

[Manuals.plus](#) /

› [Oral-B](#) /

› Oral-B iO3 Electric Toothbrush User Manual

Oral-B iOG3.1B6.0

Oral-B iO3 Electric Toothbrush User Manual

Model: iOG3.1B6.0

INTRODUCTION

Welcome to the Oral-B iO3 Electric Toothbrush user manual. This guide provides essential information for the proper setup, operation, and maintenance of your new electric toothbrush. The Oral-B iO3 combines a dentist-inspired round brush head with micro-vibrating bristles to deliver a professional purifying clean.



Image: Oral-B iO3 Electric Toothbrush and its packaging, highlighting its design and key features.

WHAT'S INCLUDED

Your Oral-B iO3 Electric Toothbrush package typically includes the following components:

- (1) Oral-B iO3 Electric Toothbrush Handle
- (1) Oral-B Ultimate Clean Brush Head
- (1) Charging Station

SETUP

Initial Charging

Before first use, fully charge your Oral-B iO3 Electric Toothbrush. This ensures optimal battery performance.

1. Plug the charging station into a suitable electrical outlet.

2. Place the toothbrush handle onto the charging station. The light ring on the handle will illuminate to indicate charging.
3. Allow the toothbrush to charge completely. A full charge typically takes several hours.

Attaching the Brush Head

To attach a brush head:

1. Align the brush head with the metal shaft on the toothbrush handle.
2. Push the brush head down firmly until it clicks into place.

To remove the brush head, simply pull it straight off the handle.

OPERATING YOUR TOOTHBRUSH

Brushing Procedure

1. Wet the brush head and apply a small amount of toothpaste.
2. Press the power button to turn on the toothbrush.
3. Place the brush head against your teeth at a slight angle towards the gumline.
4. Move the brush head slowly from tooth to tooth, allowing the micro-vibrating bristles to clean each surface.
5. Brush for the dentist-recommended 2 minutes. The built-in timer will signal every 30 seconds to prompt you to move to a different quadrant of your mouth.
6. After brushing, rinse your mouth and the brush head thoroughly.

Smart Pressure Sensor

The Oral-B iO3 features a Smart Pressure Sensor to protect your gums and ensure optimal cleaning. The light ring on the handle changes color based on the pressure applied:

- **Green Light:** Indicates optimal brushing pressure.
- **Red Light:** Signals that you are brushing too hard. Reduce pressure immediately.
- **White Light:** Indicates too little pressure or that the toothbrush is not actively brushing.

Protect your gums

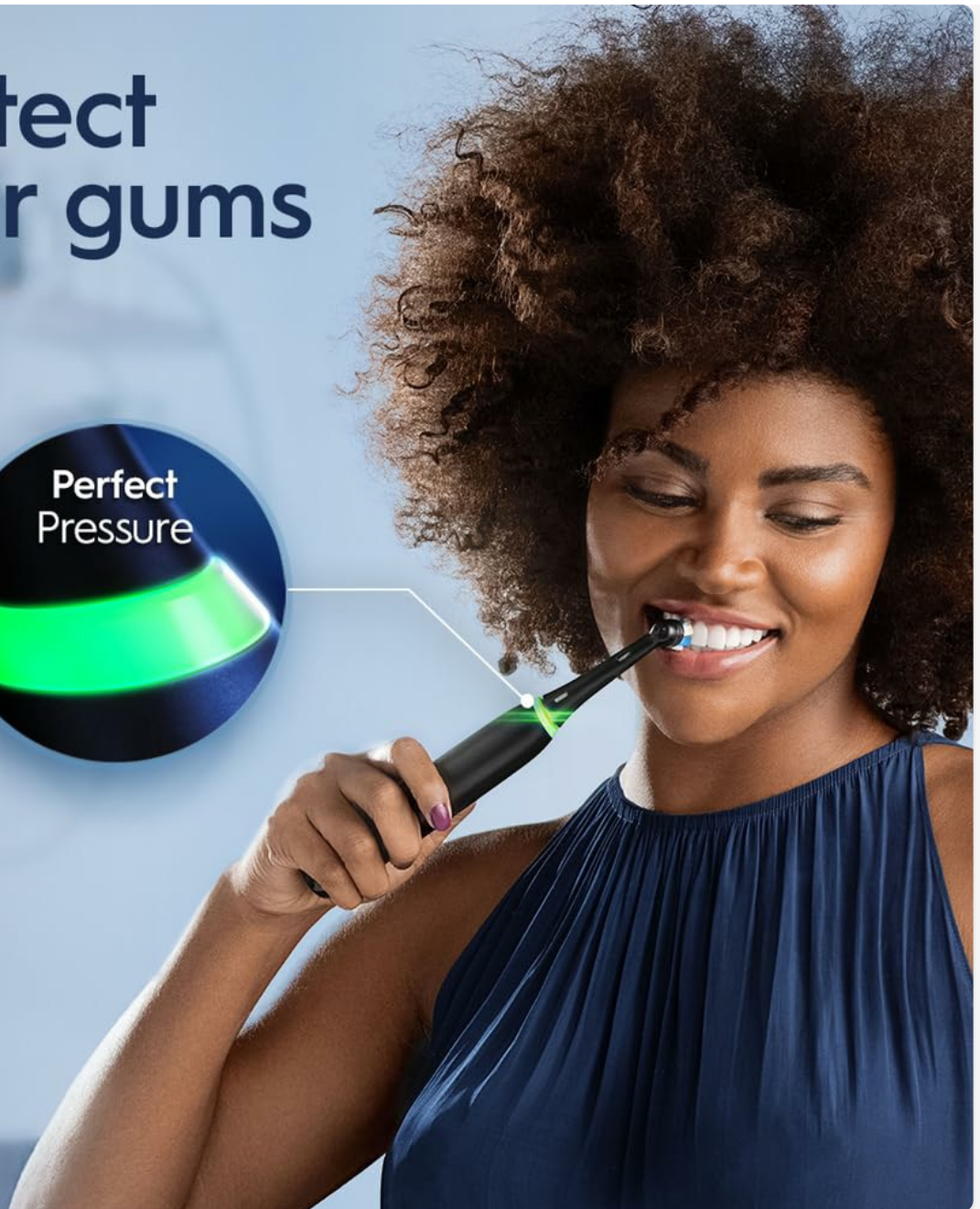


Image: Illustration of the Smart Pressure Sensor, showing a green light for optimal pressure.

Cleaning Modes

The Oral-B iO3 offers 3 cleaning modes to personalize your brushing experience:

- **Daily Clean:** For everyday cleaning.
- **Sensitive:** Gentle cleaning for sensitive areas.
- **Whitening:** Designed to polish and whiten teeth.

To switch between modes, press the mode button located below the power button.

Smart features to maximize your clean



Image: Diagram illustrating the smart features of the Oral-B iO3, including cleaning modes, 2-minute timer, and brush head change alert.

The Oral-B iO3 is clinically proven to remove 400% more plaque along the gumline compared to a regular manual toothbrush, contributing to healthier gums in just one week.

Remove 100% more plaque in 1 week*

with Oral-B's most advanced technology



Image: Visual comparison demonstrating superior plaque removal with the Oral-B iO3 electric toothbrush versus a manual toothbrush.

MAINTENANCE

Cleaning Your Toothbrush

- After each use, rinse the brush head thoroughly under running water with the handle switched on.
- Remove the brush head and rinse both the brush head and the handle separately.
- Wipe the handle with a damp cloth. Do not immerse the handle in water.
- Clean the charging station periodically with a damp cloth. Ensure it is unplugged before cleaning.

Brush Head Replacement

For optimal performance and hygiene, replace your brush head every 3 months, or sooner if the bristles become worn or frayed. The toothbrush may provide an alert when it's time to replace the brush head.

TROUBLESHOOTING

Problem	Possible Cause / Solution
Toothbrush not turning on or charging.	Ensure the charging station is properly plugged in and the toothbrush handle is correctly seated. Check the power outlet. If the issue persists, the battery may be depleted or faulty; contact customer support.
Red light constantly illuminates during brushing.	This indicates you are applying too much pressure. Reduce the pressure on your teeth and gums. The sensor is designed to protect your oral health.
Battery life seems shorter than expected.	Ensure the toothbrush is fully charged before use. Battery performance can vary with usage patterns and over time. If significantly shorter, ensure proper charging or contact customer support.
Brush head feels loose or falls off.	Ensure the brush head is pushed down firmly until it clicks into place. If it continues to be loose, the brush head or handle may be damaged and require replacement.

SPECIFICATIONS

Feature	Detail
Model Number	iOG3.1B6.0
Power Source	Battery Powered (1 Lithium Ion battery included)
Special Feature	Pressure Sensor
Item Firmness Description	Medium
Manufacturer	Procter & Gamble
UPC	069055138045
Package Dimensions	10.12 x 7.09 x 4.02 inches
Item Weight	15.52 ounces

WARRANTY AND SUPPORT

Specific warranty details for the Oral-B iO3 Electric Toothbrush are not provided in this manual. For comprehensive warranty information, product registration, and customer support, please visit the official Oral-B website or contact their customer service directly.

You can often find support resources, FAQs, and contact information on the manufacturer's official website: www.oralb.com

