

ANLAN SD-8811

ANLAN RF EMS Facial and Neck Device SD-8811 Instruction Manual

Model: SD-8811

INTRODUCTION

This manual provides detailed instructions for the safe and effective use of your ANLAN RF EMS Facial and Neck Device, model SD-8811. This device is designed to support skin care routines through Radio Frequency (RF), Electrical Muscle Stimulation (EMS), and multi-frequency vibration technologies, aiming to improve skin elasticity and reduce the appearance of wrinkles.

IMPORTANT SAFETY INFORMATION

- **Metallic Implants:** If you have metallic implants in your mouth (e.g., metal braces), it is not recommended to use the EMS mode of this product. The EMS microcurrent may form a circuit with the metallic object, potentially causing discomfort to the gums and affecting the device's effectiveness.
- **Invisible Aligners:** If you wear invisible aligners that do not contain metal, the microcurrent beauty device can generally be used normally.
- **RF Mode:** The RF mode can be used normally regardless of metallic implants.
- **Usage Frequency:** It is recommended to use the RF mode 2 to 3 times per week. The EMS mode can be used daily. Each session should last approximately 10 minutes.
- Consult a healthcare professional if you have any pre-existing medical conditions or concerns before using this device.

PRODUCT OVERVIEW

The ANLAN RF EMS Facial and Neck Device is a handheld beauty instrument designed for comprehensive skin care. It integrates Radio Frequency (RF), Electrical Muscle Stimulation (EMS), multi-frequency vibration, and light therapy (Red and White Light) to target various skin concerns.



Image: The ANLAN RF EMS Facial and Neck Device shown with its retail packaging.

Device Components

Effectively solve various skin problems



Image: Detailed diagram illustrating the device's components, including the massage area, light therapy area, level/vibration button, power/mode button, level indicators, and charging port.

- **Massage Area / Light Therapy:** The head of the device, which comes into contact with the skin, delivering RF, EMS, vibration, and light therapy.
- **Level/Vibration Button:** Used to adjust intensity levels and activate vibration.
- **Power/Mode Button:** Used to turn the device on/off and switch between operating modes.
- **Level Indicators:** Lights indicating the selected intensity level.
- **Charging Port:** USB-C port for charging the device.

KEY FEATURES AND FUNCTIONS



Image: Overview of the device's six primary functions, including RF, EMS Micro Current, White Light, Vibration, Red Light, and 3 adjustable intensity levels.

1. Radio Frequency (RF)

RF energy penetrates the dermal layer of the skin, stimulating collagen regeneration and promoting cellular rejuvenation. This process helps to reduce the appearance of wrinkles and tighten facial contours.

RF: Acts on collagen in the dermis, promotes collagen synthesis by heating, and tightens the skin

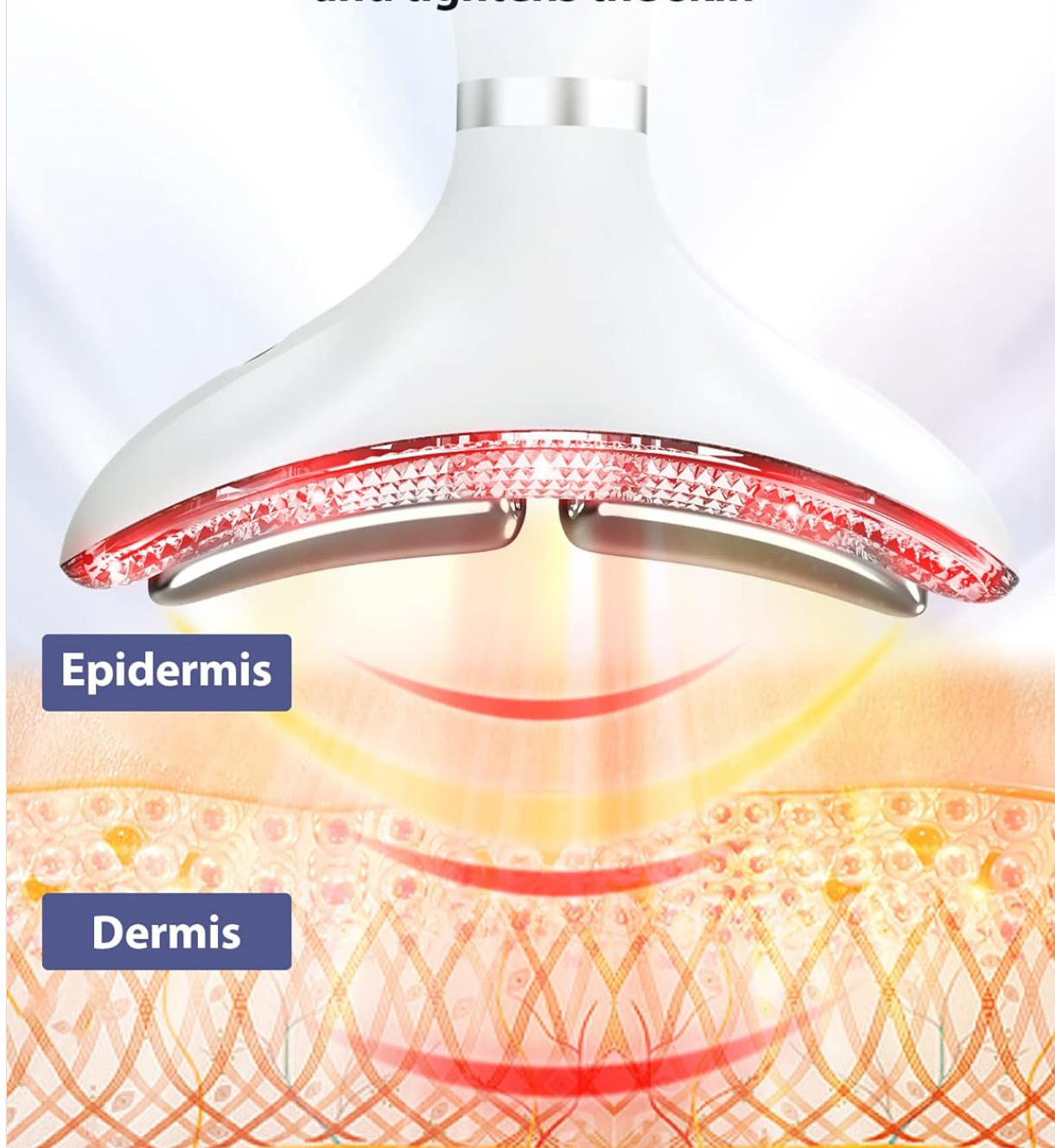


Image: Visual representation of RF energy acting on collagen in the dermis, promoting collagen synthesis and skin tightening.

2. Electrical Muscle Stimulation (EMS) Microcurrent

EMS microcurrent stimulates and massages skin cells, increasing cellular activity. This promotes better absorption of skincare products and cellular proliferation, contributing to plumper, firmer skin and restored elasticity.

EMS: stimulates and massages skin cells, increases cell activity, and promotes absorption

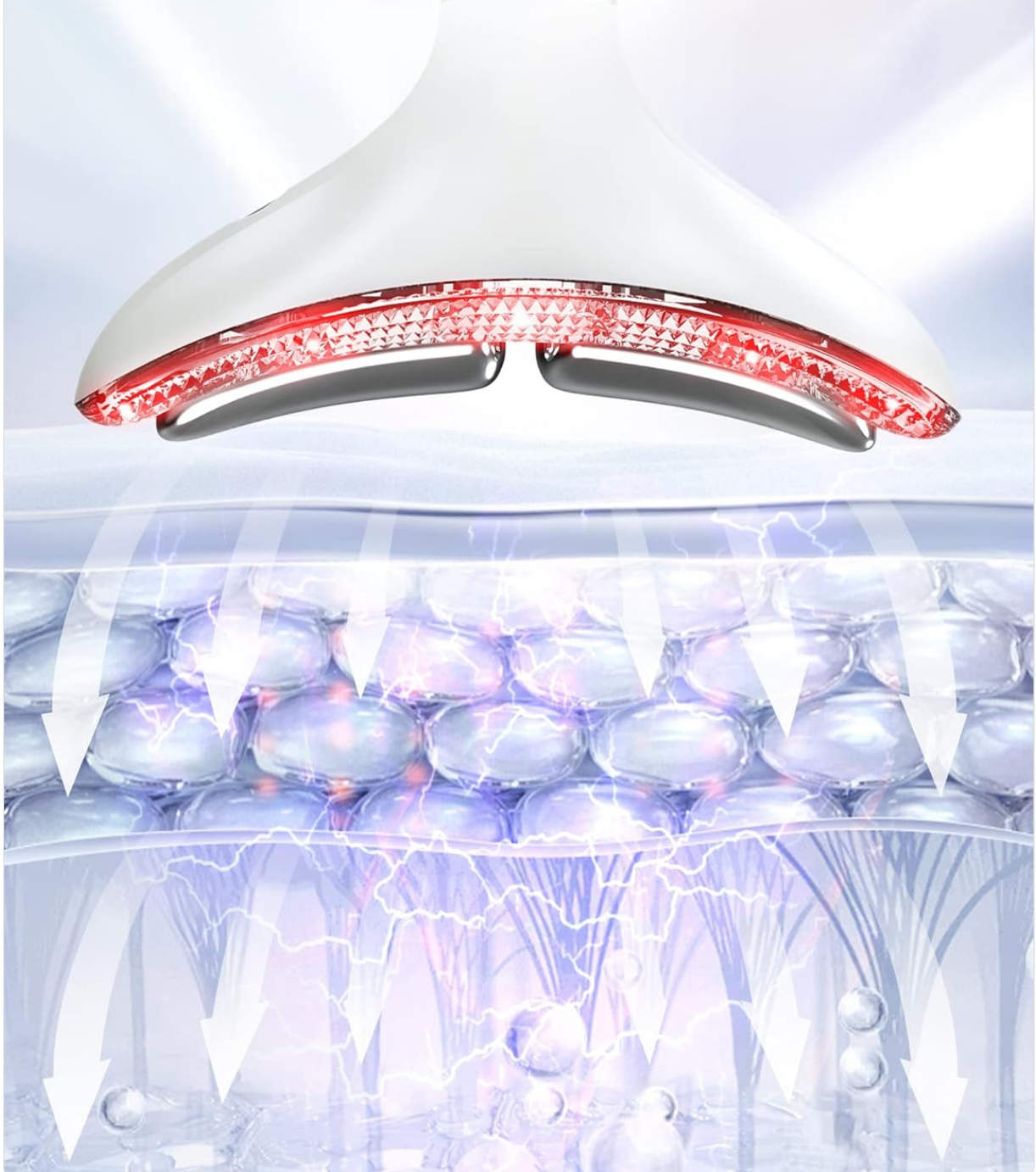


Image: Diagram showing EMS stimulating and massaging skin cells, increasing cell activity and promoting absorption.

3. Multi-Frequency Vibration

The device features adjustable vibration speeds suitable for various skin types and ages. Vibration helps to relax facial muscles and enhance the absorption of skincare products.

High frequency vibration Relax muscles and promote absorption

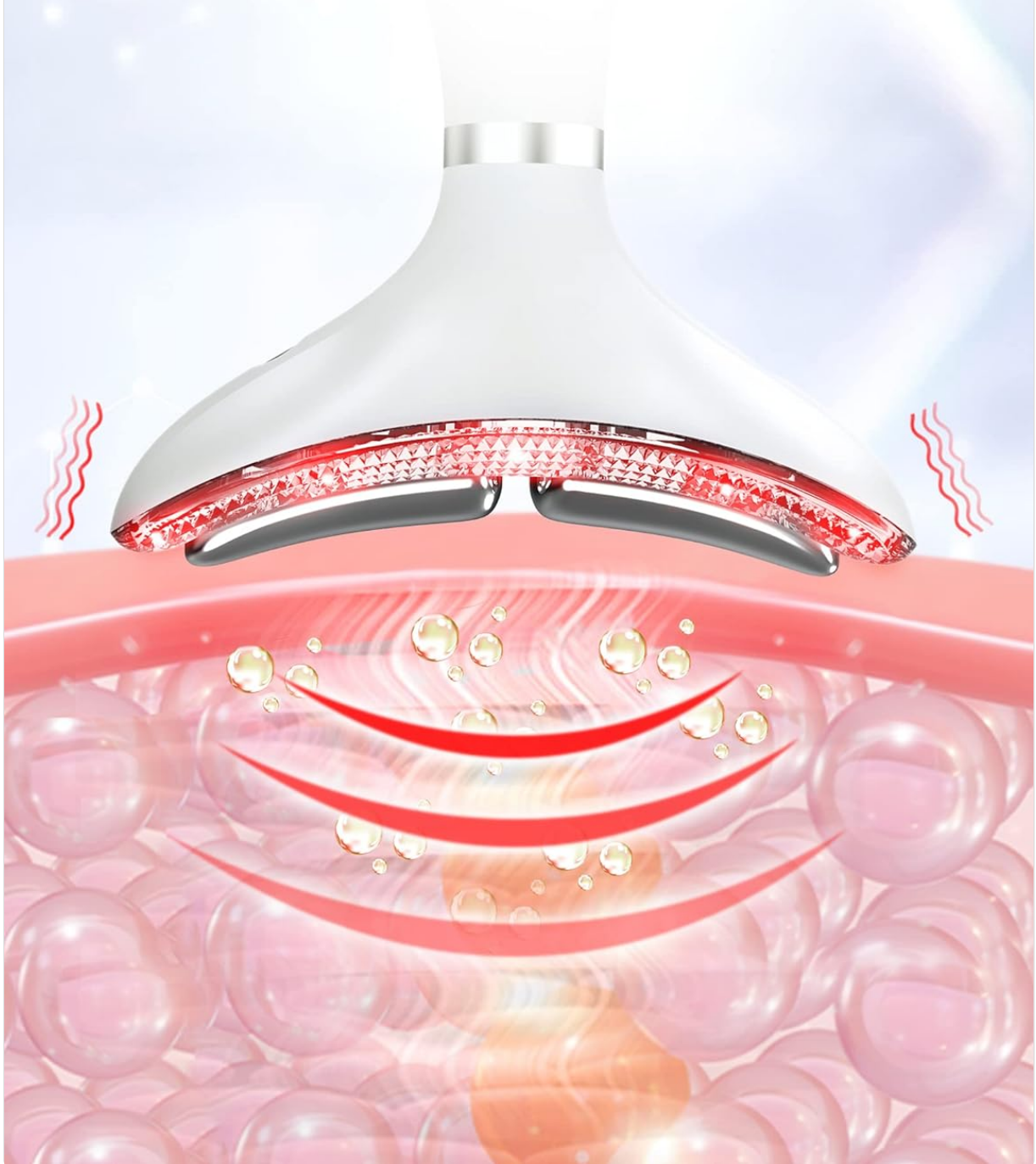


Image: Visual representation of high-frequency vibration relaxing muscles and promoting absorption within the skin.

4. Red Light Therapy (640nm)

Red light at 640nm wavelength helps to rejuvenate collagen and improve skin elasticity.

5. White Light Therapy (490nm)

White light at 490nm wavelength assists in lightening pigmentation and breaking down dark spots, contributing to a more even skin tone.



Image: Illustrates the effects of Red Light (640nm) for collagen rejuvenation and White Light (490nm) for lightening pigment and breaking down dark spots.

SETUP

1. Charging the Device

Before first use, fully charge the device. Connect the provided power adapter to the charging port on the device and plug it into a power outlet. The level indicators will show charging status and turn solid when fully charged.

2. Skin Preparation

Cleanse your face and neck thoroughly. Apply a suitable skincare product such as a serum, gel, or moisturizer to the treatment area. This ensures smooth gliding of the device and optimal conductivity for RF and EMS functions.

OPERATING INSTRUCTIONS

1. Power On/Off

- To power on, press and hold the **Power/Mode Button** for a few seconds until the device activates.
- To power off, press and hold the **Power/Mode Button** again for a few seconds.

2. Mode Selection

Briefly press the **Power/Mode Button** to cycle through the available modes:



Image: Illustrates the two main operating modes: Mode 1 (RF + Red Light) and Mode 2 (EMS + White Light), along with instructions for button usage.

- **Mode 1: RF + Red Light (621 nm)** - Focuses on deep skin penetration, accelerating active tissue metabolism, and improving fine lines and skin laxity.

- **Mode 2: EMS + White Light (473 nm)** - Targets lightening dark spots, improving dullness, and promoting skin firmness.

3. Level and Vibration Adjustment

- Briefly press the **Level/Vibration Button** to change the intensity level (3 levels available). Start with the lowest level and gradually increase as your skin adapts.
- Press and hold the **Level/Vibration Button** to activate or deactivate the vibration function.

4. Application

- Gently glide the device over the desired areas of your face and neck. Ensure continuous movement and avoid holding the device stationary on one spot for too long.
- For optimal results, follow the recommended usage frequency: RF mode 2-3 times per week, EMS mode daily. Each session should be approximately 10 minutes.

MAINTENANCE

1. Cleaning

- After each use, turn off the device and wipe the massage head with a soft, damp cloth.
- Do not use abrasive cleaners, alcohol, or harsh chemicals, as these can damage the device.
- Ensure the device is completely dry before storing.

2. Storage

Store the device in a cool, dry place, away from direct sunlight and out of reach of children.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery or not charged.	Charge the device fully before use. Ensure the charging cable is properly connected.
No sensation during EMS/RF mode.	Insufficient contact with skin or low intensity level.	Ensure the device head is in full contact with moisturized skin. Increase the intensity level using the Level/Vibration Button.
Device stops working during use.	Battery depleted or automatic shutdown.	Recharge the device. The device may have an automatic shutdown feature after a certain period of inactivity or continuous use.
Skin feels uncomfortable or irritated.	Intensity level too high or insufficient skincare product.	Reduce the intensity level. Apply more serum or gel to ensure smooth gliding and proper conductivity. Discontinue use if irritation persists.

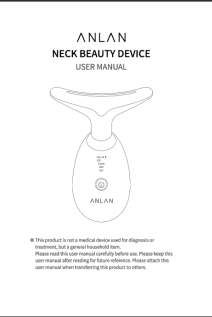

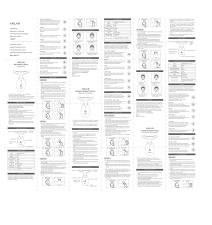
SPECIFICATIONS

Brand	ANLAN
Model Number	SD-8811
Dimensions	18 x 6 x 3.1 cm
Weight	190 grams
Material	Plastic
Color	White
Power Source	Mains operated (via adapter)
Voltage	5 Volts
Included Components	Power adapter, possibly facial tips

WARRANTY AND SUPPORT

For warranty information and customer support, please refer to the documentation included with your purchase or contact ANLAN customer service directly. Keep your proof of purchase for warranty claims.

Related Documents - SD-8811

 <p>The image shows the cover of the ANLAN Neck Beauty Device User Manual. It features a line drawing of the device, which is a handheld, oval-shaped unit with a strap. The text on the cover includes 'ANLAN NECK BEAUTY DEVICE USER MANUAL' and a disclaimer: 'This product is not a medical device used for diagnosis or treatment, but an auxiliary product for general household use. Please read the user manual carefully before use. Please keep this user manual after reading for future reference. Please attach this user manual when returning the product to us.'</p>	<p>ANLAN Neck Beauty Device User Manual: Features, Functions, and How to Use</p> <p>Comprehensive user manual for the ANLAN Neck Beauty Device (Model ES-081). Learn about its features, functions (EMS, warm, red/blue/green light), safety precautions, and how to use it for skin rejuvenation and care.</p>
 <p>The image shows the cover of the ANLAN Ultrasonic Skin Scrubber User Manual. It features a line drawing of the device, which is a handheld, rectangular unit with a strap. The text on the cover includes 'ANLAN Ultrasonic Skin Scrubber User Manual' and a disclaimer: 'This product is not a medical device for diagnosis and treatment, but an auxiliary product for general household use. Please read the user manual carefully before use. Please keep this user manual after reading for future reference. Please attach this user manual when returning the product to us.'</p>	<p>ANLAN Ultrasonic Skin Scrubber User Manual</p> <p>User manual for the ANLAN Ultrasonic Skin Scrubber (Model: JD-CP002), detailing its features, safety precautions, usage instructions, maintenance, and FAQs.</p>
 <p>The image shows the cover of the ANLAN Neck Beauty Device User Manual. It features a line drawing of the device, which is a handheld, oval-shaped unit with a strap. The text on the cover includes 'ANLAN NECK BEAUTY DEVICE USER MANUAL' and a disclaimer: 'This product is not a medical device used for diagnosis or treatment, but an auxiliary product for general household use. Please read the user manual carefully before use. Please keep this user manual after reading for future reference. Please attach this user manual when returning the product to us.'</p>	<p>ANLAN Neck Beauty Device User Manual</p> <p>User manual for the ANLAN Neck Beauty Device, detailing its features, functions, usage instructions, and warranty information. The device utilizes positive and negative ions, vibration, warm massage, and red, blue, and green light therapy for skin cleansing, EMS care, and moisturizing.</p>



Comprehensive user manual for the ANLAN Body Beauty Device (Model OFY-2102), covering safety precautions, product details, usage instructions, troubleshooting, and warranty information.