



Manuals.plus /

› Pousutong /

› Liberex FF-20S Cordless Water Flosser Instruction Manual

## Pousutong FF-20S

# Liberex FF-20S Cordless Water Flosser Instruction Manual

MODEL: FF-20S

Overview   Setup   Introduction   Operation   Safety Information   Maintenance   Troubleshooting   Product   Specifications   Support

## INTRODUCTION

Thank you for choosing the Liberex FF-20S Cordless Water Flosser. This device is designed to provide an effective and convenient way to improve your oral hygiene by removing plaque and food debris from between teeth and below the gumline. Please read this manual thoroughly before first use to ensure proper operation and maintenance.

## IMPORTANT SAFETY INFORMATION

Please observe the following safety precautions to prevent injury or damage to the device:

- Do not use the device while it is charging.
- Keep the device out of reach of children.
- Do not submerge the main body of the device in water beyond the water tank. The device is IPX7 waterproof, but the charging port should remain dry during use.
- This product is not intended for use by individuals with certain medical conditions without consulting a healthcare professional.
- Do not use the device if it appears damaged in any way.
- Always follow the instructions provided in this manual for safe and effective use.


## PRODUCT OVERVIEW

The Liberex FF-20S Cordless Water Flosser features advanced technology for superior oral cleaning.

- **Customizable Flossing Experience:** Offers 4 preset flossing modes (Soft, Medium, Strong, Pulse) and a custom mode allowing pressure adjustment from 8-110 PSI across 10 levels.
- **OLED Display:** An easy-to-read screen shows current running modes and pressure levels.
- **Effective Cleaning:** Patented 45° radian angled nozzles concentrate the water stream to 0.6mm, removing up to

99% of plaque. 360° rotatable nozzles ensure comprehensive cleaning.

- **Waterproof Design:** IPX7 waterproof rating allows for safe use in the shower.
- **Convenient Water Reservoir:** Detachable 300ml tank for easy refilling and cleaning, providing enough capacity for a full flossing session.
- **USB-C Rechargeable:** Modern USB-C charging port; a single charge lasts approximately one month.
- **Portable:** Includes a travel bag for easy transport.

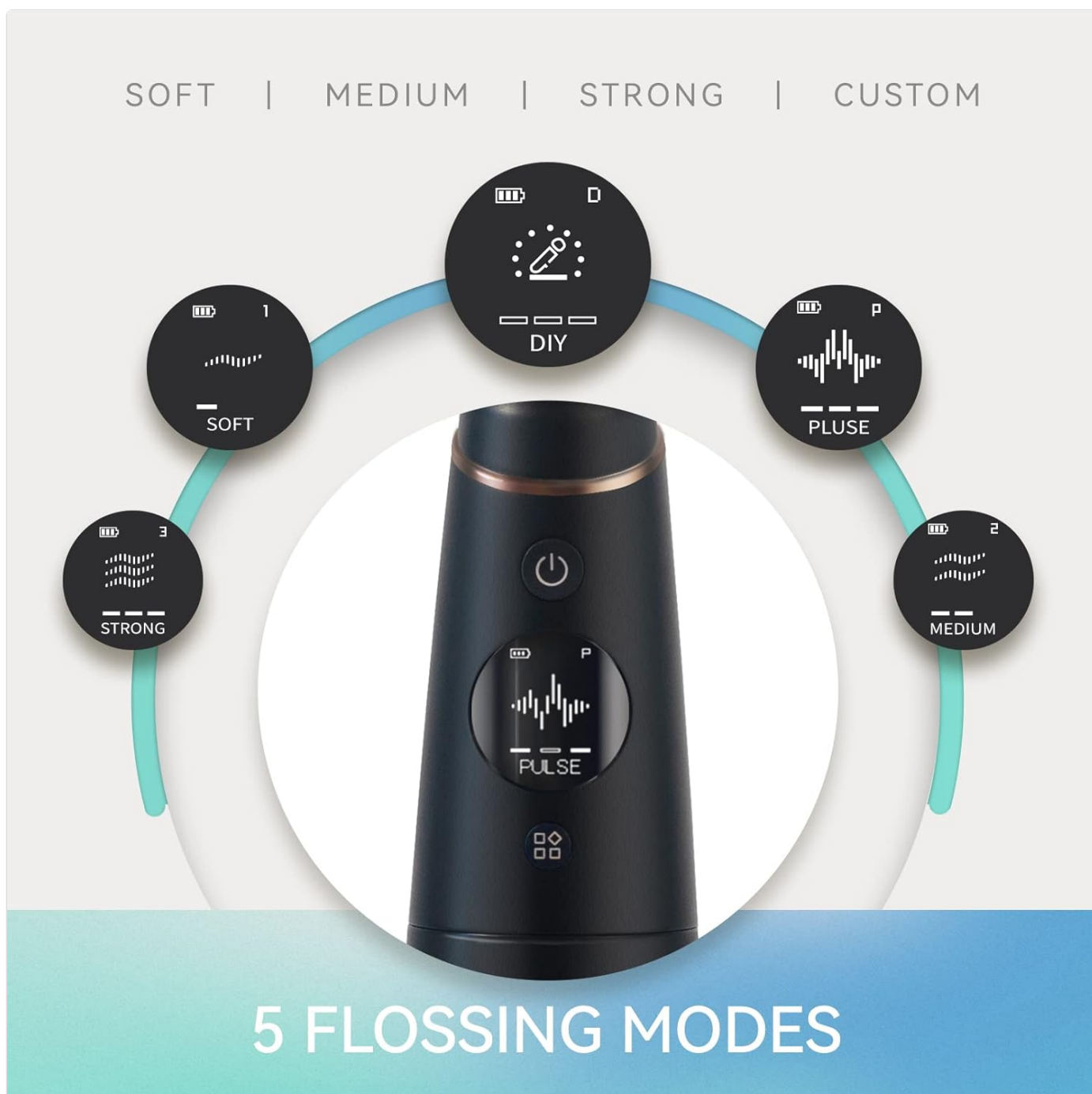

*Image: The OLED display on the Liberex FF-20S Water Flosser illustrates the five available flossing modes: Soft, Medium, Strong, Pulse, and a customizable DIY mode, each with an icon and pressure indicator.*



EVERY DAY FRESH BREATH

## IPX7 WATERPROOF

*Image: The Liberex FF-20S Water Flosser is shown being used under running water, emphasizing its IPX7 waterproof design for easy cleaning and shower use.*

# Detachable water reservoir

---



*Image: A visual representation of the Liberex FF-20S Water Flosser's detachable 300ml water reservoir, demonstrating how it can be easily removed for refilling and cleaning.*

RY DAY FRESH B



## USB-C RECHARGEABLE

*Image: A detailed view of the Liberex FF-20S Water Flosser's USB-C charging port, illustrating the connection of a charging cable.*



## FF-20S DIY OLED Cordless Water Flosser

*Image: The complete package contents of the Liberex FF-20S Water Flosser, including the device, various dental tips, charging cable, and a travel bag, are displayed alongside its retail box.*

### SETUP GUIDE

---

- 1. Charge the Device:** Before first use, fully charge the water flosser using the provided USB-C cable. Connect the cable to the USB-C port on the device and a suitable power adapter. The OLED display will indicate charging status. A full charge typically lasts for one month of regular use.
- 2. Fill the Water Tank:** Detach the 300ml water tank from the main body of the flosser. Fill it with lukewarm water. You can also fill it directly through the water port without detaching, if preferred. Ensure the tank is securely reattached to prevent leaks.
- 3. Attach a Dental Tip:** Select one of the 5 dental tips provided. Insert the tip into the opening at the top of the flosser until it clicks into place. Ensure it is firmly seated.

### OPERATING INSTRUCTIONS

---

1. **Select Mode and Pressure:** Press the power button to turn on the device. Use the mode button to cycle through Soft, Medium, Strong, and Pulse modes. For the Custom mode, press the DIY button and then use the pressure adjustment buttons to set your desired water pressure from 8-110 PSI. The OLED display will show your selection.
2. **Position the Tip:** Lean over a sink. Place the dental tip into your mouth, pointing it towards your gum line at a 90-degree angle. Keep your mouth slightly open to allow water to flow into the sink.
3. **Begin Flossing:** Turn on the device. Guide the water stream along your gum line, pausing briefly between teeth. Ensure you clean both the front and back surfaces of your teeth. The 360° rotatable nozzle allows you to reach all areas easily. Avoid directing the stream directly into gum pockets or eyes.
4. **Complete Session:** Continue until you have cleaned all areas of your mouth or the water tank is empty. The device will automatically shut off after a set period.



Image: A user demonstrates the correct way to hold and position the LibereX FF-20S Water Flosser for effective oral cleaning.

## MAINTENANCE AND CARE

---

- Clean After Each Use:** After each flossing session, remove the dental tip and rinse it thoroughly under warm water. Empty any remaining water from the tank and rinse the tank as well. Allow both to air dry.
- Deep Cleaning (Weekly):** For a more thorough cleaning, you can detach the water tank and wash it with mild soap and water. Ensure all soap residue is rinsed away. Wipe the main body of the flosser with a damp cloth.
- Store Properly:** Store the device in a clean, dry place. The included travel bag is ideal for protecting the flosser when not in use or during travel.
- Replace Dental Tips:** For optimal hygiene and performance, it is recommended to replace the dental tips every 2-3 weeks, or sooner if they show signs of wear.

## TROUBLESHOOTING

If you encounter any issues with your Liberex FF-20S Water Flosser, please refer to the table below:

Problem	Possible Cause	Solution
Device does not turn on.	Low battery or no charge.	Charge the device fully using the USB-C cable.
Weak water pressure.	Low water level in tank; clogged tip; incorrect mode/pressure setting.	Refill water tank; clean or replace tip; adjust mode/pressure setting on OLED display.
Water leaks from the device.	Water tank not properly attached; tip not securely inserted.	Ensure water tank is firmly seated; ensure dental tip clicks into place.
Battery not holding charge.	Battery degradation or faulty charging cable/port.	Try a different USB-C cable or power adapter. If issue persists, contact customer support.

## SPECIFICATIONS

Feature	Detail
Model Number	FF-20S
Water Tank Capacity	300ml
Pressure Range	8-110 PSI (10 levels customizable)
Flossing Modes	Soft, Medium, Strong, Pulse, Custom (DIY)
Charging Interface	USB-C
Waterproof Rating	IPX7
Item Weight	1.01 pounds
Package Dimensions	9.65 x 4.8 x 3.19 inches
Power Source	Battery Powered
Special Features	Multiple Operation Modes, Multiple Pressure Settings, Multiple Tips, Non-Slip Grip, Portable, Timer
Product Benefits	Gum Health, Prevents Bad Breath, Prevents Cavities, Removes Plaque

## WARRANTY AND CUSTOMER SUPPORT

---

For warranty information, technical assistance, or any questions regarding your Liberex FF-20S Cordless Water Flosser, please refer to the warranty card included with your product or contact Pousutong customer support through their official website or the retailer where the product was purchased. Please have your model number (FF-20S) and purchase date ready when contacting support.