

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [SmileDirectClub](#) /
- › [SmileDirectClub Pro Teeth Whitening System with Wireless 20-LED Light - Instruction Manual](#)

## SmileDirectClub Pro Teeth Whitening System with Wireless 20-LED Light

# SmileDirectClub Pro Teeth Whitening System with Wireless 20-LED Light - Instruction Manual

[Overview](#)   [What's Included](#)   [Setup](#)   [Operating](#)  
[Instructions](#)   [Maintenance](#)   [Troubleshooting](#)   [Specifications](#)   [Ingredients](#)   [Safety](#)  
[Information](#)   [Warranty & Support](#)

## 1. OVERVIEW

The SmileDirectClub Pro Teeth Whitening System is designed to provide professional-level teeth whitening results. This system utilizes a 3-step process involving a whitening accelerator toothpaste, an advanced whitening gel, and a wireless 20-LED light, followed by a whitening extender gel to maintain results. The system is formulated to be enamel-safe and cause low sensitivity.

### PRO WHITENING SYSTEM

**Up to 12 shades  
brighter in 1 week.\***

Get professional results in 3 easy steps to boost, whiten and extend results.

\*Maximum results when used as directed. Results may vary.



Image: The complete SmileDirectClub Pro Whitening System, showcasing the accelerator toothpaste, advanced whitening gel pens, ProLong extender gel pens, and the wireless 20-LED light.

## 2. WHAT'S INCLUDED

Your SmileDirectClub Pro Teeth Whitening System package includes the following components:

- **1 ProActivate Whitening Accelerator Toothpaste:** Preps teeth to increase brightness by up to 65%. Designed for vigorous, deep cleaning to break up tough stains and remove plaque. It is enamel-safe and

easy to use as regular toothpaste.

- **1 Advanced Whitening Gel (6mL):** Brightens teeth by up to 12 shades in 1 week. Provides up to 20 treatments, with results visible after just 1 use. Contains the same enamel-safe, low-sensitivity whitening agent used by dentists.
- **2 ProLong Whitening Extender Gels (4mL each):** Apply post-whitening or use daily to protect and prolong results. Helps prevent new stains from forming and smooths the tooth's surface.
- **Rechargeable Wireless Whitening Light:** A patented 20-LED light designed to boost whitening results. A single charge lasts a whole week, and the light recharges on its base. Features an automatic 10-minute treatment timer.



## Good things come in 3's.

- Step 1** Whitening Activator Paste preps your teeth and boosts your bright by 65%.<sup>1</sup>
- Step 2** Advanced Whitening Gel whitens your smile up to 12 shades in 1 week.<sup>2</sup>
- Step 3** ProLong™ Whitening Extender protects and prolongs your results, preventing new stains from forming.

<sup>1</sup>Maximum results when used as directed vs SmileDirectClub Premium Teeth Whitening gel alone. Results may vary. <sup>2</sup>Maximum results when used as directed. Results may vary.

Image: Close-up view of the ProActivate Whitening Accelerator Paste, Advanced Whitening Gel pens, and ProLong Whitening Extender Gel pens, which are key components of the whitening system.

## 3. SETUP

### 3.1 Charging the Wireless Whitening Light

1. Place the wireless 20-LED light onto its charging base.
2. Ensure the charging base is connected to a power source.
3. Allow the light to charge fully before first use. A full charge typically lasts for one week of treatments.

## 4. OPERATING INSTRUCTIONS

Follow these three steps daily for one week to achieve optimal whitening results. Each step is designed to work synergistically for a brighter smile.

### 4.1 Step 1: ProActivate Whitening Accelerator Toothpaste

1. Brush your teeth thoroughly with the ProActivate Whitening Accelerator Toothpaste for two minutes.
2. This toothpaste prepares your teeth for the whitening process by removing surface stains and plaque.
3. Rinse your mouth after brushing.

PRO WHITENING SYSTEM

# 3 steps to professional-level results.

## Our most advanced whitening treatment.

Brighten your smile up to 12 shades in 1 week<sup>1</sup> in just 3 easy steps for 10 minutes a day.

<sup>1</sup>Maximum results when used as directed. Results may vary.



Image: A person brushing their teeth with the ProActive Whitening Accelerator Toothpaste, demonstrating the first step of the whitening process.

### 4.2 Step 2: Advanced Whitening Gel with Wireless LED Light

1. After brushing, apply a thin layer of the Advanced Whitening Gel directly to your teeth using the pen applicator. Ensure even coverage on all visible tooth surfaces.
2. Insert the fully charged wireless 20-LED light into your mouth, ensuring it covers all teeth treated with the gel.
3. Turn on the LED light. The light has an automatic 10-minute treatment timer.
4. Keep the light in your mouth for the full 10-minute duration.
5. After 10 minutes, remove the light and rinse your mouth.

**Step 1**

# ProActivate™ Whitening Accelerator

The ProActivate Accelerator Paste helps remove tough stains and increases brightness by 65%.<sup>1</sup>

<sup>1</sup>Maximum results when used as directed. Results may vary.



Image: A person applying the Advanced Whitening Gel to their teeth using the pen applicator, preparing for the LED light treatment.

## Step 3

# ProLong™ Whitening Extender

The ProLong Extender Gel prevents new stains and extends whitening results for up to 6 months.

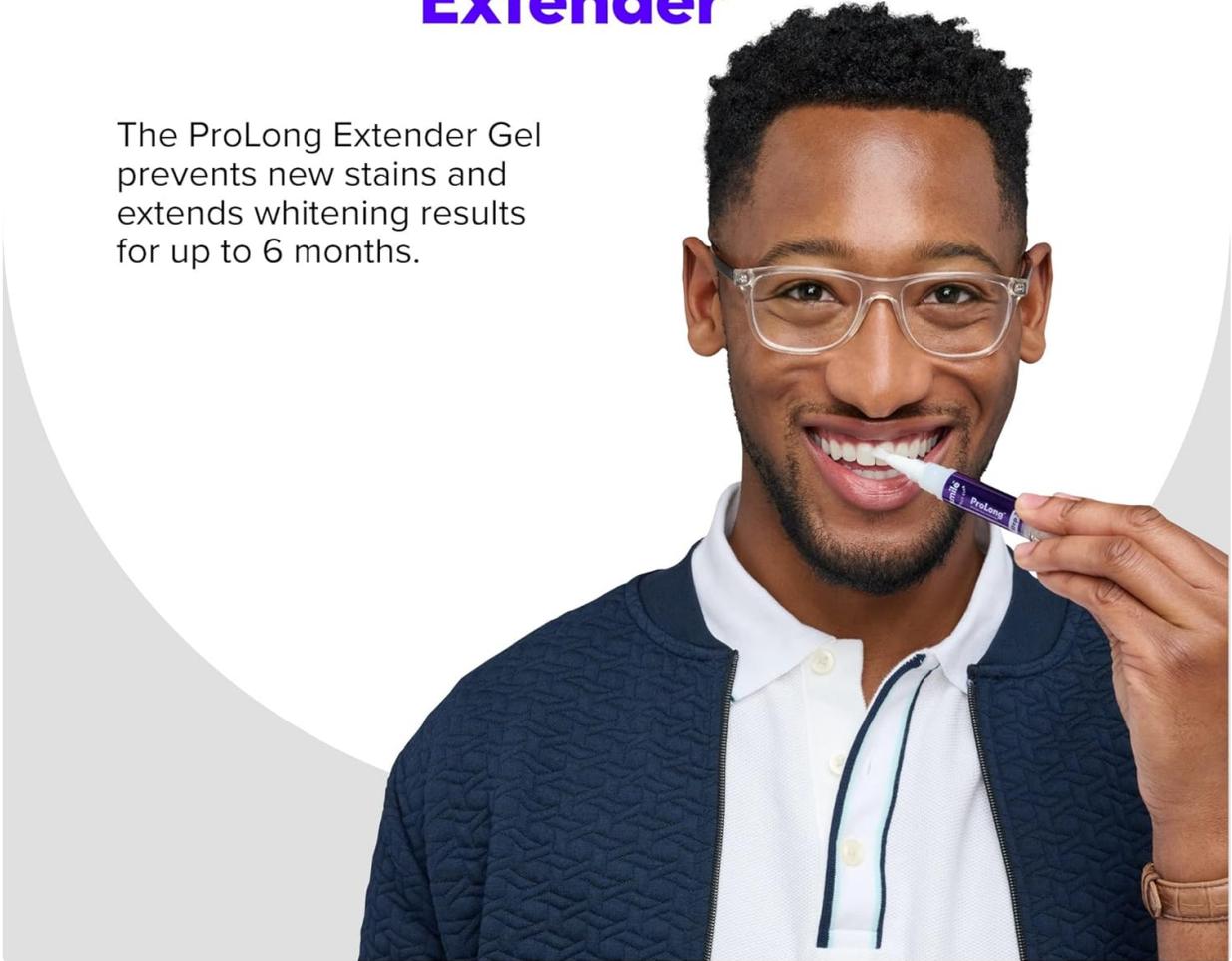


Image: A hand holding the SmileDirectClub wireless 20-LED light, ready for insertion into the mouth to activate the whitening gel.

### 4.3 Step 3: ProLong Whitening Extender Gel

1. After completing Step 2, apply the ProLong Whitening Extender Gel to your teeth.
2. This gel helps to protect your newly whitened teeth and prolong the results, preventing new stains from forming.
3. The extender gel can also be used daily as part of your regular oral hygiene routine to maintain brightness.

## Step 2

# Advanced Whitening Gel

The Advanced Whitening Gel quickly absorbs into teeth and brightens up to 12 shades in just 1 week.<sup>1</sup>

<sup>1</sup>Maximum results when used as directed. Results may vary.



Image: A person applying the ProLong Whitening Extender Gel to their teeth, the final step in the whitening system to maintain results.

## 5. MAINTENANCE

---

### 5.1 Cleaning the LED Light

- After each use, gently rinse the mouth tray portion of the LED light with cool water.
- Wipe the light clean with a soft, damp cloth.
- Do not submerge the entire LED light unit in water.
- Ensure the light is completely dry before placing it back on the charging base or storing it.

### 5.2 Storage

- Store all gels and the LED light in a cool, dry place away from direct sunlight.
- Keep out of reach of children.

## 6. TROUBLESHOOTING

---

### 6.1 LED Light Not Turning On or Malfunctioning

- **Check Charge:** Ensure the wireless LED light is fully charged by placing it on its charging base. The indicator lights on the base should confirm charging status.

- **Proper Placement:** Verify that the light is correctly seated on the charging base for proper contact.
- **Activation:** Press and hold the power button on the LED light for 2-3 seconds to activate it. A quick press may not be sufficient.
- If the issue persists after checking these points, contact SmileDirectClub customer support.

## 6.2 Tooth Sensitivity

- The system is designed for low sensitivity. However, if you experience increased sensitivity, consider reducing the frequency of use (e.g., every other day instead of daily) or the duration of the LED light treatment.
- Ensure you are not applying excessive amounts of whitening gel. A thin, even layer is sufficient.
- If sensitivity is severe or persists, discontinue use and consult your dentist.

## 6.3 Uneven Whitening Results

- Ensure the whitening gel is applied evenly to all visible tooth surfaces.
- Make sure the LED light is positioned correctly to cover all teeth being treated.
- Consistent daily use as directed is crucial for uniform results.

# 7. SPECIFICATIONS

---

- **Product Dimensions:** 2.5 x 7.5 x 8.3 inches
- **Item Weight:** 0.64 Pounds (10.24 ounces)
- **Manufacturer:** SmileDirectClub
- **ASIN:** B0BS9XVTQ2
- **UPC:** 810819032389
- **Item Form:** Gel
- **Age Range Description:** Adult
- **Number of Items:** 5 (including 3 gel pens, 1 toothpaste, 1 LED light system)
- **Unit Count:** 1.64 Fluid Ounces (total gel volume)
- **Flavor:** Unflavored

# 8. INGREDIENTS

---

## 8.1 ProActivate Whitening Accelerator Toothpaste

Water, Sorbitol, Hydrated Silica, Glycerin, Erythritol, Lithium Magnesium Sodium Silicate, Sodium Hexametaphosphate, Cellulose Gum, Xanthan Gum, Phenoxyethanol, Sodium Saccharin, Glycolipids, FD&C Blue 1, CI45430, Flavor.

## 8.2 Advanced Whitening Gel

Water, Polyethylene Glycol, Hydrogen Peroxide, Poloxamer 407, T-Butyl Alcohol, Etidronic Acid, Peppermint Oil, Sodium Stannate, EDTA.

## 8.3 ProLong Whitening Extender Gel

Water, Sorbitol, Glycerin, Erythritol, Lithium Magnesium Sodium Silicate, Hydroxyapatite, Potassium Chloride, Potassium Nitrate, Phenoxyethanol, Xanthan Gum, Tetrasodium Pyrophosphate, Tartaric Acid, Sodium Saccharin, Spearmint Oil.

*Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.*

## **9. SAFETY INFORMATION**

---

- For adult use only. Keep out of reach of children.
- Do not swallow the whitening gels or toothpaste.
- Avoid contact with eyes. In case of contact, rinse thoroughly with water.
- If irritation, sensitivity, or any adverse reaction occurs, discontinue use and consult a dental professional.
- Do not use if you are pregnant, nursing, or have gum disease, open sores in the mouth, or recent oral surgery.
- Consult your dentist before use if you have any existing dental conditions or concerns.

## **10. WARRANTY AND SUPPORT**

---

For information regarding product warranty, returns, or technical support, please refer to the packaging or contact SmileDirectClub customer service directly through their official website or the retailer where the product was purchased. Please have your product model number and purchase date available when contacting support.