

[Manuals.plus](#) /

› [MLAY](#) /

› MLAY IPL Hair Removal Device T- Series User Manual

MLAY T-

MLAY IPL Hair Removal Device T- Series User Manual

Model: T-

1. INTRODUCTION

The MLAY IPL Hair Removal Device utilizes Intense Pulsed Light (IPL) technology to target hair follicles, aiming to disrupt the hair growth cycle for permanent hair reduction. This device is designed for safe and effective home use, featuring an integrated ice cooling function for enhanced comfort during treatment.



Image 1.1: MLAY IPL Hair Removal Device and accessories.

Key Features:

- **IPL Technology:** Designed for permanent hair removal by targeting hair roots.
- **Ice Cooling Function:** Integrated plate maintains skin surface temperature between 32-41°F (0-5°C) to reduce discomfort, redness, and swelling.
- **Adjustable Power Levels:** 5 energy levels to suit various skin sensitivities and hair types.
- **Dual Flash Modes:** Manual mode for precise treatment on small areas, Auto mode for faster treatment on larger areas.
- **High Energy Output:** 12-18 joules of light energy for effective results.
- **Long Lifespan:** Up to 500,000 flashes.

2. IMPORTANT SAFETY INFORMATION

Read all instructions and warnings carefully before using this device. Failure to follow these instructions may result in injury or device damage.


2.1 Contraindications:

- Do not use on red, white, or gray hair.

- Do not use on black skin tones.
- Do not use over tattoos or dark spots (e.g., moles, freckles, birthmarks).
- Not water resistant. Avoid contact with water.

2.2 Skin and Hair Color Compatibility:

The effectiveness and safety of IPL treatment depend on your skin tone and hair color. Refer to the compatibility chart below to determine if this device is suitable for you.

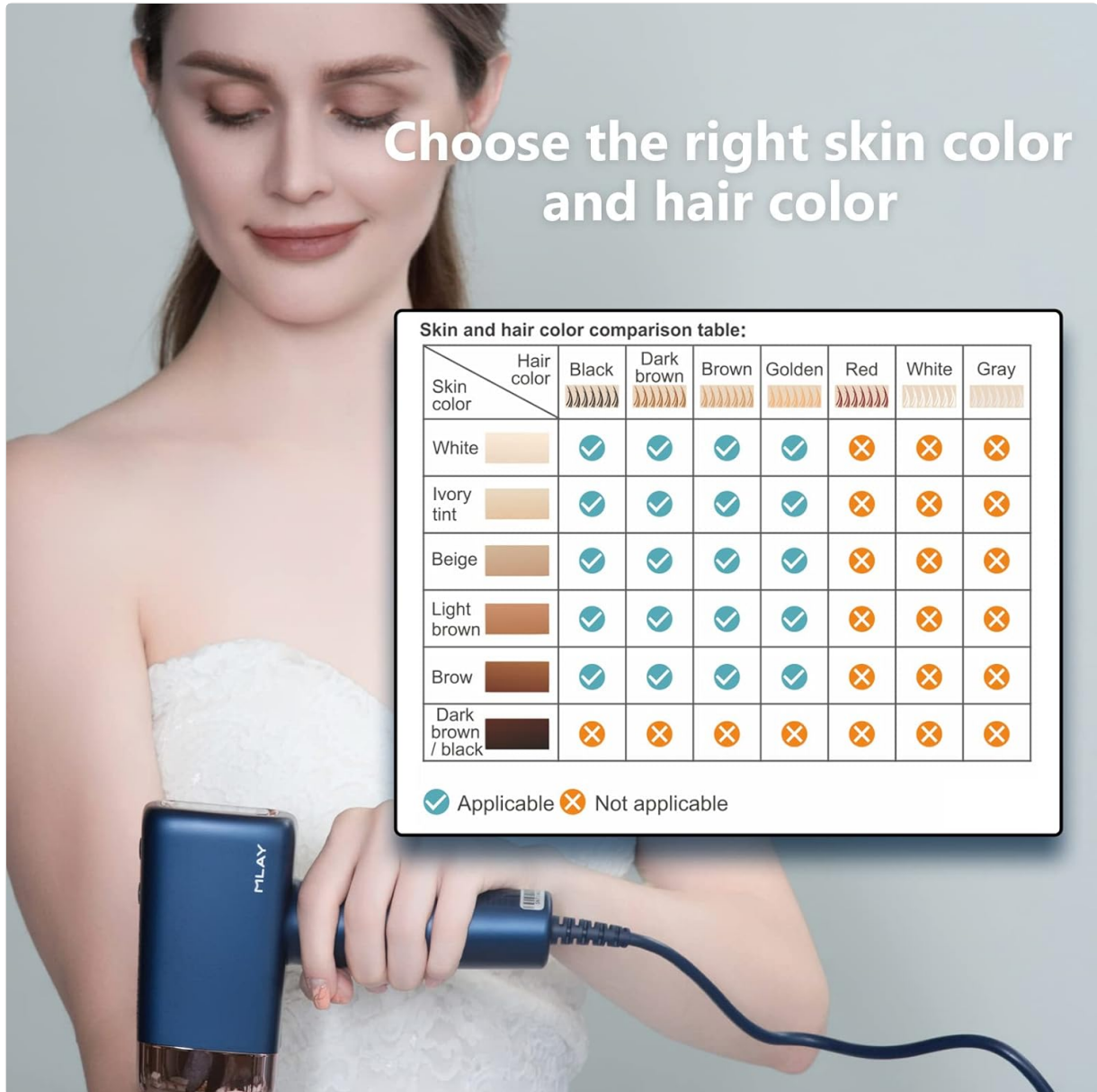

Choose the right skin color and hair color

Skin and hair color comparison table:		Black	Dark brown	Brown	Golden	Red	White	Gray
Skin color	Hair color							
White	☑	☑	☑	☑	☑	☒	☒	☒
Ivory tint	☑	☑	☑	☑	☑	☒	☒	☒
Beige	☑	☑	☑	☑	☑	☒	☒	☒
Light brown	☑	☑	☑	☑	☑	☒	☒	☒
Brown	☑	☑	☑	☑	☑	☒	☒	☒
Dark brown / black	☒	☒	☒	☒	☒	☒	☒	☒

☑ Applicable ☒ Not applicable

Image 2.1: IPL Skin and Hair Color Compatibility Chart. Green checkmark indicates applicable, red cross indicates not applicable.

2.3 General Precautions:

- Always wear the provided protective eyewear during operation to shield your eyes from intense light flashes.
- Do not hold the applicator window still over one spot on the skin. Flash once, then move to an adjacent area.
- Wait at least 20 seconds before returning to a previously treated skin spot.
- Keep out of reach of children.
- Consult a physician before use if you have any underlying skin conditions or medical concerns.

3. PACKAGE CONTENTS

Verify that all items are present in the package:



Image 3.1: Package contents including the device, adapter, power cord, eyewear, user manual, and razors.

- MLAY IPL Hair Removal Device
- Power Adapter
- Power Cord
- Protective Eyewear
- User Manual
- Disposable Razors

4. SETUP

1. **Prepare the Skin:** Shave the area to be treated. Do not wax, pluck, or use depilatory creams, as these methods remove the hair follicle, which is necessary for IPL treatment. Clean and dry the skin thoroughly.
2. **Connect Power:** Connect the power adapter to the device and plug it into a suitable electrical outlet.
3. **Wear Eyewear:** Put on the provided protective eyewear before turning on the device.

5. OPERATING INSTRUCTIONS

5.1 Device Overview:



Image 5.1: Device display indicating 500,000 flashes, 5 energy levels, mode display, and power cord connection.

5.2 Power On/Off:

Press and hold the Power button (usually located on the handle) to turn the device on or off.

5.3 Adjusting Power Levels:

The device offers 5 power levels. Press the Power button briefly to cycle through the levels. Start with Level 1 for your first treatment and gradually increase as your skin tolerates it. Higher levels provide more intense treatment.

5.4 Selecting Flash Modes:

The device has two flash modes:

- **Manual Mode:** Press the flash button for each individual flash. Ideal for small, sensitive, or curved areas like the upper lip, underarms, or bikini line.
- **Auto Mode:** Hold the flash button, and the device will automatically emit continuous flashes when in contact with the skin. Suitable for larger areas such as legs, arms, and back.

5.5 Using the Ice Cooling Function:

The ice cooling function can be activated simultaneously with the IPL treatment. This feature helps to cool

the skin surface, minimizing discomfort and protecting the skin from heat. Activate this function via a dedicated button or by a specific setting on the device (refer to the device's display for indication).



Image 5.2: Illustration of the independent cold compress technology, showing cooling to 3°C for painless, comfortable, and safe use.

5.6 Treatment Application:

1. Ensure the device's treatment window is in full contact with your skin.
2. Press the flash button (Manual Mode) or hold it (Auto Mode) to emit a flash.
3. After each flash, slowly move the device to an adjacent, untreated skin area. Do not overlap flashes or treat the same spot multiple times in quick succession.
4. Wait at least 20 seconds before treating a previously flashed area again, if necessary.

Week 1-4 New hair are thinner and begin to fal of
Week 4-8 Reasonable reduction in hair growth
Week 8-12 No more hair regrowth

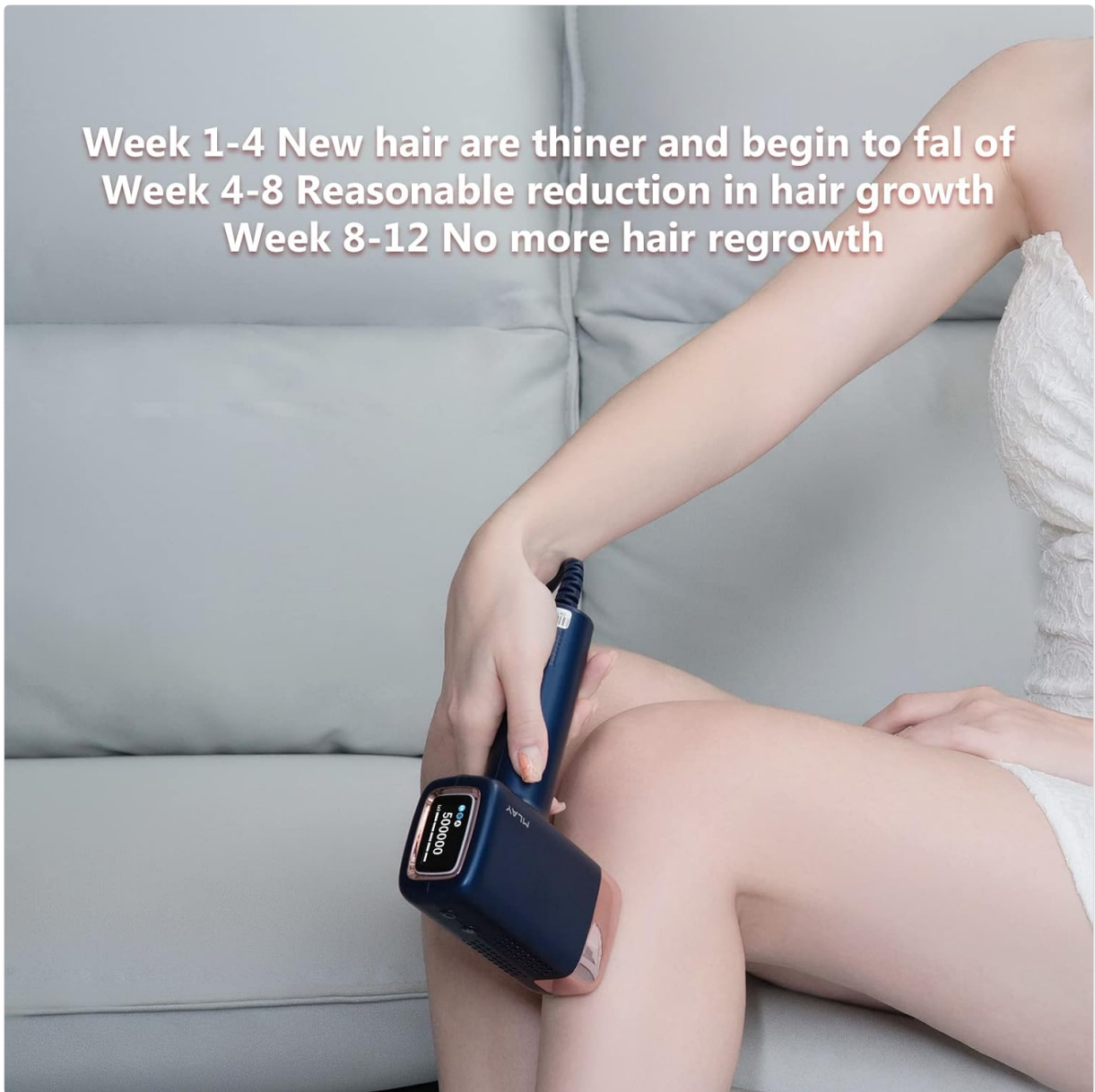


Image 5.3: Proper application of the IPL device on the leg, demonstrating how to hold and move the device.

5.7 Recommended Treatment Schedule:

Consistent use is key for optimal results. A typical treatment plan involves:

- **Weeks 1-3:** Use once per week.
- **Weeks 4-7:** Use once every 2-3 weeks.
- **After Week 8:** Use once every 1-2 months for touch-up applications as needed.

You may observe hair thinning and reduced growth within 4-8 weeks, with lasting effects after 8-12 weeks of consistent treatment.

intensity pulsed light hair removal technique

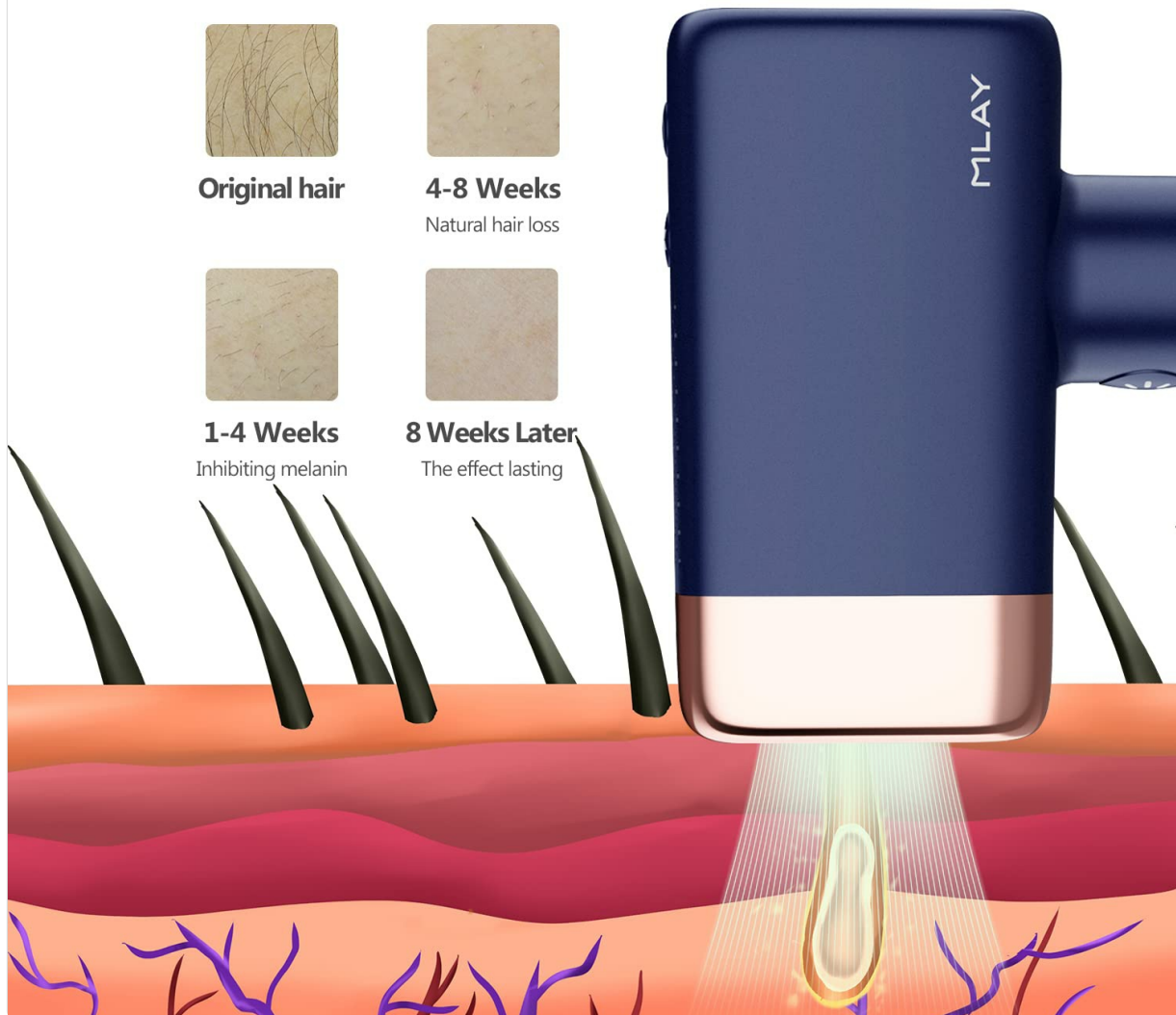


Image 5.4: Visual representation of IPL hair removal technique and expected results over time, from original hair to lasting effect after 8 weeks.

5.8 Applicable Body Areas:

The MLAY IPL Hair Removal Device is suitable for use on various body areas, including:

- Arms
- Armpits
- Legs
- Bikini Line
- Facial areas (avoiding eyes and sensitive mucous membranes)

Full body hair removal

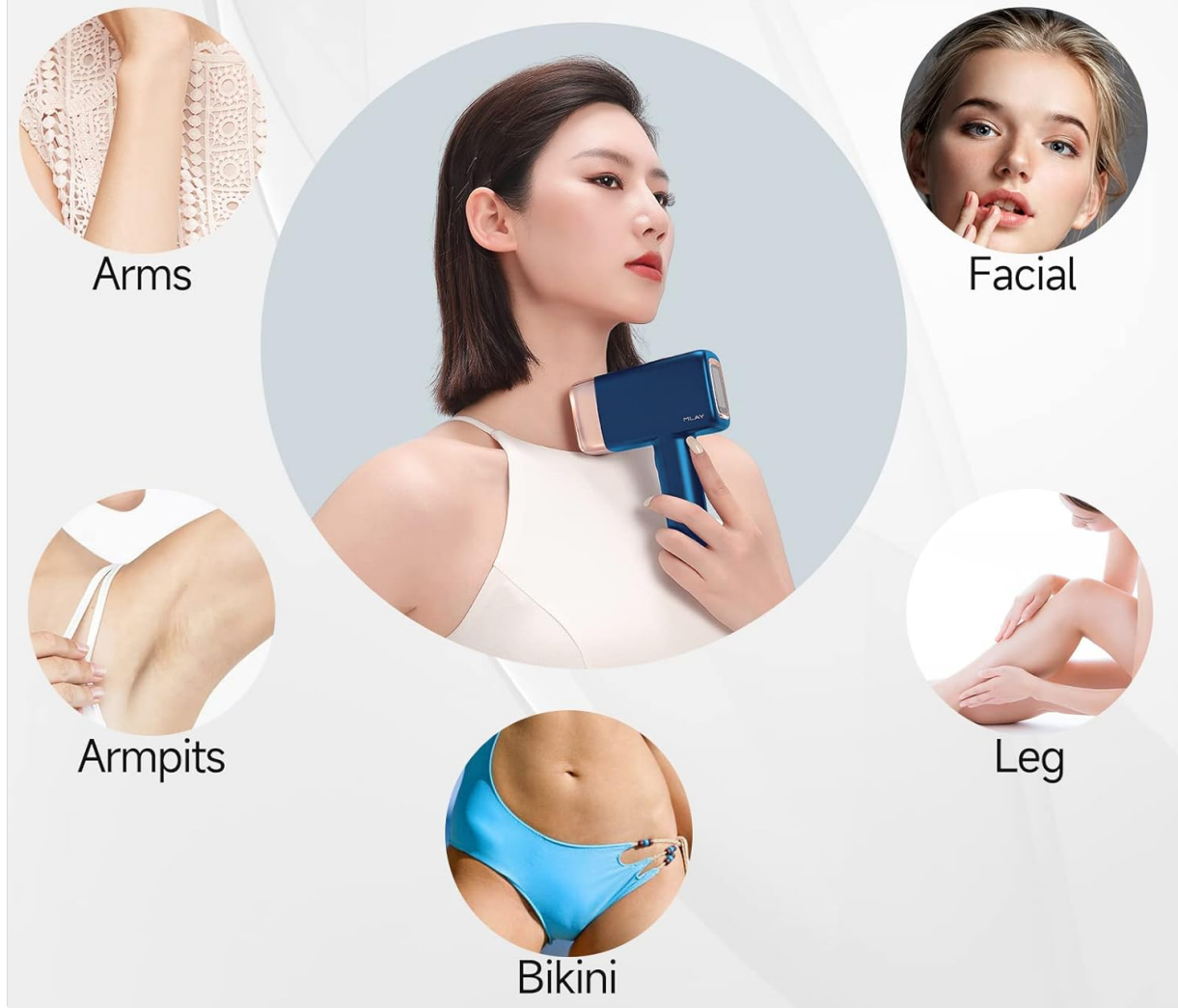


Image 5.5: Illustration of full body hair removal areas where the device can be used.

6. MAINTENANCE

Proper maintenance ensures the longevity and effectiveness of your device.

1. **Cleaning:** After each use, ensure the device is powered off and unplugged. Use a soft, dry cloth to gently wipe the treatment window and the device body. Do not use abrasive cleaners or immerse the device in water.
2. **Storage:** Store the device in a cool, dry place, away from direct sunlight and moisture. Keep it in its original packaging or a protective case to prevent damage.

7. TROUBLESHOOTING

If you encounter issues with your MLAY IPL Hair Removal Device, refer to the following common solutions:

- **Device Not Turning On:**

- Ensure the power adapter is securely connected to both the device and a working electrical outlet.
- Check if the power outlet is functional by plugging in another appliance.

- **No Flash Emitted:**

- Ensure the treatment window is in full contact with the skin. The device has a skin contact sensor that prevents flashing if not properly applied.
- Check the remaining flash count on the display. If it's near zero, the lamp may need replacement or the device has reached its end of life.
- Verify that the device is not in standby mode. Press the power button briefly to wake it up.

- **Skin Discomfort/Redness:**

- Reduce the power level.
- Ensure the ice cooling function is activated.
- Avoid treating the same area too frequently or overlapping flashes.
- If irritation persists, discontinue use and consult a healthcare professional.

For further assistance, please contact MLAY customer support.

8. SPECIFICATIONS

Feature	Detail
Model Number	T-
Brand	MLAY
Technology	IPL (Intense Pulsed Light)
Energy Levels	5
Flash Modes	Manual, Auto
Flashes Lifespan	Up to 500,000
Ice Cooling Temperature	32 ~ 41°F (0 ~ 5°C)
Light Energy	12-18 Joules
Power Source	Wired (not wireless)
Water Resistance	Not Water Resistant
Package Dimensions	9.88 x 7.56 x 3.94 inches
Item Weight	3.06 Pounds

9. WARRANTY AND SUPPORT

For information regarding warranty coverage, product registration, or technical support, please refer to the official MLAY website or contact their customer service directly. Contact details are typically found in the packaging or on the manufacturer's official online channels.

MLAY Official Store: [Visit MLAY Store on Amazon](#)

