

COMFIER CF-2405

COMFIER Back Massager with Heat

INSTRUCTION MANUAL - MODEL CF-2405

1. INTRODUCTION

Thank you for choosing the COMFIER Back Massager with Heat, Model CF-2405. This product is designed to provide a relaxing and therapeutic massage experience for your full back, upper back, and lower back, along with a vibrating seat massage. Featuring advanced shiatsu kneading and soothing heat, this massager is an ideal solution for daily relaxation and pain relief. Please read this manual thoroughly before use to ensure proper operation and safety.

2. SAFETY INSTRUCTIONS

- Always unplug the massager from the power outlet when not in use or before cleaning.
- Do not use the massager in wet or damp environments.
- Keep the massager away from heat sources and direct sunlight.
- Do not use if the power cord or plug is damaged.
- This appliance is equipped with a **15-minute automatic shut-off function** for safety. Do not override this feature.
- Consult a physician before using if you have any medical conditions, are pregnant, or have a pacemaker.
- Do not use on swollen, inflamed, or injured skin areas.
- Keep out of reach of children.

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. This product is for personal, non-medical use.

3. PRODUCT OVERVIEW AND COMPONENTS

The COMFIER Back Massager CF-2405 consists of a back massage cushion with integrated shiatsu nodes and a vibrating seat pad, controlled by a wired remote. It features a detachable flap for intensity adjustment and secure straps for attachment to various chairs.

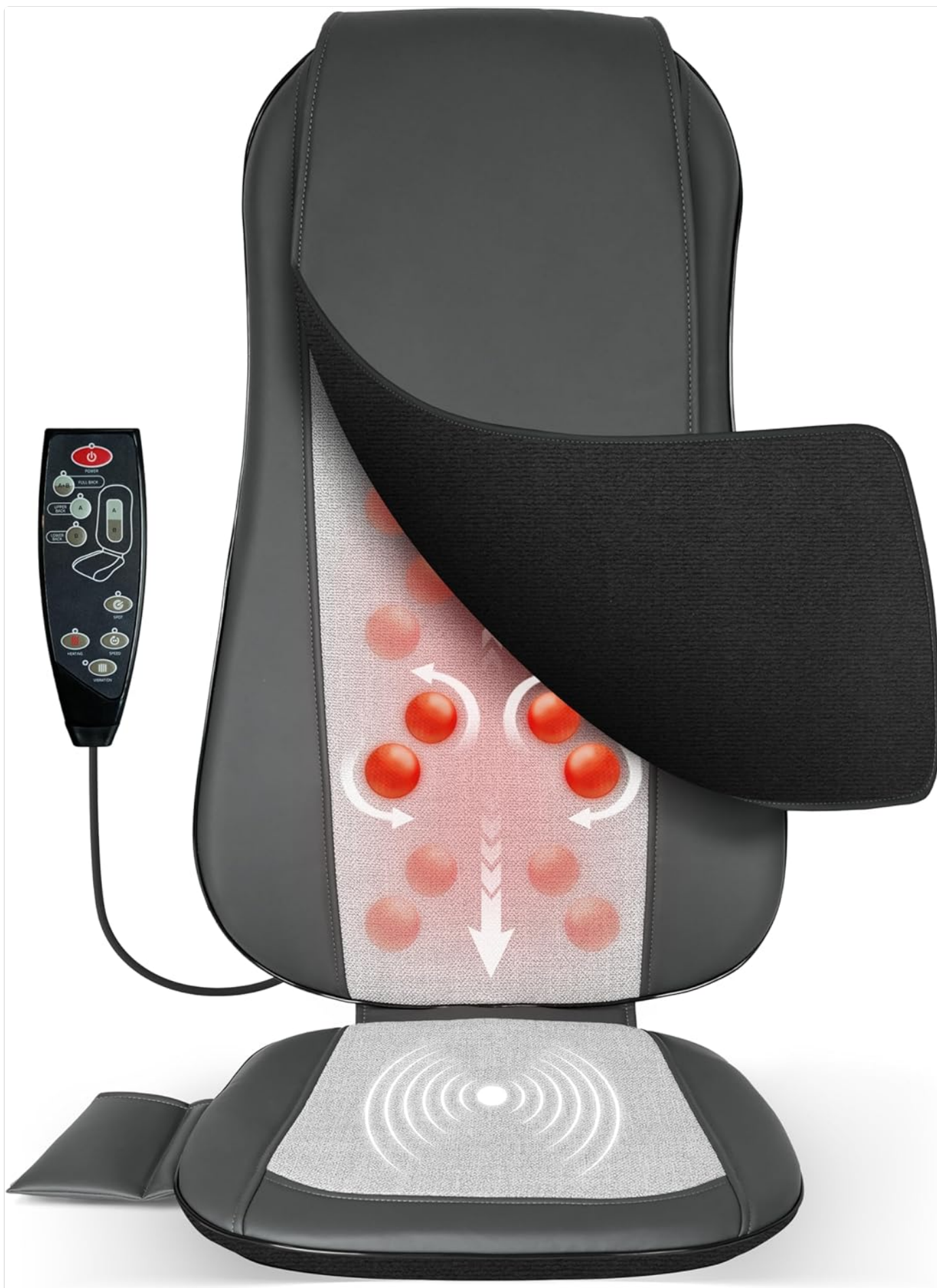


Figure 3.1: COMFIER Back Massager with Heat and its wired remote control.

Easy-to-use remote control



Figure 3.2: Key features include a detachable flap for massage intensity, breathable material for comfort, secure straps, and a convenient remote control kit.

4. SETUP

Setting up your COMFIER Back Massager is quick and easy, allowing you to enjoy a relaxing massage almost anywhere.

4.1 Attaching to a Chair

1. Place the massager on your desired chair, ensuring it is stable and upright.
2. Use the integrated straps on the back of the massager to secure it firmly to the chair. This prevents slipping during use.
3. The massager is designed to fit various chair types, including office chairs, sofa chairs, wooden chairs, and rocking chairs.



Figure 4.1: The massager is versatile and can be used on a variety of chairs for convenience.

4.2 Power Connection


- Connect the power adapter to the massager's input port.
- Plug the power adapter into a suitable electrical outlet.

5. OPERATING INSTRUCTIONS

Your COMFIER Back Massager is controlled by a user-friendly wired remote. Familiarize yourself with its functions for optimal use.

SPOT Massage

Lets you freely choose to focus your massage on one area for precise point-to-point relaxation

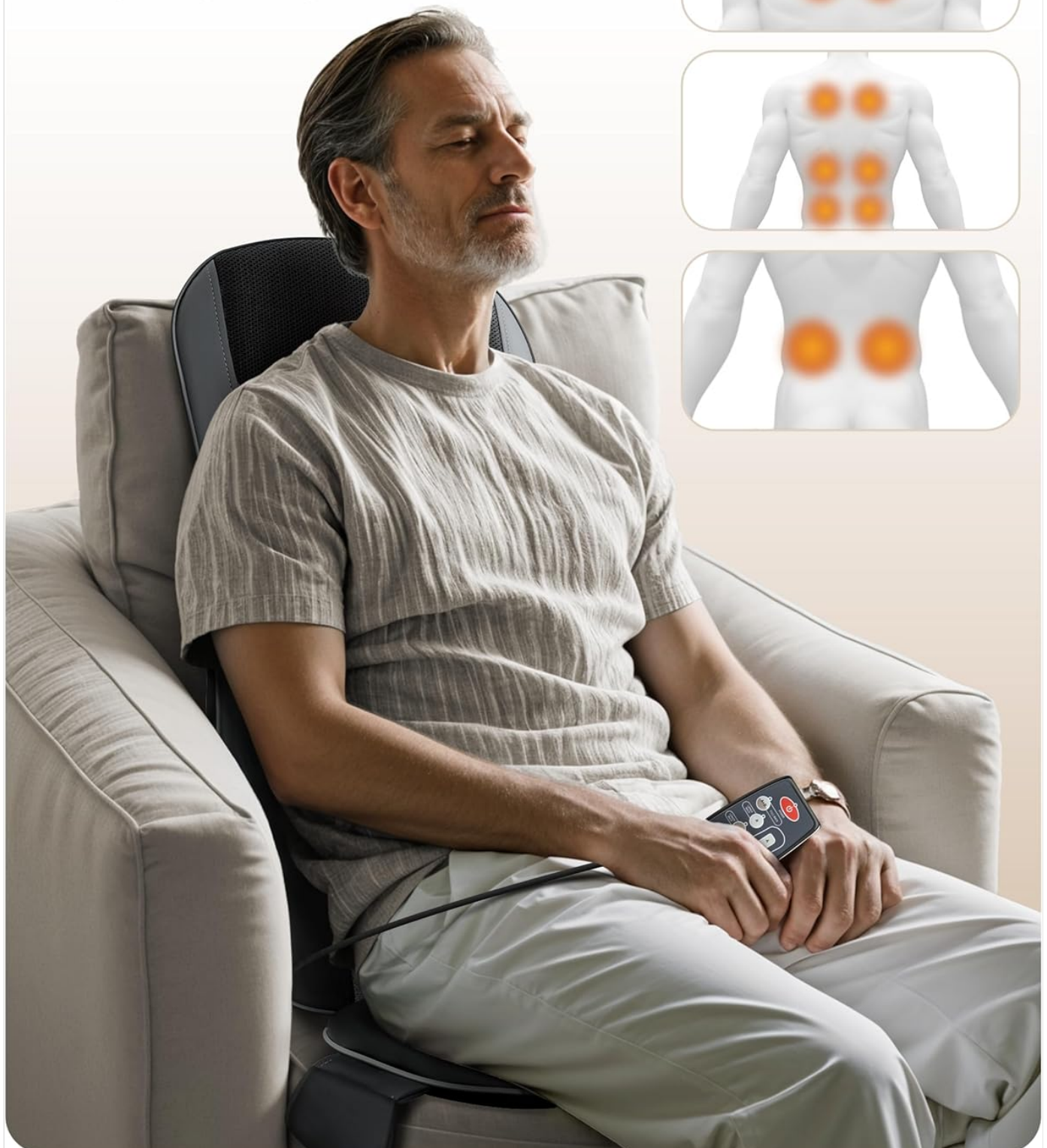

5.1 Remote Control Functions

- **Power Button:** Press to turn the massager ON/OFF.
- **Shiatsu Area Buttons (Full Back, Upper Back, Lower Back):** Select the desired massage zone.
- **SPOT Button:** Activates or deactivates the spot massage function, allowing the shiatsu nodes to focus on a specific area.
- **Shiatsu Speed Buttons:** Adjust between two adjustable kneading speeds for customized relief.
- **Heat Button:** Turn the soothing heat function ON/OFF.
- **Vibration Button:** Turn the vibrating massage function for the seat ON/OFF.

5.2 Massage Modes and Features

- **Customized Shiatsu Massage:** Choose between two adjustable kneading speeds for either soothing relief or deep tissue massage.
- **Targeted Zone Massage:** Select full back, upper back, or lower back massage.
- **SPOT Massage:** Pinpoint a specific area for concentrated massage by pressing the SPOT button.
- **Rhythmic Vibration:** The seat pad provides a comfortable vibrating massage for the buttocks and thighs.
- **Soothing Heat:** Activate the heat function for added comfort and muscle relaxation. The heat and vibration functions can be controlled independently.
- **Detachable Flap:** For a gentler massage, the removable flap can be placed over the shiatsu nodes to reduce intensity.

Relaxtion in Every Touth



2 Levels Shiatsu Speed



Rhythmic Vibration

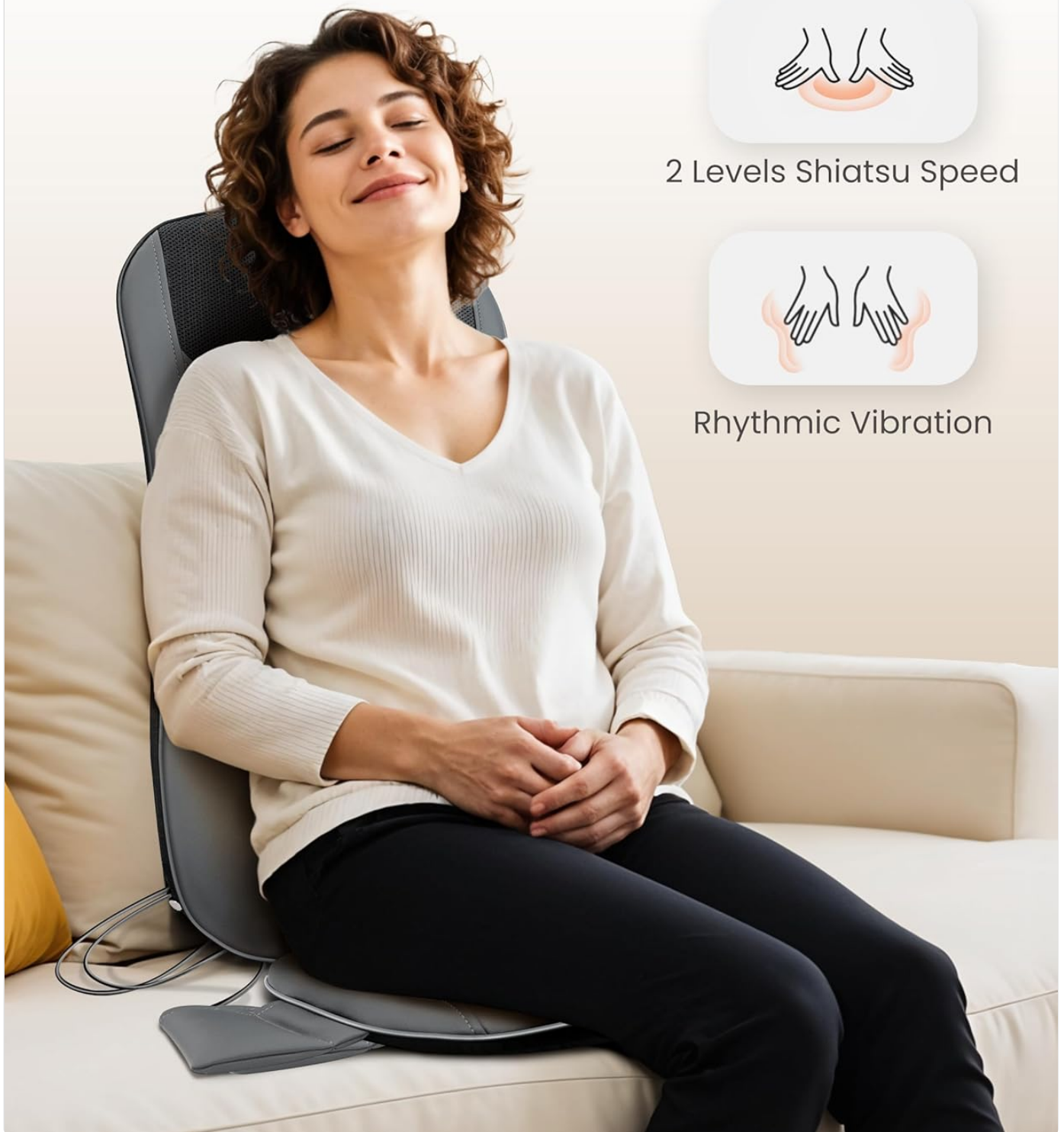


Figure 5.2: Experience personalized relaxation with adjustable shiatsu speeds and rhythmic vibration.

Intelligent Timing

Experience Safe Relaxation with COMFIER Shiatsu Massage Cushion with Heat



15 min Auto-off

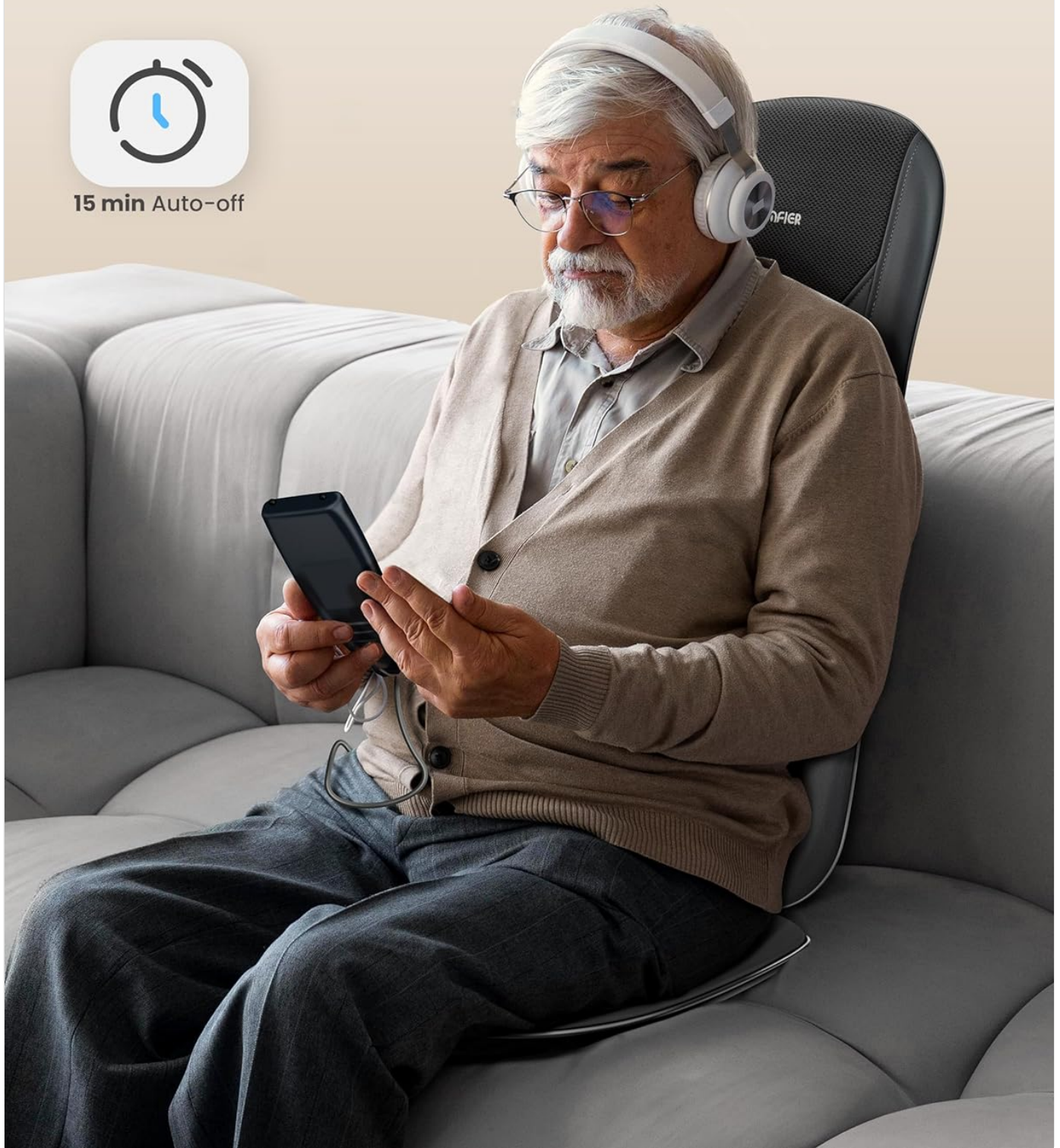


Figure 5.3: The SPOT massage function allows you to target specific areas for precise relief.

5.3 Automatic Shut-off

The massager is equipped with a **15-minute automatic shut-off timer** to ensure safety and prevent overheating. After 15 minutes of continuous use, the device will automatically turn off. If you wish to continue your massage, simply restart the unit.



Figure 5.4: The intelligent 15-minute auto-off timer ensures safe and worry-free operation.

6. MAINTENANCE

6.1 Cleaning

- Always unplug the massager before cleaning.

- Wipe the surface of the massager with a soft, slightly damp cloth.
- Do not use abrasive cleaners, brushes, or strong chemicals.
- Allow the massager to dry completely before storing or using again.

6.2 Storage

- Store the massager in a cool, dry place, away from direct sunlight.
- Avoid placing heavy objects on the massager.
- The lightweight design (6.1 lbs) makes it easy to store and transport.

7. TROUBLESHOOTING

If you encounter any issues with your COMFIER Back Massager, please refer to the following common troubleshooting tips:

Problem	Possible Cause	Solution
Massager does not turn on.	Not properly plugged in; Power outlet issue; Damaged cord/adapter.	Ensure all connections are secure. Try a different outlet. If cord/adapter is damaged, discontinue use and contact support.
Massage feels too intense.	No flap used; High intensity setting.	Place the detachable flap over the shiatsu nodes. Adjust the intensity setting on the remote.
Massager stops unexpectedly.	Automatic 15-minute shut-off activated.	This is a safety feature. Simply restart the massager using the power button on the remote.

If the problem persists after trying these solutions, please contact COMFIER customer support for assistance.





8. SPECIFICATIONS



Feature	Detail
Model Number	CF-2405
Product Dimensions	16.54 x 3.94 x 42.13 inches
Item Weight	6.1 Pounds
Material	Leather
Power Source	Corded Electric
Manufacturer	COMFIER
Country of Origin	China

9. WARRANTY AND SUPPORT

Please retain your proof of purchase for any warranty-related inquiries.

Related Documents - CF-2405

 <p>CONFIER FUSSMASSAGEBETRIEB DER FOOT MASSAGER 50 MASSEUR DE PIEDS 50 MASSAGGIATORE PIEDI 50 MASSAGEUR DE PIEDS 50</p> <p>CE User Manual CF-5913</p>	<p>Comfier CF-5913 Shiatsu Foot Massager User Manual</p> <p>User manual for the Comfier CF-5913 Shiatsu Foot Massager Machine. Learn how to use this kneading foot massager with heat for pain relief, plantar fasciitis, and general relaxation. Includes operation, safety, maintenance, and warranty details.</p>
 <p>CONFIER Shiatsu Neck and Back Massager with Heat</p> <p>CE User Manual CF-2303-KR</p>	<p>COMFIER CF-2303-KR Shiatsu Neck and Back Massager with Heat User Manual</p> <p>Comprehensive user manual for the COMFIER CF-2303-KR Shiatsu Neck and Back Massager with Heat, covering safety instructions, features, technical data, setup, operation, controller functions, care, cleaning, troubleshooting, and warranty information.</p>
 <p>CONFIER Shiatsu massagecushion mit wärmeleitung 50 Shiatsu massage cushion with heat 50 Coussin de massage Shiatsu avec chaleur 50 Cuscino per massage Shiatsu riscaldato 50 Cushion massage Shiatsu con calefacción 50</p> <p>CE User Manual CF-2606MC</p>	<p>COMFIER CF-2606MC Shiatsu Massage Cushion with Heat User Manual</p> <p>Comprehensive user manual for the COMFIER CF-2606MC Shiatsu massage cushion with heat. Includes safety instructions, operating guide, features, maintenance, and warranty information.</p>
 <p>CONFIER Shiatsu massage pillow with heat</p> <p>CE User Manual CF-1228</p>	<p>Comfier CF-1228 Shiatsu Massage Pillow with Heat User Manual</p> <p>User manual for the Comfier CF-1228 Shiatsu massage pillow with heat, detailing safety instructions, operating procedures, features, maintenance, disposal, technical specifications, and warranty information.</p>

<div><div>COMFIER</div><div>Shiatsu Neck and Back Massager</div><div></div><div>User Manual CF-2309A</div></div>	<p>Comfier Shiatsu Neck and Back Massager CF-2309A User Manual</p> <p>User manual for the Comfier Shiatsu Neck and Back Massager, model CF-2309A. Includes safety instructions, features, technical data, setup and operation, controller instructions, care and cleaning, troubleshooting, and warranty information.</p>
<div><div>COMFIER</div><div>Shiatsu neck-and-back-massager DE Shiatsu neck and back massager GB Massager de cou et de dos Shiatsu FR Masajeador Shiatsu por cuello y espalda IT Masajeador Shiatsu de para cuello y espalda ES Shiatsu neck- en rugmassagereenlat. NL</div><div></div><div>CE User Manual CF-2309A</div></div>	<p>Comfier CF-2309A Shiatsu Neck and Back Massager User Manual</p> <p>Comprehensive user manual for the Comfier CF-2309A Shiatsu Neck and Back Massager, covering safety instructions, features, operation, maintenance, technical specifications, and warranty information.</p>