

Ulike UI06 PR

Ulike Air 3 Ice-Cooling IPL Device User Manual

Model: UI06 PR | Brand: Ulike

ABOUT THE ULIKE AIR 3 IPL DEVICE

The Ulike Air 3 IPL Hair Removal Device offers a nearly painless solution for long-lasting hair reduction from the comfort of your home. Utilizing advanced IPL (Intense Pulsed Light) technology with a patented ice-cooling touch, it maintains a skin temperature of 67°F (19.4°C) during treatment, ensuring comfort and minimizing discomfort. This device is designed for both women and men, providing noticeable hair reduction within approximately 3 weeks of consistent use.

Key features include:

- **Nearly Painless Treatment:** Patented Ice-Cooling technology maintains a cool sensation on the skin.
- **Fast Results:** Visible hair reduction can be observed in as little as 3 weeks.
- **Multiple Modes:** Features Soft Mode for delicate areas, Body Mode for larger areas, and Power Mode for stubborn hair.
- **Auto-Flashing:** For quick and efficient full-body hair removal.
- **Safe IPL Technology:** Designed for safety and effectiveness across most suitable skin and hair types.

IMPORTANT SAFETY INFORMATION

Always read the user manual carefully before using the device, including all contraindications, warnings, and safety information.

- Not suitable for tattooed and over-tanned skin, wounded areas, or pregnant individuals.
- Perform a skin safety test before full treatment.
- Always wear the provided IPL goggles to protect your eyes from the intense pulsed light during the entire treatment process.
- This is a plug-in device; ensure it is always connected to power during treatment.
- Not water resistant. Do not use in wet environments.

PACKAGE CONTENTS

Your Ulike Air 3 IPL Device package includes the following items:

- Ulike Air 3 IPL Hair Removal Handset
- Power Adapter and Cord
- IPL Goggles for eye protection
- Razor for pre-treatment hair removal
- User Manual (PDF version also available online)



Image: Ulike Air 3 IPL Device, razor, and protective goggles.

SETUP INSTRUCTIONS

1. **Prepare the Skin:** Before using the device, ensure the treatment area is clean, dry, and free of any lotions or

creams. Shave any visible hair from the area using the provided razor or your preferred method. The IPL device works most effectively on hair follicles, not surface hair.

2. **Connect to Power:** Plug the power cord into the Ulike Air 3 handset and then into a power outlet. The device is designed to be used while plugged in.
3. **Wear Eye Protection:** Put on the provided IPL goggles to protect your eyes from the intense pulsed light during the entire treatment process.
4. **Turn On the Device:** Long press the power button on the device to activate it. The indicator lights will illuminate.

Check if the Air 3 right for you?

*Please notice those warm tips before using:

- 1.If you are not sure of your skin type, please refer to the skin tone table before using.
- 2.Pregnant are not suitable to use this hair removal device.
- 3.Not suitable for tattooed, allergic to light, or over-tanned skin.

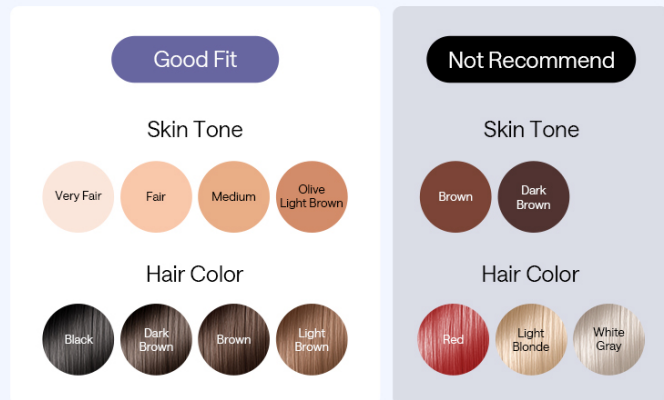


Image: Visual guide for connecting the device to power and turning it on.

OPERATING MODES AND USAGE

The Ulike Air 3 offers customized treatment modes for different body areas and hair types:

- **Soft Mode:** Ideal for sensitive areas like the face and upper lip.
- **Body Mode:** Suitable for larger areas such as legs and arms.
- **Power Mode:** Designed for coarser hair and stubborn areas like underarms and bikini line.

To select a mode, short press the power button to cycle through the intensity levels. The indicator lights will show the current mode. The device also features an auto-glide mode for continuous flashing, making full-body treatments faster.

How to Use:

1. Place the light window vertically against your skin, ensuring full contact.
2. Press the treatment button to emit a light pulse. For auto-glide mode, hold the button down and slide the device across the skin.
3. Move the device to the next area, ensuring slight overlap to cover all hair follicles.
4. After treatment, you may apply a non-scented lotion if desired.



Image: Customized modes for different body areas.



Image: Auto-Glide and Manual Mode operation.

USAGE GUIDELINES AND SKIN COMPATIBILITY

IPL technology requires a contrast between hair melanin and skin melanin for effective and safe treatment. Refer to the skin tone and hair color chart below to determine if the Ulike Air 3 is suitable for you.



Image: Skin tone and hair color compatibility chart for Ulike Air 3.

Recommended Treatment Frequency:

- **Weeks 1-3:** 3 times per week.
- **Weeks 4-6:** 2 times per week.
- **After 6 Weeks:** 1-2 times per month for maintenance.

Press Treatment Button To Start Hair Removal

Put the light window vertically to your skin, and press treatment button to start hair removal.



Image: Recommended treatment frequency for optimal results.

MAINTENANCE

Proper maintenance ensures the longevity and effectiveness of your Ulike Air 3 device:

- **Cleaning:** After each use, wipe the light window with a soft, dry cloth to remove any residue. Do not use abrasive cleaners or solvents.
- **Storage:** Store the device in its original packaging or a clean, dry place away from direct sunlight and moisture.
- **Handle with Care:** Avoid dropping the device or subjecting it to strong impacts.

TROUBLESHOOTING

If you encounter any issues with your Ulike Air 3 device, refer to the following common troubleshooting tips:

- **Device Not Turning On:** Ensure the device is properly plugged into a working power outlet. Check the power cord for any damage.
- **No Light Pulse:** Verify that the light window is in full contact with your skin. The device has a safety mechanism that prevents flashing if not properly placed.
- **Skin Irritation:** Ensure you are using the correct intensity level for your skin tone. Perform a patch test before full treatment. If irritation persists, discontinue use and consult a dermatologist.
- **Hair Not Reducing:** Ensure consistent use as per the recommended treatment frequency. Shave the area thoroughly before each treatment.

For more detailed troubleshooting or persistent issues, please refer to the full user manual or contact Ulike customer support.

SPECIFICATIONS


Feature	Detail
Brand	Ulike
Model Number	UI06 PR
Product Dimensions	3.66 x 0.9 x 0.04 inches
Item Weight	4 Pounds
Color	Purple
Power Source	Corded Electric
Water Resistance Level	Not Water Resistant
Manufacturer	Shenzhen Ulike Smart Electronics Co.,Ltd
First Available Date	March 7, 2023

WARRANTY AND SUPPORT

For information regarding product warranty, returns, or technical support, please refer to the official Ulike website or contact their customer service directly. Details can typically be found in the full user manual provided with your device or on the product's official online listing.

You can also find additional resources and contact information by visiting the [Ulike Store on Amazon](#).

Related Documents - UI06 PR

	<p>Ulike IPL Hair Removal Device Quick Start Guide</p> <p>A quick start guide for the Ulike IPL Hair Removal Device, covering safety precautions, product information, operation, and troubleshooting.</p>
---	--

 The image shows the cover of the Ulike IPL Hair Removal Device Quick Start Guide. It features the Ulike logo at the top, followed by the title 'QUICK START GUIDE' and 'IPL Hair Removal Device'. Below this, it mentions '90-Day Hassle-Free Return' and '2-Year Warranty'. The central image is a white, rectangular device with a screen and buttons. At the bottom, it says 'The hair removal options Ulike'.	<p>Ulike IPL Hair Removal Device Quick Start Guide</p> <p>This guide provides instructions for using the Ulike IPL Hair Removal Device, including safety precautions, operation principles, application ranges, treatment methods, maintenance, and troubleshooting.</p>
 The image shows the cover of the Ulike IPL Hair Removal Device Quick Start Guide. It features the Ulike logo at the top, followed by the title 'QUICK START GUIDE' and 'IPL Hair Removal Device'. Below this, it mentions '90-Day Hassle-Free Return' and '2-Year Warranty'. The central image is a white, rectangular device with a screen and buttons. At the bottom, it says 'The hair removal options Ulike'.	<p>Ulike IPL Hair Removal Device Quick Start Guide</p> <p>A quick start guide for the Ulike IPL Hair Removal Device, detailing its features, usage instructions, and safety precautions. Learn about IPL technology, how to use the device for effective hair removal, and what to expect over a four-week period.</p>
 The image shows the cover of the Ulike IPL Hair Removal Device Quick Start Guide. It features the Ulike logo at the top, followed by the title 'QUICK START GUIDE' and 'IPL Hair Removal Device'. Below this, it mentions '90-Day Hassle-Free Return' and '2-Year Warranty'. The central image is a white, rectangular device with a screen and buttons. At the bottom, it says 'The hair removal options Ulike'.	<p>Ulike IPL Hair Removal Device Quick Start Guide</p> <p>Comprehensive guide to using the Ulike IPL Hair Removal Device, detailing product features, step-by-step instructions, safety warnings, and expected results for effective hair removal.</p>
 The image shows the cover of the Ulike Air IPL Hair Removal Device: User Manual, Safety, and How-To Guide. It features the Ulike logo at the top, followed by the title 'PRODUCT DESCRIPTION' and 'IPL Hair Removal Device'. The central image is a white, rectangular device with a screen and buttons. At the bottom, it says 'AIR'.	<p>Ulike Air IPL Hair Removal Device: User Manual, Safety, and How-To Guide</p> <p>Comprehensive guide to the Ulike Air IPL Hair Removal Device, covering safety precautions, usage instructions, suitability, maintenance, and troubleshooting for effective and safe hair removal.</p>
 The image shows the cover of the Ulike IPL Hair Removal Guide: Treatment Areas, Schedule & Intensity Settings. It features the Ulike logo at the top, followed by the title 'PRODUCT DESCRIPTION' and 'IPL Hair Removal Device'. The central image is a white, rectangular device with a screen and buttons. At the bottom, it says 'AIR'.	<p>Ulike IPL Hair Removal Guide: Treatment Areas, Schedule & Intensity Settings</p> <p>A concise guide to using the Ulike IPL hair removal device, covering preparation, treatment areas (armpit, face, arm, leg, Brazilian wax area), recommended treatment schedules, and intensity settings. Includes important usage precautions.</p>

