Nuvita NU-IBKT0041

Nuvita 1147 Newborn Baby Care Beauty Set Instruction Manual

Model: NU-IBKT0041

[Instructions](#) [Maintenance](#) [Product Overview](#) [Setup](#) [Operating](#) [Troubleshooting](#) [Specifications](#) [Warranty & Support](#)

1. PRODUCT OVERVIEW

The Nuvita 1147 Newborn Baby Care Beauty Set is a comprehensive kit designed for the daily hygiene and care of infants. This set includes 9 essential accessories, all neatly organized within a stylish and functional travel case.

Included Components:

- Digital Thermometer with Flexible Tip
- Natural Wool Brush
- Comb
- Baby Scissors
- Baby Nail Clippers
- 4 Nail Files
- Nasal Aspirator with Silicone Tip
- Eco-leather Travel Case

Figure 1: The complete Nuvita 1147 Baby Care Beauty Set, showcasing all included items within and outside the travel case.

2. SETUP

Before first use, please inspect all components for any damage. Clean all items according to the 'Maintenance' section. Ensure the digital thermometer has a working battery.

Unpacking and Initial Inspection:

1. Carefully remove all items from the travel case.
2. Verify that all 9 components listed in the 'Product Overview' are present.
3. Dispose of all packaging materials responsibly.

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Figure 2: The Nuvita Baby Care Beauty Set with its travel case opened, showing the organized compartments for each accessory.

3. OPERATING INSTRUCTIONS

3.1 Digital Thermometer with Flexible Tip

This thermometer provides accurate temperature measurements in approximately 30 seconds. It is suitable for oral, axillary (underarm), or rectal use.

1. Press the ON/OFF button to activate the thermometer.
2. Place the flexible tip in the desired measurement area (oral, axillary, or rectal).
3. Wait for approximately 30 seconds until the reading stabilizes and a beep sounds.
4. Read the temperature displayed on the screen.
5. Press the ON/OFF button again to turn off the thermometer, or it will automatically shut off after a period of inactivity.

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Figure 3: The digital thermometer with its flexible tip, designed for precise and comfortable temperature measurement.

3.2 Natural Wool Brush and Comb

These tools are designed for gentle grooming of your baby's hair and scalp.

- **Natural Wool Brush:** Features soft bristles ideal for delicate baby hair and stimulating the scalp.
- **Comb:** Designed with an ergonomic shape and rounded tips to prevent scratching the scalp. It is particularly useful for gently removing cradle cap after applying oil.

Figure 4: The natural wool brush being used to gently comb a baby's hair.

Figure 5: Close-up of the comb and natural wool brush, highlighting their ergonomic design and soft materials.

3.3 Baby Nail Scissors, Nail Clippers, and Nail Files

These tools are designed for safe and gentle nail care for infants.

- **Baby Nail Scissors:** Feature rounded tips for safety and are made from hypoallergenic, nickel-free material. Use with caution to trim nails straight across.
- **Baby Nail Clippers:** Made from hypoallergenic material. Ensure a firm grip and carefully trim nails to avoid cutting the skin.
- **Nail Files (4 pieces):** Use these cardboard files to gently smooth any sharp edges after trimming, preventing scratches.

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4 limette in
cartone

Forbici con punta
arrotondata

Figure 6: The baby nail care tools: clippers, scissors with rounded tips, and nail files.

3.4 Nasal Aspirator with Silicone Tip

This aspirator helps clear your baby's nasal passages gently and effectively.

1. Gently squeeze the bulb of the aspirator to expel air.
2. Carefully insert the soft silicone tip into one of your baby's nostrils.
3. Slowly release the bulb to create suction, drawing mucus into the aspirator.
4. Remove the aspirator and squeeze the bulb to expel the collected mucus into a tissue.
5. Repeat for the other nostril if necessary.

Figure 7: The nasal aspirator, featuring a soft and flexible silicone tip for comfortable use.

4. MAINTENANCE

Proper cleaning and maintenance will ensure the longevity and hygiene of your Nuvita Baby Care Beauty Set.

General Cleaning Instructions:

- **Thermometer:** Wipe the tip with an alcohol swab or a soft cloth dampened with mild soap and water. Do not immerse the entire thermometer in water.
- **Brush and Comb:** Wash with mild soap and warm water. Rinse thoroughly and allow to air dry completely.
- **Nail Scissors, Clippers, and Files:** Clean with a damp cloth after each use. Ensure they are dry before storing to prevent rust (for metal parts). The cardboard files are disposable and should be replaced when worn.
- **Nasal Aspirator:** Disassemble the aspirator (if possible) and wash all parts with warm soapy water. Rinse thoroughly to remove all soap residue and allow to air dry completely before reassembling.
- **Travel Case:** Wipe the eco-leather exterior with a damp cloth. Do not machine wash.

Always store all components in the provided travel case in a clean, dry place when not in use.

5. TROUBLESHOOTING

This section addresses common issues you might encounter with the Nuvita 1147 Beauty Set.

Problem	Possible Cause	Solution
---------	----------------	----------

Problem	Possible Cause	Solution
Digital thermometer not turning on.	Battery is low or depleted.	Replace the battery (refer to thermometer's specific instructions for battery type and replacement).
Inaccurate thermometer readings.	Improper placement or insufficient measurement time.	Ensure the tip is correctly placed and allow the full 30 seconds for measurement. Avoid movement during measurement.
Nasal aspirator not creating suction.	Tip not sealed against nostril or aspirator is clogged.	Ensure a tight seal around the nostril. Disassemble and clean the aspirator thoroughly as per maintenance instructions.
Brush bristles or comb teeth are dirty.	Accumulation of hair, skin flakes, or product residue.	Clean the brush and comb with mild soap and water, then rinse and air dry.

6. SPECIFICATIONS

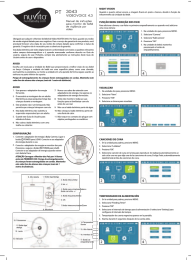
Feature	Detail
Product Model Number	NU-IBKT0041
Brand	Nuvita
Product Dimensions	0.17 x 0.19 x 0.04 cm (approximate case dimensions)
Item Weight	290 grams
Batteries Required	No (except for thermometer, which typically uses a button cell battery, not specified in JSON if included or type)
Materials	Natural wool (brush), hypoallergenic/nickel-free (scissors/clippers), silicone (aspirator tip), eco-leather (case)

7. WARRANTY & SUPPORT

For information regarding warranty coverage, product support, or to purchase replacement parts, please contact Nuvita customer service directly. Details can typically be found on the product packaging or the official Nuvita website.

Keep your purchase receipt as proof of purchase for any warranty claims.

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