

[Manuals.plus](#) /

› [JOVS](#) /

› JOVS IPL Hair Removal Device JR5-E User Manual

JOVS JR5-E

JOVS IPL Hair Removal Device JR5-E User Manual

Model: JR5-E | Brand: JOVS

1. INTRODUCTION

This manual provides essential instructions for the safe and effective operation of your JOVS IPL Hair Removal Device, Model JR5-E. This device utilizes Intense Pulsed Light (IPL) technology for hair removal and skin rejuvenation, featuring up to 24J energy output, a 37°F Sapphire Ice-Cooling system for comfort, and 6 specialized modes for various body areas.

Please read this manual thoroughly before first use and retain it for future reference.



Image: The JOVS IPL Hair Removal Device JR5-E, showcasing its ergonomic design and control panel.

2. SAFETY INFORMATION AND PRECAUTIONS

Your safety is paramount. Please adhere to the following guidelines:

- **Eye Protection:** Always wear the provided safety goggles during operation to protect your eyes from intense pulsed light. Do not look directly into the light emission window.
- **Skin Preparation:** Ensure skin is clean, dry, and shaved before treatment. Avoid waxing, plucking, or chemical depilatories for at least 4 weeks prior to treatment, as the device targets hair follicles.
- **Skin Sensitivity Test:** Before full treatment, perform a patch test on a small skin area. Wait 24 hours to check for adverse reactions.
- **Skin Tone and Hair Color:** This device is effective on specific skin tones and hair colors. Refer to the "Skin Tone and Hair Color Suitability" section to ensure compatibility. Not suitable for very dark skin tones or very light/red/grey hair.

- **Avoid Certain Areas:** Do not use on tattoos, permanent makeup, moles, freckles, birthmarks, or areas with skin disorders, infections, burns, or recent sun exposure. Avoid using on mucous membranes (e.g., lips, nostrils, genitals).
- **Medical Conditions:** Consult a physician before use if you have any medical conditions, are pregnant, breastfeeding, or taking photosensitizing medications.
- **Children and Pets:** Keep the device out of reach of children and pets. This device is not intended for use by individuals under 18 years of age.
- **Power Supply:** Use only the original power adapter provided with the device.
- **Water Exposure:** Do not use the device near water or immerse it in water.

3. PACKAGE CONTENTS

Verify that all items are present in your package:



Image: Diagram showing the contents of the JOVS IPL Hair Removal Device package.

- 1 x JOVS Venus Pro II Device
- 1 x SR Head (Skin Rejuvenation)
- 1 x Guide Card
- 1 x Instructions Manual
- 1 x Safety Goggles
- 1 x Razor
- 1 x Power Adapter
- 1 x Power Cable

4. SETUP AND DEVICE OVERVIEW

4.1 Initial Setup

1. **Prepare Skin:** Shave the surface hair from the area you intend to treat. Do not wax or pluck.
2. **Connect Power:** Connect the power adapter to the device's DC socket, then plug it into a power outlet.
3. **Wear Goggles:** Put on the provided safety goggles.

4.2 Device Controls

Packing list


① x Packing box


② x Venus Pro II device
③ x SR head


④ x Guide card


⑤ x Instructions


⑦ x Goggles
⑥ x Razor


⑨ x Power Adapter
⑧ x Power Cable

Packaged & Specification

Product name	JOVS Venus pro II hair removal		
Weight	318g	Energy output	2-6 J/cm²
Cool sapphire area	12.1cm²	Flash interval	0.7s
Rated power	60 W	Wavelength	590 ~ 1200 nm
Number of flashes	Infinite	System	OPT technology
Input	100-240V~50/60Hz, 1.5A		
Output	24V=2.5A		
Function	Hair Removal, Skincare		
After-sales service	24 months		

Image: Detailed diagram of the JOVS IPL device controls and their functions.

- **Energy Level Switching Ring:** Rotate this ring (0-6) to adjust the energy intensity. Level 0 is OFF, Levels 1-6 are ON with increasing intensity.
- **Six Modes Status:** Indicators for the selected treatment mode.
- **Auto Mode Touch Button:** Touch to turn auto mode ON/OFF.
- **Six Modes Switching Button:** Click to switch between the six available treatment modes.
- **Control Button:** Click once for a single IPL treatment flash. Long press for 3 seconds to switch to Auto Mode.
- **DC Socket:** For connecting the power adapter.

5. OPERATING INSTRUCTIONS

5.1 Treatment Steps

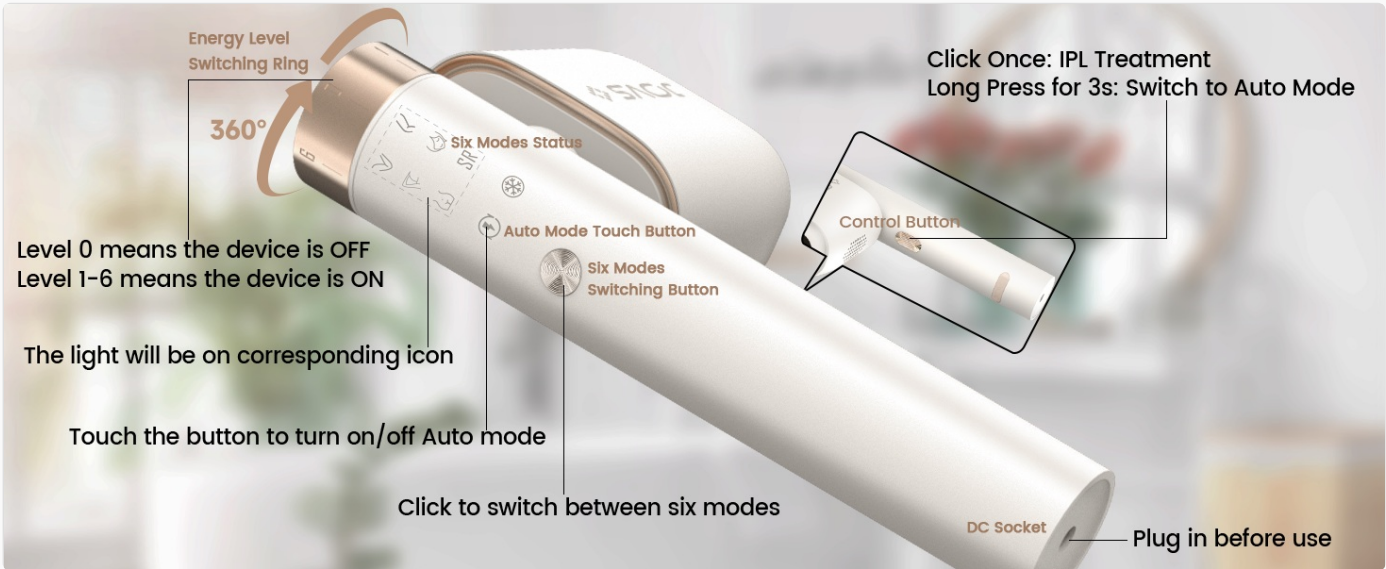


Image: Step-by-step guide on how to use the JOVS IPL device for hair removal.

1. **Shave:** Ensure the treatment area is shaved clean.
2. **Connect & Power On:** Connect the device to power and turn it on using the Energy Level Switching Ring.

3. **Protect Eyes:** Wear the safety goggles.
4. **Select Mode & Intensity:** Use the Six Modes Switching Button to choose the appropriate mode for the body area (e.g., Arm, Leg, Face, Bikini, Armpit, SR). Adjust the energy level (1-6) using the Energy Level Switching Ring based on your skin sensitivity and hair type.
5. **Apply to Skin:** Place the device's light emission window vertically and firmly against your skin. The device will only flash when full contact is detected.
6. **Initiate Flash:**
 - For single flashes: Click the Control Button once.
 - For continuous flashes (Auto Mode): Long press the Control Button for 3 seconds. The device will flash automatically as you glide it over your skin.
7. **Post-Treatment Care:** After treatment, you may apply a soothing skincare product to the treated area.

5.2 Treatment Modes and Areas

The JOVS IPL device offers 6 specialized modes, each optimized for different body areas and hair characteristics:



Image: Chart illustrating the 6 different light wave modes, their wavelengths, and target hair depths/diameters for precise hair removal.

- **SR Mode (Skin Rejuvenation):** 640-1200nm wavelength.
- **Face Mode:** 530-1200nm wavelength.
- **Arm Mode:** 550-1200nm wavelength.
- **Leg Mode:** 560-1200nm wavelength.
- **Armpit Mode:** 580-1200nm wavelength.
- **Bikini Mode:** 620-1200nm wavelength.

Each mode adjusts the pulsed light's power and wavelength to suit the specific characteristics of hair and skin in that area.

5.3 Treatment Frequency

How to use ?

Perfect Salon Level Effect



Image: Recommended treatment schedule for the JOVS IPL device, showing frequency over several weeks.

- **Weeks 1-2:** 3-4 times per week.
- **Weeks 2-3:** 2-3 times per week.
- **After 4 Weeks:** 2 times monthly or once monthly, as needed for maintenance.

Consistency is key for optimal results. Hair growth cycles vary, so continued use is necessary to target all hair follicles.

6. SKIN TONE AND HAIR COLOR SUITABILITY

The effectiveness and safety of IPL hair removal depend on your skin tone and hair color. IPL technology targets the melanin in hair follicles. Therefore, it is most effective on darker hair colors and lighter skin tones.



Image: Chart indicating suitable skin tones and hair colors for JOVS IPL device use.

- **Suitable Skin Tones:** Light to medium skin tones.
- **Unsuitable Skin Tones:** Very dark skin tones (Fitzpatrick skin types V and VI) due to the risk of skin damage.
- **Suitable Hair Colors:** Black, dark brown, brown, light brown hair.
- **Unsuitable Hair Colors:** White, grey, red, or very light blonde hair, as these colors contain insufficient melanin for IPL to be effective.

Always consult the skin tone chart in the image to determine if this device is suitable for you.

7. MAINTENANCE AND STORAGE

- **Cleaning:** After each use, ensure the device is unplugged and cooled down. Gently wipe the light emission window with a soft, dry, lint-free cloth. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the device in its original packaging or a protective case in a cool, dry place, away from direct sunlight and moisture. Keep it out of reach of children.
- **Cord Care:** Do not wrap the power cord tightly around the device. Avoid kinking or stressing the cord.

8. TROUBLESHOOTING

If you encounter any issues with your JOVS IPL device, please consider the following general troubleshooting steps:

- **Device Not Turning On:** Ensure the power adapter is securely connected to both the device and a working power outlet. Check if the Energy Level Switching Ring is set to a level other than 0.
- **No Flash:** Verify that the light emission window is in full and firm contact with your skin. The device has a skin contact sensor for safety. Ensure the device is not in standby mode.
- **Discomfort During Use:** Reduce the energy intensity level. Ensure your skin is properly prepared (shaved, clean, dry). The Sapphire Ice-Cooling system should minimize discomfort, but individual sensitivity varies.
- **Inconsistent Results:** Ensure you are following the recommended treatment schedule and applying the device consistently to all areas. Verify your skin tone and hair color are suitable for IPL treatment.

For more detailed troubleshooting or persistent issues, please refer to the full user manual (PDF link below) or contact JOVS customer support.

9. SPECIFICATIONS

Feature	Detail
Product Name	JOVS Venus Pro II Hair Removal
Model Number	JR5-E
Weight	318g
Energy Output	2-6 J/cm²
Cool Sapphire Area	12.1 cm²
Flash Interval	0.7s (fastest)
Rated Power	80 W
Wavelength	590 - 1200 nm
Number of Flashes	Infinite
System Technology	OPT technology

Feature	Detail
Input	100 - 240V~50/60Hz, 1.5A
Output	24V = 2.5A
Function	Hair Removal, Skincare

10. WARRANTY AND SUPPORT

The JOVS IPL Hair Removal Device Model JR5-E comes with a**24-month warranty** from the date of purchase. Please retain your proof of purchase for warranty claims.

For technical support, warranty service, or any questions regarding the use of your device, please contact JOVS customer support through the official JOVS website or your retailer.

You can also refer to the official PDF user manual for more detailed information[Download PDF User Manual](#)

© 2025 JOVS. All rights reserved.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page’s manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question

