

Waterpik WF-17CD017-1

Waterpik Cordless Slide Professional Water Flosser (Model WF-17CD017-1) Instruction Manual

Model: WF-17CD017-1 | Brand: Waterpik

1. INTRODUCTION

The Waterpik Cordless Slide Professional Water Flosser is designed for convenient and effective oral hygiene, offering a compact and portable solution for home use and travel. Its collapsible design allows for easy storage, reducing its size by half. This model features a removable 6-ounce reservoir, 3 pressure settings, and is waterproof for shower use. It is equipped with a rechargeable lithium-ion battery providing up to 4 weeks of use on a single charge.

The Waterpik Cordless Slide Professional is clinically proven to remove up to 99.9% of plaque bacteria from treated areas and is accepted by the American Dental Association (ADA) for safety and effectiveness.

2. WHAT'S INCLUDED

Your Waterpik Cordless Slide Professional Water Flosser package includes the following items:

- Waterpik Cordless Slide Professional Water Flosser Unit
- 4 Water Flossing Tips (various types for different dental needs)
- USB-C to USB-A Charging Cable
- Microfiber Travel Bag

Included in the box



Figure 1: Waterpik Cordless Slide Professional Water Flosser with 4 tips, USB charging cable, and travel bag.

3. SETUP

3.1 Initial Charging

Before first use, charge the unit for 4 hours. The rechargeable lithium-ion battery provides up to 4 weeks of use from a single charge. Use the included USB-C to USB-A cable to charge via a laptop, car, or wall adapter (adapter not included).

Rechargeable battery

Up to 4 weeks of use
LED recharge indicator



Figure 2: The Waterpik Cordless Slide charging via USB-C cable.

3.2 Attaching a Tip

Select the desired water flossing tip. Insert the tip into the center of the handle until it clicks into place. To remove a tip, press the eject button on the side of the handle.

3.3 Filling the Reservoir

The unit features a removable 6-ounce reservoir. To fill, open the flip-top cap on the reservoir and fill with warm water. The reservoir is top-rack dishwasher safe for easy cleaning.

waterpik® waterflosser™

CORDLESS SLIDE

- Convenient collapsible design saves space for storage and/or travel
- Waterproof, can be used in the shower
- Ultra quiet technology

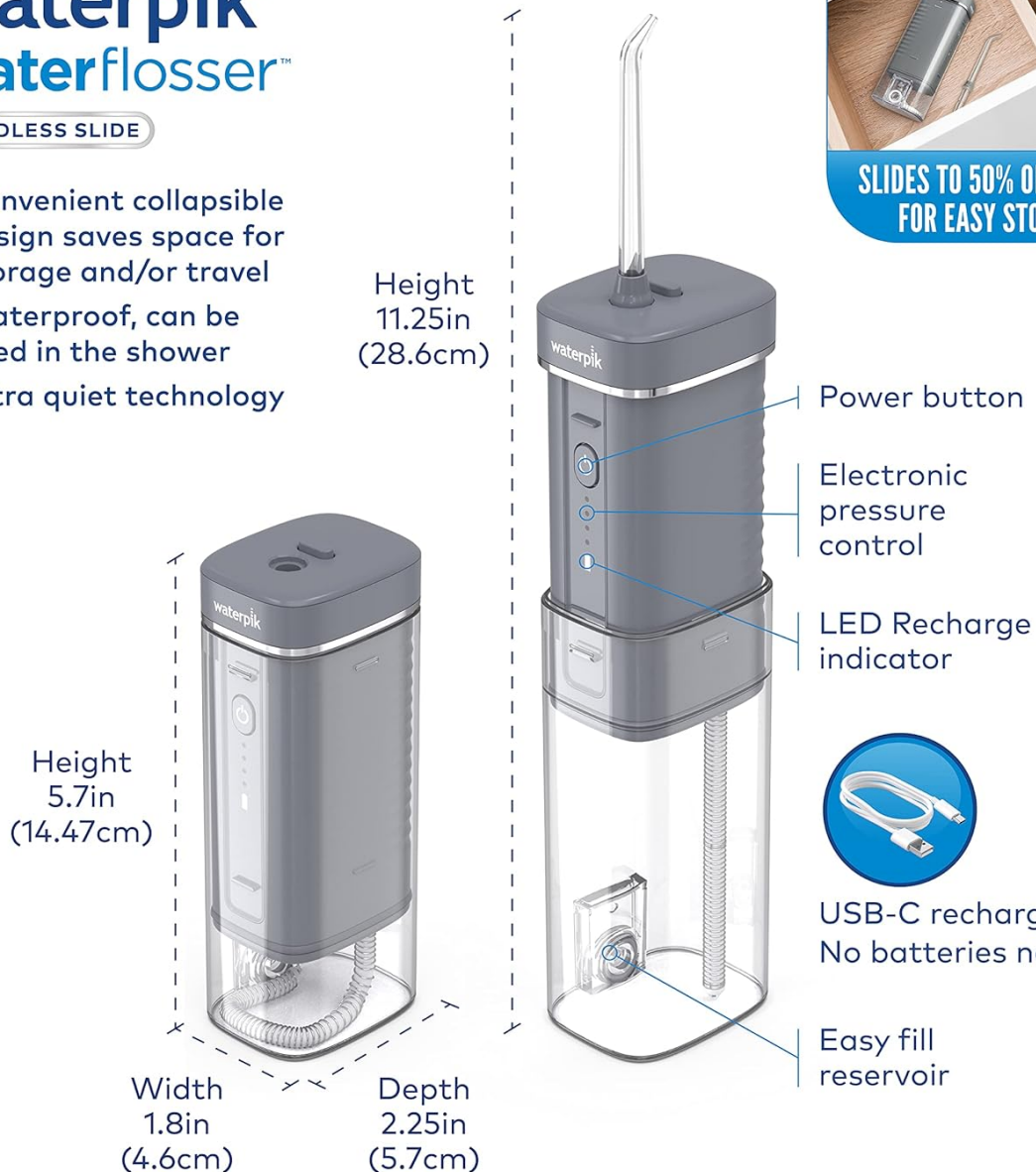


Figure 3: Overview of the Waterpik Cordless Slide, highlighting the reservoir and controls.

4. OPERATING INSTRUCTIONS

4.1 General Use Directions

1. Fill the reservoir with warm water and click the tip into the handle.
2. Lean over a sink and place the tip in your mouth.
3. Turn the unit ON. For first-time users, start at the lowest pressure setting and gradually increase over time.
4. Aim the tip along the gumline at a 90-degree angle, pausing briefly between teeth.
5. Close your lips enough to prevent splashing while allowing water to flow from your mouth into the sink.
6. Floss inside and outside of upper and lower teeth.
7. When finished, empty any remaining water from the reservoir.

Your browser does not support the video tag.

Video 1: Official Waterpik guide on how to use the Cordless Slide Water Flosser.

4.2 Pressure Settings

The Waterpik Cordless Slide Professional offers 3 pressure settings (50-80 PSI) to customize your cleaning experience. Adjust the pressure according to your comfort and dental needs.

3 pressure settings

For gentle to maximum clean



Figure 4: The control panel indicating the three available pressure settings.

5. MAINTENANCE

5.1 Cleaning the Reservoir

The 6-ounce reservoir is removable and top-rack dishwasher safe for easy and thorough cleaning. Ensure the reservoir is completely dry before reattaching.

5.2 General Cleaning

Wipe the main unit with a soft, damp cloth. Do not immerse the main unit in water. Periodically clean the tips with warm, soapy water.

6. TROUBLESHOOTING

- If you encounter any issues with your Waterpik Cordless Slide Professional Water Flosser, please refer to the following general advice:
- **No Power/Not Charging:** Ensure the USB-C cable is securely connected to both the unit and a power source. Verify the power source is active. If the unit has been exposed to water, ensure it is completely dry before attempting to charge.
 - **Weak Water Stream:** Check the water level in the reservoir. Ensure the tip is correctly inserted and not clogged. Clean the reservoir and tips thoroughly.
 - **Leaking:** Ensure the reservoir is properly seated and the flip-top cap is securely closed. Check for any cracks or damage to the reservoir or tips.

For persistent issues, please contact Waterpik customer support.

7. SPECIFICATIONS

Feature	Detail
Product Dimensions	2.25 x 1.8 x 5.7 inches; 8 ounces
Item Model Number	WF-17CD017-1
Batteries	1 Lithium Ion batteries required. (included)
UPC	073950321211
Manufacturer	Waterpik
Power Source	Battery Powered
Special Feature	Rechargeable and Portable, Travel Accessories, Waterproof for Shower Use
Product Benefits	ADA Accepted for Safety and Effectiveness, Gum Health, Removes Plaque
Unit Count	1 Count
GTIN	073950321211

8. SAFETY INFORMATION

Please observe the following safety guidelines:

- Do not immerse the main unit in water.
- Keep away from children.
- Do not use with certain medical conditions or while pregnant without consulting a healthcare professional.
- Use appropriate voltage adapters when traveling internationally.
- Always follow the instructions provided in this user manual for safe and proper use.

9. WARRANTY AND SUPPORT

The Waterpik Cordless Slide Professional Water Flosser is backed by a 2-year limited manufacturer's warranty. For detailed warranty information, please refer to the product manual included in the packaging.

For any product questions or support needs, Waterpik's Colorado-based support team is available to assist you.



© 2025 Waterpik. All rights reserved.

Related Documents - WF-17CD017-1

	<p>Waterpik Cordless Slide Water Flosser WF-16/WF-17 Quick Start Guide</p> <p>A quick start guide for the Waterpik Cordless Slide Water Flosser, models WF-16 and WF-17. Includes instructions on charging, use, cleaning, and troubleshooting.</p>
	<p>Waterpik Water Flosser Models WF-16 WF-17 User Manual and Guide</p> <p>Comprehensive guide for Waterpik Water Flosser models WF-16 and WF-17, covering setup, usage, cleaning, maintenance, safety, and warranty information.</p>
	<p>Waterpik Water Flosser WF-11 & WF-12 User Manual</p> <p>Comprehensive user manual for the Waterpik Water Flosser models WF-11 and WF-12, providing detailed instructions on setup, usage, cleaning, troubleshooting, and warranty information. Enhance your oral hygiene with Waterpik.</p>
	<p>Waterpik WF-02 Water Flosser: User Manual, Features, and Safety Guide</p> <p>Comprehensive user manual and guide for the Waterpik WF-02 Water Flosser. Learn about product features, usage instructions, cleaning, maintenance, battery care, and warranty information.</p>
	<p>Waterpik Water Flosser WF-02 User Manual and Guide</p> <p>Official user manual for the Waterpik Water Flosser Model WF-02. Includes setup, operation, cleaning, troubleshooting, and warranty information for optimal oral hygiene. Learn how to use and maintain your Waterpik water flosser for effective dental care.</p>
	<p>Waterpik Kids Water Flosser: Quick Start Guide and Registration</p> <p>Welcome to the Waterpik Kids Water Flosser. This guide provides essential steps for getting started, priming the unit, water flossing techniques, and details on registration, reviews, and customer support.</p>