

## YteaTh Y01

# YteaTh Shcall Y01 Sonic Electric Toothbrush User Manual

Model: Y01

## IMPORTANT SAFETY WARNINGS

---

Please read all instructions carefully before using the YteaTh Shcall Y01 Sonic Electric Toothbrush. Keep this manual for future reference.

- This product is not suitable for children under 3 years of age.
- Children should use this product under adult supervision.
- Replace the brush head every 3 months, or sooner if the bristles become worn or discolored (e.g., blue indicator bristles fade).
- Do not immerse the charging base in water. The toothbrush handle is IPX7 waterproof, but the charging port cover must be securely closed during use and cleaning.
- Consult a dentist before use if you have severe gum disease, oral surgery, or other oral health concerns.

## PRODUCT OVERVIEW

---

The YteaTh Shcall Y01 Sonic Electric Toothbrush is designed to provide an effective and gentle cleaning experience. It features advanced sonic technology, multiple brushing modes, and a long-lasting battery.



Image: The YteaTh Shcall Y01 Sonic Electric Toothbrush in light green, accompanied by six replacement brush heads. The toothbrush handle is shown with water splashing around it, highlighting its waterproof design.

### Key Features:

- **Advanced Sonic Motor:** Delivers 42,000 vibrations per minute for efficient plaque removal.
- **Long Battery Life:** Up to 180 days of use on a single full charge (based on normal usage).
- **5 Brushing Modes:** Customizable cleaning for various oral care needs.
- **Smart Timer:** 2-minute auto-timer with 30-second quadrant reminders.
- **IPX7 Waterproof:** Suitable for use in the shower.
- **USB Type-C Charging:** Convenient and fast charging.
- **6 Replacement Brush Heads:** Includes high-quality DuPont bristles for extended use.

### SETUP AND FIRST USE

## 1. Charging the Toothbrush:

1. Before first use, fully charge the toothbrush.
2. Locate the USB Type-C charging port at the bottom of the toothbrush handle.
3. Open the charging port cover.
4. Connect the provided USB Type-C cable to the toothbrush and a compatible USB power adapter (not included).
5. The indicator light will show charging status. A full charge can take several hours.
6. Once fully charged, securely close the charging port cover to maintain its IPX7 waterproof rating.



Image: Close-up of the YteaTh Shcall Y01 Sonic Electric Toothbrush's USB Type-C charging port with its protective cover open, showing a charging cable connected. Text indicates 1800mAh battery and 180 days of use.

## 2. Attaching the Brush Head:

- Align the brush head with the metal shaft on the toothbrush handle.
- Push the brush head down firmly until it clicks into place. Ensure it is securely attached.

# Whiten Your Teeth & Improve Gum Health

## Clean

Clean your teeth daily

## White

1 times per week

## Sensitive

Beneficial for fresh user or sensitive teeth

## Polish

Quick whitening in an emergency

## Massage

Improve gum health



Image: The YteaTh Shcall Y01 Sonic Electric Toothbrush handle with a brush head attached, ready for use.

## OPERATING INSTRUCTIONS

### 1. Selecting a Brushing Mode:

The YteaTh Shcall Y01 offers 5 distinct brushing modes to suit various oral care needs. Press the power button repeatedly to cycle through the modes:

- **Clean:** For daily dental care and thorough cleaning.
- **White:** Recommended for use once per week to help whiten teeth.

- **Sensitive:** Gentle cleaning for new users or individuals with sensitive teeth and gums.
- **Polish:** Designed for quick whitening in specific situations.
- **Massage:** Improves gum health through gentle stimulation.



Image: Diagram illustrating the five brushing modes (Clean, White, Sensitive, Polish, Massage) available on the YteaTh Shcall Y01 Sonic Electric Toothbrush, with descriptions of their intended use.

## 2. Brushing Your Teeth:

1. Wet the brush head and apply a small amount of toothpaste.
2. Place the brush head against your teeth at a slight angle towards the gum line.
3. Press the power button to turn on the toothbrush and select your desired mode.
4. Gently move the brush head slowly across your teeth, allowing the sonic vibrations to do the work. Do not scrub aggressively.

- The toothbrush features a **Smart Timer** that pauses briefly every 30 seconds to indicate when to move to a different quadrant of your mouth (upper right, upper left, lower right, lower left).
- After 2 minutes, the toothbrush will automatically turn off, signaling the completion of the recommended brushing time.

## Newly Upgraded Functional Configuration



VS



Electric Motor  
**More stable and stronger**

Electric Motor  
**Average**

Battery Life  
**180 days**

Battery Life  
**50 days**

Craft  
**Water flow curve  
texture design:  
comfortable  
and non-slip**

Craft  
**Run of the mill**

Decibel  
**<55dB**

Decibel  
**65dB**

Image: The YteaTh Shcall Y01 Sonic Electric Toothbrush with an animated circular arrow around it, representing the smart vibration timer that notifies users every 30 seconds to change brushing quadrants.

### 3. Turning Off the Toothbrush:

To turn off the toothbrush before the 2-minute cycle is complete, press and hold the power button for at least three seconds.

## MAINTENANCE AND CARE

## 1. Cleaning the Toothbrush:

- After each use, rinse the brush head thoroughly under running water.
- Remove the brush head and rinse the metal shaft area on the handle.
- Wipe the handle with a damp cloth. Ensure the charging port cover is closed during cleaning.
- Allow all parts to air dry before storing.

## 2. Replacing Brush Heads:

For optimal hygiene and brushing performance, replace the brush head every 3 months, or sooner if the bristles appear frayed or the blue indicator bristles fade halfway. The YteaTh Shcall Y01 comes with 6 replacement brush heads, providing an 18-month supply.

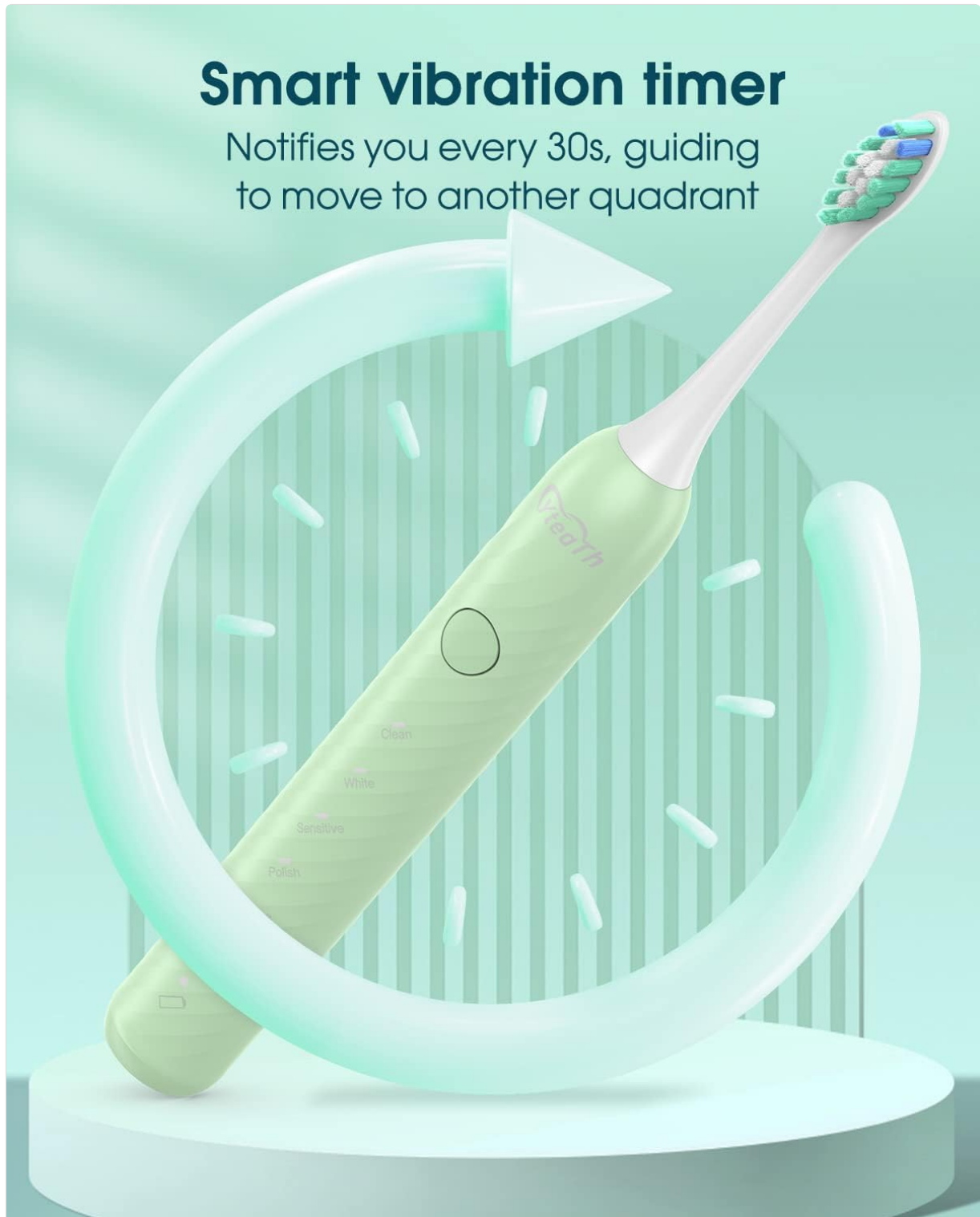


Image: A collection of six replacement brush heads for the YteaTh Shcall Y01 Sonic Electric Toothbrush. An inset shows a close-up of bristles with blue indicator bristles, and text advises replacement every three months or when blue bristles wear out.

## TROUBLESHOOTING

| Problem                                 | Possible Cause                                   | Solution   |
|---|--|--|
| Toothbrush does not turn on.            | Battery is depleted.                             | Charge the toothbrush fully using the USB Type-C cable.  |
| Weak vibrations or reduced performance. | Low battery or worn brush head.                  | Charge the toothbrush. Replace the brush head if bristles are worn.  |
| Toothbrush stops before 2 minutes.      | Power button held down too long, or battery low. | Ensure you are not holding the power button for more than 3 seconds. Charge the toothbrush.                |
| Charging port not working.              | Debris in port or cable issue.                   | Inspect the charging port for debris and clean gently. Try a different USB Type-C cable and power adapter. |

## SPECIFICATIONS

- **Model:** Y01
- **Brand:** YteaTh
- **Product Name:** Sonic Electric Toothbrush
- **Vibration Frequency:** 42,000 VPM (Vibrations Per Minute)
- **Brushing Modes:** 5 (Clean, White, Sensitive, Polish, Massage)
- **Battery Type:** Lithium-ion
- **Battery Life:** Up to 180 days (approx. 6 months) on a single charge
- **Charging:** USB Type-C
- **Waterproof Rating:** IPX7
- **Timer:** 2-minute auto-timer with 30-second quadrant reminder
- **Bristle Material:** DuPont
- **Dimensions:** 20 x 8.7 x 3.2 cm (Package)
- **Weight:** 250 grams (Package)
- **Country of Origin:** China

## WARRANTY AND SUPPORT

For warranty information or technical support, please refer to the retailer's return policy or contact YteaTh customer service through the platform where the product was purchased. Please retain your proof of purchase for any warranty claims.

**Note:** The product description mentions a 30-day return policy through Amazon. Specific manufacturer warranty details are not provided in the product data.