

AMINZER SKB2008

AMINZER SKB2008 IPL Hair Removal Device User Manual

Model: SKB2008

1. INTRODUCTION

The AMINZER SKB2008 IPL Hair Removal Device utilizes advanced Intense Pulsed Light (IPL) technology to provide a safe and effective solution for permanent hair reduction at home. This device incorporates a unique double ice cooling system and 600nm red light wave technology to target hair follicles, preventing regrowth and ensuring a comfortable, painless experience. With 5 adjustable energy levels and multiple operating modes, it is designed for various skin tones and body areas.

2. SAFETY INFORMATION AND WARNINGS

Before using the AMINZER SKB2008 device, please read all safety instructions carefully. Improper use can lead to skin irritation or injury. This device is intended for hair removal only.

2.1 Important Safety Precautions


- Always perform a patch test on a small skin area 24 hours before full treatment to check for adverse reactions.
- Do not use on irritated, broken, or infected skin.
- Avoid direct eye exposure to the pulsed light. Always wear protective eyewear (not included, but recommended) or ensure the device is pressed firmly against the skin to activate the safety sensor.
- Keep the device dry. Do not use it in wet environments or near water.
- Keep out of reach of children.
- Consult a doctor before use if you have any underlying medical conditions or are taking medications that may increase skin sensitivity.

2.2 Contraindications and Skin Suitability


The AMINZER SKB2008 device is not suitable for everyone. Refer to the chart below to determine if this device is appropriate for your skin tone and hair color, and to identify conditions where use is contraindicated.

Check if this Device is Right for you?


Skin tone & body hair type chart.




Tattoos
















Skin Disease



Sunburn



Pregnant Women

HAIR COLOUR		SKIN TONE							
									
White		✓	✓	✓	✓	✗	✗	✗	
Beige		✓	✓	✓	✓	✗	✗	✗	
Cream Color		✓	✓	✓	✓	✗	✗	✗	
Light Brown		✓	✓	✓	✓	✗	✗	✗	
Dark brown		✗	✗	✗	✗	✗	✗	✗	
Brownishblack		✗	✗	✗	✗	✗	✗	✗	

✓ It's match! ✗ Try a different hair removal method

Image: Skin Tone and Hair Color Suitability Chart. This chart helps determine if the device is suitable for your specific skin tone and hair color, and lists conditions such as tattoos, skin disease, sunburn, and pregnancy as contraindications for use.

- **Do not use if:** You have tattoos in the treatment area, skin diseases, sunburn, or are pregnant.
- The device is most effective on lighter skin tones with darker hair. It may not be effective on very dark skin tones or very light (blonde, red, grey) hair.

3. PACKAGE CONTENTS

Upon unpacking, please ensure all items are present and in good condition:

- AMINZER SKB2008 IPL Hair Removal Device
- Power Adapter
- User Manual

4. DEVICE OVERVIEW

The AMINZER SKB2008 device features a sleek design with intuitive controls. Key components include the light emission window, skin contact sensor, power button, mode selection button, and energy level indicators.

Innovative Dual Ice Cooling Function

Lower Skin Surface Temperature for Painless Hair Removal

10° ❄️

The image illustrates the device's cooling mechanism. The top part shows a cutaway of the device with a fan and a cooling system. Below this, two diagrams compare skin temperatures. The left diagram, labeled 'Other', shows a device emitting a red light beam, with the skin temperature indicated as 60°C. The right diagram, labeled 'Our', shows the device emitting a blue light beam with snowflake icons, indicating a cooling effect that reduces the skin temperature to 10°C - 20°C.

For the non-ice-cooling device, the temperature of skin is 60°C - 100°C

Our device with cooling effect will cool the skin temperature to 10°C - 20°C

Image: AMINZER SKB2008 IPL Hair Removal Device. This image displays the overall design of the device, highlighting its ergonomic shape and color scheme.

5. SETUP

5.1 Initial Preparation

1. **Unpack:** Carefully remove the device and all accessories from the packaging.
2. **Inspect:** Check the device for any visible damage. Do not use if damaged.
3. **Connect Power:** Plug the power adapter into the device and then into a suitable electrical outlet.

5.2 Pre-Treatment Skin Preparation

- **Shave:** Shave the area to be treated 12-24 hours before your session. Do not wax, pluck, or use

depilatory creams, as these remove the hair root which is necessary for IPL to be effective.

- **Cleanse:** Ensure the skin is clean, dry, and free from any lotions, creams, deodorants, or makeup.

6. OPERATING INSTRUCTIONS

6.1 Powering On and Off

- **Power On:** Press and hold the power button for a few seconds until the indicator lights illuminate.
- **Power Off:** Press and hold the power button again until the lights turn off. The device will also automatically power off after a period of inactivity.

6.2 Selecting Energy Levels and Modes

The device offers 5 energy levels and 3 operating modes to customize your treatment.

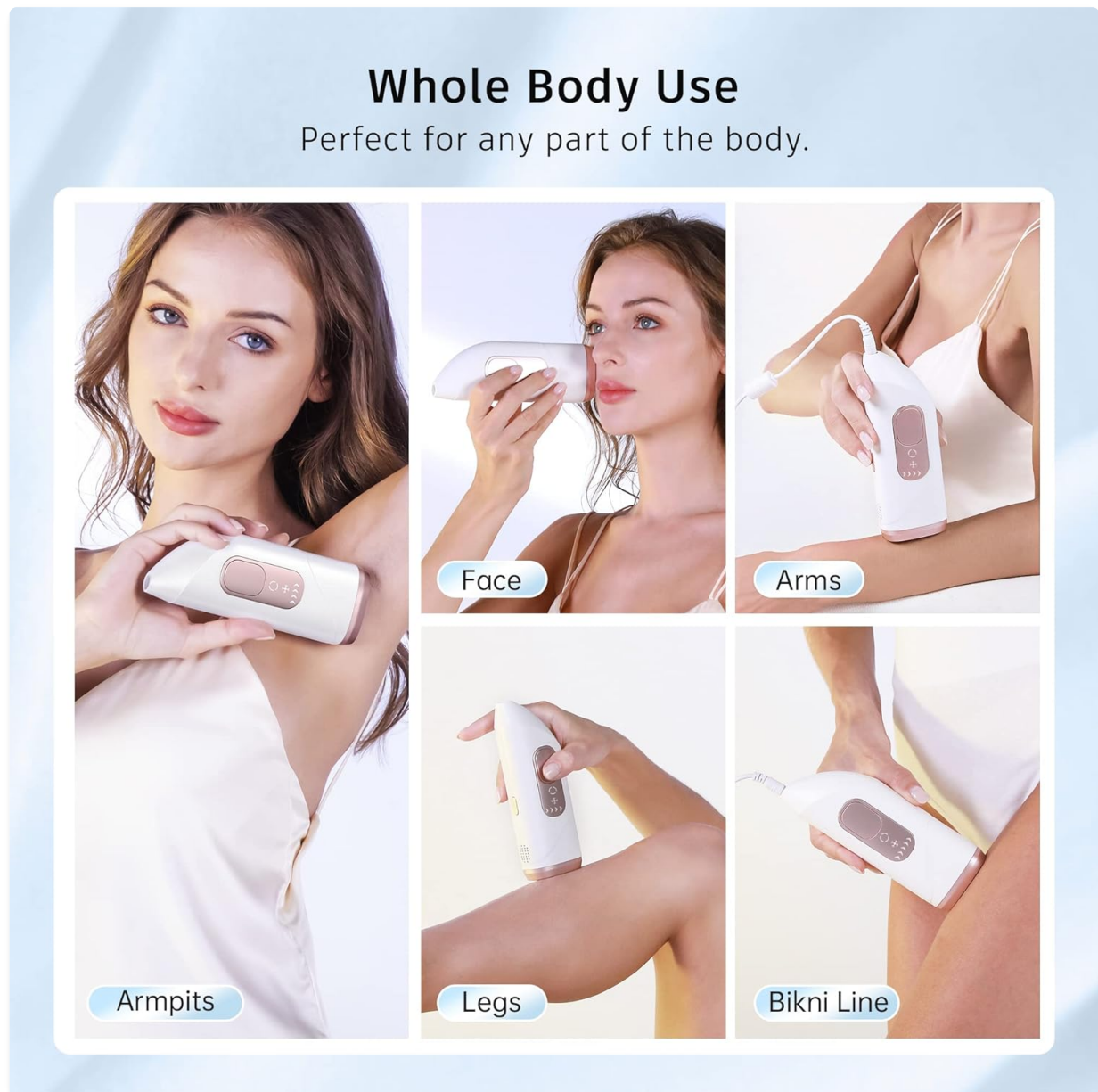


Image: Energy Levels and Operating Modes. This image illustrates the five energy intensity levels (Level 1 to Level 5) and the two primary operating modes: Auto Mode for large areas and Manual Mode for precise treatment.

- **Energy Levels:** Press the power button briefly to cycle through the 5 energy levels. Start with Level 1 and gradually increase as your skin tolerates it. Higher levels provide more effective results.
- **Auto Mode:** Ideal for large areas like legs, arms, and back. Press and hold the flash button for

continuous flashes when the device is in full contact with the skin.

- **Manual Mode:** Suitable for smaller, sensitive, or precise areas such as the bikini line, armpits, and face. Press the flash button once for each individual flash.
- **Ice Applying Mode:** This mode activates the double ice cooling system, which cools the skin to approximately 10°C, providing a comfortable and painless experience during and after treatment. This function is active during IPL flashes.



Image: Innovative Dual Ice Cooling Function. This diagram demonstrates how the device's cooling system reduces skin surface temperature to 10°C, ensuring a painless hair removal experience compared to devices without cooling.

6.3 Application Technique

1. Ensure the device's light emission window is in full contact with your skin. The device will not flash unless proper contact is made.
2. For Manual Mode, press the flash button once for each area. Move to the next adjacent area, ensuring slight overlap to avoid missed spots.
3. For Auto Mode, hold the flash button down and glide the device slowly over the skin, maintaining continuous contact.
4. Treat each area 1-2 times per session. Avoid treating the same spot multiple times in one session to prevent overheating.

6.4 Treatment Schedule

Consistent use according to the recommended schedule is crucial for optimal results.

Whole Body Use

Perfect for any part of the body.



The image shows two side-by-side photographs of a woman's armpit. The left photo, labeled 'Before', shows the armpit with visible hair. The right photo, labeled 'After', shows the same armpit with significantly reduced hair. Below each photo is a diagram illustrating the light wavelength used for treatment.

Before

Others
530nm light waves-unable to remove hair completely.

530nm

After

Ours
600nm red light waves- Effectively removes hair from the root.

600nm

Image: IPL Treatment Schedule. This illustration outlines the recommended frequency of use for the first 8 weeks and beyond to achieve effective hair reduction.

- **Weeks 1-4:** Use 2 times per week.
- **Weeks 5-8:** Use 1 time per week.
- **After 8 Weeks:** Use 1-2 times per month for maintenance, or as needed.

6.5 Areas of Use

The AMINZER SKB2008 is designed for use on various body parts.

Choose the Suitable Mode & Level as You Like



2 Modes

Auto Mode

Suitable for large body areas



Manual Mode

Suitable for smaller and sensitive areas



Image: Whole Body Use and Effectiveness. This image demonstrates the device's suitability for various body areas and illustrates the superior effectiveness of 600nm red light waves in removing hair from the root compared to 530nm light waves.

- **Face:** Upper lip, chin, sideburns (avoiding eyes).
- **Body:** Arms, armpits, legs, bikini line, back, chest.
- Always ensure the skin is clean and dry before treatment.

7. MAINTENANCE

7.1 Cleaning the Device

- Always unplug the device before cleaning.
- Use a soft, dry cloth to wipe the light emission window and the device body.
- Do not use abrasive cleaners, alcohol, or harsh chemicals.
- Ensure no water or liquids enter the device.

7.2 Storage

- Store the device in a cool, dry place, away from direct sunlight and moisture.

- Keep the device in its original packaging or a protective case to prevent damage.

8. TROUBLESHOOTING

If you encounter any issues with your AMINZER SKB2008 device, refer to the following common problems and solutions:

Problem	Possible Cause	Solution
Device does not turn on	Not properly plugged in; Power button not pressed correctly; Power adapter fault.	Ensure power adapter is securely connected to both device and outlet. Press and hold the power button for a few seconds. Try a different outlet.
Device does not flash	Light emission window not in full contact with skin; Skin tone sensor not detecting suitable skin; Device malfunction.	Ensure the light emission window is fully pressed against the skin. Check skin tone suitability (refer to Section 2.2). If problem persists, contact customer support.
Skin irritation or discomfort	Energy level too high; Skin not properly prepared; Device used on contraindicated skin.	Reduce energy level. Ensure skin is clean and shaved. Review contraindications in Section 2.2. Use Ice Applying Mode for comfort. Apply a soothing gel after treatment if needed.
No noticeable hair reduction	Inconsistent use; Energy level too low; Unsuitable hair/skin type.	Follow the recommended treatment schedule (Section 6.4). Increase energy level if comfortable. Review skin and hair suitability (Section 2.2).

If the problem persists after attempting these solutions, please contact AMINZER customer support for further assistance.

9. SPECIFICATIONS

- **Model:** SKB2008
- **Brand:** AMINZER
- **Technology:** IPL (Intense Pulsed Light)
- **Light Wave:** 600nm Red Light Wave
- **Flashes:** 999,990+
- **Energy Levels:** 5
- **Cooling System:** Double Ice Cooling
- **Modes:** Auto Mode, Manual Mode, Ice Applying Mode
- **Units:** 1.0 Count

10. WARRANTY AND SUPPORT

Your AMINZER SKB2008 IPL Hair Removal Device comes with a manufacturer's warranty. Please refer to the warranty card included in your product packaging for specific terms and conditions, including warranty

duration and coverage.

For technical support, troubleshooting assistance, or warranty claims, please contact AMINZER customer service through the contact information provided in your product documentation or on the official AMINZER website. Please have your model number (SKB2008) and purchase details ready when contacting support.