

Ulike UI06

Ulike Air 3 Ice-Cooling IPL Hair Removal Device Instruction Manual

Model: UI06 | Brand: Ulike

INTRODUCTION

The Ulike Air 3 Ice-Cooling IPL Hair Removal Device is designed to provide a nearly painless and long-lasting hair removal solution from the comfort of your home. Utilizing advanced IPL technology with a patented ice-cooling touch, it maintains a skin temperature of approximately 67°F (19.4°C) during treatment, ensuring comfort while effectively targeting hair follicles. This device is suitable for full-body use and offers noticeable hair reduction within weeks.

Key features include:

- **Nearly Painless Treatment:** Patented Ice-Cooling technology keeps skin cool during use.
- **Visible Results:** Users may experience noticeable hair reduction within 3-5 weeks.
- **Versatile Use:** Safe and effective for various body areas including legs, arms, bikini area, and face.
- **Multiple Modes:** Features Soft, Body, and Power modes for tailored treatment.
- **Auto-Glide Mode:** For fast and efficient full-body treatments.

SETUP

Before your first use, ensure all components are present and familiarize yourself with the device.

What's in the Box

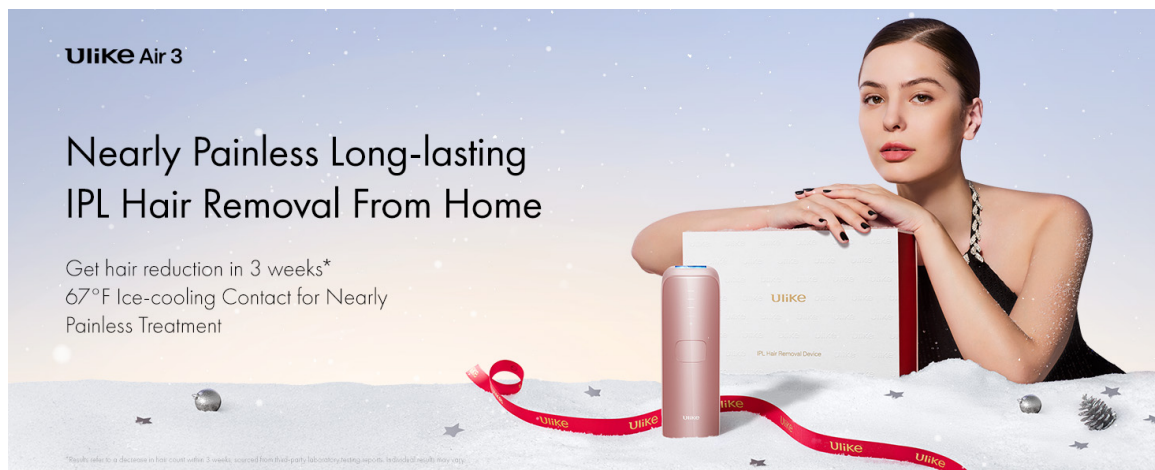


Image: The Ulike Air 3 device, product packaging, power cord, UV protective glasses, razor, and adapter.

Initial Preparation

1. **Unpack:** Carefully remove all items from the packaging. Verify that all parts listed in the 'What's in the Box' section are present.
2. **Connect Power:** Plug the power cord into the adapter, then connect the adapter to the device's power input port. Finally, plug the adapter into a power outlet. The device must remain connected to power during treatment.
3. **Wear Safety Glasses:** Always wear the provided UV protective glasses before turning on the device to protect your eyes from the intense light.
4. **Prepare Skin:** Ensure the treatment area is clean, dry, and freshly shaved. Do not wax, pluck, or use depilatory creams, as these methods remove the hair follicle, which is necessary for IPL treatment to be effective.



Image: Steps to connect the device to power and activate it.

OPERATING INSTRUCTIONS

Follow these steps for effective and safe hair removal.

Power On and Mode Selection

1. **Turn On:** Long-press the power button to activate the device.
2. **Select Mode:** The Ulike Air 3 offers three comfort modes:

- **Soft Mode:** Ideal for sensitive areas like the face and upper lip.
- **Body Mode:** Suitable for larger areas such as legs and arms.
- **Power Mode:** For stubborn hair in areas like underarms and the bikini line.

Press the mode button to cycle through the settings and choose the appropriate one for your treatment area.

4-7 Weeks*

After 7 weeks of consistent bi-weekly treatments, the hair follicles start to shrink and the hair root begins to fall out.



Image: Illustration of the three modes: Soft Mode for face/upper lip, Body Mode for legs/arms, and Power Mode for bikini/armspits.

Treatment Application

1. **Position Device:** Place the treatment window flat against your skin at a 90-degree angle. The device will detect skin contact.
2. **Auto-Glide Mode:** For continuous flashing, simply glide the device slowly over the treatment area. The device will automatically emit light pulses at regular intervals.
3. **Manual Mode:** For precise treatment, press the lighting button to emit a single light pulse. Lift the device, move to the next area, and repeat.
4. **Cooling Sensation:** The integrated Sapphire Ice-Cooling technology will provide a cool sensation on your skin throughout the treatment.

7-10 Weeks*

Most people will have achieved desired results by this stage, with skin becoming smooth and nearly hairless.



Image: Demonstrates Auto-Glide Mode for large areas and Manual Mode for targeted treatment.

Recommended Treatment Schedule

- **Weeks 1-3:** 3 treatments per week.
- **Weeks 4-6:** 2 treatments per week.
- **Weeks 7-10:** 1-2 times for long-lasting results.
- **Maintenance:** After 10 weeks, use as needed to maintain smooth skin.

3 Modes for Full Body Hair Removal

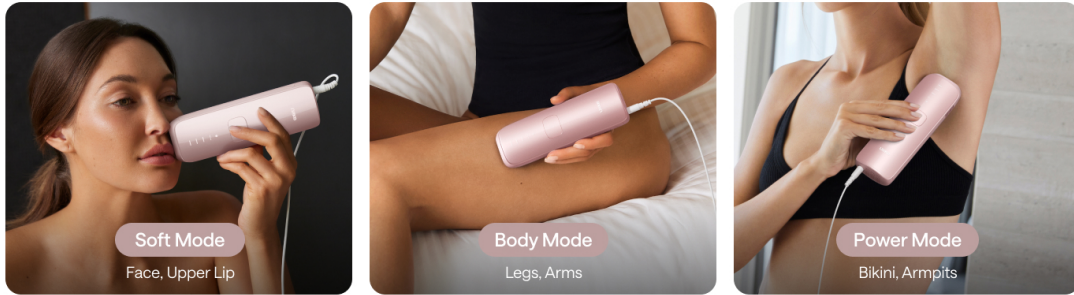


Image: Visual representation of hair reduction progress after 1-3 weeks.



Image: Visual representation of hair reduction progress after 4-7 weeks.



Image: Visual representation of hair reduction progress after 7-10 weeks.

SAFETY INFORMATION

Your safety is paramount. Please read and adhere to all safety guidelines.

- Before using the device, please read the user manual carefully, including all contraindications, warnings, and safety information.
- Not suitable for tattooed and over-tanned skin, wounded areas, and pregnant individuals.

- The device is a plug-in unit; ensure it is always connected to power during treatment.
- Always wear the provided UV protective glasses during operation.

Skin and Hair Tone Compatibility

The Ulike Air 3 IPL device is effective for most skin and hair types. Refer to the chart below to determine if it is suitable for you.

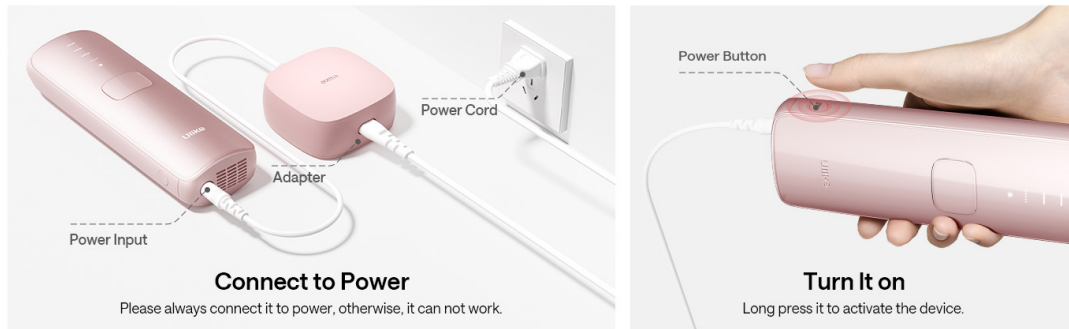


Image: Skin tone and hair color chart for IPL device compatibility.

MAINTENANCE

Proper maintenance ensures the longevity and optimal performance of your Ulike Air 3 device.

- **Cleaning:** After each use, ensure the device is unplugged and cool. Gently wipe the treatment window with a soft, dry cloth. Do not use abrasive cleaners or harsh chemicals.
- **Storage:** Store the device and its accessories in a cool, dry place, away from direct sunlight and moisture. Keep it in its original packaging or a protective case to prevent damage.

TROUBLESHOOTING

If you encounter any issues with your Ulike Air 3 device, refer to the following common troubleshooting tips:

- **Device Not Turning On:** Ensure the device is properly plugged into a power outlet and the power button is long-pressed. Check if the power adapter and cord are securely connected.
- **No Flash:** Verify that the treatment window is placed flat against your skin at a 90-degree angle. The device has a skin contact sensor that prevents flashing if not properly positioned.
- **Discomfort During Treatment:** Ensure you are using the appropriate mode for your skin sensitivity and body area. The ice-cooling feature should minimize discomfort. If discomfort persists, reduce the intensity level or consult a dermatologist.
- **Hair Regrowth Not Reducing:** Ensure you are following the recommended treatment schedule consistently. IPL requires regular use to be effective. Also, confirm your skin and hair tone are compatible with IPL treatment (refer to the compatibility chart).

For further assistance or issues not covered here, please refer to the full user manual or contact Ulike

customer support.

SPECIFICATIONS

Feature	Detail
Model Number	UI06
Product Dimensions	7.05 x 2.29 x 1.34 inches
Item Weight	9.6 ounces
UPC	749749575614
Manufacturer	Shenzhen Ulike Smart Electronics Co.,Ltd

WARRANTY AND SUPPORT

Ulike is committed to customer satisfaction and product quality.

- **Money-Back Guarantee:** Enjoy a 100-day money-back guarantee.
- **Sales Service:** The product comes with 5 years of sales service.
- **Contact Support:** For any inquiries or support needs, please refer to the contact information provided in your full user manual or visit the official Ulike website.

OFFICIAL PRODUCT VIDEO

Video: Official Ulike Air 3 IPL Hair Removal demonstration, showcasing its features and ease of use.