

VOYOR WF200

VOYOR Water Flosser WF200 Oral Irrigator Instruction Manual

Model: WF200

1. INTRODUCTION

Thank you for choosing the VOYOR Water Flosser WF200. This cordless, rechargeable dental water jet flosser is designed to provide powerful and effective oral cleansing, promoting healthy teeth and gums. Please read this manual thoroughly before use to ensure proper operation and maintenance.



Image: VOYOR Water Flosser WF200, showing the device, its packaging, and the included nozzle tips.

2. PRODUCT OVERVIEW

The VOYOR Water Flosser WF200 is a compact and portable oral irrigator featuring advanced technology for superior dental hygiene. Key components and design elements include:

- **On/Off Button:** Controls power to the device.
- **Mode Selection Button:** Cycles through available cleaning modes.
- **Working Modes:** Normal, Soft, and Pulse settings for varied intensity.
- **Charging Indicator Light:** Displays charging status.
- **Nozzle Release Button:** Allows for easy attachment and detachment of nozzles.
- **Interchangeable Nozzles:** Includes various tips for different cleaning needs.
- **Detachable Water Tank:** 250ml capacity for easy filling and cleaning.
- **Charging Port:** Located on the device for convenient recharging.

ELEGANT BODY & ERGONOMIC HANDLE DESIGN

Easy to handle, a great gift choice for your friends and loved ones.



Image: Detailed diagram of the VOYOR Water Flosser WF200, highlighting its ergonomic design and labeled functional parts.


3. FEATURES

The VOYOR Water Flosser WF200 is equipped with several features to enhance your oral hygiene routine:


- **Powerful & Effective Cleansing:** 1800t/m water pulse frequency removes 99.99% of food residue.
- **Multiple Working Modes:** 3 cleansing modes (Normal, Soft, Pulse) to suit different oral needs.
- **360° Rotating Nozzle:** Allows for easy cleaning at any angle, ideal for braces.
- **IPX7 Waterproof:** Safe for use in the shower or bathroom.
- **Rechargeable Battery:** 1500mAh lithium-ion battery provides up to 1 month of use on a 4-hour charge.
- **Compact & Portable:** Lightweight (260 grams) and slim design for easy handling and travel.

8 CORE ADVANTAGES

Offer You Professional Oral Cleansing To
Maintain a Healthy Gums & Brighter Teeth In Daily Use.

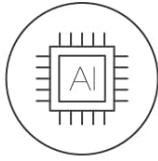

3 WORKING
MODES

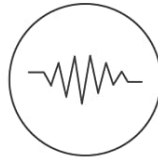

HIGH IMPULSE
FREQUENCY




2 MINUTES AUTO
TIMER



0.6MM PRECISE
WATER JET



360 DEGREE
CLEANSING

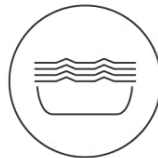

250ML BIG
WATERTANK



CORDLESS &
RECHARGEABLE



IPX7
WATERPROOF



Image: Infographic detailing the 8 core advantages of the VOYOR Water Flosser WF200, emphasizing its professional oral cleansing capabilities.

WITH 1800 TIMES/MINUTE HIGH FREQUENCY WATER PULSE

Deep cleanse every corner precisely

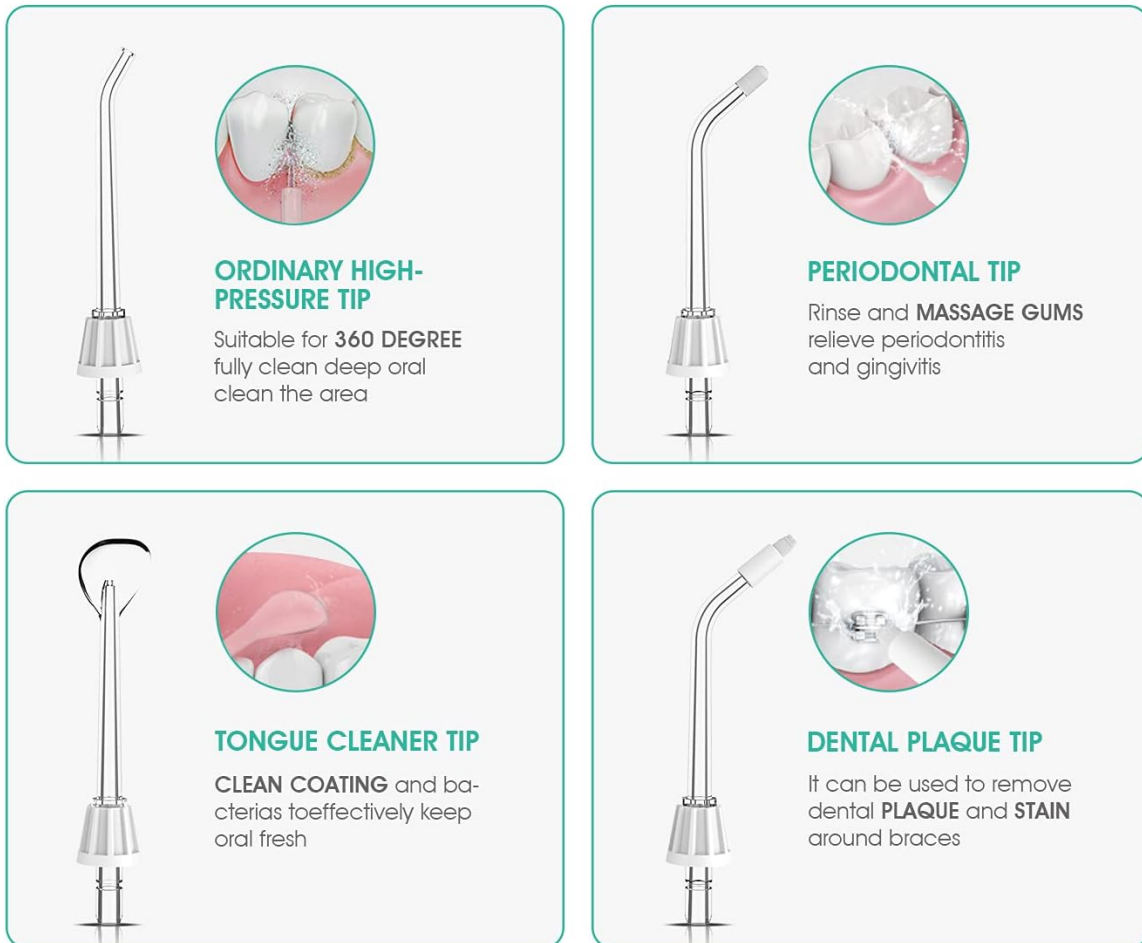


Image: Visual representation of the four specialized nozzle tips included with the VOYOR Water Flosser WF200, each designed for specific oral care needs.

4. SETUP

4.1 Charging the Device

Before first use, fully charge the water flosser. Connect the provided USB charging cable to the charging port on the device and to a compatible USB power adapter (not included). The charging indicator light will illuminate during charging and turn off or change color when fully charged. A full charge takes approximately 4 hours and provides up to 1 month of use.

4.2 Attaching the Nozzle

1. Select the desired nozzle tip for your oral care needs.
2. Insert the nozzle into the opening at the top of the device until you hear a click, ensuring it is securely in place.

4.3 Filling the Water Tank

1. To detach the water tank, twist the bottom part of the device counter-clockwise to the unlock

position and pull down.

2. Fill the tank with water (temperature lower than 40°C). You may also add mouthwash if desired.
3. Reattach the water tank by aligning it with the main body and twisting clockwise to the lock position. Alternatively, you can open the small flap on the side of the tank and fill directly without detaching.



Image: Instructions for filling the water tank, illustrating both detachment and direct filling methods.

5. OPERATING INSTRUCTIONS

1. Ensure the water tank is filled and a nozzle is securely attached.
2. Place the nozzle in your mouth before turning on the device to prevent water splash.
3. Press the **On/Off** button to start the water flosser.
4. Press the **Mode** button to select your preferred cleaning mode: Normal, Soft, or Pulse.
5. Direct the water jet flow at a 90-degree angle against your teeth and gums. Move the nozzle slowly along the gumline and between teeth for effective cleaning.
6. After use, press the **On/Off** button to turn off the device.

VOYOR
Enjoy A Healthy And Beautiful Life



Image: A user demonstrating the correct way to hold and position the VOYOR Water Flosser WF200 for oral cleaning.

6. MAINTENANCE

6.1 Cleaning the Water Tank

The detachable water tank allows for thorough cleaning to prevent limescale and bacterial buildup. After each use, empty any remaining water and rinse the tank. For deeper cleaning, detach the tank and wash it with mild soap and water. Ensure it is completely dry before reattaching.

6.2 Cleaning Nozzles

Remove the nozzle after each use by pressing the nozzle release button. Rinse the nozzle under running water. For optimal hygiene, replace nozzles every 3-6 months or as recommended by your dental professional.

6.3 Device Exterior

Wipe the exterior of the device with a soft, damp cloth. Do not immerse the main unit in water, despite its IPX7 rating, to protect internal components, especially the charging port when uncovered.



DETACHABLE WATER TANK

Helps to clean the limescale and bacteria inside of the water tank.

THE STORAGE BOX

Allows you to store the nozzles more hygienically and carry it with you.



IPX7 WATERPROOF

Completely washable and worry free to use in your shower.

Image: Visual guide to the detachable water tank, nozzle storage, and the IPX7 waterproof feature of the VOYOR Water Flosser WF200.

7. TROUBLESHOOTING

- **No Power:** Ensure the device is fully charged. If the charging indicator does not light up, check the charging cable and adapter connection.
- **Weak Water Pressure:** Check if the water tank is sufficiently filled. Ensure the nozzle is not clogged and is securely attached.
- **Water Leakage:** Verify that the water tank is correctly attached and locked into place. Ensure the charging port cover is properly sealed.
- **Device Stops During Use:** The device has a 2-minute auto-timer. If it stops before 2 minutes, check the battery level.

8. SPECIFICATIONS

- **Product Dimensions:** 4.02 x 3.43 x 8.78 inches
- **Weight:** 8.8 ounces
- **Item Model Number:** WF200
- **Batteries:** 1 Lithium Ion battery required (included)
- **Power Source:** Battery Powered
- **Water Pulse Frequency:** 1800t/m
- **Water Tank Capacity:** 250ml
- **Waterproof Rating:** IPX7

9. SAFETY INFORMATION

- Remember to put the nozzle in your mouth before turning on the device to prevent water splash.
- Direct the water jet flow at a 90-degree angle against your teeth and gums for effective oral cleansing.
- Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

10. WARRANTY & SUPPORT

For warranty information or product support, please refer to the contact details provided with your purchase or visit the official VOYOR website. Keep your purchase receipt for any warranty claims.