



[Manuals.plus](#) /

> [Breco](#) /

> Breco N5 Mini Neck Massager User Manual

Breco Neck N5 MINI

Breco N5 Mini Neck Massager User Manual

Model: Neck N5 MINI | Brand: Breco

1. INTRODUCTION

Thank you for choosing the Breco N5 Mini Neck Massager with Heat. This manual provides essential information for the safe and effective use of your device. Please read it thoroughly before operation and retain it for future reference.



Figure 1: Breo N5 Mini Neck Massager

2. SAFETY INFORMATION

To ensure safe operation, please observe the following warnings:

- Do not use near water or in a damp environment.
- Keep away from children.
- Do not use if you have a pacemaker or other medical implant. Consult a physician before use if you have any pre-existing medical conditions.
- Do not cover the device during use to prevent overheating.
- Do not use for more than 15 minutes at a time. Allow the device to cool down between uses.
- Do not bend or twist the power cord excessively.
- Do not force the massage heads against your body.
- Use only the included power adapter for charging.

3. PACKAGE CONTENTS

Verify that all items are present in the package:

- Breo N5 Mini Neck Massager
- USB Charging Cable
- User Manual (this document)



Figure 2: Package Contents

4. PRODUCT FEATURES

- **Patented 5D Massage Heads:** Features a "4+2" 5D massage head configuration designed to cover larger areas from the neck to the trapezius and shoulder, providing a comprehensive massage.
- **Ergonomic Design:** Engineered to mimic human hand massage techniques, targeting neck and shoulder muscles for effective relief.
- **Wearable Backpack-wearing Design:** Allows for hands-free operation, enabling use during daily activities.
- **Integrated Heat Function:** Provides constant heat compress at approximately 107.6°F (42°C) for enhanced muscle relaxation.
- **Long-lasting Battery:** Equipped with a high-quality lithium battery, offering extended use on a single charge.
- **Multiple Massage Modes:** Offers two distinct massage modes: Deep Massage and Soothing Massage, with adjustable intensities.

6 SOFT SILICONE MASSAGE HEADS

Gentle touch to cervical & shoulders

Upper Massage Heads *2

for neck muscle relaxation

Mid Massage Heads *2 +

Bottom Massage Heads*2

for shoulder muscle relaxation



Figure 3: 5D Massage Head Configuration

CONSTANT HEAT COMPRESS

107.6°F heat for fatigue relief



Figure 4: Constant Heat Compress

FREE YOUR ARMS WITH BACKPACK-WEARING DESIGN

Non-Limiting, Hands Free, Enjoyable Massage

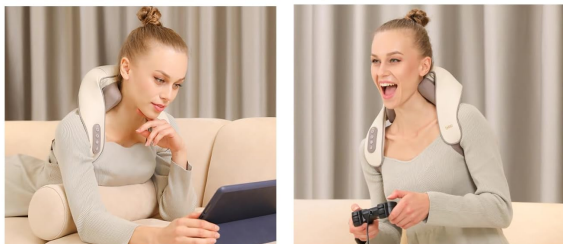


Figure 5: Hands-Free Design

5. SETUP

5.1 Initial Charging

1. Before first use, fully charge the massager. Connect the provided USB charging cable to the device's charging port and to a suitable USB power source.
2. The charging indicator light will show the charging status. Once fully charged, the indicator will change (refer to the device's specific light behavior in the quick start guide if available).
3. A full charge can provide approximately one week of 10-minute daily use.

5.2 Positioning the Massager

1. Place the massager around your neck and shoulders.
2. Adjust the position of the massage heads to ensure they comfortably target the desired areas of your neck and trapezius.
3. Secure the massager using the adjustable straps, wearing it like a backpack to keep your hands free during use.

6. OPERATING INSTRUCTIONS

6.1 Power On/Off

- Press and hold the power button (usually marked with a power symbol) to turn the massager on or off.

6.2 Selecting Massage Mode and Intensity

- The Breo N5 Mini offers two massage modes: **Deep Massage** and **Soothing Massage**.
- Press the mode button (refer to control panel icons) to cycle through the available massage modes.
- Press the intensity button (refer to control panel icons) to adjust the massage intensity to your preference.

2 MODES & 2 INTENSITIES

Unique massaging technique



Mode 1: Deep Massage



Mode 2: Soothing Massage



Figure 6: Control Panel and Modes

6.3 Activating Heat Function

- Press the heat button (usually marked with a heat or wave symbol) to turn the heat function on or off.
- The heat function provides a constant temperature of approximately 107.6°F (42°C) for added comfort and muscle relaxation.

6.4 Usage Duration

- The device is designed for sessions of up to 15 minutes. It may have an automatic shut-off feature after this period.
- Allow the device to cool down for a period before starting another session.

7. MAINTENANCE

7.1 Cleaning

- Ensure the massager is turned off and unplugged before cleaning.
- The massage heads can be detached for cleaning. Gently remove them and wash separately with mild soap and water.
- Wipe the main body of the massager with a soft, damp cloth. Do not use abrasive cleaners or immerse the

device in water.

- Allow all parts to dry completely before reassembling or storing.

7.2 Storage

- Store the massager in a cool, dry place, away from direct sunlight and extreme temperatures.
- Keep it away from sharp objects that could damage the faux leather material.

8. TROUBLESHOOTING

If you encounter issues with your Breo N5 Mini Neck Massager, please refer to the following common solutions:

Problem	Possible Cause	Solution
Device does not turn on	Low battery or not charged	Charge the device fully using the provided USB cable.
Massage intensity feels weak	Incorrect mode/intensity setting	Adjust the massage mode and intensity using the control panel buttons.
Heat function not working	Heat function not activated	Press the heat button to activate the heat function. Allow a few moments for it to warm up.
Device stops unexpectedly	Automatic shut-off or low battery	The device may have reached its automatic shut-off time (e.g., 15 minutes). Recharge if the battery is low.

If the problem persists after trying these solutions, please contact Breo customer service for further assistance.

9. SPECIFICATIONS

Feature	Detail
Product Dimensions	9.25 x 6.3 x 6.3 inches
Item Weight	2.75 Pounds (device only), 2.8 Pounds (shipping)
Model Number	Neck N5 MINI
Power Source	Corded Electric (for charging)
Material	Faux Leather
Color	Beige
Manufacturer	Breo

10. WARRANTY AND SUPPORT

Breo is committed to your satisfaction. The Breo N5 Mini Neck Massager comes with a **365-Day warranty**. We also offer lifetime top customer service to address any questions or concerns you may have.

For detailed warranty information, product support, or to download the official user manual PDF, please visit the

Breo official website or refer to the provided PDF document: [Breo N5 Mini User Manual \(PDF\)](#).
You can also visit the official Breo Store on Amazon for more products and support: [Breo Store](#).