

YA-MAN YJFC0B

YA-MAN Blue Green Shot IPL LED Facial Device User Manual

Model: YJFC0B

Brand: YA-MAN

1. IMPORTANT SAFETY INFORMATION

Please read all safety warnings and instructions carefully before using the device. Keep this manual for future reference.

Warning:

The face line is particularly prone to skin entering the irradiation port, so please be careful not to press too hard during flash irradiation.

2. PRODUCT OVERVIEW

The YA-MAN Blue Green Shot is an innovative IPL LED facial device designed to address various skin concerns using a combination of green LED light and IPL (Intense Pulsed Light) technology. It features a smart focus sensor that detects skin dullness and automatically adjusts irradiation power for targeted care.





Figure 2.1: The YA-MAN Blue Green Shot IPL LED Facial Device.

Key Features:

- **Green LED x IPL Technology:** Japan's first device combining green LED (505nm wavelength) and IPL for comprehensive skin care.
- **Focus Sensor:** Quickly detects skin dullness and automatically switches to high-power irradiation for targeted treatment.
- **Two Modes:** Daily Mode for regular maintenance and Special Mode for intensive weekly care.
- **Minimal Downtime:** Designed to be gentle with approximately 95% of users reporting no pain during use.
- **Compact and Handy:** Ergonomic design for comfortable home use.

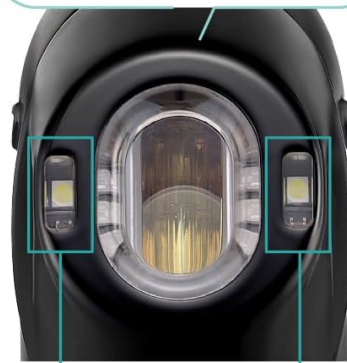


Figure 2.2: The device highlighting its core "Green" technology.

フォーカスセンサーが くすみ*をキャッチ

フォーカスセンサーが肌表面のくすみ*をすばやく検知。
パワーをスマートに切りかえ、濃い肌色の部位には自動でハイパワーで照射。

くすみ*に自動で
ハイパワー照射



フォーカスセンサー



薄い肌色



濃い肌色

*くすみ…うるおい不足による

Figure 2.3: Illustration of the focus sensor detecting skin tone and adjusting power.

3. SETUP

1. **Unpack the Device:** Carefully remove the Blue Green Shot device, AC adapter, and user manual from the packaging.
2. **Inspect for Damage:** Before first use, inspect the device for any visible damage. Do not use if damaged.
3. **Connect Power:** Connect the AC adapter to the device and plug it into a suitable power outlet. The device is wired for stable high-power output.
4. **Initial Patch Test:** Before using on your face, perform a patch test on a small, inconspicuous area of your skin (e.g., inner arm) to check for any adverse reactions. Follow the operating instructions for the patch test.

4. OPERATING INSTRUCTIONS

Follow these steps for effective use of your YA-MAN Blue Green Shot device.

1. **Preparation:** Remove all makeup and thoroughly wash your face before use. Ensure your skin

is clean and dry.

2. **Power On:** Press and hold the power button for 3 seconds to turn on the device. The cooling fan will start.
3. **Select Mode:** Press the "MODE" button to switch between Daily Mode and Special Mode.
4. **Select Level:** Press the "LEVEL" button to adjust the intensity level. Start with the lowest level and gradually increase as your skin tolerates.
5. **Begin Treatment:** Place the irradiation head firmly against your skin. The device will automatically detect skin contact and emit flashes.

4.1. Daily Mode

Daily Mode is designed for regular, quick facial care. It provides 5 continuous flashes at short intervals, allowing for speedy and easy light facial treatment.

- **Frequency:** Use daily for approximately 10 minutes.
- **Application:** Gently glide the device over your skin, ensuring full contact. The device will flash automatically.
- **Recommended Use:** Apply after applying your regular lotion or skincare product to reduce friction.



毎日の照射で
透明印象肌*を **キープ**

秒間隔の5回連続照射で、スピーディー
かつ手軽に光フェイシャル。



*...美白有効成分アルブチン配合の化粧品との併用において


Figure 4.1: Using the device in Daily Mode for continuous flashes.

4.2. Special Mode

Special Mode uses powerful IPL irradiation for intensive care, targeting specific areas affected by

daily stressors like UV rays and dryness.

- **Frequency:** Use once a week. Do not combine with Daily Mode on the same day.
- **Application:** Place the irradiation head on the desired area and press firmly to activate a single flash, similar to stamping. Move to the next area and repeat. Close your eyes during flashes.
- **Caution:** Avoid using on areas with hair growth or around the eyes (eyelids, under-eye bags).

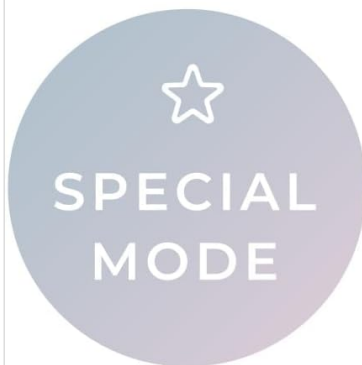

IPLの強力照射で 気になる箇所を **集中ケア**

毎日の生活に潜む紫外線や乾燥などの
ストレス。週の終わりに強力照射でス
ペシャルケアを。



Figure 4.2: Using the device in Special Mode for targeted care.

4.3. Power Off

Press and hold the power button for 3 seconds to turn off the device. The cooling fan may continue to run for up to 10 minutes after shutdown. You may unplug the power cord once the fan stops or if you need to store the device immediately.

5. MAINTENANCE

Proper maintenance ensures the longevity and optimal performance of your device.

- **Cleaning the Irradiation Head:** After each use, wipe the irradiation head with a soft, dry cloth to remove any residue. Do not use abrasive cleaners or solvents.
- **Device Body:** Wipe the main body of the device with a clean, damp cloth.
- **Storage:** Store the device in a cool, dry place away from direct sunlight and extreme temperatures. Keep out of reach of children.

- **Water Resistance:** The device is designed to be water-resistant, allowing for use with lotions. However, do not submerge the device in water.

6. TROUBLESHOOTING

If you encounter any issues with your device, please refer to the following common troubleshooting steps:

Problem	Possible Cause	Solution
Device does not turn on.	Not properly plugged in; power button not held long enough.	Ensure AC adapter is securely connected and plugged into a working outlet. Press and hold the power button for 3 seconds.
No flashes during treatment.	Irradiation head not in full contact with skin; incorrect mode/level selected.	Ensure the irradiation head is firmly and completely pressed against the skin. Check that the correct mode and level are selected.
Device feels warm.	Normal operation; fan running.	This is normal. The device generates heat during operation, and the fan helps dissipate it. If it feels excessively hot, discontinue use and contact support.
Fan noise is loud.	Normal operation.	The fan noise is part of normal operation to cool the device. It may be noticeable but does not indicate a malfunction.

If the problem persists after trying these steps, please contact customer support.

7. SPECIFICATIONS

Feature	Detail
Model Number	YJFC0B
Dimensions (W x D x H)	Approx. 43mm x 214mm x 65mm (1.7 x 8.4 x 2.6 inches)
Weight	Approx. 290g (0.64 lbs) (main unit only)
Power Rating	DC12V 5.0A
Power Consumption	Standby: Approx. 5.0W (Instantaneous Max: Approx. 55W)
Materials	Polycarbonate resin, Acrylic resin
Color	Black
Country of Manufacture	Japan

8. WARRANTY AND SUPPORT




This product comes with a manufacturer's warranty. Please refer to the warranty card included in your package for specific terms and conditions, including the warranty period and coverage.

For technical support, troubleshooting assistance, or warranty claims, please contact YA-MAN customer service. Contact information can typically be found on the official YA-MAN website or on the warranty card.

You may also find additional support resources and FAQs on the [YA-MAN Brand Store on Amazon](#).

© 2023 YA-MAN. All rights reserved.

Related Documents - YJFC0B

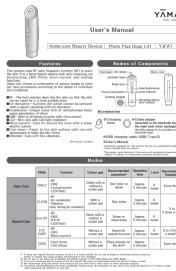	<p>YA-MAN YJMF4 Blue Green Mask Lift User Manual</p> <p>Comprehensive user manual for the YA-MAN YJMF4 Blue Green Mask Lift, providing detailed instructions on usage, safety precautions, maintenance, and troubleshooting for this home-care beauty device.</p>
	<p>YA-MAN YJMF4 Blue Green Mask Lift User Manual</p> <p>Comprehensive user manual for the YA-MAN YJMF4 Blue Green Mask Lift home-care beauty device, covering setup, usage instructions, safety precautions, troubleshooting, and product specifications.</p>
	<p>YA-MAN Liftlogy SP User Manual</p> <p>User manual for the YA-MAN Liftlogy SP home-care beauty device, providing instructions on usage, safety precautions, maintenance, and troubleshooting.</p>
	<p>YA-MAN Blue Green Mask YJMFO User Manual</p> <p>Comprehensive user manual for the YA-MAN Blue Green Mask (Model YJMFO), detailing safety precautions, usage instructions, maintenance, and troubleshooting for this home-care beauty device.</p>
	<p>Yaman YJMF4 Blue Green Mask Lift User Manual</p> <p>This manual provides instructions for using the Yaman YJMF4 Blue Green Mask Lift, including setup, operation, safety precautions, and maintenance. Learn how to achieve optimal results with your skincare device.</p>
	<p>YA-MAN Photo Plus Deep Lift YJFA1 User Manual - Home Care Beauty Device</p> <p>Comprehensive user manual for the YA-MAN Photo Plus Deep Lift YJFA1, a home-care beauty device featuring RF, EMS, Iontophoresis, LED, and Micro-current technologies for facial skin treatment.</p>