

Fxtiaa Ls-001/Ls-002

Fxtiaa 7-in-1 LED Light Therapy Facial Tool Instruction Manual

Model: Ls-001/Ls-002

1. INTRODUCTION

This Fxtiaa 7-in-1 LED Light Therapy Facial Tool is designed to provide various skin benefits through different light wavelengths and microcurrent vibration. It is suitable for most skin types and aims to support skin rejuvenation and overall skin health.

2. PRODUCT FEATURES

The Fxtiaa LED Light Therapy Facial Tool incorporates advanced technology to offer multiple benefits:

- **7 LED Light Colors:** Each color targets specific skin concerns.
- **High-Frequency Vibration:** Promotes tissue metabolism and enhances product absorption.
- **EMS Microcurrent:** Helps to soothe and firm the skin.
- **100% UV-Free:** Safe for regular use without harmful UV exposure.
- **Non-Invasive:** Provides gentle skincare without invasive procedures.
- **Adjustable Head:** The massage head can rotate 90 degrees for comfortable application on various facial contours.

LED Light Color Benefits:

- **Red Light (630nm):** Promotes blood circulation and helps relieve allergic skin conditions.
- **Blue Light (465nm):** Calms allergic skin and inhibits wound infection.
- **Green Light (525nm):** Calms sensitive skin, promoting a sense of well-being.
- **Yellow Light (610nm):** Relieves skin erythema (redness) and swelling.
- **Purple Light (390nm):** Aids in whitening and brightening the skin.
- **Cyan Light (490nm):** Promotes skin metabolism and strengthens skin tissue.
- **White Light (MIX):** Helps to tighten the skin, improve fine lines, and address loose skin.



Figure 1: Diagram illustrating how various LED light colors penetrate different layers of the skin, targeting concerns such as smoothing, skin tightening, pore reduction, oil control, anti-aging, and freckle removal.

7 Colors Led Light Therapy

— Suitable for most skin types



Figure 2: A comprehensive chart outlining the 7 available LED light colors and their respective wavelengths and benefits.

3. SETUP

Charging the Device:


- Before first use, fully charge the device using the provided Type-C USB charging cable.
- Connect the Type-C end of the cable to the charging port on the device and the USB-A end to a compatible power adapter (not included).
- The indicator light will show charging status and turn solid when fully charged.

Adjusting the Head:

- The massage head can be rotated 90 degrees to better adapt to different areas of the face and body. Gently pivot the head to the desired angle.

4. OPERATING INSTRUCTIONS

4 Easy Steps

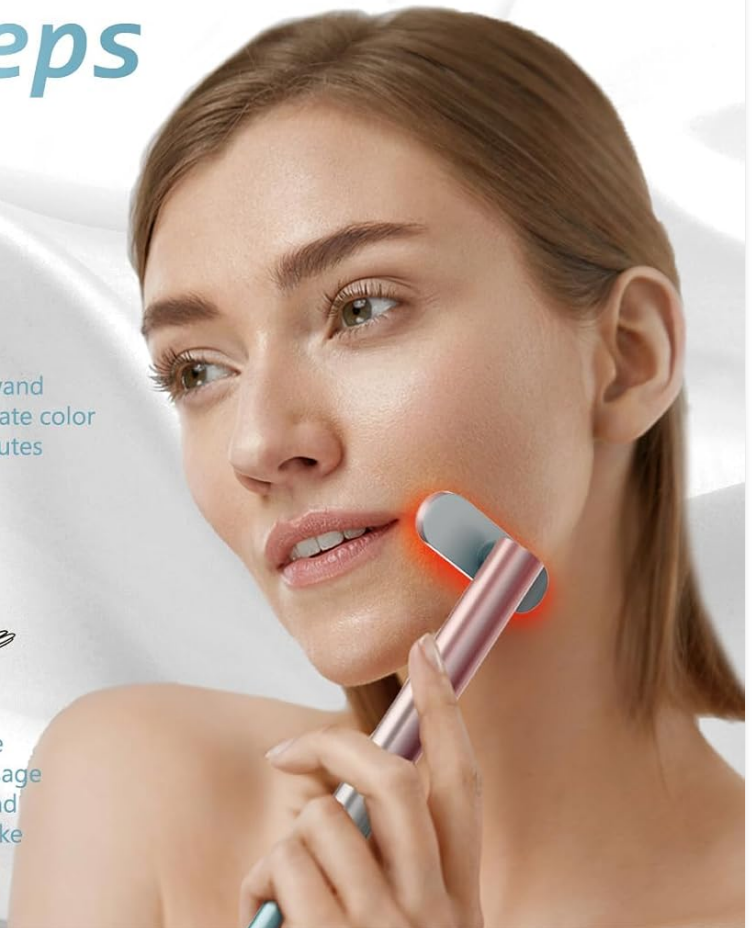

Cleansing



Facial mask



Water cream
essence



Face cream



It can be used with skin care products

Figure 3: This image illustrates the four easy steps for using the facial wand.

Preparation:

- Ensure your skin is clean and dry before use.
- Apply your preferred skincare products, such as serum or cream, to the treatment area. This helps the device glide smoothly and enhances product absorption.

Power On/Off:

- To turn on the device, press and hold the bottom power button for a few seconds until the lights illuminate.
- To turn off the device, press and hold the bottom power button again.

Changing Light Modes:

- Once the device is on, press the top button repeatedly to cycle through the 7 different LED light colors. Select the color that best suits your skincare needs.

Application:

- Gently press the massage head against your skin. The device will automatically activate its vibration and microcurrent functions upon contact with skin.
- Slowly glide the device over the desired treatment areas, such as the face, neck, hands, and around the eyes.
- For optimal results, focus on each area for approximately 10 seconds, repeating the motion 2-3 times.
- Recommended usage is 5-10 minutes per session, 3-5 times per week.



Figure 4: A visual representation of the LED light therapy wand being used on the face, highlighting its ergonomic design and ease of application.

5. PRODUCT CARE AND MAINTENANCE

- After each use, wipe the massage head with a soft, damp cloth to remove any skincare product residue.
- Do not immerse the device in water or use harsh chemicals for cleaning.
- Store the device in a cool, dry place away from direct sunlight and extreme temperatures.
- Keep out of reach of children.

6. TROUBLESHOOTING

- **Device not turning on:** Ensure the device is fully charged. If the issue persists, try using a different charging cable or adapter.
- **No vibration/microcurrent:** Ensure the massage head is in full contact with your skin. The device is designed to activate these functions upon skin contact.
- **Lights not changing:** Press the top button firmly to cycle through the modes. If unresponsive, try restarting the device.

If you encounter any other issues, please refer to the manufacturer's support information.

7. SPECIFICATIONS

Feature	Detail
Product Dimensions	8.11 x 3.38 x 1.29 inches
Weight	4.66 ounces
Item Model Number	Ls-001, Ls-002
Manufacturer	Shenzhen Jingxin Electronics Co., Ltd

Power Source

Type-C USB Rechargeable

8. PRODUCT DEMONSTRATION

Your browser does not support the video tag.

Video 1: This official seller video provides a detailed demonstration of the Fxtiaa LED Light Therapy Wand, showcasing its features, how to switch between the 7 light colors, and various application techniques for the face and eye area.



Figure 5: This banner showcases the LED face light therapy wand, emphasizing its 7 LED colors, 100% UV-free operation, non-invasive nature, and high-intensity LEDs for professional-level treatment.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question