

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

› [Hyperice](#) /

› [Hyperice Normatec 3 Recovery System User Manual \(Standard Size Legs and Backpack\)](#)

Hyperice Normatec 3

Hyperice Normatec 3 Recovery System User Manual

Standard Size Legs and Backpack

1. INTRODUCTION

The Hyperice Normatec 3 Recovery System utilizes dynamic air compression technology to enhance wellness, accelerate recovery, and optimize performance. This system is designed to provide a restorative massage experience for your legs, aiding in circulation and muscle recovery. This manual provides essential information for the proper setup, operation, and maintenance of your Normatec 3 system with standard size leg attachments and the accompanying backpack.



Image 1.1: Hyperice Normatec 3 Legs and Backpack. This image displays the two leg compression attachments connected to the control unit, alongside the dedicated Hyperice backpack for transport and storage.

2. SAFETY INFORMATION

Before using the Normatec 3 Recovery System, please read and understand all safety instructions. Improper use can lead to injury or device malfunction.

- **Body Positioning:** Ensure proper body positioning during use. Consult with a healthcare professional if you have any pre-existing medical conditions or contraindications for compression therapy.
- **Electrical Safety:** Follow all electrical safety guidelines for powered massagers. Do not use the device near water or in excessively humid environments. Ensure the charging port and cables are dry before connecting.
- **Medical Advice:** This device is intended for personal recovery and wellness. It is not a substitute for professional medical treatment. If you experience pain or discomfort, discontinue use immediately and consult a healthcare professional.
- **Children and Pets:** Keep the device out of reach of children and pets.

3. SETUP

3.1 Unpacking and Inspection

1. Carefully remove all components from the packaging.
2. Inspect the control unit, leg attachments, and hoses for any signs of damage. If any components are damaged, do not use the system and contact customer support.

3.2 Sizing and Fit

The Normatec 3 Legs are available in a standard size on Amazon. To ensure proper fit and effectiveness, verify your measurements against the following guidelines:

- **Standard Size:** Inseam 31" - 35"
- **Maximum Thigh Circumference:** 30"

To measure your inseam, measure from just below your crotch to your ankle. For thigh circumference, measure around the fullest part of your mid-thigh.

Normatec size chart

Our size chart helps ensure the standard Normatec fits you comfortably. Use your inseam, height, and thigh circumference to confirm compatibility.

Size

Standard

Inseam

31" - 35"

Height

5'4" - 6'3"

Thigh circumference

30" max

How to measure

Measure from just below your crotch to your ankle for inseam. For thigh, measure around the fullest part, typically mid-thigh.




Image 3.1: Normatec Size Chart. This chart illustrates how to measure inseam and thigh circumference to ensure the standard Normatec leg attachments fit comfortably.

3.3 Connecting the System

1. Place the control unit on a stable, flat surface.
2. Connect the air hoses from the leg attachments to the hose connection port on the control unit. Ensure a secure connection.
3. Plug the power adapter into the control unit's charging port and then into a wall outlet to charge the device before first use, or to operate while charging.

4. OPERATING INSTRUCTIONS

4.1 Control Unit Overview

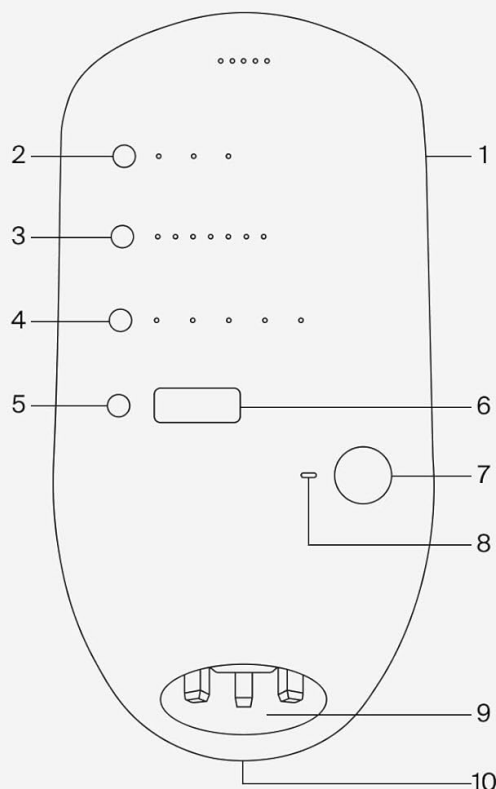

Control unit

- 1 Power button (ON/OFF)¹
- 2 Attachment selection button
- 3 Pressure level button
- 4 ZoneBoost™ button
- 5 Time adjustment
- 6 Display screen
- 7 Start/Stop button
- 8 Bluetooth® status indicator
- 9 Hose connection
- 10 Charging port

Footnotes

- 1. The power button is located on the side of the control unit. Not to be confused with the Start/Stop button on the front of the device (lower-right).

Image 4.1: Normatec 3 Control Unit Diagram. This diagram labels the various buttons and indicators on the control unit, including the power button, attachment selection, pressure level, ZoneBoost, time adjustment, display screen, start/stop button, Bluetooth status indicator, hose connection, and charging port.

- 1. **Power Button (ON/OFF):** Located on the side of the device. Press and hold to turn on or off.
- 2. **Attachment Selection Button:** Used to select the type of attachment (e.g., legs, hips, arms).
- 3. **Pressure Level Button:** Adjusts the compression intensity (7 levels available).
- 4. **ZoneBoost™ Button:** Activates targeted compression in specific zones.
- 5. **Time Adjustment:** Sets the duration of the session.
- 6. **Display Screen:** Shows current settings and session time.
- 7. **Start/Stop Button:** Initiates or pauses the compression session.
- 8. **Bluetooth® Status Indicator:** Shows Bluetooth connectivity status for the Hyperice App.
- 9. **Hose Connection:** Port for connecting the air hoses from the attachments.
- 10. **Charging Port:** For connecting the power adapter.

4.2 Starting a Session

- 1. Put on the leg attachments. Ensure they are snug but not overly tight before inflation.
- 2. Press and hold the **Power Button** (located on the side) to turn on the control unit.
- 3. Use the **Pressure Level Button** to select your desired compression intensity (levels 1-7).
- 4. Adjust the session time using the **Time Adjustment** button.
- 5. Optionally, use the **ZoneBoost™ Button** to apply extra compression to specific areas.
- 6. Press the **Start/Stop Button** to begin the compression session. The display screen will show the remaining time.



The ultimate massage experience

Warm-up and recover like never before. The Normatec 3 uses dynamic air compression to create a restorative massage that helps you feel refreshed faster.

Image 4.2: User operating Normatec 3 Legs. This image shows an individual comfortably seated, using the Normatec 3 leg attachments and interacting with the control unit.

4.3 Using the Hyperice App (Bluetooth Connectivity)

The Normatec 3 features Bluetooth connectivity, allowing you to control and customize your sessions via the Hyperice App.

1. Download the Hyperice App from your device's app store.
2. Ensure Bluetooth is enabled on your mobile device.
3. Open the Hyperice App and follow the on-screen instructions to pair with your Normatec 3 control unit. The **Bluetooth® Status Indicator** on the control unit will confirm connection.
4. Through the app, you can access additional features, customize settings, and track your recovery sessions.



Intuitive, refreshed design

Normatec 3's refreshed system design features an intuitive interface, reduced weight, and the ability to unlock next-level features and customization through the Bluetooth connected Hyperice App.

Image 4.3: Hyperice App interface for Normatec 3. This image shows a user interacting with the Hyperice mobile application to control their Normatec 3 session, demonstrating the intuitive, refreshed design.

5. MAINTENANCE

5.1 Cleaning

- **Control Unit:** Wipe the control unit with a soft, damp cloth. Do not use abrasive cleaners or submerge the unit in water.
- **Leg Attachments:** The nylon material of the leg attachments can be wiped clean with a damp cloth and mild soap. Ensure they are completely dry before storage.
- **Hoses:** Wipe hoses clean with a damp cloth. Ensure no moisture enters the connectors.

5.2 Storage

When not in use, store the Normatec 3 system in a cool, dry place. The included backpack is designed for convenient and protective storage and transport.



Image 5.1: Hyperice Backpack rear view. This image shows the rear of the Hyperice backpack, highlighting its ventilated shoulder straps and ergonomic design.



Image 5.2: Hyperice Backpack open view. This image displays the interior of the Hyperice backpack, showing its various compartments, including a padded lining, compression-molded sunglass case, elastic side pockets, and a laptop compartment, designed to hold the Normatec 3 system.

6. TROUBLESHOOTING

If you encounter issues with your Normatec 3 system, refer to the following common troubleshooting steps. For more detailed assistance, please consult the official Troubleshooting Guide PDF available on the product page or contact Hyperice customer support.

Problem	Possible Cause	Solution
---------	----------------	----------

Problem	Possible Cause	Solution
Device does not power on	Low battery; Power adapter not connected; Device malfunction	Charge the device; Ensure power adapter is securely connected; Contact support if issue persists.
No compression or weak compression	Hoses not securely connected; Leak in attachment or hose; Incorrect pressure setting	Check all hose connections; Inspect attachments and hoses for damage; Increase pressure level.
Bluetooth connection issues	Bluetooth off on mobile device; Device not in pairing mode; Interference	Ensure Bluetooth is on; Restart control unit and app; Move closer to the device.
Attachments feel uncomfortable	Incorrect sizing; Too high pressure setting; Improper positioning	Verify sizing (refer to Section 3.2); Reduce pressure; Adjust attachment position.

7. SPECIFICATIONS

- **Product Dimensions:** 5.71 x 4.43 x 1.34 inches (Control Unit)
- **Item Weight:** 12.4 Pounds (Total System)
- **Model Number:** REJ6 63085 001-00, 67070 001-00, 61020 001
- **Batteries:** 1 Lithium Ion battery required (included)
- **Power Source:** Battery Powered
- **Material:** Nylon (Leg Attachments)
- **Use For:** Legs
- **Manufacturer:** Hyperice
- **First Available:** December 7, 2022



Image 7.1: Normatec 3 Legs and Backpack with dimensions. This image provides visual representation of the approximate dimensions for the Normatec 3 system components.

8. WARRANTY AND SUPPORT

For information regarding your product's warranty, please refer to the documentation included with your purchase or visit the official Hyperice website. For technical support, replacement parts, or any other inquiries, please contact Hyperice customer service directly.


Online Resources:

- Safety Information (PDF): [View PDF](#)
- Troubleshooting Guide (PDF): [View PDF](#)
- User Manual (PDF): [View PDF](#)

Normatec 3 EN Operating Instructions EN Instructions d'usage	Normatec 3 Operating Instructions and Safety Manual Comprehensive guide to the Normatec 3 air compression recovery system by Hyperice, covering setup, operation, safety, and maintenance.
Normatec 3 EN Operating Instructions EN Instructions d'usage	Normatec 3 Operating Instructions and User Manual This document provides comprehensive operating instructions, safety guidelines, technical specifications, and troubleshooting tips for the Hyperice Normatec 3 recovery system.
Normatec 3 EN Operating Instructions EN Instructions d'usage	Normatec 3 User Manual & Safety Guide Hyperice Official user manual and safety guide for the Hyperice Normatec 3 dynamic air compression massager. Learn about setup, operation, features like ZoneBoost, troubleshooting, technical specifications, and warranty.
Normatec 3 EN Operating Instructions EN Instructions d'usage	Normatec 3 Operating Instructions and Safety Information Comprehensive operating instructions, safety guidelines, technical specifications, and warranty information for the Hyperice Normatec 3 air pressure massager.
Let's get started 	Normatec 3 Recovery System User Guide A comprehensive guide to the Normatec 3 recovery system, including setup, usage, and care instructions. Learn how to maximize your recovery with this advanced compression device.

Let's get started 	Normatec 3 Recovery System Quick Start Guide Hyperice A comprehensive guide to setting up and using the Hyperice Normatec 3 dynamic compression recovery system, including component descriptions, setup instructions, usage tips, and care information.
---	--

Documents - Hyperice – Normatec 3

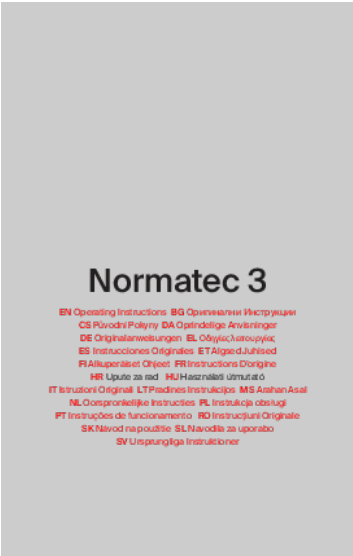


[pdf] User Manual Instructions

Untitled Normatec 3 jui 2023 — IMPORTANT SAFETY INSTRUCTIONS ORIGINAL Read the entire instruction manual before using System not used in accordance with hyperice normatec 30 manuals coolblue be a3 |||

Normatec 3 EN Operating Instructions BG CS Pvodn Pokyny DA Oprindelige Anvisninger DE Originalanwe ... nd/or pains, and to temporarily increase circulation to the treated area. RISKS AND BENEFITS OF THE **NORMATEC 3** SYSTEM The risks and benefits of using the **Normatec 3** System are similar to having a mass...

lang:en score:63 filesize: 2.02 M page_count: 264 document date: 2023-07-03

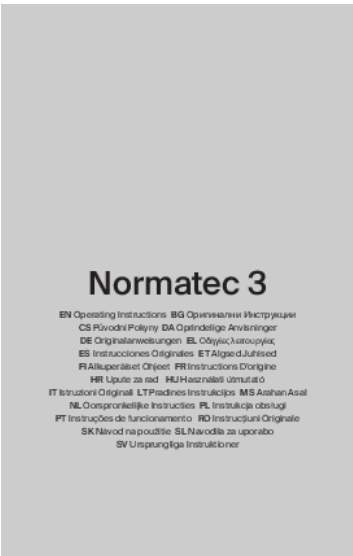

[pdf] User Manual Instructions Declaration of Conformity

Normatec 3 IMPORTANT SAFETY INSTRUCTIONS ORIGINAL Read the entire instruction manual before using System WARNING No modification of this hyperice normatec 30 jalkojen palautumisjarjestelma 63010 006 03 gigantti fi document dv doc D1800012671106 529523 |||

Normatec 3 EN Operating Instructions BG CS Pvodn Pokyny DA Oprindelige Anvisninger DE Originalanwe ... nd/or pains, and to temporarily increase circulation to the treated area. RISKS AND BENEFITS OF THE **NORMATEC 3** SYSTEM The risks and benefits of using the **Normatec 3** System are similar to having a mass...

lang:en score:57 filesize: 27.11 M page_count: 264 document date: 2022-03-24

Let's
get
started



[\[pdf\]](#) Quick Start Guide Guide Warranty

Normatec 3 Quick Start Guide Hyperice Arm Attachment Garment No Control Unit quick start guide
hyperice normatec arm attachment vitalitymedical assets 1 |||

Let s get started Meet your new **Normatec 3** Control unit 2 3 4 5 1. Power button 2.
Attachment selection button 1 3. Pressure level button 4. ZoneBoost™ button 5.
Time adjustment button 6 6. Display screen 7 7. Start/Stop button 8 8. Bluetooth
status indicator 9. Air outlet and power inlet 9 Ho...

lang:en score:56 filesize: 76.59 K page_count: 8 document date: 2024-12-05

Let's
get
started



[\[pdf\]](#) Quick Start Guide Guide Warranty

Quick Start Guide Hyperice Normatec 3 Hip Attachment Garment No Control Unit Vitality Medical quick
start guide hyperice normatec hip attachment vitalitymedical assets 1 |||

Let s get started Meet your new **Normatec 3** Control unit 2 3 4 5 1. Power button 2.
Attachment selection button 1 3. Pressure level button 4. ZoneBoost™ button 5.
Time adjustment button 6 6. Display screen 7 7. Start/Stop button 8 8. Bluetooth
status indicator 9. Air outlet and power inlet 9 Ho...

lang:en score:55 filesize: 280.84 K page_count: 8 document date: 2024-12-09

Normatec 3

EN Operating Instructions FR Instructions D'origine

[Normatec 3 Operating Instructions and User Manual](#)

This document provides comprehensive operating instructions, safety guidelines, technical specifications, and troubleshooting tips for the Hyperice Normatec 3 recovery system.

lang:en score:51 filesize: 87.28 K page_count: 13 document date: 2024-12-05

Normatec 3

EN Operating Instructions FR Instructions D'origine

[Normatec 3 Operating Instructions and Safety Information](#)

Comprehensive operating instructions, safety guidelines, technical specifications, and warranty information for the Hyperice Normatec 3 air pressure massager.

lang:en score:50 filesize: 1.25 M page_count: 13 document date: 2024-12-09

Let's
get
started

[Normatec 3 Recovery System Quick Start Guide | Hyperice](#)

A comprehensive guide to setting up and using the Hyperice Normatec 3 dynamic compression recovery system, including component descriptions, setup instructions, usage tips, and care information.

lang:en score:50 filesize: 283.89 K page_count: 8 document date: 2023-07-07



Let's get
started



[Normatec 3 Recovery System User Guide](#)


A comprehensive guide to the Normatec 3 recovery system, including setup, usage, and care instructions. Learn how to maximize your recovery with this advanced compression device.

lang:en score:50 filesize: 280.97 K page_count: 8 document date: 2024-12-05

Normatec 3 Pickleball CentralNe pas laver ni sécher en machine Ne faire nettoyer à sec MAINTENANCE
DE L APPAREIL unité de commande le tuyau chargeur Normatec3 UserManualpickleballcentral content
Normatec3 UserManual srsItd AfmBOorWTXYbpYrOzGsHvu n8WwO53h79Jz5tpyg2ErDYPn LKyPaOPy
CentralNormatec chargeur La fonction ZoneBoost est destinée être utilisée lorsque vous voulez
concentrer massage sur une zone spécifique La ajoute 60 secondespickleballcentral
LKyPaOPypickleballcentral AfmBOOpEjJQ2sPy4dcyG0Kfs9TsQmCd9f8rBwx5Fh QA5Jp 2zISVtM
secondesIf you do not understand these operating instructions contact Hyperice at 1 949 565 4994 SET
UP THE SYSTEM Step Plug the charger into an electrical pickleballcentral 2zISVtMpickleballcentral
AfmBOoDPIUT4DJf7j yMVGKGN3rSLKIk0ksE6QAfzywS6VwFA0A7oMO Recovery System with
Patented Dynamic Compression Massage Technology Standard Size Legs and Backpack
BackpackNormatec CentralREJ6 63085 001 00 67070 61020 User Manual B0BPDZLVTH B0B94WYYKP
B0BPTJF75D B0B72QBWHC B0C8WB35YJNe electrical B1 KeMhBKrLm media amazon images I B1
KeMhBKrLmpickleballcentral electrical This equipment complies FCC ISSED radiation exposure limits set
forth for uncontrolled environment meets radio frequency RF Exposure B1

[illegible]

Normatec 3 EN Operating Instructions FR Instructions D origine EN TO REDUCE RISKS OF ELECTRIC SHO ... and/or pains and to temporarily increase circulation to the treated area. RISKS AND BENEFITS OF THE **NORMATEC 3** SYSTEM The risks and benefits of using the **Normatec 3** System are similar to having a mass...
lang:en score:45 filesize: 2.54 M page count: 28 document date: 2022-02-18

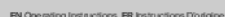




[\[pdf\]](#) User Manual Instructions

Normatec 3 2023 10 13 — CAUTION Do not attempt to take apart the system The has no user serviceable parts When service or repairs are required please contact customer Manual files bbystatic WOUAQ72LrplOL1MlyQtzJg |||


Normatec 3 EN Operating Instructions FR Instructions D origine EN TO REDUCE RISKS OF ELECTRIC SHO ... and/or pains and to temporarily increase circulation to the treated area. RISKS AND BENEFITS OF THE **NORMATEC 3** SYSTEM The risks and benefits of using the **Normatec 3** System are similar to having a mass... lang:en score:38 filesize: 2.54 M page_count: 28 document date: 2022-02-18

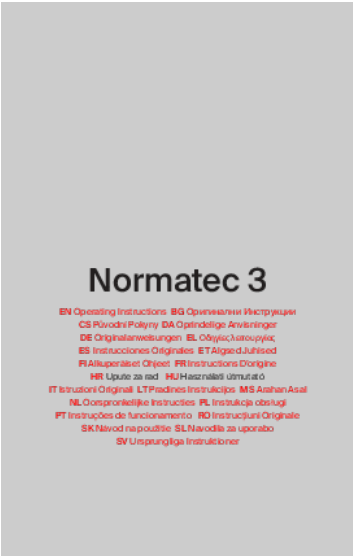

[\[pdf\]](#) User Manual Instructions

Untitled Normatec 3 CoolblueЗа деактивиране на ZoneBoost натискайте бутона докато вече няма светещи светодиодни индикатори Зоните приставката са номерирани във възходящ v1manuals coolblue nl a3 hyperice normatec 30 |||

Normatec 3 EN Operating Instructions BG CS Pvodn Pokyny DA Oprindelige Anvisninger DE Originalanwe ... nd/or pains, and to temporarily increase circulation to the treated area. RISKS AND BENEFITS OF THE **NORMATEC 3** SYSTEM The risks and benefits of using the **Normatec 3** System are similar to having a mass... lang:en score:35 filesize: 2.06 M page_count: 264 document date: 2023-07-03

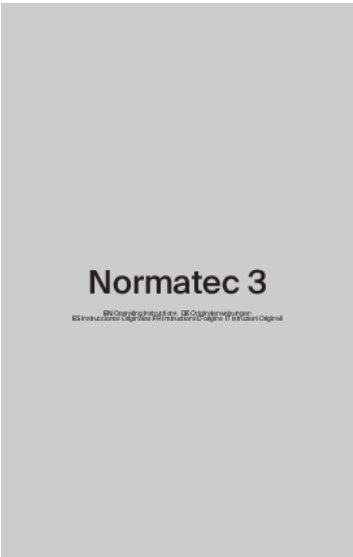

[\[pdf\]](#) User Manual Instructions

Normatec 3 sont garantis par Hyperice contre les défauts de fabrication et pendant une période an à compter la date d achat chez C1R2CC8QPNL m media amazon images I |||

Normatec 3 EN Operating Instructions DE Originalanweisungen ES Instrucciones Originales FR Instructi ... nd/or pains, and to temporarily increase circulation to the treated area. RISKS AND BENEFITS OF THE **NORMATEC 3** SYSTEM The risks and benefits of using the **Normatec 3** System are similar to having a mass... lang:en score:35 filesize: 7.3 M page_count: 64 document date: 2024-10-04




[\[pdf\]](#) Troubleshooting Guide Guide

Hyperice Normatec 3 Recovery System with Patented Dynamic Compression Massage Technology
Standard Size Legs FSA HSA Approved Troubleshooting Guide B0BPDZLVTH B0B94WYYKP
B0BPTJF75D B0B72QBWHC B0C8WB35YJ 9183FLKiSJL m media amazon images I |||

Normatec 3 Dynamic air compression for rapid recovery What s it for How do I use it
Normatec 3 is the next evolution of the iconic Normatec series, using dynamic air compression to advance your wellness, recover faster, improve your training, and maximize your performance. Normatec s patented pr...


lang:en score:34 filesize: 915.24 K page_count: 2 document date: 2023-09-28

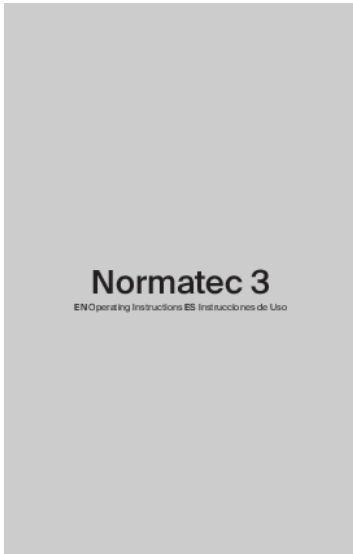

[Normatec 3 User Manual & Safety Guide | Hyperice](#)

Official user manual and safety guide for the Hyperice Normatec 3 dynamic air compression massager. Learn about setup, operation, features like ZoneBoost, troubleshooting, technical specifications, and warranty.

lang:en score:32 filesize: 2.57 M page_count: 28 document date: 2021-11-30

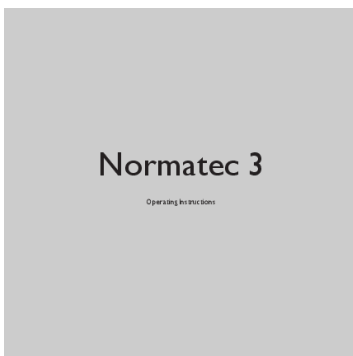


[\[pdf\]](#) User Manual Instructions

Johnny Northrup Manual NormaTec Industries LP NT3 2AY3Y 2AY3YNT3 nt3

Normatec 3 Operating Instructions EN WARNINGS Read the entire instruction manual before using the ... d/or pains and to temporarily increase circulation to the treated areas. RISKS AND BENEFITS OF THE **NORMATEC 3** SYSTEM The risks and benefits of using the **Normatec 3** System are the same as having a mas...

lang:en score:32 filesize: 401.21 K page_count: 23 document date: 2021-06-25

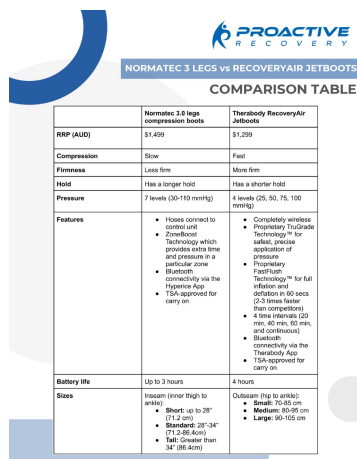

[\[pdf\]](#) Guide

Ice Bathing Beginner s Guide Proactive Recovery DAGGwYOcRmM BAFf5EQNycU Jhay Arianne


Lombog Get the copy of this tableCompression Boots Normatec 3 Legs vs RecoveryAir JetBoots —

RecoveryNormatec by Hyperice and Therabody comparison tablecdn shopify files 1 0709 1352 3994

table v 1718081186 |||

NORMATEC 3 LEGS vs RECOVERYAIR JETBOOTS COMPARISON TABLE ...

lang:en score:30 filesize: 571.7 K page_count: 1 document date: 2024-06-11

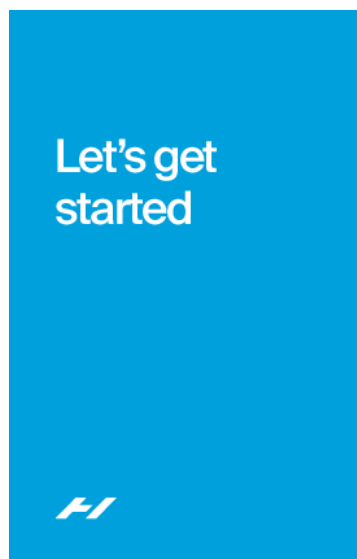

[\[pdf\]](#) User Manual Guide Warranty

Let s get started Pickleball Central 2 giorni fa — Bluetooth® status indicator one attachment please see

the user manual for more information on customizing your settings Normatec3 QuickStartGuide

pickleballcentral content |||

Let's get started Meet your new **Normatec 3** Control unit 2 3 4 5 1. Power button 2.

Attachment selection button 1 3. Pressure level button 4. ZoneBoost™ button 5.

Time adjustment button 6 6. Display screen 7 7. Start/Stop button 8 8. Bluetooth

status indicator 9. Air outlet and power inlet 9 Ho...

lang:en score:30 filesize: 284.04 K page_count: 8 document date: 2021-09-28



NORMATEC

DECLARATION OF CONFORMITY

1. Manufacturer: Normatec Industries, Inc.
 Address: 4600 Pleasant Street, Suite A205, Watertown, Massachusetts, United States of America, 02472
 SIN: US-PA-000000018

2. Authorized Representative: Energo Europe B.V.
 Address: Pijnackerstraat 10, 2514 AP Hap, The Netherlands, +31 20 365 5075
 SIN: NL-NL-000000116

3. Product(s):

Product Name (Model Number)	Basic UDI-DI
Normatec 1 Leg Abutment – Short (30060-201-00)	00010000041410
Normatec 1 Leg Abutment – Standard (30391-001-00)	00010000041720
Normatec 1 Leg Abutment – Tall (30390-201-00)	00010000041730
Normatec 1 Leg Abutment – Short Power (30720-001-00)	00010000041881
Normatec 1 Leg Abutment – Standard Power (31040-001-00)	00010000042191
Normatec 1 Leg Abutment – Tall Power (30760-001-00)	00010000042192
Normatec 3 Axis Abutment – Left – Single (30910-001-00)	00010000042605
Normatec 3 Axis Abutment – Right – Single (30920-001-00)	00010000042606
Normatec 3 Axis Abutment – Single (30901-001-00)	00010000042608

4. The product(s) described above is/are in conformity with:

Title	Document No.
Medical Device Regulation	(EU) 2017/745
RoHS Directive	2011/65/EU

These products are fully compliant and do not contain the restricted substances above listed noted in EU/US Directive 2017/745/65.

5. Intended Purpose: Intended for use with Normatec 3 Control Device.

6. Additional Information:

Declaration statement according to:	Case in accordance with
Classification:	Class II in accordance with Annex VII, table I of Regulation MDR (EU) 2017/745

The declaration of conformity is issued under the sole responsibility of the manufacturer, Normatec Industries, LP.


 Steve Heidmann
 Director of Quality & Regulatory Systems

2023 January 3 | Watertown, MA USA
 Date | Location

[\[pdf\]](#) Decleration of Conformity

EU Declaration of Conformity Normatec 3 Attachments Hyperice hyperice zendesk hc en us article attachments 24899244052251 |||

DECLARATION OF CONFORMITY 1 Manufacturer: NormaTec Industries, LP
Address: 480 Pleasant Street, Su ... Netherlands, 31 70 345 8570 SRN: NL-AR-
000000116 3 Product s : Product Trade Name Model Number **Normatec 3** Leg
Attachment Short 63085 001-00 **Normatec 3** Leg Attachment Standard 63091 001-00
N...

lang:en score:29 filesize: 163.42 K page_count: 1 document date: 0000-00-00

Normatec 3

Normatec 3 Operating Instructions and Safety Manual

Comprehensive guide to the Normatec 3 air compression recovery system by Hyperice, covering setup, operation, safety, and maintenance.

lang:en score:28 filesize: 2.22 M page_count: 285 document date: 2025-10-02



[pdf]

Commençons Le Normatec 3 est très résistant et a été conçu pour une usure normale système ne nécessite aucun entretien ni aucune maintenance Nettoyage de votre 81z7xMV71cL ref dp product quick view m media amazon images I |||

Commençons Découvrez votre nouveau Unit de commande 2 3 4 5 1. Bouton d'alimentation 2. Bouton de ... nts faciles sur hyperice.com/ register-product.

Chargement de votre appareil Branchez votre nouveau **Normatec 3** une prise murale avec le chargeur Hyperice fourni. Vous aurez besoin de charger la batte...

lang:fr **score:28** filesize: 286.05 K page_count: 8 document date: 2023-10-27

 DECLARATION OF CONFORMITY	
1) Manufacturer: NormaTec Industries, LP Address: 480 Pleasant Street, Suite A200, Watertown, Massachusetts, United States of America, 02472 SRN: NL-AR-000000116 815	
2) Authorized Representative: Emurge Europe B.V. Address: Pissinweg 20, 2514 AP The Hague, The Netherlands, +31 70 345 8570 SRN: NL-AR-00000116	
3) Product(s):	
Product Trade Name (Model Number)	Basic UDI-DI
Normatec 3 Hose (81500-001-00)	00810050201713
4) The product(s) described above is/are in conformity with:	
Title	Document No.
Medical Device Regulation	(EU) 2017/745
Notified Body	001155633
These products are fully compliant and do not contain the restricted substances above listed in RoHS Directive 2011/65/EU	
5) Intended Purpose: Air tubing for connecting Normatec 3 Control Device to Normatec 3 Attachments.	
6) Additional Information:	
Technical Documentation compiled according to:	Annex I and II of Regulation MDR (EU) 2017/745
Classification	Class I in accordance with Annex VII, rule 1 of Regulation MDR (EU) 2017/745
This declaration of conformity is issued under the sole responsibility of the manufacturer, NormaTec Industries, LP.	
 Elise Hendriksen Director of Quality & Regulatory Systems	2023 January 31 Watertown MA USA Date / Location
36-WJ228 REV 1	

[pdf] Documentation Declaration of Conformity

EU Declaration of Conformity Normatec 3 Hose Hyperice hyperice zendesk hc en us article attachments

24899244063003 |||

DECLARATION OF CONFORMITY 1 Manufacturer: NormaTec Industries, LP

Address: 480 Pleasant Street, Su ... Netherlands, 31 70 345 8570 SRN: NL-AR-

000000116 3 Product s : Product Trade Name Model Number **Normatec 3** Hose

63080 001-00 Basic UDI-DI 00810050281713 4 The product s described above is/ a...

lang:en **score:28** filesize: 162.27 K page_count: 1 document date: 0000-00-00

 DECLARATION OF CONFORMITY	
NormaTec Industries, LP CE Certificate Number: CE621634 480 Pleasant Street, Suite A200 Watertown, Massachusetts United States of America, 02472 This declaration of conformity is issued under the sole responsibility of the manufacturer. Product Trade Name (Model Number): Normatec 3 REJ6 The object of the declaration described above is in conformity with: • Council Directive 90/269/EEC concerning medical devices • Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment (RoHS) The manufacturer has been verified by the Notified Body listed below to ISO 13485 and complies with Annex I including section 4 of the Medical Device Directive. • Classification: Class I in accordance with Annex II, see 9 of the Medical Device Directive (MDD/CEC) Notified Body: 001155633 800 George The Netherlands BV 1066 CP Amsterdam The Netherlands +31 20 345 8570 Note that the Notified Body number does not apply to the RoHS Directive. Authorized Representative: Emurge Europe Pissinweg 20 2514 AP The Hague The Netherlands +31 70 345 8570	
The Quality Management Representative is Elise Hendriksen Normatec 3, as indicated above, is a system that includes charger, control device, hose, and attachments (Legs, Arms or RPS). Compliance of the device system with the CE mark is the subject of the table below. They are packaged as a system to be used with relevant user instructions from the manufacturer. All activity is governed by internal control and inspection by the manufacturer. Products described above and labeled with the CE Mark are in conformity with the current versions of: Standard No. Title	
ANSI Z535.7-2017	Medical electrical equipment and system electromagnetic immunity test for exposure to radio frequency electromagnetic fields
EN 301 488-1	Electromagnetic Compatibility (EMC) standard for radio equipment and services; Part 1: Common technical requirements
EN 301 488-17	Electromagnetic Compatibility (EMC) standard for radio equipment and services; Part 17: Specific conditions for broadcast and mobile systems
EN 60601-1	Medical electrical equipment Part 1: General requirements for basic safety and essential performance
EN 60601-1-2	Collateral standard: Electromagnetic Compatibility - Requirements and Tests
EN 60601-1-11	Requirements for medical electrical equipment and medical electrical systems used in the home healthcare environment
EN 60601-1-6	Collateral standard: Usability including IEC 62366: Application of usability engineering to medical devices
EN 62366-1	Medical devices - Application of usability engineering to medical devices
 Elise Hendriksen Director of Quality Systems & Regulatory Systems	27 December 2022 Date

[pdf] Declaration of Conformity

EU Declaration of Conformity Normatec 3 Hyperice hyperice zendesk hc en us article attachments

24899211996955 |||

DECLARATION OF CONFORMITY 2797 NormaTec Industries, LP CE Certificate

Number CE621634 480 Pleasant ... y is issued under the sole responsibility of the

manufacturer. Product Trade Name Model Number : **Normatec 3** REJ6 The object of


the declaration described above is in conformity with: Council Dir...

lang:en **score:28** filesize: 142.9 K page_count: 1 document date: 0000-00-00



[pdf] User Manual Instructions

User Manual Hyperice Inc VENOMGO 2AWQY 2AWQYVENOMGO venomgo
Venom Go EN Operating Instructions BG CS Nvod k použit DA Betjeningsvejledning
DE Bedienungsanlei ... oole. TEAVE SISEMISE AKU KOHTA Seadmel Venom Go on
laetav liitium-polmeeraku. Sisemine aku vimaldab **Normatec 3** ssteemi kasutada
kikjal isegi kui seinapistikupesa pole keprast. Enne esmakordset kasuta...
lang:en score:27 filesize: 2.28 M page_count: 164 document date: 2022-06-07

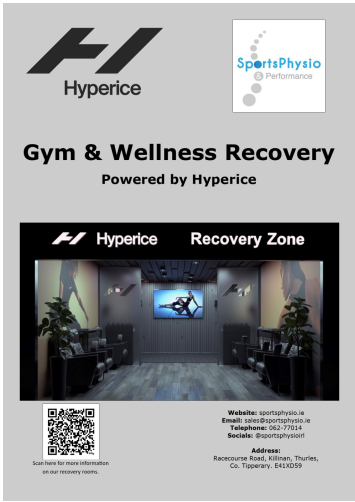


[pdf] Catalog

padraig CatalogsRR Gym Recovery Powered by Hypericesportsphysio ie pub media s RR Hyperice |||
Gym Wellness Recovery Powered by Hyperice Scan here for more information on our
recovery rooms. ... 4 VAT Package Cost with wall box 2,775 VAT Standard Recovery
Room 4 x Zero Gravity Chairs 4 x **Normatec 3** Leg System 1 x Percussion Wall Box
These are suggested expansions of a recovery. These a...
lang:en score:27 filesize: 7.25 M page_count: 16 document date: 2025-04-25



[pdf] Accessories Borchure

Recovery Room pricing 3—Venom Go 2 Device 3—Normatec Lower Leg 1—Hypervolt Pelican Case
MSRP value €5 604 Talk with your local Hyperice account lead for partner hype recoveryroom tierpricing
euros trifoldbrochure 001 ricatto bleed pub mediabox storage rxweb prd exhibitor document exh
2ec1568b 1893 4610 914a e315e2900665 d80b7863 53fb 4d1d b4e0 492145dcd929 |||
Travel Recovery Room Accessories and additional products Recovery Room pricing
Included 6--Hyperv ... Venom 2 Range Normatec Lower Legs Venom 2 Leg Normatec
DreamSeat Recline Normatec Zero Gravity **Normatec 3** Cup Holder Chair Questions
or to order: Jason Hussain Director of Fitness - EMEA jhussai...
lang:en score:27 filesize: 1.26 M page_count: 2 document date: 2024-04-05

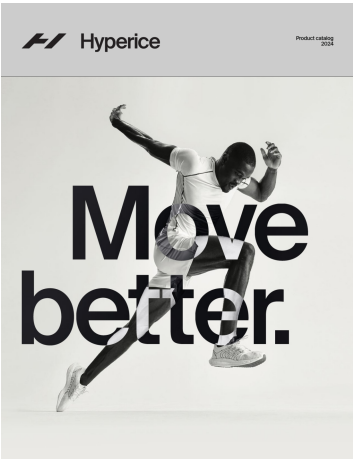

[pdf] Accessories Catalog

Product catalog 2024 3 hour battery life Travel friendly TSA approved carry on What s included Normatec units 2 5V dual USB charger charging cables Page 15 12 HYPE24 INTL ProductCatalog No Price nobleeds pub mediabox storage rxweb prd exhibitor document exh 2ec1568b 1893 4610 914a e315e2900665 921dd6fc ed2c 45de adab eb21ae922af3 |||

Product catalog 2024 Table of contents Hypervolt Normatec Venom Vyper Hypersphere Hyperice Apps ... n or sore muscles, it is a proven way to increase circulation, and reduce pain and soreness. 7 7 **Normatec 3** Legs Dynamic air compression system Features Control unit Bluetooth connected to ...

lang:en score:23 filesize: 3.94 M page_count: 36 document date: 2024-02-16

HYPERICE Normatec 3.0 Leg Control Unit Σύστημα Αποκατάστασης Δυναμικής Συμπίεσης



SKU: IT59148

€399.96

Specifications

Τεχνολογία

- Τεχνολογία Pulse: παρέχει ακριβή θεραπεία σε 5 επεκτασιμότητες ζώνες
- Τεχνολογία ZoneBoost: στοχεύει συγκεκριμένες περιοχές του ποδιού

Τεχνικά Χαρακτηριστικά

- Τύπος: Δυναμικής συμπίεσης
- Εύρος πίεσης: 30-110 mmHg
- Αποσπώμενη χρήση
- Διάρκεια ζωής μπαταρίας: 3 ώρες
- Εγκεκριμένο από TSA

[pdf] Specifications

title HYPERICE Normatec 3 0 Leg Control Unit Σύστημα Αποκατάστασης Διάρκεια ζωής μπαταρίας ώρες Εγκεκριμένο από TSA Page 2 7 επίπεδα συμπίεσης Συνδεσιμότητα Bluetooth μέσω Hyperice app Γενικά Χαρακτηριστικά Μονάδα electroline cy products small appliances personal care massage devices hyperice normatec leg control unit for recovery dynamic compression system |||

HYPERICE **Normatec 3.0** Leg Control Unit SKU: IT59148 399.96 Specifications Pulse: 5 ZoneBoost: : : 30-110 mmHg : 3 TSA 7 Bluetooth Hyperice app. Premium Soft Touch : , , ZoneBoostTM 3 : , , Custom : 1.45 kg ...

lang:en score:23 filesize: 59.48 K page_count: 2 document date: 2025-07-12

HYPERICE Normatec 3.0 Leg Recovery System



SKU: IT59147

€799.00

Specifications

Technical Specifications

- Type: Leg Recovery System
- Compression levels: 7
- Connection: Bluetooth
- ZoneBoost™ Technology: Provides extra time and pressure in a specific zone
- Battery: Up to 3 hours

General Specifications

- Compact, lightweight control unit
- All attachments feature high-quality materials and overlapping zones for gap-free compression
- Intuitive, refreshed design
- The ultimate massage experience

[pdf] Specifications

title HYPERICE Normatec 3 0 Leg Recovery System SKU IT59147 €999 00 Specifications Technical Type Compression levels 7 electroline cy en products small appliances personal care massage devices hyperice normatec leg recovery system |||

HYPERICE **Normatec 3.0** Leg Recovery System SKU: IT59147 799.00 Specifications Technical Specifications Type: Leg Recovery System Compression levels: 7 Connection: Bluetooth ZoneBoostTM Technology: Provides extra time and pressure in a specific zone Battery: Up to 3 hours General Specifications Compac...

lang:en score:22 filesize: 119.75 K page_count: 2 document date: 2023-12-10

HYPERICE Normatec 3.0 Legs
Αποθεραπεία και Αποκατάσταση Με
Συμπίεση Αέρα για τα Πόδια



SKU: IT59147
€999.00

Specifications

Τεχνικά Χαρακτηριστικά

- Τύπος: Αποθεραπεία και αποκατάσταση με συμπίεση αέρα για τα πόδια
- Επίπεδα συμπίεσης: 7
- Σύνδεση: Bluetooth
- Τεχνολογία ZoneBoost™: Παρέχει επιπλέον χρόνο και πίεση σε μια συγκεκριμένη ζώνη
- Μπαταρία: Έως 3 ώρες

Γενικά Χαρακτηριστικά

- Συμμετοχή, ελαφρύ μονάδα ελέγχου
- Όλα τα προσάρτηματα περιλαμβάνουν υλικά υψηλής ποιότητας και

[\[pdf\]](#) Specifications

title HYPERICE Normatec 3 0 Legs Αποθεραπεία και Αποκατάσταση Με Τεχνικά Χαρακτηριστικά Τύπος αποκατάσταση με συμπίεση αέρα για τα πόδια Επίπεδα συμπίεσης 7 Σύνδεση Bluetooth electroline cy products small appliances personal care massage devices hyperice normatec leg recovery system ||| HYPERICE **Normatec 3.0** Legs SKU: IT59147 999.00 Specifications : : 7 : Bluetooth ZoneBoost™: : 3 , , // ...

lang:en **score:21** filesize: 127.23 K page_count: 2 document date: 2025-07-12



February 25, 2022

NormaTec Industries, LP
 N. Prithul Bom
 Responsible Third Party Official
 Regulatory Technology Services, LLC
 1000 Westgate Drive,
 Suite 510k
 Saint Paul, Minnesota 55114

Re: K220217
 Trade/Device Name: Normatec 3
 Regulation Number: 21 CFR 890.5650
 Regulation Name: Powered inflatable tube massager
 Regulatory Class: Class II
 Product Code: BFR
 Dated: January 25, 2022
 Received: January 26, 2022

Dear Prithul Bom:

We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use listed in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. Although this letter refers to your product as a device, please be aware that some cleared products may instead be combination products. The 510(k) Premarket Notification Database located at <https://www.accessdata.fda.gov/cdrh/cdrhinfo/510knewdevices.cfm> identifies combination product submissions. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration. Please note: CDRH does not evaluate information related to contract liability warranties. We remind you, however, that device labeling must be truthful and not misleading.

If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the Federal Register.

U.S. Food & Drug Administration
 1080 Locust Street, Room 5600
 Philadelphia, PA 19106
 (215) 855-1200
 www.fda.gov

Doc ID: 3245712.0.0

[\[pdf\]](#)

K220217 Letter SE Safety Information Hyperice Normatec 3 Recovery System with Patented Dynamic Compression Massage Technology Standard Size Legs and Backpack Health Household 81pK0dJv7kL m media amazon images I ||| February 25, 2022 NormaTec Industries, LP Prithul Bom Responsible Third Party Official Regulatory Technology Services, LLC 1000 Westgate Drive, Suite 510k Saint Paul, Minnesota 55114 Re: K220217 Trade/Device Name: **Normatec 3** Regulation Number: 21 CFR 890.5650 Regulation Name: Powered inflatable tub... lang:en **score:21** filesize: 258.12 K page_count: 2 document date: 2022-02-25




Commercial Recovery Room




RECOVERY ROOM
 Recovery of athletes

Website: sportsphysio.ie
 Email: sale@sportsphysio.ie
 Telephone: 062-77014
 Social: @sportsphysio

Address:
 Racecourse Road, Killin, Thurles, Co. Tipperary, E41XD59

[\[pdf\]](#) Catalog


padraig CatalogsRR Commercial Recovery Roomsportsphysio ie pub media s RR Room ||| Commercial Recovery Room Website: sportsphysio.ie Email: sale@sportsphysio.ie Telephone: 062-77014 S ... 4 VAT Package Cost with wall box 2,775 VAT Standard Recovery Room 4 x Zero Gravity Chairs 4 x **Normatec 3** Leg System 1 x Percussion Wall Box These are suggested expansions of a recovery. These a... lang:en **score:20** filesize: 12.45 M page_count: 23 document date: 2025-04-25

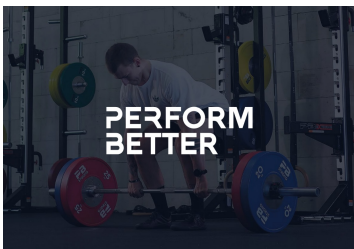

PERFORM BETTER

[\[pdf\]](#) User Manual

PB Image Broschuere 2023 web perform better de media c5 50 cd 1686644668 ||| 2 Perform Better verbindet Expertise und Leidenschaft, um Menschen dabei zu helfen, fitter, strke ... als ein herkmmlisches Wrmekissen und sorgt fr sofortige, wohltuende Linderung bei Muskelkater. 39 **NORMATEC 3** BEINE NormaTec 3 ist die nchste Generation der kultigen NormaTec-Serie und nutzt dynamisc... lang:de **score:18** filesize: 8.7 M page_count: 96 document date: 2023-05-24

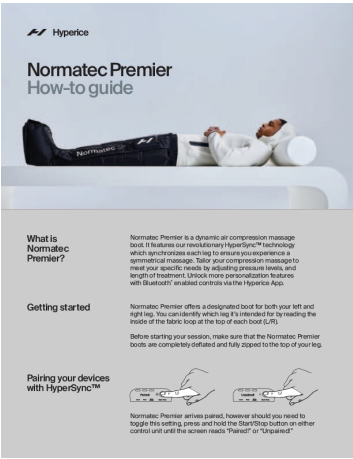

[\[pdf\]](#) Quick Start Guide Guide

Handbuch EN Recovery Boots NormaTec 3 Premier motion sports Normatec v 1729176546 cdn shopify s files 1 0823 4896 6217 |||

Normatec Premier How-to guide What is Normatec Premier Getting started Normatec Premier is a dynamic air compression massage boot. It features our revolutionary HyperSync™ technology which synchronizes each leg to ensure you experience a symmetrical massage. Tailor your compression massage to me...

lang:en score:16 filesize: 437.35 K page_count: 2 document date: 2024-09-20



[\[pdf\]](#) User Manual

docId 84AB2644D56BF68D340E7D340E5B97F1 Handbuch DE Recovery Boots NormaTec 3 Premier motion sports Normatec v 1729176524 cdn shopify s files 1 0823 4896 6217 |||

Normatec Premier Gebrauchsanweisung Was ist Normatec Premier Der Normatec Premier ist ein dynamischer LuftdruckMassageschuh. Er verfügt über unsere revolutionäre HyperSync™-Technologie, die jedes Bein synchronisiert, um eine symmetrische Massage zu gewährleisten. Passen Sie Ihre Kompressionsmassage ...

lang:de score:16 filesize: 609.98 K page_count: 2 document date: 2024-10-08



[\[pdf\]](#)

PB Broschuere ENG web perform better de media 8b a9 82 1712564683 |||

Let visions become reality: Experience meets passion Following the philosophy „Turning visions int ... ining Rope Medizinball Soft Medizinball Jam Ball SAQ Hrden 6er Set Speed Rope SAQ Ladder 89 **Normatec 3** Leg Recovery System Normatec Lower Legs Hypervolt 2 Pro Hypervolt Go 2 Hypervolt Wrme...

lang:en score:14 filesize: 9.41 M page_count: 96 document date: 2024-04-08



[\[pdf\]](#)

PB Broschuere ENG web perform better de media ef 57 02 1724078327 |||

Let visions become reality: Experience meets passion Following the philosophy „Turning visions int ... ining Rope Medizinball Soft Medizinball Jam Ball SAQ Hrden 6er Set Speed Rope SAQ Ladder 93 **Normatec 3** Leg Recovery System Normatec Lower Legs Hypervolt 2 Pro Hypervolt Go 2 Hypervolt Wrme...

lang:en score:14 filesize: 10.48 M page_count: 100 document date: 2024-08-19



[\[pdf\]](#) User Manual

PB Broschuere D web perform better de media d0 c3 52 1712564675 |||

Visionen Wirklichkeit werden lassen: Erfahrung trifft auf Leidenschaft Unter dem Motto „Visionen W ... ining Rope Medizinball Soft Medizinball Jam Ball SAQ Hrden 6er Set Speed Rope SAQ Ladder 89 **Normatec 3** Leg Recovery System Normatec Lower Legs Hypervolt 2 Pro Hypervolt Go 2 Hypervolt Wrme...

lang:de **score:14** filesize: 9.4 M page_count: 96 document date: 2024-04-08




[\[pdf\]](#) User Manual

PB Broschuere D web perform better de media fd 9c 6b 1724078279 |||

Visionen Wirklichkeit werden lassen: Erfahrung trifft auf Leidenschaft Unter dem Motto „Visionen W ... ining Rope Medizinball Soft Medizinball Jam Ball SAQ Hrden 6er Set Speed Rope SAQ Ladder 93 **Normatec 3** Leg Recovery System Normatec Lower Legs Hypervolt 2 Pro Hypervolt Go 2 Hypervolt Wrme...

lang:de **score:14** filesize: 10.48 M page_count: 100 document date: 2024-08-19

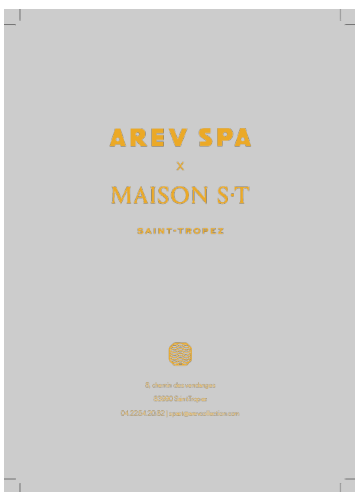

[\[pdf\]](#) Catalog

FR CatalogN°7 Summer feel in the Citadelle Maison STAREV SPA MENU A5 FRcdn shopify s files 1 0554 7540 4995 AREV v 1754434210 |||

8, chemin des vendanges 83990 Saint-Tropez 04.22.54.20.82 spast arevcollection.com 2 AREV SPA ... essotherapie. Fruit de l expertise accumule, elles visent optimiser vos performances sportives. Les **NORMATEC 3** acclrent la rcupration grce leur technologie d impulsion brevete, stimulant la circulati...

lang:fr **score:11** filesize: 1.62 M page_count: 28 document date: 2025-07-07