

[Manuals.plus](#) /

› [Healthkeep](#) /

› Healthkeep BG260W Digital Personal Scale Instruction Manual

Healthkeep BG260W

Healthkeep BG260W Digital Personal Scale

Instruction Manual

1. INTRODUCTION

This manual provides instructions for the proper use and maintenance of your Healthkeep BG260W Digital Personal Scale. Please read this manual thoroughly before use to ensure accurate measurements and safe operation. This scale is designed for personal weight measurement and features high-precision sensors and step-on technology.

2. SAFETY INFORMATION

- Place the scale on a hard, flat surface. Using it on soft surfaces like carpet may affect accuracy.
- The scale's surface is made of tempered safety glass. Avoid dropping heavy objects on it or subjecting it to strong impacts, as this may cause damage.
- Keep the scale dry. Avoid contact with water or other liquids.
- Do not use the scale if the glass surface is cracked or broken.
- Do not stand on the edge of the scale; stand centrally to prevent tipping.
- Keep out of reach of children.
- The maximum weight capacity is 180 kg (400 lbs). Do not exceed this limit.

3. PACKAGE CONTENTS

The Healthkeep BG260W Digital Personal Scale package includes:

- 1 x Healthkeep BG260W Digital Personal Scale
- 2 x 1.5V AAA Batteries (pre-installed)
- 1 x Instruction Manual (this document)

4. PRODUCT FEATURES

- **High-Precision Sensors:** Provides accurate measurements up to 180 kg (400 lbs) in 0.1 kg (0.2 lb) increments.
- **Step-On Technology:** Automatically turns on and displays weight upon stepping on the platform.
- **LED Display:** Bright, backlit display ensures clear visibility even in low-light conditions.
- **Auto On/Off Function:** Conserves battery life by automatically turning off after use.
- **Durable Design:** Constructed with tempered safety glass for robustness.
- **Anti-Skid Padding:** Non-slip feet ensure stability and safety during use.
- **Unit Conversion:** Easily switch between pounds (lb) and kilograms (kg).



Image: An exploded view of the Healthkeep BG260W scale, highlighting its internal high-precision sensors for accurate measurements.

5. SETUP

1. **Unpack the Scale:** Carefully remove the scale from its packaging.
2. **Battery Installation:** The scale comes with 2x AAA batteries pre-installed. If the display does not light up, check the battery compartment on the underside of the scale. Ensure the batteries are correctly inserted with the correct polarity. Remove the plastic insulation tab if present.
3. **Placement:** Place the scale on a hard, flat, and stable surface. Avoid carpets or uneven flooring for accurate readings.
4. **Initial Calibration (First Use or After Moving):** For the first use or after moving the scale, step on it briefly until numbers appear, then step off. The scale will display "0.0" or "CAL" and then turn off. This calibrates the scale to its surface.



Image: Underside of the scale showing the battery compartment and the unit conversion button. The scale uses 2x AAA batteries, which are pre-installed.

6. OPERATING INSTRUCTIONS

1. **Step-On Measurement:** Simply step onto the center of the scale with bare, dry feet. The scale will automatically turn on and display your weight. Remain still until the reading stabilizes.
2. **Reading the Display:** Your weight will be shown on the bright LED display. The display is designed for clear visibility even in dim lighting.
3. **Unit Conversion (lb/kg):** To switch between pounds (lb) and kilograms (kg), locate the "UNIT" button on the underside of the scale (near the battery compartment). Press this button to toggle between the units.
4. **Automatic Shut-off:** The scale will automatically turn off after a few seconds of inactivity to conserve battery life.

Effortless Readability at a Glance

Experience effortless readability with our large, bright, and clear numbers on the backlit display.

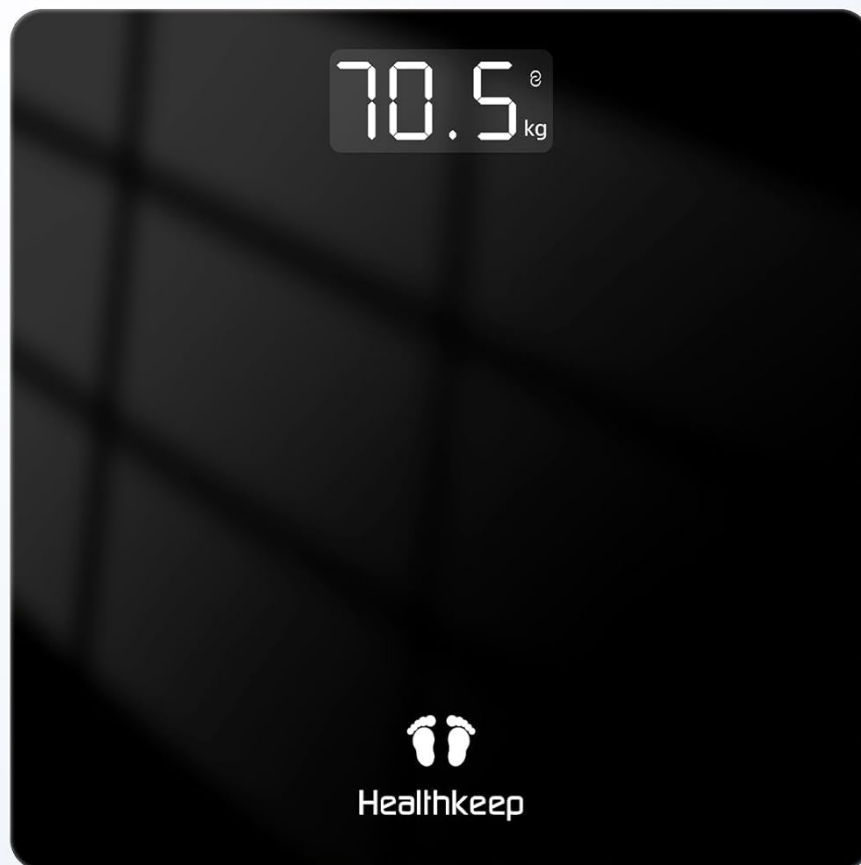


Image: The Healthkeep BG260W scale showing a clear weight reading of 70.5 kg on its bright LED display.

Effortless Measurement with Just One Step

Step onto the scale and watch as your weight appears instantly on the display.



Image: A person stepping onto the Healthkeep BG260W scale, demonstrating the instant weight measurement feature. The display shows 48.2 kg.

7. MAINTENANCE

- **Cleaning:** Wipe the scale surface with a damp cloth and mild detergent. Do not use abrasive cleaners or immerse the scale in water.
- **Storage:** Store the scale in a cool, dry place. Avoid extreme temperatures.
- **Battery Replacement:** If the display shows "Lo" or becomes dim, replace the batteries. Open the battery compartment on the underside, remove old batteries, and insert new 2x AAA batteries, ensuring correct polarity.

8. TROUBLESHOOTING

Problem	Possible Cause	Solution
---------	----------------	----------

Problem	Possible Cause	Solution
No display / "Lo" on display	Batteries are low or incorrectly installed.	Replace batteries with new 2x AAA batteries, ensuring correct polarity. Check if the plastic insulation tab has been removed.
Inaccurate readings	Scale is on an uneven or soft surface (e.g., carpet). Scale needs re-calibration. Weight exceeds maximum capacity.	Place the scale on a hard, flat surface. Step on the scale briefly, then step off to re-calibrate. Ensure weight is within 180 kg (400 lbs) limit.
"Err" on display	Overload.	Step off the scale immediately. The weight capacity has been exceeded.

9. SPECIFICATIONS

- **Model:** BG260W
- **Brand:** Healthkeep
- **Material:** Tempered Glass
- **Color:** Black
- **Dimensions:** Approximately 26 cm x 26 cm (10.24 inches x 10.24 inches)
- **Weight Capacity:** 180 kg / 400 lbs
- **Accuracy:** 0.1 kg / 0.2 lb
- **Display Type:** LED
- **Power Source:** 2 x AAA Batteries
- **Item Weight:** 1.15 kg

Unwavering Precision You Can Trust

Leverage our **decade-long** expertise and the unwavering trust of over **10 million** customers. Experience the pinnacle of precision with our industry-leading scales.

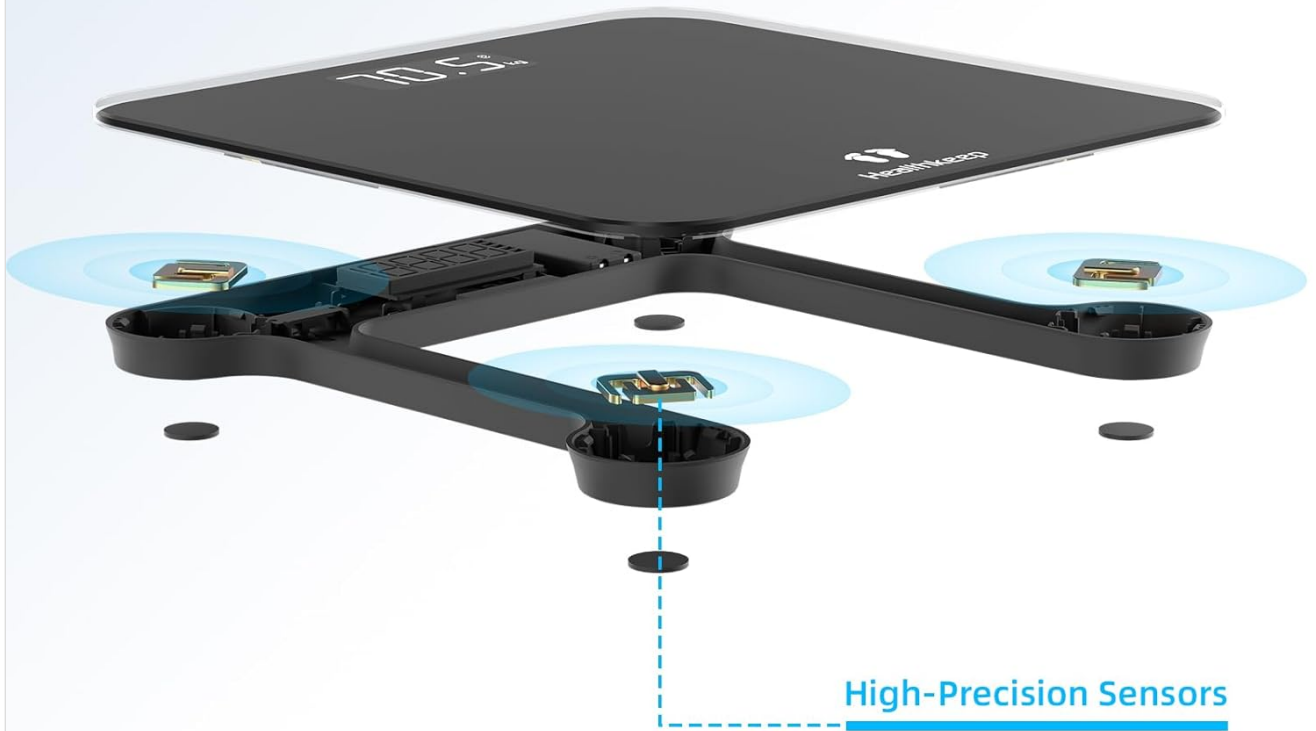


Image: The Healthkeep BG260W scale highlighting its compact dimensions of 10.24 inches by 10.24 inches.

10. WARRANTY AND SUPPORT

Healthkeep products are manufactured with high-quality standards. For warranty information or technical support, please refer to the contact details provided with your purchase documentation or visit the official Healthkeep website. Please retain your proof of purchase for warranty claims.