

DGYAO UV-2212

DGYAO 311nm UV Light Therapy Narrow Band User Manual

Model: UV-2212

1. INTRODUCTION

The DGYAO 311nm UV Light Therapy Narrow Band device is a handheld phototherapy unit designed for home use. It utilizes a high-quality PHILIP PL-S 9W/01 light bulb to emit narrow-band UVB light at 311-312nm, which is commonly used in the treatment of various skin conditions. This manual provides detailed instructions for the safe and effective operation, maintenance, and troubleshooting of your device.

Please read this manual thoroughly before using the device and retain it for future reference.

2. IMPORTANT SAFETY INFORMATION

WARNING: UV radiation can be harmful. Always follow safety guidelines.

- **Consult a Physician:** Before beginning any UV light therapy, consult with a qualified medical professional to determine if this treatment is appropriate for your condition and to establish a suitable treatment plan, including exposure times and frequency.
- **Eye Protection:** Always wear the provided UV-protective eyewear during treatment to prevent eye damage. Do not look directly at the UV light source.
- **Skin Protection:** Protect areas of skin not intended for treatment by covering them with clothing or sunscreen. Avoid overexposure, which can lead to burns, premature skin aging, and increased risk of skin cancer.
- **Children and Pets:** Keep the device out of reach of children and pets. This device is not intended for use by individuals with reduced physical, sensory, or mental capabilities unless supervised.
- **Indoor Use Only:** This device is designed for indoor use only. Do not expose it to water or high humidity.
- **Power Source:** Use only the provided power adapter. Ensure the power outlet matches the voltage requirements of the device. Disconnect from power when not in use or during cleaning.
- **Damaged Device:** Do not use the device if the power cord, plug, or device itself is damaged. Contact customer support for assistance.

3. PACKAGE CONTENTS

Carefully unpack your DGYAO 311nm UV Light Therapy device and check for the following items:



Image: The DGYAO 311nm UV Light Therapy device shown with all its accessories, including the power adapter, protective eyewear, various shaping panels, cleaning cloths, a comb attachment, and a hook.

- 1 x DGYAO 311nm UV Light Therapy Handheld Device (with pre-installed PHILIP PL-S 9W/01 Light Bulb)
- 1 x Power Adapter
- 1 x UV Protective Eyewear
- 4 x Shaping Panels (Ellipse, Round, Rectangle, Long Ellipse)
- 1 x Comb Attachment
- 2 x Cleaning Cloths
- 1 x Hook
- 1 x Stand/Rack

4. SETUP

1. **Unpack the Device:** Remove all components from the packaging.
2. **Connect Power:** Insert the power adapter cable into the device's power port and then plug the adapter into a suitable electrical outlet.
3. **Attach Comb or Shaping Panels (Optional):**
The device comes with a comb attachment and various shaping panels. These can be used to target specific areas or to ensure a consistent distance from the skin.

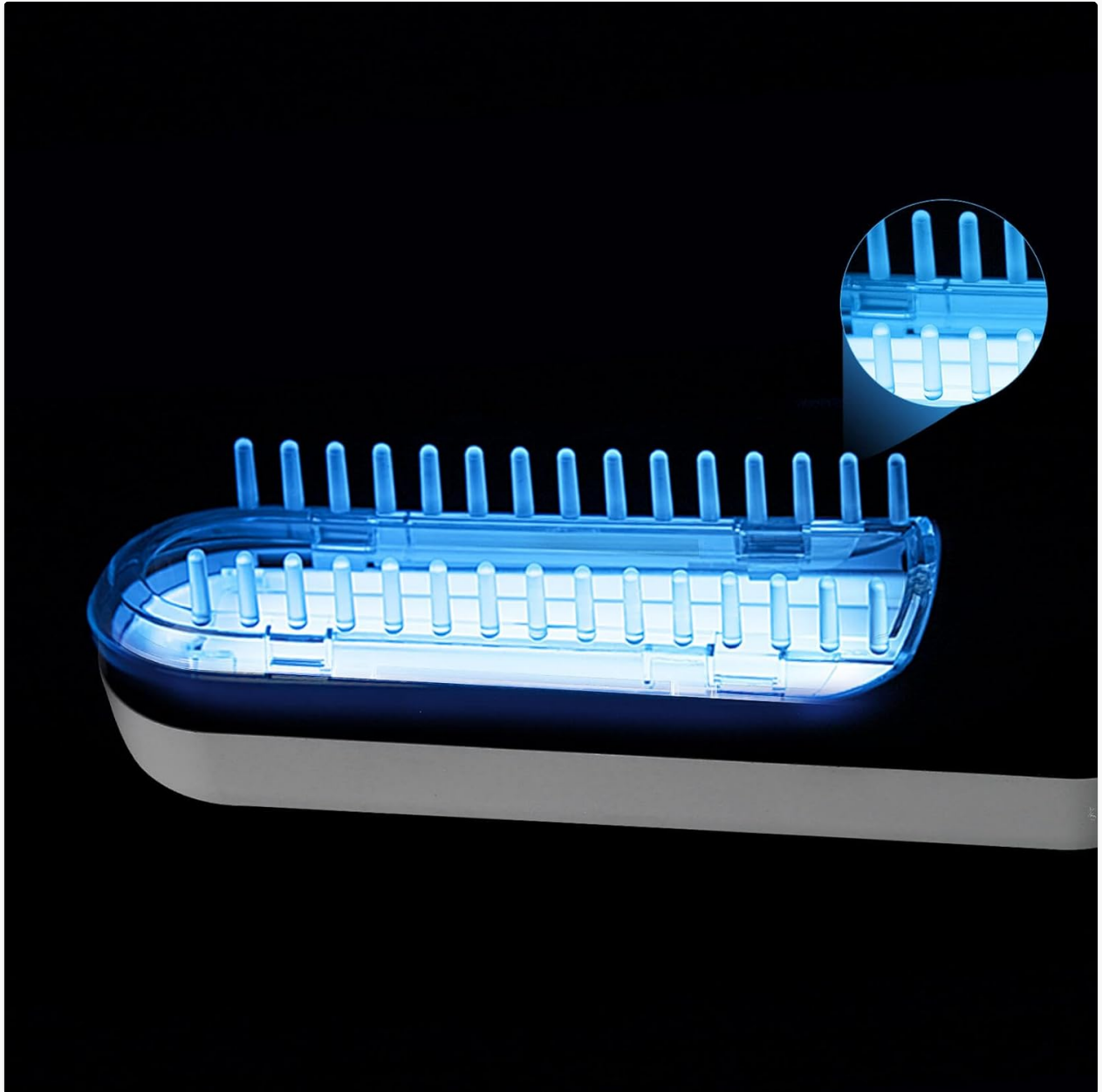


Image: A close-up view of the device showing the comb attachment, which helps maintain a consistent distance from the skin during treatment.

To attach the comb, slide it onto the light head until it clicks into place. To use a shaping panel, simply install the desired panel above the comb teeth. The panels are designed to shape the light for specific skin areas.

Use Block Panels

Simply install the block panel above to the comb teeth. There are four shapes of block panels: Long Ellipse, Ellipse, Rectangle, Round. The block panel is only used for some certain shape of skin area.

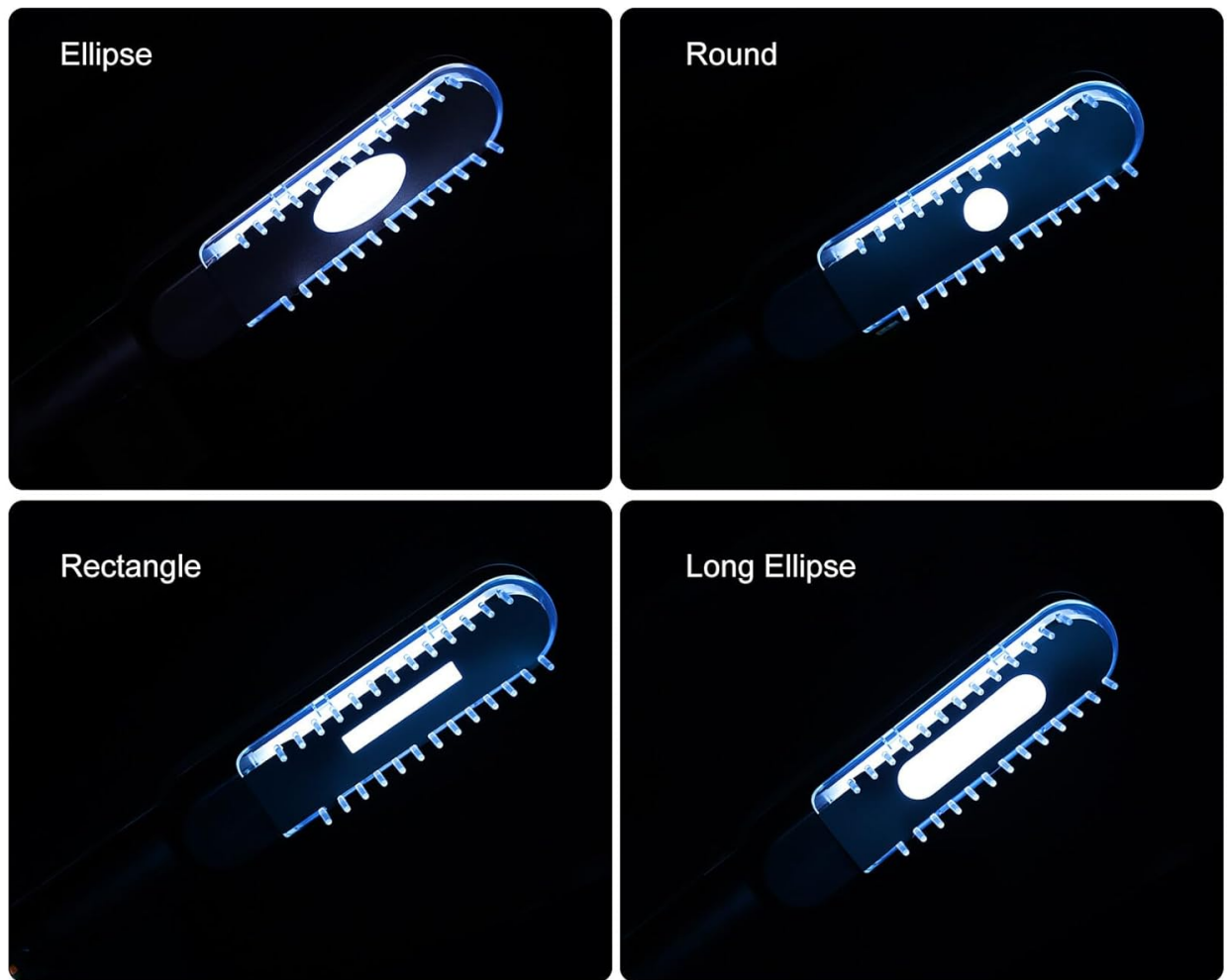


Image: A visual guide demonstrating the four available shaping panels (Ellipse, Round, Rectangle, Long Ellipse) and how they are installed on the device to target different skin areas.

5. OPERATING INSTRUCTIONS

Always wear the provided UV protective eyewear before turning on the device.

1. **Power On:** Press the power button on the device. The display will illuminate.
2. **Set Treatment Time:** Use the 'TIME' button and the '+' / '-' buttons to adjust the desired treatment duration as recommended by your physician. The device features a digital display to show the set time.
3. **Begin Treatment:** Place the light head of the device directly over the affected skin area. Maintain a consistent distance, typically 1-3 cm, or use the comb attachment to ensure proper spacing.



Image: The DGYAO UV therapy device in active use, demonstrating its handheld application on a person's leg, highlighting its ease of use for targeted treatment.

4. **Treat Various Body Areas:** The device is suitable for use on various parts of the body, including the scalp, face, neck, hands, fingers, arms, legs, and the front and back of the body.

Suitable For The Whole Body



Image: An illustrative diagram indicating different body areas such as the scalp, face, neck, hands, fingers, arms, legs, and torso, where the DGYAO UV therapy device can be effectively applied.

5. **Monitor Treatment:** The device will automatically turn off once the set time has elapsed.
6. **Power Off:** After treatment, press the power button to turn off the device. Disconnect the power adapter from the wall outlet.

Note: Always follow your physician's recommendations for treatment duration and frequency. Do not exceed recommended exposure times.

6. MAINTENANCE

- **Cleaning:** Ensure the device is unplugged before cleaning. Use a soft, dry cloth to wipe the exterior of the device. For the light surface and comb, use the provided cleaning cloths. Do not use abrasive cleaners or immerse the device in water.
- **Bulb Replacement:** The PHILIP PL-S 9W/01 bulb is designed for long life. If replacement is necessary, contact DGYAO customer support for guidance or professional service. Do not attempt to replace the bulb yourself unless specifically instructed and qualified.
- **Storage:** Store the device in a cool, dry place, away from direct sunlight and moisture. Keep it in its original packaging or a protective case when not in use.

7. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not power on.	Not plugged in correctly; Power outlet issue; Damaged power cord/adapter.	Ensure power cord is securely plugged into device and wall outlet. Test outlet with another appliance. Inspect power cord and adapter for damage; if damaged, do not use and contact support.
Light does not emit.	Device not powered on; Bulb failure.	Ensure device is powered on. If powered on and light still does not emit, the bulb may need replacement. Contact customer support.
Display is not working.	Device not powered on; Internal malfunction.	Ensure device is powered on. If the device powers on but the display remains blank or faulty, contact customer support.






8. SPECIFICATIONS

Feature	Detail
Brand	DGYAO
Model Number	UV-2212
Light Type	Fluorescent (PHILIP PL-S 9W/01)
Wavelength	311-312nm (Narrow Band UVB)
Wattage	9 Watts
Bulb Base	E26 Medium
Power Source	Electrical
Controller Type	Push Button
Indoor/Outdoor Usage	Indoor
Item Weight	1 kg 110 g

9. WARRANTY AND SUPPORT

The DGYAO 311nm UV Light Therapy Narrow Band device comes with a **1-Year warranty** on the entire system from the date of purchase. This warranty covers manufacturing defects and malfunctions under normal use. For warranty claims, technical support, or any questions regarding your device, please contact DGYAO customer service through the retailer where you purchased the product or refer to the contact information provided in your product packaging.

Related Documents - UV-2212

	<p>DGYAO Red & Infrared Light Therapy Slipper User Manual and Benefits</p> <p>Comprehensive guide to the DGYAO Red Light Therapy Slipper (Model RLT-S21-FTP), covering package contents, usage instructions, benefits like pain relief and improved circulation, warranty, FAQs, and customer reviews.</p>
	<p>DGYAO RLT-S22-WT Red & Infrared Light Therapy Waist & Back User Manual</p> <p>User manual for the DGYAO RLT-S22-WT Red & Infrared Light Therapy Waist & Back device, detailing package contents, introduction, getting started, specifications, warranty, and FAQs for pain relief and circulation enhancement.</p>
	<p>DGYAO RLT-WB01 Red & Infrared Light Therapy Device User Manual</p> <p>Comprehensive user manual for the DGYAO RLT-WB01 Red & Infrared Light Therapy Device. Learn about features, benefits, usage instructions, warranty, and FAQs for pain relief and recovery.</p>
	<p>dermalight® 80 psoracomb UV-B Therapy Device - Operating Instructions</p> <p>Comprehensive operating instructions for the dermalight® 80 psoracomb UV-B narrow-band phototherapy device by Dr. Höhle Medizintechnik. Learn about safety, usage, therapy, and maintenance.</p>
	<p>Daavlin Home Phototherapy Operating Instructions and General Information</p> <p>This document provides comprehensive operating instructions and general information for Daavlin home phototherapy devices, including the UV Series, 7 Series, 2 Series, M Series, and 1 Series. It covers setup, usage, maintenance, and warranty information.</p>