

[Manuals.plus](#) /

> [Elitecran](#) /

> Elitecran 5-in-1 Wireless Full Body Massager with 6 Modes User Manual

Elitecran 1

Elitecran 5-in-1 Wireless Full Body Massager with 6 Modes User Manual

Model: 1

IMPORTANT SAFETY INFORMATION

Please read all safety warnings and instructions carefully before using this device. Failure to follow these instructions may result in electric shock, fire, or serious injury.

- Do not touch the electrodes while the device is operating.
- Do not use the device on the heart, head, internal carotid artery, mouth, intimate areas, surgical areas, infected areas, or damaged skin.
- Do not use the device on individuals who cannot move or communicate, such as children, the elderly, and mentally disabled persons.
- Keep adhesive electrodes away from the skin when the device is still running.
- Do not use the device on a table or unstable surface.
- This massager should not be used by pregnant women.
- If you have any medical conditions, consult a doctor before use.

PRODUCT OVERVIEW

The Elitecran 5-in-1 Wireless Full Body Massager is a portable electric mini massager designed to provide muscle relief and relaxation. It utilizes micro-current technology and offers multiple massage modes and intensity levels.

- **Micro-current Technology:** Made of high-quality, odorless, skin-friendly plastic. Ideal for reducing arm flaccidity, lymphatic drainage, and improving blood circulation.
- **6 Massage Modes & 10 Intensity Levels:** Features 6 high-tech comfortable massage modes

with 10 adjustable intensity levels, catering to different individual needs.

- **Rechargeable Battery:** Can be used for a long period after a full charge via USB. Includes a battery level indicator.
- **Wide Application:** Suitable for various body parts including the neck, shoulders, legs, waist, and arms. Effectively relaxes the body and relieves fatigue.
- **Easy and Safe Use:** Activated by a single button press for intelligent massage. Recommended to start on the back due to varying body sensitivities.

PACKAGE CONTENTS

- 1 x Elitecran 5-in-1 Wireless Full Body Massager Unit
- 1 x USB Charging Cable
- Adhesive Pads (quantity may vary)
- User Manual

SETUP

Charging the Device

1. Connect the provided USB charging cable to the massager's charging port.
2. Plug the other end of the USB cable into a compatible USB power adapter (not included) or a computer's USB port.
3. The battery indicator will show charging status. Charge until the indicator shows a full charge.
4. A full charge allows for extended use.

Attaching the Adhesive Pads


1. Ensure your skin is clean and dry before applying the pads.
2. Remove the protective film from the adhesive pads.
3. Attach the pads firmly to the desired area of your body. Ensure good contact with the skin for effective operation.
4. After use, reattach the protective film to the pads to preserve their adhesive quality.

OPERATING INSTRUCTIONS

Power On/Off and Intensity Adjustment

6 Massage Models

Each mode automatically switches off after 12 minutes of operation

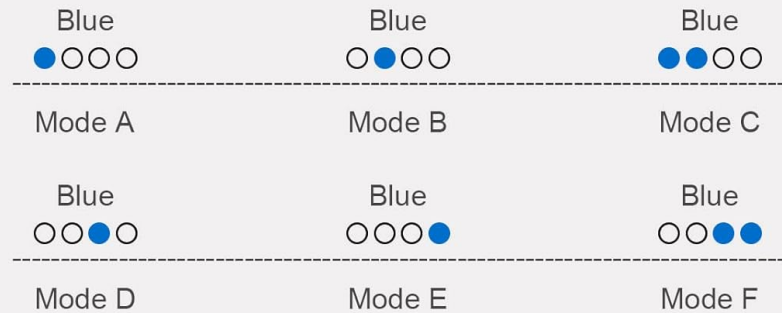

Close-up of the massager's control panel, showing the 'ON/INC' button for power and intensity increase, the 'OFF/DEC' button for power off and intensity decrease, and the 'PROGRAM' button for mode selection.

1. **Power On:** Press the **ON/INC** button once to turn on the device.
2. **Increase Intensity:** After powering on, press the **ON/INC** button repeatedly to increase the massage intensity by one level with each press. There are 10 intensity levels.
3. **Decrease Intensity:** Press the **OFF/DEC** button repeatedly to decrease the massage intensity by one level with each press.
4. **Power Off:** Press and hold the **OFF/DEC** button until the device turns off, or decrease the intensity to level 0.

Selecting Massage Modes

Widely Applicable




Diagram illustrating the six distinct massage modes (A, B, C, D, E, F) available on the Elitecran massager, each indicated by a different pattern of blue lights. Each mode automatically switches off after 12 minutes.

1. Press the **PROGRAM** button to cycle through the 6 available massage modes (Mode A to Mode F).
2. Each mode offers a different massage pattern. Select the mode that best suits your preference.
3. The device is designed to automatically switch off after 12 minutes of operation in any mode.

Application Areas and Benefits

TENS Microelectric Pulse Massage


Simulates massage techniques through bionic microcurrent pulses to relieve muscle discomfort. Continuous, gentle stimulation of neuromuscles and cells to activate muscle vitality and achieve a soothing effect.

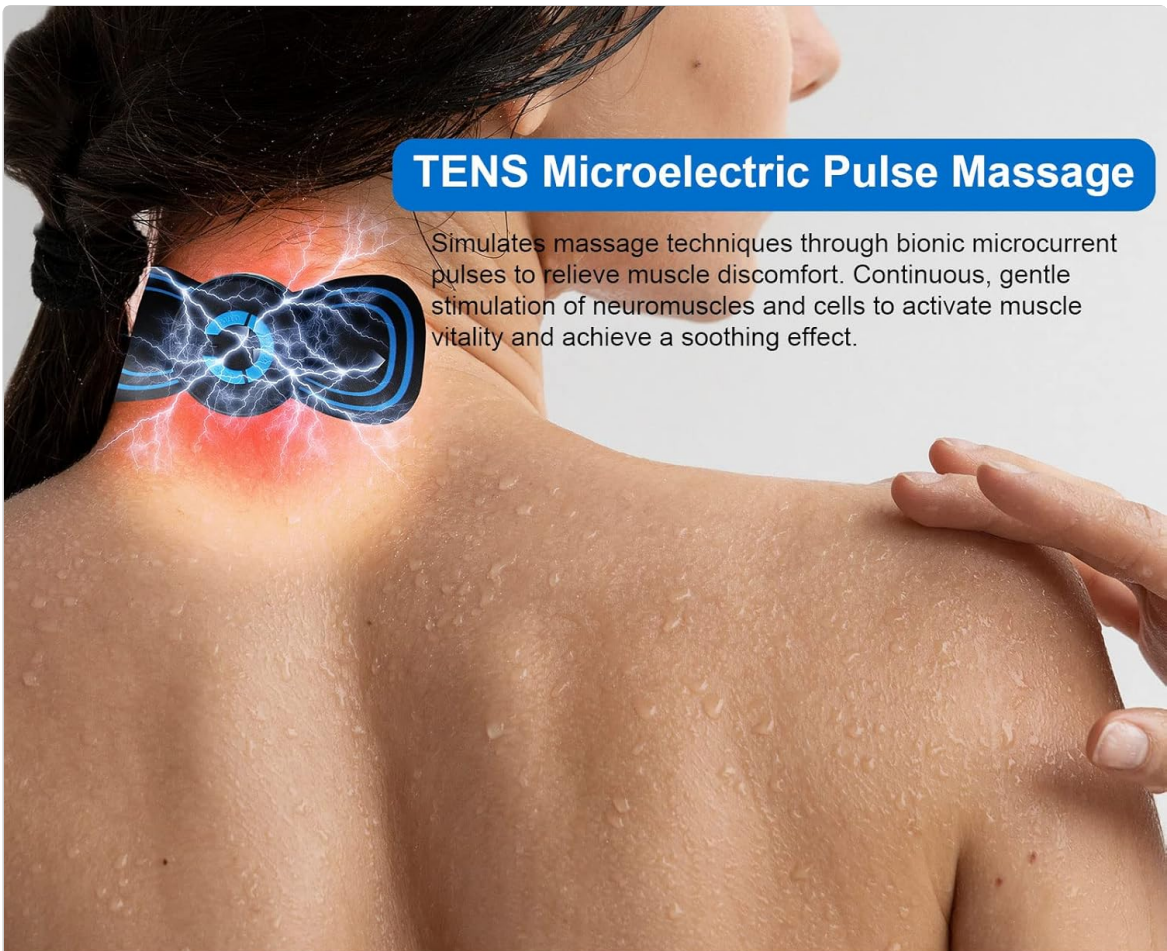

Stiff muscles



TENS Pulse
Massage

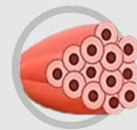

Activation
recovery

Images demonstrating the versatile application of the Elitecran massager on various body parts, including the neck, legs, arms, and waist, highlighting its wide applicability for targeted relief.

Effective Pain Relief



An image showing the Elitecran massager applied to the upper back/neck area, with visual cues indicating pain relief and relaxation. Surrounding images depict individuals experiencing headaches, neck pain, and restful sleep, suggesting the massager's benefits.

Easy to Use

ON/INC

ON: Power on INC: Increase
Click on to switch on and increase the intensity by one level with each press



PROGRAM

Mode switch key
click to switch modes

OFF/DEC

OFF: Power on DEC: Decrease
Each press decreases the intensity by one level until off

An illustration explaining the TENS microelectric pulse massage technology used by the Elitecran device. It shows how bionic microcurrent pulses stimulate neuromuscular cells to activate muscles and provide a soothing effect, transitioning from stiff muscles to activation recovery.

The massager can be applied to various body parts to target specific areas of discomfort:

- **Neck and Shoulders:** For tension and stiffness.
- **Back and Waist:** For general muscle soreness and fatigue.
- **Arms and Legs:** For muscle relaxation and improved circulation.

The TENS microelectric pulse massage simulates bionic microcurrents to relieve muscle discomfort, gently stimulating neuromuscles and cells to activate muscle vitality and achieve a soothing effect, aiding in activation recovery from stiff muscles.

MAINTENANCE AND CARE

- **Cleaning:** Wipe the device with a soft, damp cloth. Do not use abrasive cleaners or immerse the device in water.
- **Adhesive Pads:** Clean the adhesive surface of the pads gently with a damp cloth if they lose

stickiness. Replace pads when they no longer adhere effectively.

- **Storage:** Store the massager in a cool, dry place, away from direct sunlight and extreme temperatures. Ensure the protective film is reattached to the pads before storage.
- **Battery:** To prolong battery life, charge the device regularly, even if not in frequent use.

TROUBLESHOOTING

- **Device not turning on:** Ensure the device is fully charged. Press and hold the ON/INC button firmly.
- **Weak or no massage sensation:** Check if the adhesive pads are firmly attached to clean, dry skin. Increase the intensity level using the ON/INC button. Ensure the battery is charged.
- **Pads not sticking:** Clean the adhesive surface of the pads with a damp cloth. If they are worn out, replace them with new ones.
- **Device turns off unexpectedly:** The device has an automatic shut-off feature after 12 minutes of operation. If it turns off before 12 minutes, ensure the battery is charged.

SPECIFICATIONS

| Feature | Specification |
|------------------------|--------------------------------|
| Brand | Elitecran |
| Model Number | 1 |
| Dimensions (L x W x H) | 15 x 7 x 6 cm |
| Weight | 120 grams |
| Material | Plastic |
| Color | Black |
| Power Source | Battery Powered (Rechargeable) |
| Massage Modes | 6 |
| Intensity Levels | 10 |
| Targeted Body Part | Full Body |
| Massage Technique | Vibrations (Micro-current) |
| Water Resistance | Not Waterproof |
| Wireless | Yes |

WARRANTY AND SUPPORT

Specific warranty information for the Elitecran 5-in-1 Wireless Full Body Massager is not provided in the product data. For details regarding warranty coverage, returns, or technical support, please refer to the product packaging or contact the manufacturer directly.

Keep your purchase receipt as proof of purchase for any warranty claims.