

VRAIKO AURORA-1

VRAIKO Aurora LED Gua Sha Device User Manual

Model: AURORA-1 | Brand: VRAIKO

1. INTRODUCTION

This manual provides detailed instructions for the safe and effective use of your VRAIKO Aurora LED Gua Sha Device. This device combines traditional Gua Sha techniques with modern technology, including heat, vibration, and LED light therapy, to support skin rejuvenation and facial massage.

Please read this manual thoroughly before first use and retain it for future reference.

2. IMPORTANT SAFETY INFORMATION

WARNING: To reduce the risk of injury, electric shock, fire, or damage to the device, always follow these basic safety precautions:

- Avoid using on broken or irritated skin.
- Do not use excessively.
- Keep away from water and do not immerse the device in liquids.
- Dispose of batteries properly according to local regulations.
- Keep out of reach of children.
- Consult a physician before use if you have any pre-existing medical conditions, especially skin conditions, or if you are pregnant.

3. PACKAGE CONTENTS

Verify that all items are present in the package:

- VRAIKO Aurora LED Gua Sha Device
- USB-C Charge Cable
- Silk-like Storage Bag
- User Manual (this document)

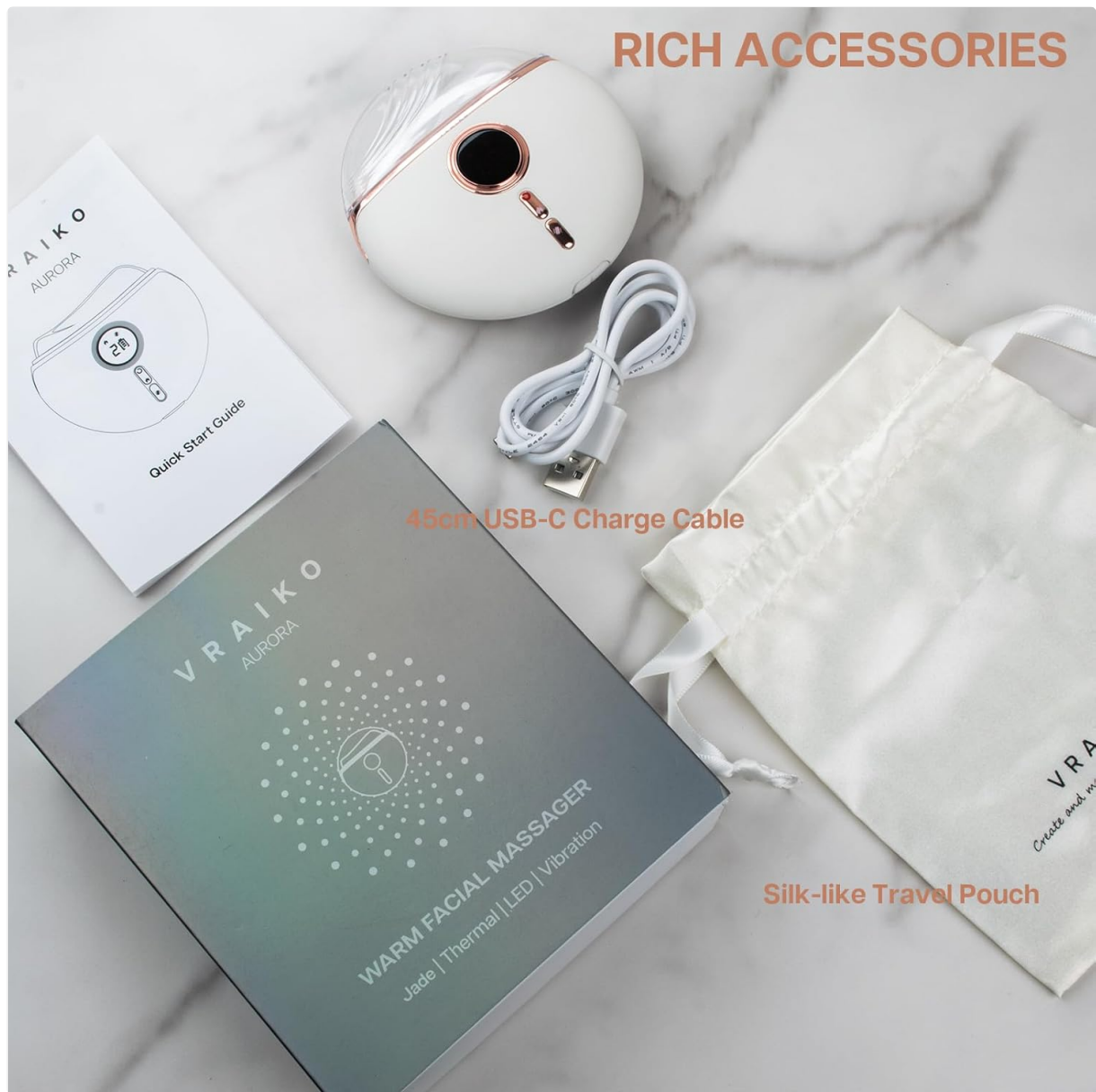


Figure 1: Package Contents. The image displays the VRAIKO Aurora device, a USB-C charging cable, a silk-like storage pouch, and the quick start guide.

4. PRODUCT FEATURES

The VRAIKO Aurora LED Gua Sha Device is designed with multiple functions to enhance your skincare routine:

- **Natural White Jade Material:** Crafted from 100% natural white jade, known for its smooth texture and ability to retain temperature. Each piece is finely polished and carefully selected.
- **Heated Therapy:** The device provides quick and efficient heating, offering a comforting sensation that promotes blood circulation and aids in lymphatic drainage. The temperature can be adjusted between 43~53°C. It takes approximately 30 seconds for the stone to fully warm.
- **Vibration Function:** A gentle vibration function can be activated to provide a relaxing massage, helping to alleviate muscle soreness and enhance product absorption.
- **630nm Far-Infrared LED Light:** This red LED light therapy is integrated to support skin rejuvenation and promote collagen growth. The light can be switched on or off by long-pressing the vibration button.

LED & HEATED MASSAGE



Red-LED gives aids in oil reduction & rejuvenating skin for more radiant and glowing complexion.
The temperature can be easily adjusted between 43~53°C
note: it would take about 30secs to get the stone fully warmed

Figure 2: LED & Heated Massage Feature. The image highlights the device's red LED light and the digital display indicating temperature.



Figure 3: 100% Natural White Jade. The device is shown surrounded by raw white jade, illustrating its primary material.

5. SETUP

5.1 Charging the Device

1. Locate the USB-C charging port on the device.
2. Connect the provided USB-C cable to the device and plug the other end into a compatible USB power adapter (not included) or a computer USB port.
3. The battery indicator on the device's display will show charging status.
4. Charge the device fully before first use. A full charge typically takes approximately 2-3 hours.



Figure 4: Charging the Device. The image shows the device plugged in with its USB-C cable, indicating it is charging.

6. OPERATING INSTRUCTIONS

Follow these steps for optimal use of your VRAIKO Aurora LED Gua Sha Device:

1. **Prepare Your Skin:** After cleansing your face, apply your preferred serum or facial oil to ensure the device glides smoothly over your skin.
2. **Power On and Select Mode:** Press and hold the power button to turn on the device. Use the mode selection buttons to choose your desired combination of heat, vibration, and LED light.
3. **Apply to Face and Neck:** Use a consistent, gentle scraping massage motion. Start from the center of your face and move outwards and upwards. For the neck, move upwards from the collarbone towards the jawline.
4. **Target Specific Areas:** Apply slightly more pressure to areas where you feel muscle soreness or tension, such as the jawline or temples.
5. **Duration:** Use for approximately 5-10 minutes per session, 3-5 times a week.



Figure 5: How to Use the Device. A visual guide illustrating the five-step process for using the VRAIKO Aurora LED Gua Sha Device on the face and neck.

7. MAINTENANCE

7.1 Cleaning

After each use, gently wipe the massage stone clean with a soft, damp cloth. Do not use abrasive cleaners or immerse the device in water.

7.2 Storage

Store the device in its silk-like storage bag in a dry, cool place, away from direct sunlight and extreme temperatures.

8. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery or no charge.	Charge the device fully using the provided USB-C cable.
No heat or vibration.	Incorrect mode selected or device malfunction.	Ensure the correct mode is selected. If the issue persists after charging, contact customer support.
LED light not working.	LED function is off.	Long-press the vibration button to toggle the LED light on/off.

9. SPECIFICATIONS

Feature	Detail
Product Dimensions	1.97 x 1.97 x 1.97 inches
Item Weight	9.17 ounces (0.26 Kilograms)
Item Model Number	AURORA-1
Batteries	1 Lithium Ion battery (included)
Power Source	Battery Powered
Material	White Jade
Manufacturer	VRAIKO



Figure 6: Device Dimensions and Power. This image illustrates the compact size and battery capacity of the VRAIKO Aurora device.

10. WARRANTY AND SUPPORT

10.1 Warranty Information

The VRAIKO Aurora LED Gua Sha Device comes with a worry-free 12-month warranty from the date of purchase. This warranty covers manufacturing defects and malfunctions under normal use.

10.2 Customer Service

For any questions, concerns, or warranty claims, please feel free to contact our friendly customer service team. Refer to the contact information provided on the product packaging or the official VRAIKO website.

	<p>SKIN XTC™ Charis User Manual: Microcurrent & Warm Heat Skincare Device</p> <p>Comprehensive user guide for the SKIN XTC™ Charis device, covering its features, benefits, technical specifications, usage instructions, and important safety contraindications for effective microcurrent and warm heat facial treatments.</p>
	<p>PALSAR7 Galvanická žehlička na tvár 5V1: Návod na použitie a funkcie</p> <p>Objavte PALSAR7 Galvanická žehlička na tvár 5V1. Tento multifunkčný prístroj hĺbkovo čistí, osviežuje a vyhladzuje pokožku pomocou štyroch režimov: teplé čistenie, fototerapia, skin lifting a starostlivosť o oči. Získajte žiarivú a mladistvú pleť.</p>
	<p>7 in 1 H2 O2 Hydra Facial User Manual - UK Beauty Machines</p> <p>Comprehensive user manual for the 7 in 1 H2 O2 Hydra Facial machine by UK Beauty Machines, detailing its principles, functions, installation, operation, and maintenance procedures.</p>
	<p>HECO Aurora XT Loudspeaker Owner's Manual</p> <p>Comprehensive owner's manual for the HECO Aurora XT series of loudspeakers, including models XT 1000, XT 700, XT 300, and XT CENTER 30. This guide provides detailed instructions on setup, connection, bi-wiring, bi-amping, placement, and maintenance for optimal audio performance.</p>
	<p>AURORA DA-CN-0353 Dry Aging Cabinet - Instruction Manual</p> <p>Comprehensive instruction manual for the AURORA DA-CN-0353 Dry Aging Cabinet, covering installation, operation, safety, troubleshooting, and technical specifications. Learn how to properly use and maintain your dry aging unit.</p>
	<p>AS161 6-in-1 Hydro Dermabrasion Facial Care Machine User Manual</p> <p>User manual for the AS161 6-in-1 Hydro Dermabrasion Facial Care Machine, detailing safety precautions, contraindications, operation guide, technical specifications, maintenance, and troubleshooting for skin rejuvenation.</p>