

FOREO BEAR 2 go

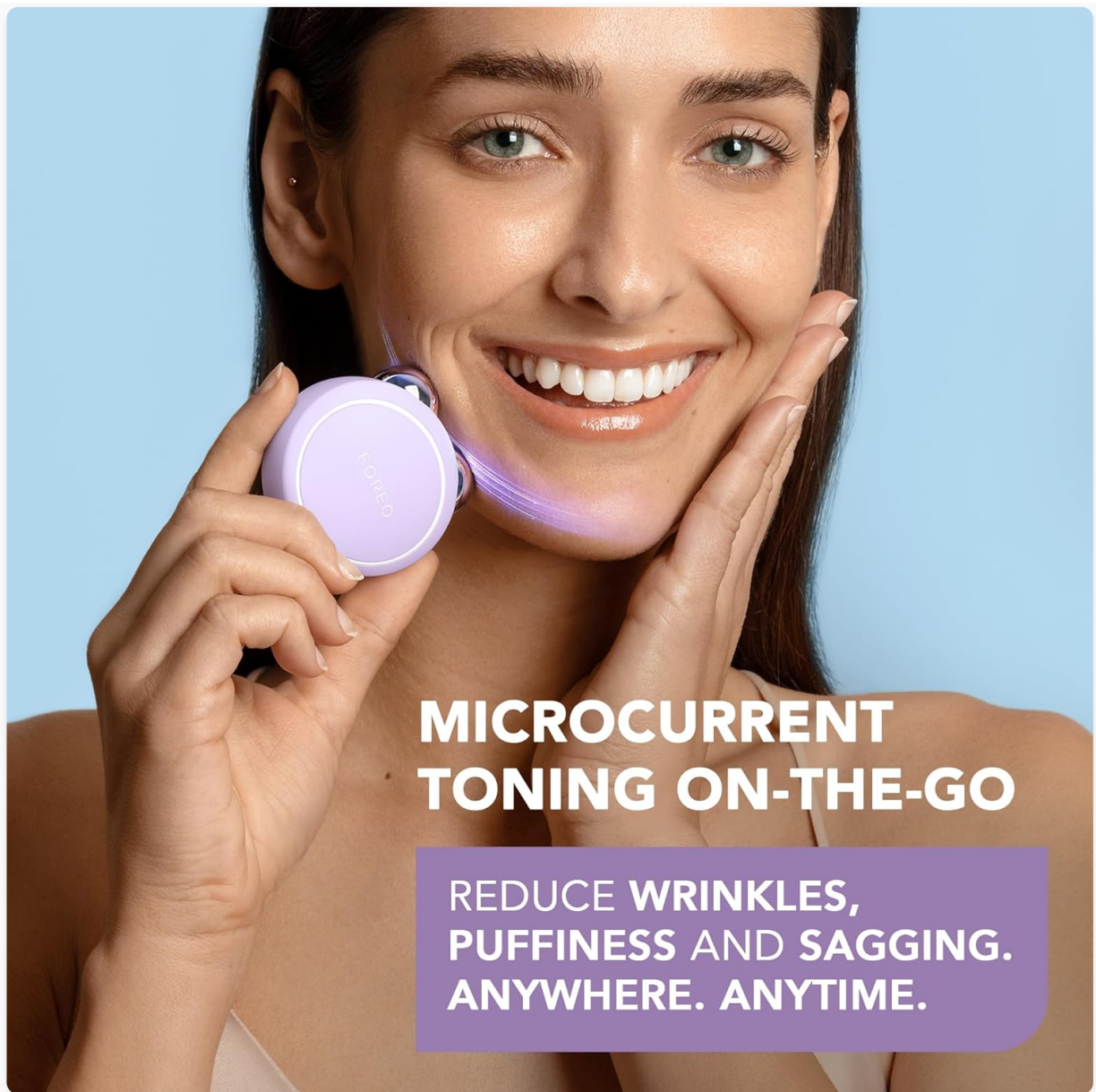
FOREO BEAR 2 go Instruction Manual

Model: BEAR 2 go (F1825)

Brand: FOREO

INTRODUCTION

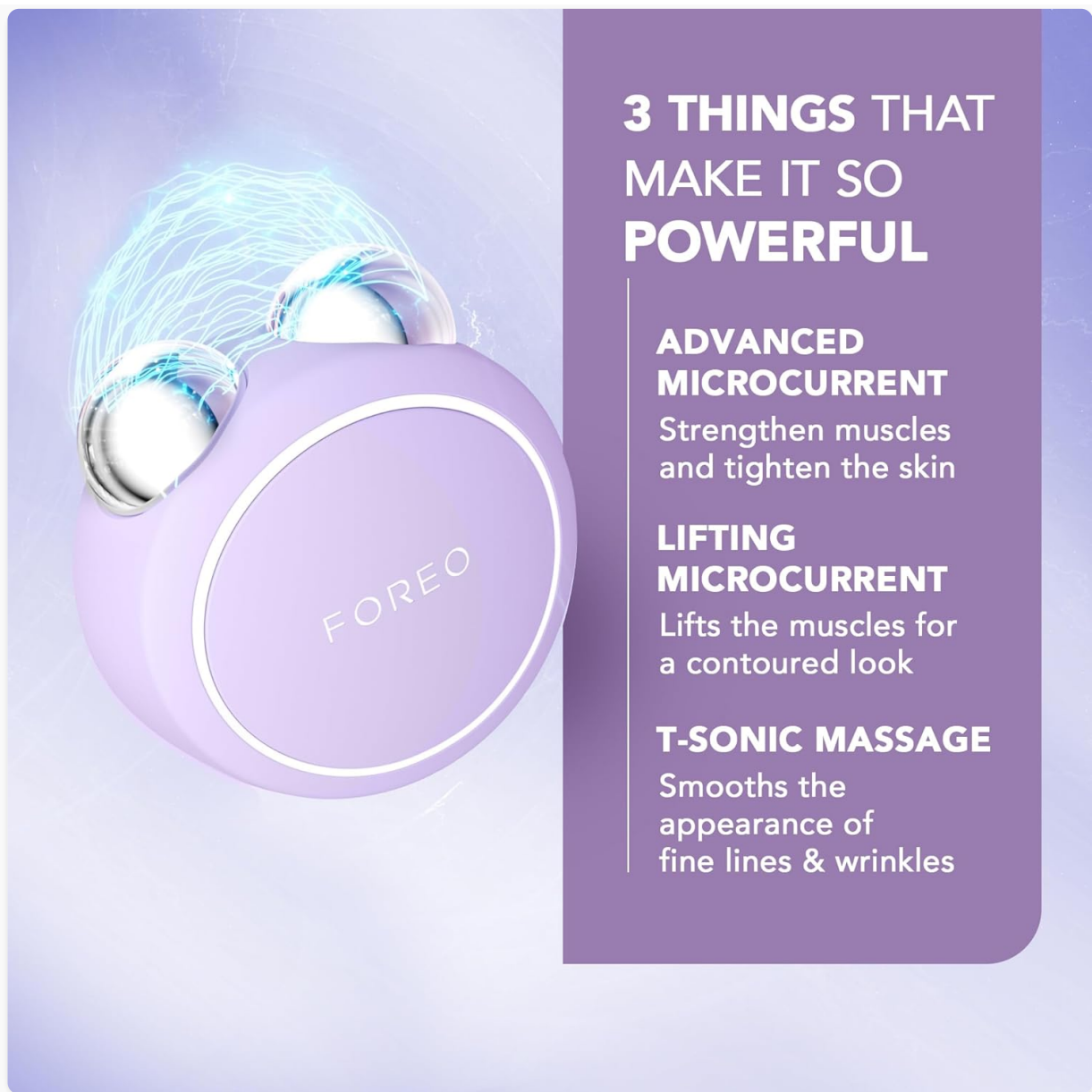
The FOREO BEAR 2 go is a travel-friendly microcurrent facial device designed to tone and firm the skin. This compact and lightweight device is clinically proven to improve the appearance of fine lines and skin firmness, offering a non-invasive solution for facial sculpting and anti-aging care.



MICROCURRENT TONING ON-THE-GO

**REDUCE WRINKLES,
PUFFINESS AND SAGGING.
ANYWHERE. ANYTIME.**

The FOREO BEAR 2 go device in use, highlighting its microcurrent toning capabilities for reducing wrinkles, puffiness, and sagging.



Key technologies of the device include Advanced Microcurrent for muscle strengthening, Lifting Microcurrent for contoured appearance, and T-Sonic Massage for smoothing fine lines.

WHAT'S IN THE BOX?

Your FOREO BEAR 2 go package includes the following items:

- FOREO BEAR 2 go Device (Lavander)
- Quick Start Guide
- General Manual
- USB Charging Cable



✓ YOUNGER-LOOKING FACE
& LIFTED CHEEKBONES

✓ HEALTHIER, BRIGHTER,
PLUMPER SKIN

✓ LESS PUFFINESS
& SAGGING

✓ MORE CONTOURED
FACIAL FEATURES

The complete package contents of the FOREO BEAR 2 go, showing the device, quick start guide, general manual, and USB charging cable.

SETUP

Charging the Device

Before first use, ensure your FOREO BEAR 2 go device is fully charged. Connect the provided USB charging cable to the device and a suitable power source. The device will indicate charging status, typically taking up to 2 hours for a full charge.

Unlocking and Pairing Your Device

To unlock and register your device, follow these simple steps:

1. Download the FOREO For You mobile app on your mobile phone. You can find the app on the App Store or Google Play.
2. Log in to your account or sign up for a new one.
3. Add device (located at the top of your screen within the app).
4. Choose the device series.
5. Press and hold the Power button on your device to connect it to the mobile app.

6. Fill in the purchase information as prompted.

Once these steps are completed, your device is ready for use.


For app download, visit: foreo.com/app

OPERATING INSTRUCTIONS

The FOREO BEAR 2 go is designed for ease of use, integrating seamlessly into your skincare routine. Always ensure your face is clean before use.

Application of Conductive Serum

For optimal conductivity and to ensure a smooth glide, it is essential to apply a suitable serum or conductive gel to your face and neck before using the device. FOREO recommends their Serum Serum Serum, which is lightweight, hydrating, and specifically formulated to enhance microcurrent transfer.



The image shows a gold-colored cylindrical bottle of FOREO Supercharged Serum 2.0 and a purple FOREO BEAR 2 go device. The bottle is labeled 'FOREO SUPERCHARGED SERUM SERUM SERUM - 2.0 - 30 ml e 1 fl.oz.'. The device is purple with two silver electrodes on top. A pink diagonal banner in the top left corner reads 'OUR RECOMMENDATION'. The background is a light purple with faint, glowing blue lines.

BEAR go's MUST-HAVE SIDEKICK

SUPERCHARGED SERUM 2.0

**Improves microcurrent transfer from the device to the skin,
ensuring you get the most out of each use.**

The Supercharged Serum 2.0 is recommended for use with the BEAR 2 go to improve microcurrent transfer and maximize treatment benefits.

Usage Steps

1. Clean your face thoroughly.
2. Apply a suitable serum or conductive gel to enhance conductivity.

3. Turn on the device and select the desired intensity level.
4. Gently glide the device over your face and neck, focusing on areas with fine lines and wrinkles.
5. Use in a circular motion for about 5 minutes.
6. Rinse your face and follow with your regular skincare routine.

The device offers adjustable intensity levels, allowing you to customize the treatment to your comfort and skin needs. The accompanying FOREO For You app provides guided routines and allows for personalized settings, acting as a personal beauty coach.

MAINTENANCE

Cleaning the Device

Cleaning your FOREO BEAR 2 go is simple due to its durable and waterproof material. After each use, wipe the device down with a damp cloth. For thorough sanitization, you may occasionally spritz it with a sanitizing spray. This ensures hygiene and longevity of your device.

TROUBLESHOOTING

Precautions to be taken in the event of changes in the performance of FOREO devices:

- **If your FOREO device is not activated when pressing the center button:**
 - Battery is empty. Recharge using USB charging cable for up to 2 hours until fully charged and then restart your device by holding down the power button.
 - If device cannot be switched off and/or interface buttons do not respond: Microprocessor is temporarily malfunctioning. Press and hold the power button to restart the device.
- **If your FOREO device won't sync to the FOREO app:**
 - Ensure that your device is fully charged and that your Bluetooth is turned on.
 - Switch your Bluetooth off and then on again to try reconnecting.
 - Close the FOREO app and then reopen it to start the process over.
 - Delete and reinstall the application - perhaps the application simply needs to be updated.
 - After deleting and reinstalling the application, make sure to restart your phone.
 - Delete your phone's cookies and cache memory.
 - Ensure that your mobile's operating system has been updated to the latest version.

SPECIFICATIONS

Feature	Detail
Product Dimensions	1.46 x 3.78 x 3.43 inches
Item Weight	2.12 ounces (0.06 Kilograms)
Item Model Number	F1825
Department	Unisex
Batteries	1 Lithium Ion batteries required (included)

Manufacturer	FOREO
Color	Lavander
Skin Type	All
Model Name	BEAR 2 go

WARRANTY & SUPPORT

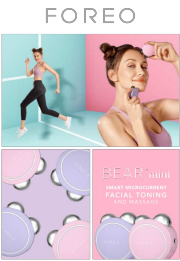
The FOREO BEAR 2 go device comes with a **2-year warranty** from the date of purchase. For any product support, warranty claims, or additional information, please visit the official FOREO support page.

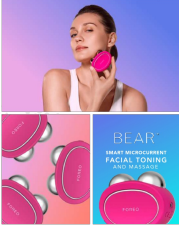

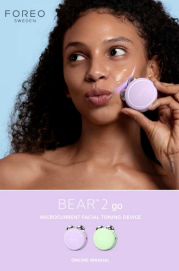
Official Support: foreo.com/support

No specific safety warnings are required for this product.

© 2024 FOREO. All rights reserved.

Related Documents - BEAR 2 go

	<p>FOREO BEAR™ 2 Advanced Microcurrent Facial Toning Device User Manual</p> <p>Comprehensive user manual for the FOREO BEAR™ 2, an advanced FDA-cleared microcurrent facial toning device. This guide covers setup, usage, types of microcurrent, safety warnings, cleaning, troubleshooting, and warranty information.</p>
	<p>FOREO BEAR™ mini: Smart Microcurrent Facial Toning Device User Manual</p> <p>User manual for the FOREO BEAR™ mini, a smart microcurrent facial toning and massage device. Learn about its features, usage, safety, and benefits for firmer, younger-looking skin.</p>
	<p>FOREO BEAR Smart Microcurrent Facial Toning Device User Manual</p> <p>Comprehensive user manual for the FOREO BEAR smart microcurrent facial toning device, covering getting started, features, usage, safety warnings, troubleshooting, warranty, and specifications.</p>

<p>FOREO</p>  The image shows the cover of the FOREO BEAR Smart Microcurrent Facial Toning and Massage User Manual. It features a woman's face in the top left, a close-up of the pink device in the top right, and a collection of the device in various colors (pink, red, orange, yellow) in the bottom left. The text 'FOREO' is at the top, and 'BEAR SMART MICROCURRENT FACIAL TONING AND MASSAGE' is in the top right.	<p>FOREO BEAR Smart Microcurrent Facial Toning and Massage User Manual</p> <p>Comprehensive user manual for the FOREO BEAR smart microcurrent facial toning and massage device, covering setup, usage, safety warnings, troubleshooting, and warranty information.</p>
<p>FOREO</p>  The image shows the cover of the FOREO BEAR mini User Manual. It features a woman's face in the top right, a close-up of the purple device in the top left, and a collection of the device in various colors (purple, blue, green, yellow) in the bottom left. The text 'FOREO' is at the top, and 'BEAR mini SMART MICROCURRENT FACIAL TONING' is in the top right.	<p>FOREO BEAR™ mini User Manual: Smart Microcurrent Facial Toning Device</p> <p>Comprehensive user manual for the FOREO BEAR™ mini, detailing its features, usage instructions, safety warnings, troubleshooting, and warranty information for smart microcurrent facial toning.</p>
<p>FOREO</p>  The image shows the cover of the FOREO BEAR 2 go Microcurrent Facial Toning Device User Manual. It features a woman's face in the top left, a close-up of the green device in the top right, and a collection of the device in various colors (green, yellow, orange, red) in the bottom left. The text 'FOREO' is at the top, and 'BEAR 2 go SMART MICROCURRENT FACIAL TONING' is in the top right.	<p>FOREO BEAR 2 go Microcurrent Facial Toning Device User Manual</p> <p>Comprehensive user manual for the FOREO BEAR 2 go, a microcurrent facial toning device. Learn about features, usage, safety warnings, cleaning, troubleshooting, warranty, and disposal.</p>