

FOREO BEAR 2 go

FOREO Bear 2 go Microcurrent Facial Toning Device Instruction Manual

Model: BEAR 2 go (F1849)

1. INTRODUCTION

The FOREO Bear 2 go is a compact, travel-friendly microcurrent facial toning device designed to smooth, firm, and lift the skin. This device utilizes advanced microcurrent technology and T-Sonic pulsations to target facial muscles, contributing to a more toned and contoured appearance.

It is clinically proven to improve lines and skin firmness within a week of consistent use. This manual provides essential information for the safe and effective operation, maintenance, and troubleshooting of your device.



Image: The FOREO Bear 2 go device, presented in its transparent retail packaging.

2. WHAT'S IN THE BOX

Upon opening your FOREO Bear 2 go package, you should find the following items:

- FOREO Bear 2 go Device (Pistachio color)
- Quick Start Guide
- General Manual
- USB Charging Cable



Image: All components included in the FOREO Bear 2 go product box.

3. SETUP AND UNLOCKING YOUR DEVICE

Before first use, you must download the FOREO For You mobile application to unlock and register your device. Follow these simple steps:

1. Download the FOREO For You app on your mobile phone.
2. Log in to your account or sign up for a new one.
3. Add your device (locate the "Add device" option at the top of your screen within the app).
4. Choose the correct device series.
5. Press and hold the Power button on your device to connect it to the mobile app via Bluetooth.
6. Fill in the required purchase information.

Your device is now ready for use.

For app download and support, visit: [FOREO App](#) and [FOREO Support](#).

UNLOCK YOUR DEVICE

Before first use, download the FOREO For You mobile app to unlock and register your device. Follow these simple steps;

1. Download the FOREO For You app on your mobile phone
2. Log in to your account or sign up for a new one
3. Add device (on the top of your screen)
4. Choose device series
5. Press and hold the Power button to connect your device to the mobile app
6. Fill in the purchase information

And your device is ready to use!

TROUBLESHOOTING

PRECAUTIONS TO BE TAKEN IN THE EVENT OF CHANGES IN THE PERFORMANCE OF FOREO DEVICES.

If your FOREO device is not activated when pressing the center button:

- Battery is empty. Recharge using USB charging cable for up to 2 hours until fully charged and then restart your device by holding down the power button. If device cannot be switched off and/or interface buttons do not respond:
- Microprocessor is temporarily malfunctioning. Press and hold the power button to restart the device.

If your FOREO device won't sync to the FOREO app:

- Ensure that your device is fully charged and that your Bluetooth is turned on.
- Switch your Bluetooth off and then on again to try reconnecting.
- Close the FOREO app and then reopen it to start the process over.
- Delete and reinstall the application - perhaps the application simply needs to be updated.
- After deleting and reinstalling the application, make sure to restart your phone.
- Delete your phone's cookies and cache memory
- Ensure that your mobile's operating system has been updated to the latest version.

Image: Detailed steps for initial device setup and unlocking via the mobile application.

4. OPERATING INSTRUCTIONS

The FOREO Bear 2 go utilizes microcurrent technology and T-Sonic pulsations to provide a facial toning experience. For optimal results, always use the device with a conductive serum.

Key Technologies:

- **Advanced Microcurrent:** Designed to strengthen and tighten facial muscles.
- **Lifting Microcurrent:** Aims to lift facial muscles for a more contoured appearance.
- **T-Sonic Massage:** Helps to smooth the appearance of fine lines and wrinkles.



Image: Visual representation of the three core technologies integrated into the device.

Usage Steps:

1. **Cleanse:** Ensure your face and neck are clean and dry.
2. **Apply Serum:** Apply a thin layer of a conductive serum (e.g., FOREO Supercharged Serum 2.0 or a similar water-based conductive gel) to the areas you intend to treat. This is crucial for the microcurrent to work effectively and comfortably.
3. **Activate Device:** Turn on your FOREO Bear 2 go device.
4. **Treatment:** Gently glide the device across your skin in upward and outward motions, following the contours of your face and neck. Refer to the FOREO For You app for guided routines and specific treatment areas.
5. **Duration:** Follow the recommended treatment duration as guided by the app, typically a few minutes per session.
6. **After Use:** Turn off the device and clean it.

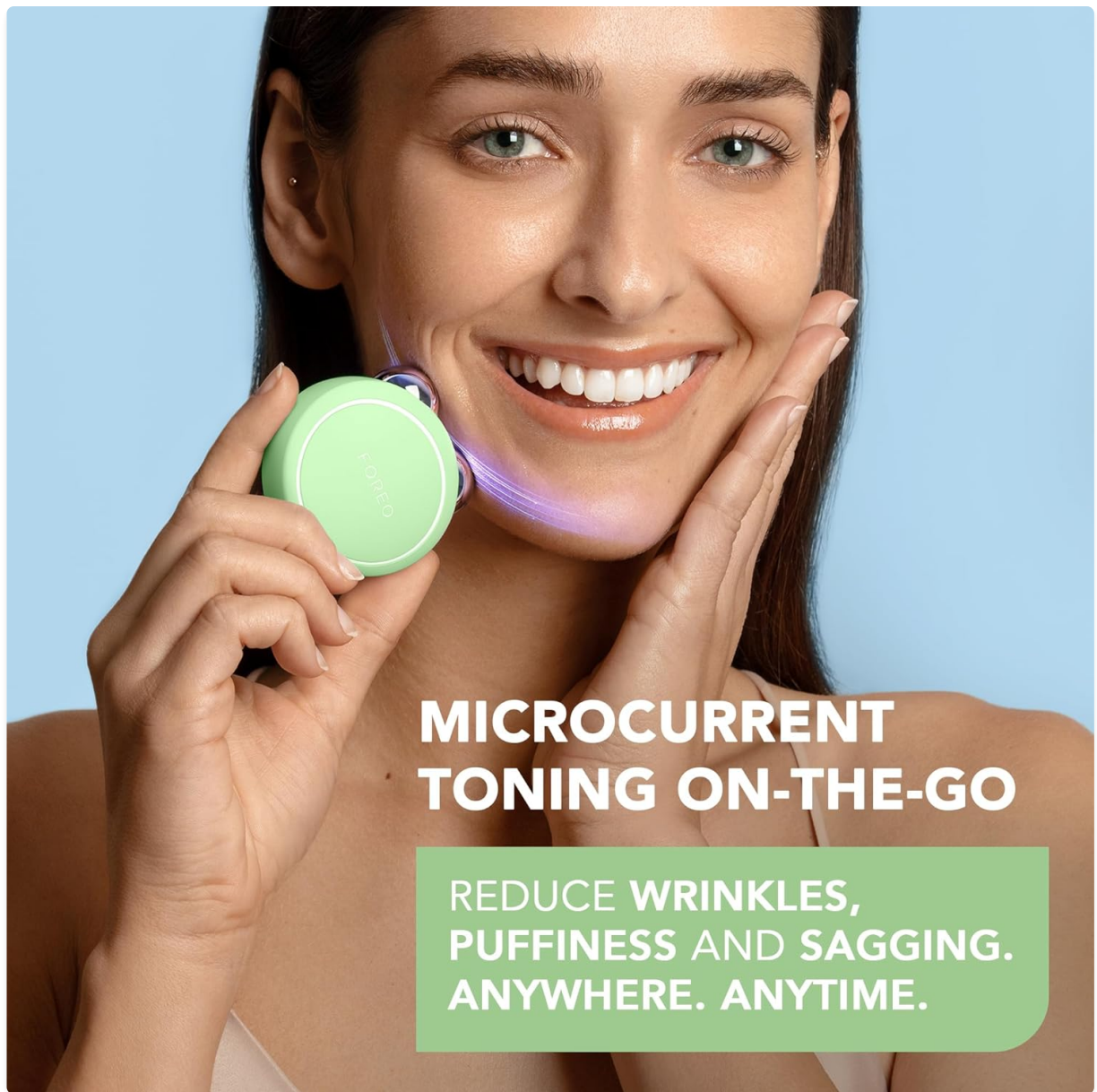


Image: A user demonstrating the application of the FOREO Bear 2 go device on the face.



BEAR go's MUST-HAVE SIDEKICK

SUPERCHARGED SERUM 2.0

Improves microcurrent transfer from the device to the skin,
ensuring you get the most out of each use.

Image: The FOREO Bear 2 go device paired with the recommended Supercharged Serum 2.0 for optimal conductivity.

WHAT'S IN THE BOX?

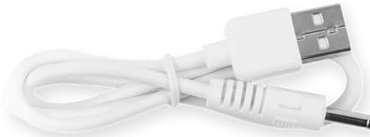
BEAR 2 Pistachio



Quick Start Guide



General Manual



USB Charging Cable

Image: Visual comparison of facial appearance before and after using the device, highlighting improvements in contour and reduction of puffiness.

5. MAINTENANCE

Proper maintenance ensures the longevity and hygiene of your FOREO Bear 2 go device.

- **Cleaning:** After each use, gently wipe the device, especially the microcurrent spheres, with a damp cloth and a mild, non-abrasive cleanser.
- **Drying:** Pat dry with a soft, lint-free cloth.
- **Storage:** Store the device in a clean, dry place away from direct sunlight and extreme temperatures.
- **Avoid Harsh Chemicals:** Do not use alcohol, petrol, or acetone-based cleaning products, as these may damage the silicone material.

6. TROUBLESHOOTING

If you experience issues with your FOREO Bear 2 go, refer to the following troubleshooting steps:

Device Not Activating (No response when pressing the center button):

- **Battery Empty:** The battery may be depleted. Recharge the device using the provided USB charging cable for up to 2 hours until fully charged. Then, attempt to restart the device by holding down the power button.
- **Device Malfunction:** If the device cannot be switched off and/or interface buttons do not respond, the microprocessor might be temporarily malfunctioning. Press and hold the power button to restart the device.

Device Not Syncing to FOREO App:

- **Bluetooth:** Ensure your device is fully charged and that Bluetooth is enabled on your mobile phone.
- **Reconnect:** Switch your phone's Bluetooth off and then on again to attempt reconnection.
- **Restart App:** Close the FOREO app completely and then reopen it to restart the connection process.
- **Reinstall App:** Delete and reinstall the application. This can resolve issues if the app is outdated or corrupted.
- **Restart Phone:** After deleting and reinstalling the application, restart your mobile phone.
- **Clear Cache:** Delete your phone's cookies and cache memory.
- **Software Update:** Ensure that your mobile's operating system has been updated to the latest version.

7. SPECIFICATIONS

Brand	FOREO
Model Name	BEAR 2 go
Item Model Number	F1849
Color	Pistachio
Material	Plastic
Product Dimensions	1.46 x 3.78 x 3.43 inches
Item Weight	2.12 ounces (0.06 Kilograms)
Batteries	1 AA battery required (included)
Manufacturer	FOREO

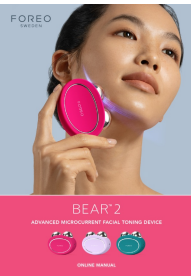
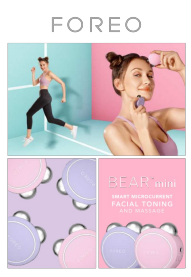
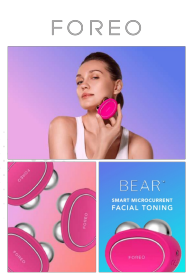
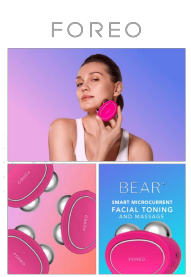
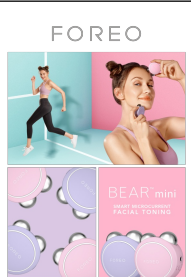
8. WARRANTY AND SUPPORT


For detailed warranty information, product registration, and further support, please refer to the General Manual included in your package or visit the official FOREO support website.

Official Support Website: www.foreo.com/support

Safety Information: No specific warnings are required for this product beyond general electronic device safety. Always follow the instructions provided in this manual and the FOREO For You app.

Related Documents - BEAR 2 go

	<p>FOREO BEAR™ 2 Advanced Microcurrent Facial Toning Device User Manual</p> <p>Comprehensive user manual for the FOREO BEAR™ 2, an advanced FDA-cleared microcurrent facial toning device. This guide covers setup, usage, types of microcurrent, safety warnings, cleaning, troubleshooting, and warranty information.</p>
	<p>FOREO BEAR™ mini: Smart Microcurrent Facial Toning Device User Manual</p> <p>User manual for the FOREO BEAR™ mini, a smart microcurrent facial toning and massage device. Learn about its features, usage, safety, and benefits for firmer, younger-looking skin.</p>
	<p>FOREO BEAR Smart Microcurrent Facial Toning Device User Manual</p> <p>Comprehensive user manual for the FOREO BEAR smart microcurrent facial toning device, covering getting started, features, usage, safety warnings, troubleshooting, warranty, and specifications.</p>
	<p>FOREO BEAR Smart Microcurrent Facial Toning and Massage User Manual</p> <p>Comprehensive user manual for the FOREO BEAR smart microcurrent facial toning and massage device, covering setup, usage, safety warnings, troubleshooting, and warranty information.</p>
	<p>FOREO BEAR™ mini User Manual: Smart Microcurrent Facial Toning Device</p> <p>Comprehensive user manual for the FOREO BEAR™ mini, detailing its features, usage instructions, safety warnings, troubleshooting, and warranty information for smart microcurrent facial toning.</p>

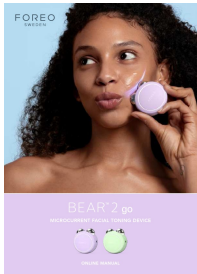

[FOREO BEAR 2 go Microcurrent Facial Toning Device User Manual](#)

Comprehensive user manual for the FOREO BEAR 2 go, a microcurrent facial toning device. Learn about features, usage, safety warnings, cleaning, troubleshooting, warranty, and disposal.