

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [Sejoy](#) /
- › [Sejoy Cordless Water Flosser Instruction Manual CYQ-F5037-BLA](#)

## Sejoy CYQ-F5037-BLA

# Sejoy Cordless Water Flosser Instruction Manual

Model: CYQ-F5037-BLA

## INTRODUCTION

---

Thank you for choosing the Sejoy Cordless Water Flosser. This device is designed to provide effective oral hygiene by removing plaque and food debris from between teeth and below the gumline, promoting healthier gums and fresher breath. Please read this manual thoroughly before first use to ensure safe and optimal operation.

## IMPORTANT SAFETY PRECAUTIONS

---

- Do not immerse the device in water for extended periods. The device is IPX7 waterproof, suitable for showering, but not for prolonged submersion.
- Keep out of reach of children.
- Do not direct water spray under the tongue, into the ear, nose, or other sensitive areas.
- Use only water or dentist-recommended mouthwash in the reservoir. Do not use corrosive or abrasive cleaning agents.
- Ensure the charging port cover is securely closed before use to maintain water resistance.
- If you have any medical conditions or concerns, consult your dentist before using this product.

## PACKAGE CONTENTS

---

Verify that all items are present in your package:

- 1 x Sejoy Cordless Water Flosser Unit
- 1 x Instruction Manual
- 1 x DC Charging Cable
- 5 x Standard Nozzles
- 1 x Orthodontic Nozzle
- 1 x Tongue Scraper
- 1 x Periodontal Pocket Nozzle


# Package Include :

1\* Instruction manual  
1\* Orthodontic Nozzle  
1\* Water Flosser

1\* DC Charging Cable  
1\* Tongue Scraper  
5\* standard nozzles  
1\* Periodontal Pocket Nozzle




Image: Package contents of the Sejoy Cordless Water Flosser, including the main unit, manual, charging cable, and various nozzles.

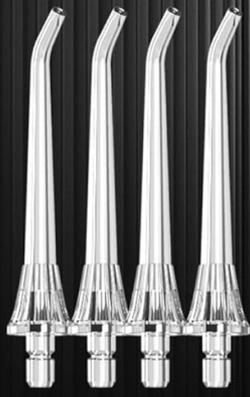

**Standard Nozzle \*5**  
Normal Cleaning



**Tongue Scraper \*1**  
Removes bacteria build up  
for fresher breath

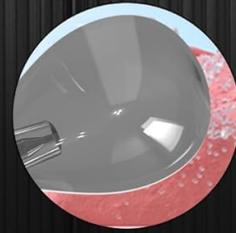

## Comes with **6** specialized Nozzles

Standard nozzles are recommended to be replaced every 6 months  
and the remaining nozzles are recommended  
to be replaced every 3 months



**Orthodontic Nozzle \*1**  
Braces Normal Cleaning



**Periodontal Pocket Nozzle \*1**  
Gum cleaning



Image: A visual representation of the 8 specialized nozzles provided with the water flosser, designed for different oral care needs.

### PRODUCT COMPONENTS

---

Familiarize yourself with the parts of your Sejoy Cordless Water Flosser:



Image: Labeled diagram illustrating the various components of the water flosser, such as the nozzle, control buttons, indicators, and charging port.

1. **Nozzle:** Directs the water stream.
2. **On/Off Button:** Powers the device on or off.
3. **Mode Button:** Cycles through different cleaning modes.
4. **Charging Indicator:** Shows battery charging status.
5. **Mode Indicator:** Displays the currently selected cleaning mode.
6. **DC Charging Port:** For connecting the charging cable.
7. **Nozzle Release Button:** Press to detach the nozzle.
8. **Water Tank:** Holds water for flossing (300ml capacity).

## SETUP

---

### 1. Charging the Device

Before first use, fully charge the water flosser. Connect the DC charging cable to the charging port on the device and plug the other end into a suitable USB power adapter (not included). The charging indicator will show the charging status. A full charge typically takes several hours.

## 2. Attaching a Nozzle

Select the desired nozzle. Insert the nozzle into the opening at the top of the flosser until it clicks into place. Ensure it is firmly seated.

## 3. Filling the Water Tank

Open the water tank cover located at the bottom or side of the device. Fill the 300ml tank with lukewarm water. Close the cover securely to prevent leakage.



Image: The water flosser with its 300ml water tank, demonstrating its capacity.

## OPERATING INSTRUCTIONS

### 1. Selecting a Mode

Press the Mode button to cycle through the available cleaning modes. The mode indicator will show the selected mode. The device offers the following modes:

- **Soft (70-110psi):** For gentle cleaning of sensitive gums.
- **Normal (25-60psi):** For regular cleaning and gum massage.
- **Shoot (20-110psi):** For targeted cleaning and gum massage.

- **Custom (60-110psi):** Allows customization of force cycle to find the right water pressure for sensitive gums or first-time users.



Image: Visual display of the four distinct operating modes and their respective pressure ranges.

## 2. Using the Water Flosser

1. Lean over a sink and place the nozzle tip into your mouth, aiming it at your gumline.
2. Press the On/Off button to start the water flow.
3. Guide the water stream along your gumline, pausing briefly between teeth.
4. Allow the water to flow into the sink.
5. Continue until you have cleaned all areas of your mouth.
6. Press the On/Off button again to turn off the device.

## MAINTENANCE

### 1. Cleaning the Water Tank

After each use, empty any remaining water from the tank. The tank can be detached and rinsed thoroughly with clean water. For deeper cleaning, use a soft brush and mild soap, then rinse completely.

# 10 Easy to Clean oz(300ML)

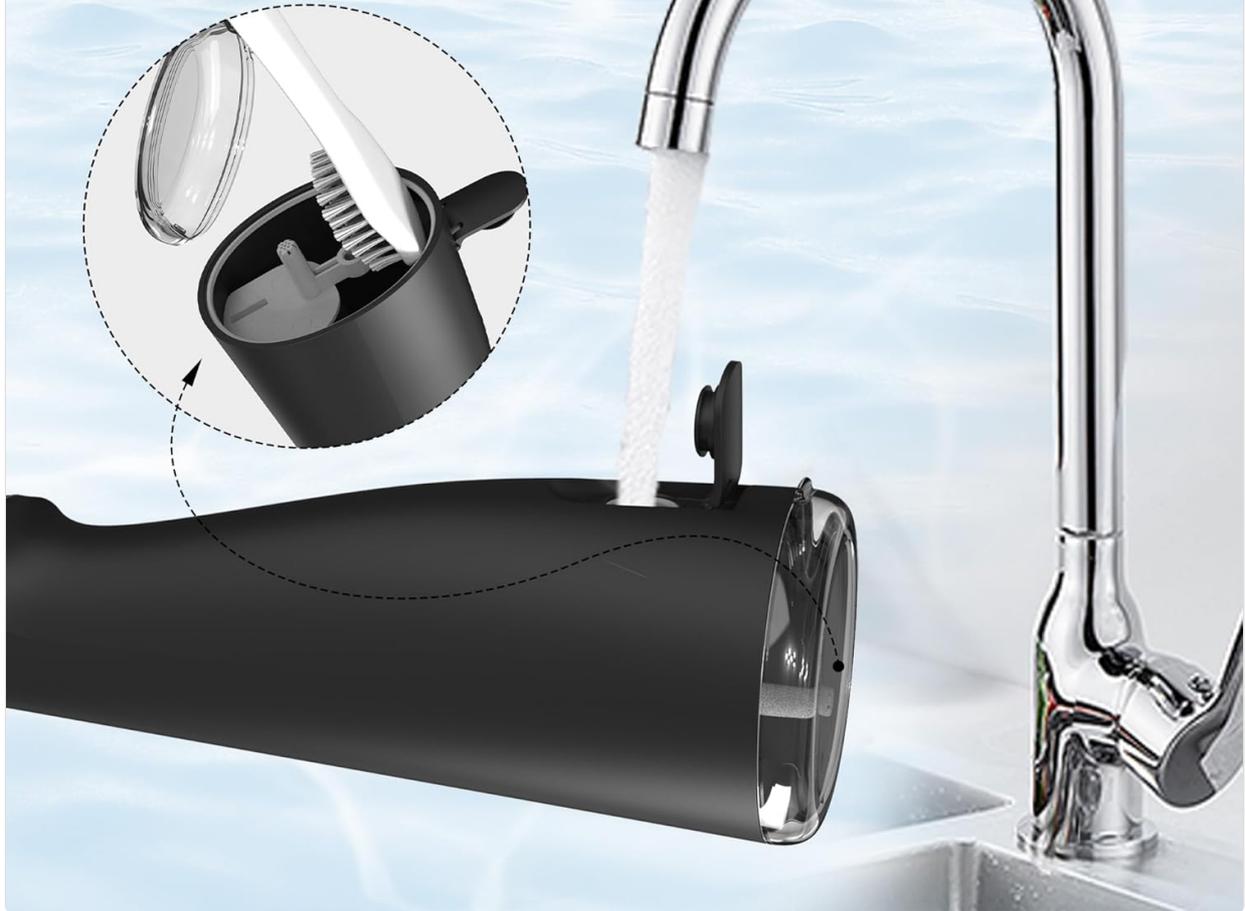


Image: Instructions for easy cleaning of the 300ml water tank under a faucet.

## 2. Cleaning the Nozzles

Rinse nozzles under running water after each use. For thorough cleaning, soak them in a solution of 2 parts water to 1 part white vinegar for 5-10 minutes, then rinse well.

## 3. Nozzle Replacement

Standard nozzles are recommended to be replaced every 6 months. Specialized nozzles (orthodontic, tongue scraper, periodontal pocket) are recommended to be replaced every 3 months for optimal hygiene and performance.

## 4. Storing the Device

Ensure the device is clean and dry before storing. Store it in a cool, dry place away from direct sunlight.

## TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery or no charge.	Charge the device fully.

Problem	Possible Cause	Solution
Weak water pressure.	Low battery, clogged nozzle, or incorrect mode.	Charge the device. Clean or replace the nozzle. Select a higher pressure mode.
Water leakage.	Water tank cover not closed properly or nozzle not seated correctly.	Ensure the water tank cover is securely closed. Re-insert the nozzle until it clicks.
Device stops during use.	Water tank empty or low battery.	Refill the water tank. Charge the device.

## SPECIFICATIONS

Feature	Detail
Brand	Sejoy
Model Number	CYQ-F5037-BLA
Power Source	Battery Powered
Special Feature	Portable
Product Benefits	Gum Health
Unit Count	1 Count
Color	Black
Tank Capacity	300 Milliliters
Manufacturer	JOYTECH Healthcare Co.,Ltd.
Dimensions	24.7 x 10.2 x 8.7 cm; 480 grams
GTIN (UPC)	06942284600047

## WARRANTY AND SUPPORT

For warranty information or technical support, please refer to the contact details provided on the product packaging or visit the official Sejoy website. Keep your purchase receipt as proof of purchase for any warranty claims.