

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [Aquasonic](#) /
- › [AquaSonic Proflex Replacement Brush Heads User Manual](#)

Aquasonic Proflex

AquaSonic Proflex Replacement Brush Heads User Manual

Model: Proflex | Brand: AquaSonic



INTRODUCTION

This manual provides instructions for the AquaSonic Proflex Replacement Brush Heads. These brush heads are designed to deliver optimal cleaning performance for your AquaSonic electric toothbrush, promoting whiter teeth and improved gum care. They are ADA Accepted, ensuring high quality and effectiveness.

The AquaSonic Proflex brush heads feature a unique bristle design: softer bristles around the perimeter for gentle gum care, and firmer bristles at the center to effectively polish away surface stains and remove plaque and tartar.

COMPATIBILITY

AquaSonic Proflex Replacement Brush Heads are compatible with the following AquaSonic electric toothbrush handles:

- Black Series
- Black Series Pro
- Vibe Series
- Vibe Series Pro
- HDC Pro
- Duo Series Pro
- Elite



This image illustrates the AquaSonic Proflex brush head and lists the compatible AquaSonic toothbrush models: Black Series, Black Series Pro, Vibe Series, Vibe Series Pro, HDC Pro, Duo Series Pro, and Elite.

SETUP AND INSTALLATION

To install a new AquaSonic Proflex Replacement Brush Head:

1. Ensure your AquaSonic toothbrush handle is turned off.
2. Gently pull the old brush head straight off the metal shaft of the toothbrush handle.
3. Align the new Proflex replacement brush head with the metal shaft.
4. Push the new brush head firmly onto the shaft until it clicks securely into place.



The AquaSonic Proflex Replacement Brush Heads are packaged with protective caps, ready for installation on your compatible AquaSonic toothbrush handle.

OPERATING INSTRUCTIONS

Once the brush head is securely attached, you can operate your AquaSonic toothbrush as usual. The Proflex brush heads are designed to work with the various cleaning modes of your AquaSonic toothbrush, such as Clean, Soft, Whiten, and Massage, to provide a comprehensive oral care experience.

For optimal cleaning, apply a pea-sized amount of toothpaste to the bristles. Place the brush head against your teeth at a slight angle towards the gum line. Turn on the toothbrush and guide the brush head slowly across all surfaces of your teeth, allowing the sonic vibrations to do the work. Do not scrub vigorously.

ProFlex

Replacement Brush Heads

Optimized for Whiter
Teeth & Gum Care




AquaSonic

The AquaSonic ProFlex brush heads are engineered for effective plaque removal and gum care, contributing to whiter teeth.

ProFlex

Replacement Brush Heads

soft bristles around the perimeter are gentle on gums, firm bristles at the center are tough on plaque

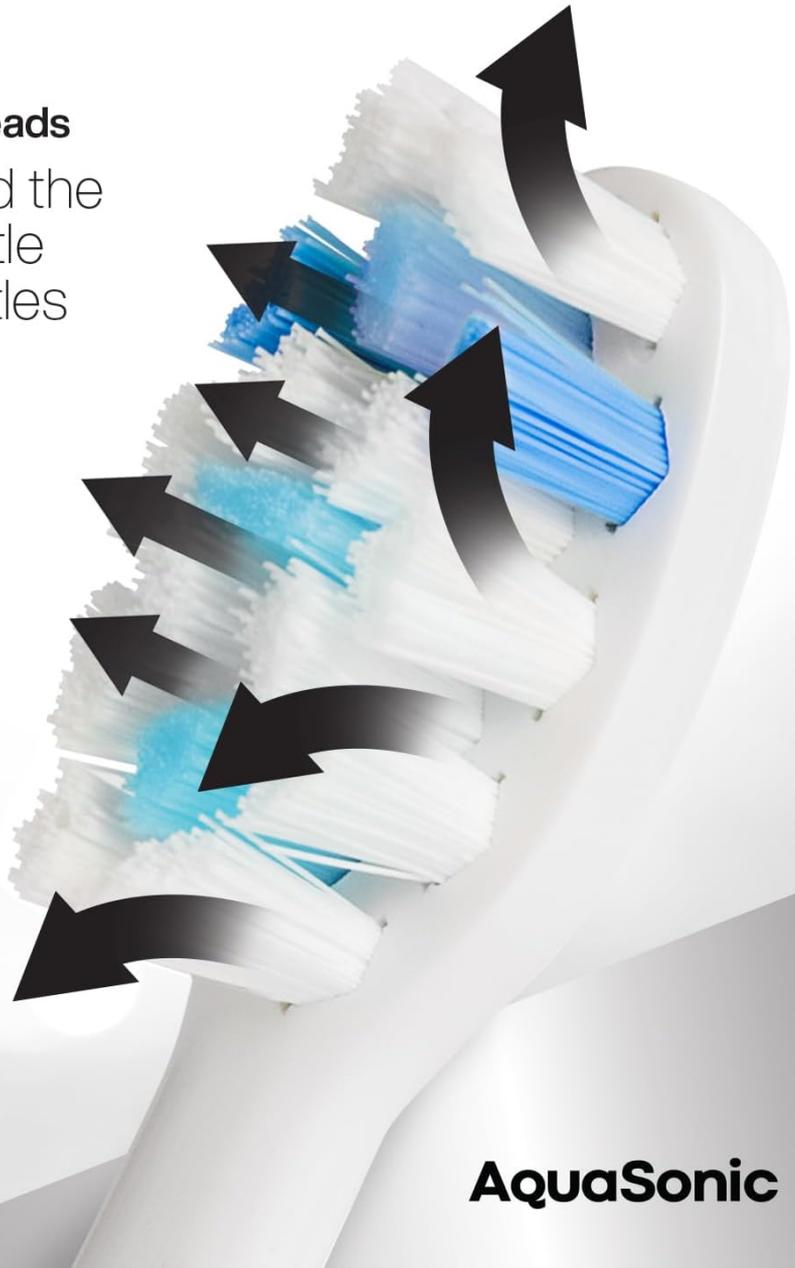

This image illustrates the dual-bristle design of the Proflex head: soft bristles on the outer edge for gum comfort and firm bristles in the center for plaque removal.

MAINTENANCE AND REPLACEMENT

For optimal brushing performance and hygiene, it is recommended to replace your AquaSonic Proflex brush head every 3 months. You should also replace the brush head sooner if the bristles become worn or frayed. After each use, rinse the brush head thoroughly under running water to remove any toothpaste residue and debris. Allow it to air dry.

Recommended

Change your brush head every **3-4 months** for healthier teeth and gums, and fresher breath



Regular replacement of your brush head, ideally every 3-4 months, is crucial for maintaining oral hygiene and effective cleaning.

SAFETY INFORMATION

Do not use these brush heads on sensitive areas or by children under 3 years of age. Always supervise children when they are brushing their teeth.

SPECIFICATIONS

Feature	Detail
Brand	Aquasonic
Model	Proflex Replacement Brush Heads
Color	White (8 Pack)
Material	Nylon (bristles)
Age Range (Description)	Adult

Package Dimensions	8.78 x 4.88 x 0.83 inches
Item Weight	3.53 ounces

WARRANTY AND SUPPORT

For any questions regarding the AquaSonic Proflex Replacement Brush Heads or for support with your AquaSonic toothbrush, please refer to the original toothbrush manual or contact AquaSonic customer service. AquaSonic is a USA-based company dedicated to customer satisfaction.

