

## Hyperice HH0001

# Hyperice Hypervolt Heated Head Attachment User Manual

Model: HH0001

## 1. INTRODUCTION

---

The Hyperice Hypervolt Heated Head Attachment is designed to enhance your percussion massage experience. This attachment is compatible with all Hypervolt models and helps to relax tight muscles and maintain range of motion by combining heat with percussion therapy.

# Details and specifications

## Weight

0.2 lbs

## Size

2.8 x 2.4 in

## Battery life

75 minutes

## Heat

3 levels of heat. Level 1, 109°F (43°C) /  
Level 2, 115°F (46°C) / Level 3, 120°F  
(49°C)



Image 1.1: The Hyperice Hypervolt Heated Head Attachment.

## 2. SAFETY INFORMATION


Please read all safety warnings and instructions before using this product. Failure to follow these instructions may result in injury or damage to the device.

- Use only as directed.
- Do not use on open wounds or broken skin.
- Consult a physician before use if you have a medical condition or are taking any medications.
- Stop use if you experience pain or discomfort.
- Keep out of reach of children. This product is not intended for use by children.
- Do not use while sleeping or in the bathtub.
- Unplug the device when not in use.
- Do not use a damaged power cord or plug.
- Do not submerge in water or expose to rain or wet conditions.
- Do not use with oils, lotions, or creams.

- 
- This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.
- FCC
- UKCA
- CE
- NOM
- SGS
- 801404
- HECHO EN: CHINA
- Model/Modelo/모델/HH-000T
- Especificaciones técnicas/특성규격: 5V=2.10A
- 제조업체: MMY
- Hyperice, Inc.  
525 Technology Dr, #100  
Irvine, CA 92618 USA  
170-00349 Rev 3
- Masajeador/加熱塗抹器  
Marca/인증 받은 자의 상호; Hyperice

### 3. WHAT'S INCLUDED

- One (1) Hyperice Hypervolt Heated Head Attachment
- One (1) USB-C Charging Cable

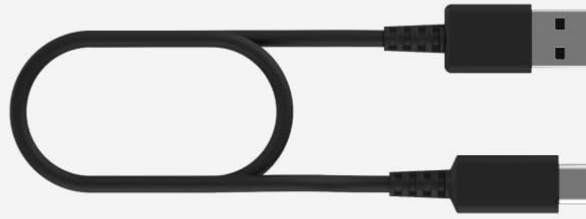

## What's included with your Heated Head Attachment

### Heated Head Attachment

The actual attachment

### Charger

USB-C charging cable to keep your device charged

Image 3.1: Contents of the package: Heated Head Attachment and USB-C cable.

## 4. SETUP AND CHARGING

### 4.1 Attaching to Your Hypervolt Device

1. Ensure your Hypervolt massager is powered off and compatible with the Heated Head Attachment.
2. Gently push and twist the Heated Head Attachment onto the Hypervolt massager's head slot until securely seated.



## Heated Head Attachment

The ultimate addition to your percussion massage experience. Compatible with all Hypervolt models, the Heated Head Attachment helps to relax tight muscles and maintain your full range of motion, providing an even more effective massage.

Image 4.1: The Heated Head Attachment being connected to a Hypervolt device.

### 4.2 Charging the Heated Head Attachment

The Heated Head Attachment has an internal battery that powers its heating function. Charge the attachment using the provided USB-C cable.

- Connect the USB-C cable to the charging port on the Heated Head Attachment.
- Connect the other end of the USB-C cable to a compatible USB power adapter (not included).
- A full charge provides approximately 75 minutes of heated operation.

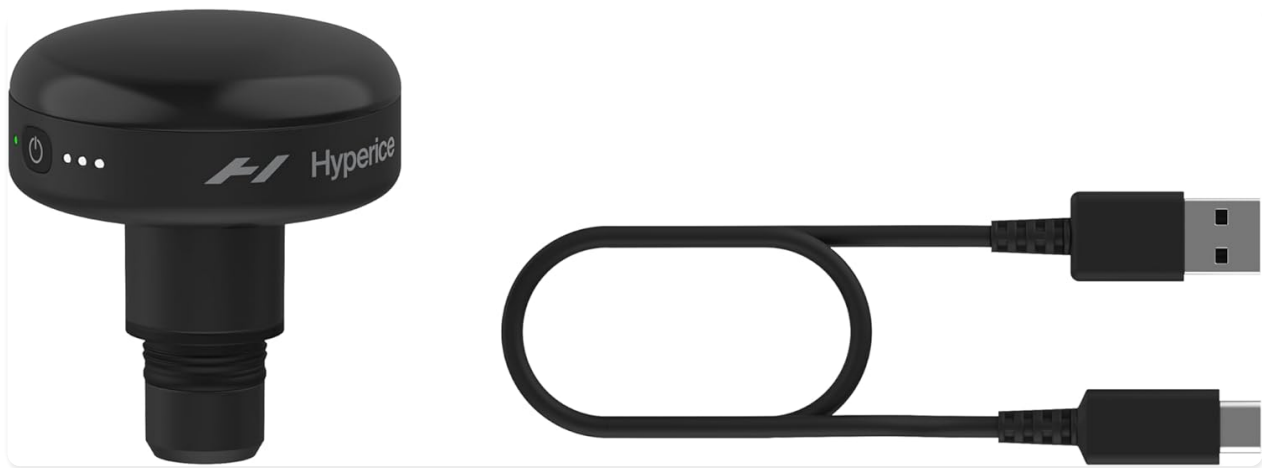


Image 4.2: The Heated Head Attachment with its USB-C charging cable.

## 5. OPERATING INSTRUCTIONS

---

Once the Heated Head Attachment is securely attached and charged, you can begin your heated percussion massage.

1. Power on your Hypervolt massager.
2. Press the single button on the Heated Head Attachment to cycle through the three heat levels:
  - **Level 1:** 109°F (43°C)
  - **Level 2:** 115°F (46°C)
  - **Level 3:** 120°F (49°C)
3. Apply the massager with the Heated Head Attachment to your desired muscle group for relaxation, pain relief, and increased blood flow.
4. To turn off the heat, cycle through the heat levels until the indicator lights are off, or power off the Hypervolt device.



## The perfect combination

The Heated Head Attachment has three levels of heat, comes with a USB-C charger and provides up to 75 minutes of heated bliss per charge.

Image 5.1: Using the Heated Head Attachment for neck massage.





## Upgrade your massage

The Heated Head Attachment has three levels of heat, comes with a USB-C charger and provides up to 75 minutes of heated bliss per charge.

Image 5.2: Applying the Heated Head Attachment to a thigh muscle.

## 6. MAINTENANCE

---

To ensure the longevity and optimal performance of your Heated Head Attachment, follow these maintenance guidelines:

- **Cleaning:** Wipe the attachment with a damp cloth after each use. Do not use harsh chemicals or abrasive cleaners. Ensure the device is completely dry before storing or recharging.
- **Storage:** Store the attachment in a cool, dry place away from direct sunlight and extreme temperatures.
- **Battery Care:** For best battery life, avoid fully discharging the battery frequently. Charge the device regularly, even if not in use for extended periods.

## 7. TROUBLESHOOTING

---

If you encounter any issues with your Hyperice Hypervolt Heated Head Attachment, consider the following:

- **No Heat:** Ensure the attachment is fully charged. Verify that the heat setting is activated by pressing the button.



- **Attachment Not Secure:** Re-attach the head, ensuring it is pushed and twisted firmly into place on your Hypervolt device.
- **Charging Issues:** Check the USB-C cable for damage and ensure it is properly connected to both the attachment and a working power source.

For further assistance, please contact Hyperice customer support.

## 8. SPECIFICATIONS

Feature	Specification
Model Number	HH0001
Product Dimensions	2.89 x 4.09 x 3.29 inches
Item Weight	0.4 Pounds (6.4 ounces)
Material	Plastic
Power Source	Battery Powered
Battery Type	1 Lithium Polymer battery (included)
Battery Life	Up to 75 minutes
Heat Levels	3 levels: 109°F (43°C), 115°F (46°C), 120°F (49°C)
Compatibility	All Hypervolt models
Manufacturer	Hyperice



Image 8.1: Detailed view of the Heated Head Attachment.






## 9. WARRANTY AND SUPPORT

For warranty information and customer support, please visit the official Hyperice website or contact their customer service department. Keep your proof of purchase for warranty claims.

**Manufacturer:** Hyperice

**Website:** [www.hyperice.com](http://www.hyperice.com)

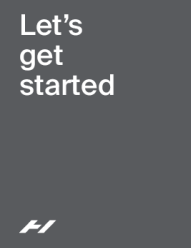Related Documents - HH0001

 <p><b>Hypervolt Go 2</b></p>	<p><a href="#">Hypervolt Go 2: User Manual &amp; Safety Guide   Hyperice</a></p> <p>Discover how to use your Hypervolt Go 2 percussion massager safely and effectively. This guide covers operation, charging, maintenance, and important safety precautions from Hyperice.</p>
 <p><b>Hypervolt 2 Pro</b></p>	<p><a href="#">Hypervolt 2 Pro Operating Instructions and Safety Information</a></p> <p>Comprehensive operating instructions, safety warnings, specifications, and warranty information for the Hyperice Hypervolt 2 Pro percussion massage device.</p>
 <p><b>Let's get started</b></p>	<p><a href="#">Hyperice Hypervolt 2 Quick Start Guide: Setup, Usage, and Care</a></p> <p>Your essential guide to the Hyperice Hypervolt 2. Learn how to set up, charge, use, and care for your percussion massage device, plus discover the benefits of the Hyperice App.</p>
 <p><b>Hypervolt Go 2</b></p>	<p><a href="#">Hypervolt Go 2 Percussion Massage Device: Operating Instructions and Safety Guide</a></p> <p>Comprehensive guide to the Hypervolt Go 2 handheld percussion massage device by Hyperice, covering operating instructions, safety warnings, charging, cleaning, specifications, and warranty information.</p>
 <p><b>Hypervolt GO</b></p>	<p><a href="#">Hypervolt GO User Manual and Instructions</a></p> <p>Discover the Hypervolt GO, a handheld percussion massage device designed to relieve muscle soreness and stiffness, accelerating warm-up and recovery post-activity. This manual provides operating instructions, safety guidelines, specifications, and care information.</p>



### [Hypervolt: High-Performance Vibration Massage Device | Hyperice](#)

Discover the Hypervolt, a state-of-the-art cordless vibration massage device by Hyperice. Learn about its features, specifications, usage instructions, and benefits for muscle recovery, athletic performance, and overall wellness.