

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [COULAX](#) /

› [COULAX Replacement Electric Toothbrush Head Instruction Manual](#)

## COULAX COU007 M304

# COULAX Replacement Electric Toothbrush Head Instruction Manual

Model: COU007 M304 Compatible Heads

## INTRODUCTION

---

This manual provides essential information for the proper use and maintenance of your COULAX Replacement Electric Toothbrush Heads. These heads are designed for compatibility with COULAX electric toothbrushes, specifically models COU007 and M304, ensuring optimal oral hygiene. Each pack contains 4 replacement brush heads, featuring U-shape bristles and soft, flexible DuPont bristles for effective cleaning and gum protection.

## IMPORTANT SAFETY INFORMATION

---

Please read all instructions carefully before using the replacement brush heads. Keep this manual for future reference.


- Ensure the brush head is securely attached to the toothbrush handle before use.
- Do not use if the brush head appears damaged or worn. Replace immediately.
- Supervise children aged 8 years and above when using the toothbrush.
- Consult a dentist if you experience excessive bleeding from the gums or if bleeding persists after 2 weeks of use.
- Do not attempt to modify the brush head.

## SETUP AND INSTALLATION

---

Follow these steps to properly install the replacement brush head onto your COULAX electric toothbrush handle:

1. Ensure your toothbrush handle is turned off.
2. Align the replacement brush head with the metal shaft on the toothbrush handle.
3. Push the brush head down firmly until it clicks into place. A secure fit is essential for proper operation.
4. To remove, pull the brush head straight up and off the handle.

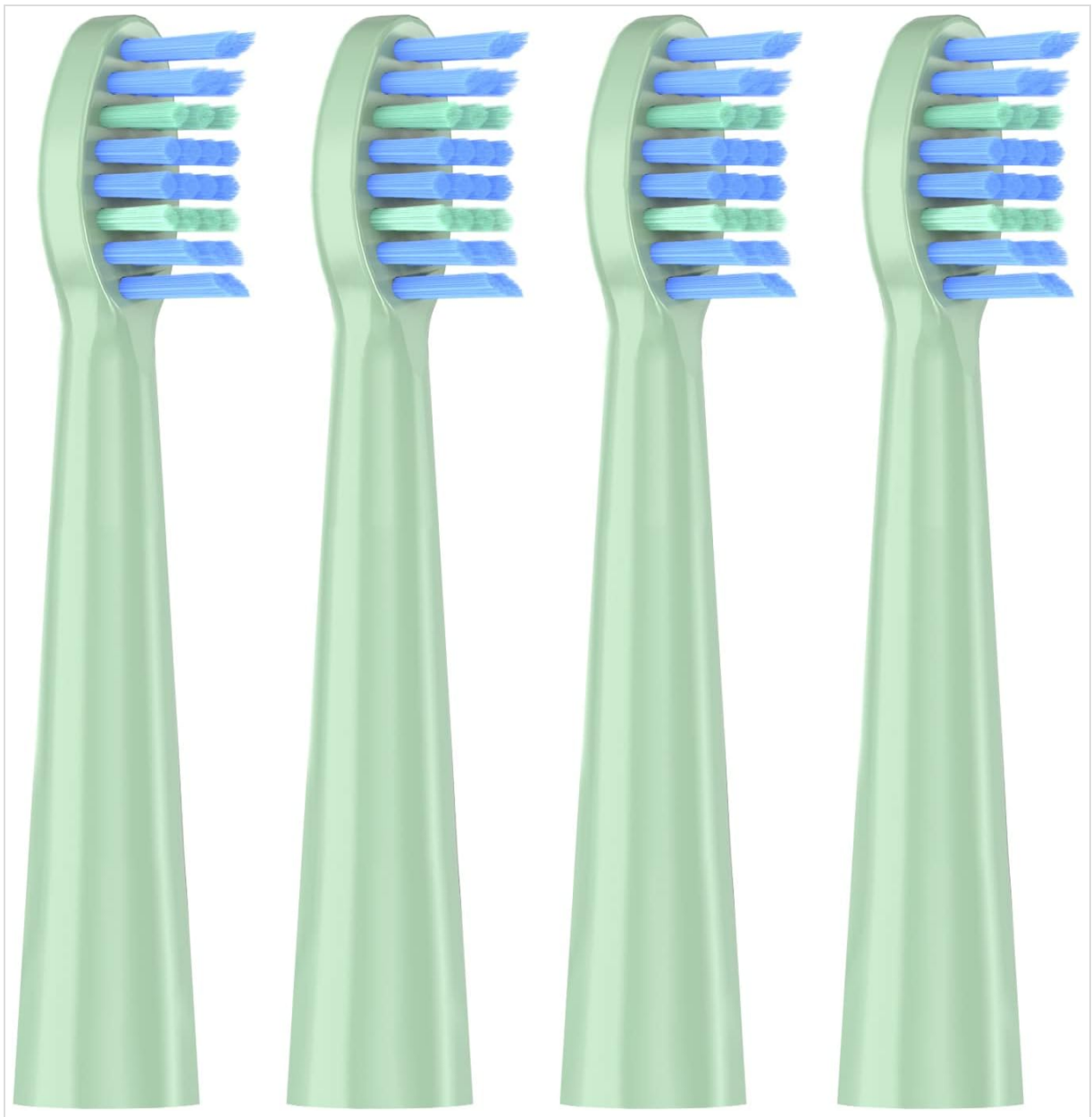

**Figure 1:** Four COULAX replacement electric toothbrush heads. Each head features a light green base and a brush head with alternating rows of blue and light green U-shaped bristles, designed for effective cleaning.

## OPERATING INSTRUCTIONS

---

Once the replacement head is installed, operate your electric toothbrush as per its original instructions. The COULAX replacement heads are designed to integrate seamlessly with your COU007 or M304 model.

- Wet the bristles and apply a small amount of toothpaste.
- Place the brush head against your teeth at a slight angle towards the gumline.
- Turn on the toothbrush and guide the brush head slowly from tooth to tooth, allowing the bristles to clean each surface.
- Do not scrub vigorously; let the electric toothbrush do the work.

## MAINTENANCE AND CARE

---

Proper maintenance extends the life of your brush heads and ensures optimal hygiene.

### Cleaning the Brush Head

- After each use, rinse the brush head thoroughly under running water while the toothbrush is still on (briefly) to remove any remaining toothpaste or debris.
- Turn off the toothbrush and remove the brush head. Rinse the brush head and the metal shaft of the handle separately.
- Allow the brush head to air dry in an upright position.

## Replacement Frequency

It is recommended to replace your brush head every 3 months, or sooner if the bristles become frayed or worn. Worn bristles are less effective at cleaning and can potentially harm gums.

## TROUBLESHOOTING

Problem	Possible Cause	Solution
Brush head does not fit securely.	Incorrect alignment or not pushed down firmly enough.	Ensure proper alignment and push until it clicks. Verify compatibility with your COULAX model.
Bristles appear worn quickly.	Excessive pressure during brushing.	Reduce pressure. Let the electric toothbrush do the work. Replace the brush head.
Cleaning effectiveness seems reduced.	Worn bristles or improper technique.	Replace the brush head. Review operating instructions for proper technique.

## SPECIFICATIONS

**Product Name:** COULAX Replacement Electric Toothbrush Head

**Model Compatibility:** COU007, M304 COULAX basic electric toothbrushes

**Quantity:** 4 pieces per pack

**Bristle Material:** DuPont Nylon

**Bristle Firmness:** Soft

**Bristle Design:** U-shape

**Color:** Light Green (base)

**Age Range:** 8 years old +

**Manufacturer:** Ningbo Jianjie Daily Chemical Co., Ltd

**Dimensions (Package):** 10 x 9 x 1.5 cm

**Weight (Package):** 16 g

## WARRANTY AND SUPPORT

For information regarding warranty coverage or technical support for your COULAX electric toothbrush or replacement heads, please refer to the documentation provided with your original toothbrush purchase or contact COULAX customer service directly.

Ensure you retain your proof of purchase for any warranty claims.

