

Karite B0CHJRPYC4

KARITE Paraffin Wax Machine Instruction Manual

Model: B0CHJRPYC4

INTRODUCTION

This manual provides detailed instructions for the safe and effective operation, maintenance, and troubleshooting of your KARITE Paraffin Wax Machine. Please read this manual thoroughly before first use and retain it for future reference. The KARITE Paraffin Wax Machine is designed for paraffin wax therapy for hands, feet, elbows, forearms, and knees. It features a 9000ml capacity, durable stainless steel inner pot, and intelligent heating modes including manual, smart, and time reservation.

IMPORTANT SAFETY PRECAUTIONS

- Always unplug the machine when not in use, before cleaning, or before adding/removing wax.
- Do not operate the machine if the power cord or plug is damaged.
- Keep the machine away from water or other liquids. Do not immerse the unit in water.
- Ensure the machine is placed on a stable, heat-resistant surface.
- Do not touch the melted wax directly with bare hands until its temperature has been checked.
- Always use the provided heat insulating grid to prevent direct contact with the hot bottom of the pot.
- Keep out of reach of children and pets.
- Consult a physician before use if you have sensitive skin, diabetes, circulatory problems, open wounds, or any skin conditions.
- Do not use the machine for purposes other than those described in this manual.

PACKAGE CONTENTS

Verify that all items are present in the package:

- KARITE Paraffin Wax Machine (1 unit)
- Lavender Paraffin Wax Blocks (16 pieces)
- Thermal Mitts (2 pieces)


- Thermal Foot Covers (2 pieces)
- Plastic Gloves for Hands (30 pieces)
- Plastic Gloves for Feet (30 pieces)
- Brush and Silicone Scraper (1 set)
- Heat Insulating Grid (1 piece)
- Instruction Manual (1 piece)



Image: Overview of the KARITE Paraffin Wax Machine and its included accessories, such as wax blocks, thermal mitts, and plastic gloves.

SETUP INSTRUCTIONS

1. **Unpack the Machine:** Carefully remove all components from the packaging. Place the paraffin wax machine on a flat, stable, and heat-resistant surface.
2. **Insert Heat Insulating Grid:** Place the heat insulating grid at the bottom of the inner pot. This prevents direct contact with the heating element.
3. **Add Paraffin Wax:** Place the desired number of paraffin wax blocks into the inner pot. The machine has a 9000ml capacity, allowing for a substantial amount of wax. Ensure the wax level does not exceed the maximum fill line.
4. **Connect Power:** Plug the power cord into a standard electrical outlet.

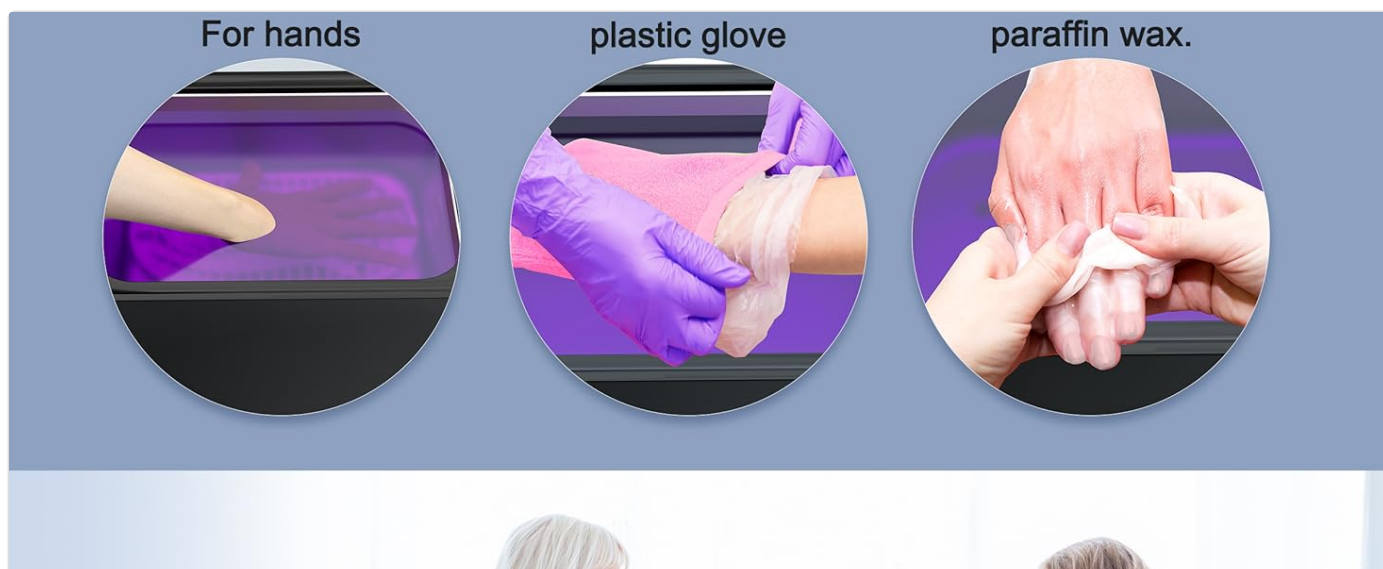

OPERATING INSTRUCTIONS

The KARITE Paraffin Wax Machine offers three intelligent heating modes:

1. Manual Heating Mode

- Press the power button to turn on the machine.
- Use the temperature control buttons (+ / -) to set the desired temperature between 86°F (30°C) and 167°F (75°C).
- The machine will begin heating the wax. Monitor the melting process.
- Once the wax is fully melted and reaches the desired temperature, it is ready for use. Always test a small amount of wax on your wrist before immersing a body part.

2. Smart Mode (One-Touch Operation)

- Press the "Smart Mode" button. The machine will automatically heat the wax to an optimal temperature for therapy.
- This mode is designed for convenience, eliminating the need for manual temperature adjustments.
- Wait for the wax to fully melt and reach the ready temperature.

3. Time Reservation Mode

- Press the "Reservation" button.
- Use the time adjustment buttons to set a desired heating start time (e.g., 3-99 hours in advance).
- The machine will automatically begin heating the wax at the scheduled time, ensuring it is ready when needed.



15 pounds of wax

How big a wax melter is needed ?

Others

Five needed



OURS

Only one is needed

9000ml/305oz



Image: Close-up of the control panel, illustrating the buttons for Smart Mode, Reservation Mode, and manual temperature adjustment.

Application Steps for Paraffin Wax Therapy:

1. **Prepare Skin:** Ensure the area to be treated (hands, feet, elbows) is clean and dry.
2. **Test Wax Temperature:** Before immersing, test the wax temperature by dipping a small area of your wrist or finger. The wax should feel warm and comfortable, not uncomfortably hot.
3. **Dip Body Part:** Slowly dip the hand, foot, or elbow into the melted wax. Remove it immediately and allow the wax to solidify for a few seconds. Repeat this process 3-5 times to build up a thick layer of wax.
4. **Cover with Plastic and Thermal Mitts:** Immediately after the final dip, cover the waxed area with a plastic glove/bag (provided) and then with a thermal mitt or foot cover (provided) to retain heat.
5. **Relax:** Allow the wax to remain on for 10-15 minutes. This allows the heat to penetrate and provide therapeutic benefits.
6. **Remove Wax:** Gently peel off the solidified wax. The wax should come off easily, taking dead skin cells with it. Dispose of used wax. Do not reuse wax that has been removed from the body.



Image: Visual guide demonstrating the application process of paraffin wax, including dipping a hand, covering with a plastic glove, and removing the solidified wax.

MAINTENANCE AND CLEANING

- **Cleaning the Machine:**
 - Always unplug the machine and allow it to cool completely before cleaning.
 - To remove residual wax, turn the machine on for a short period to slightly soften the wax, then unplug and carefully scrape out the wax using the provided silicone scraper or a wooden spatula.
 - Wipe the interior and exterior surfaces with a soft, damp cloth. Do not use abrasive cleaners or immerse the unit in water.
- **Wax Replacement:** It is recommended to replace the paraffin wax after approximately 20-30 uses, or sooner if it appears dirty or contaminated.
- **Storage:** Store the machine in a cool, dry place when not in use. Ensure the lid is securely closed to prevent dust accumulation.

TROUBLESHOOTING GUIDE

Problem	Possible Cause	Solution
Machine does not turn on.	Not plugged in; power outlet issue; damaged cord.	Ensure power cord is securely plugged in. Test outlet with another device. Inspect cord for damage; if damaged, discontinue use and contact support.
Wax not melting or melting slowly.	Insufficient heating time; too much wax; low temperature setting.	Allow more time for heating. Ensure the machine is set to an appropriate temperature (e.g., Smart Mode or higher manual setting). Reduce the amount of wax if overfilled.
Wax is too hot.	Temperature set too high.	Reduce the temperature setting. Allow wax to cool before use. Always test wax temperature before application.
Wax is too cold/hard.	Temperature set too low; not enough heating time.	Increase the temperature setting or allow more heating time.

PRODUCT SPECIFICATIONS

- **Model:** B0CHJRPYC4
- **Capacity:** 9000ml / 305oz
- **Temperature Range:** 86°F - 167°F (30°C - 75°C)
- **Inner Pot Material:** Stainless Steel



Image: Diagram illustrating the large capacity (9000ml) and dimensions of the KARITE Paraffin Wax Machine.

WARRANTY AND CUSTOMER SUPPORT

Every product purchased from KARITE comes with a 24-month warranty. If you encounter any questions or issues with your KARITE Paraffin Wax Machine, please contact our customer service for assistance.

Customer Service Contact: Please refer to the contact information provided in your product packaging or on the official KARITE website.

