

INIA D-1198

INIA Fond IPL Laser Hair Removal Device User Manual

Model: D-1198 | Brand: INIA

1. INTRODUCTION

The INIA Fond IPL Laser Hair Removal Device utilizes advanced Intense Pulsed Light (IPL) technology to target hair follicles, aiming for long-lasting hair reduction. Designed for safe and effective at-home use, this device offers a convenient alternative for managing unwanted hair on various body and facial areas. Regular use can lead to a significant reduction in hair growth, promoting smoother skin.

2. SAFETY INFORMATION

Your safety is paramount. Please read and understand all safety precautions before using the INIA Fond device.

- **Not Suitable For:** Tattoos, skin diseases, sunburn, scalp, eyebrows, or injured skin.
- **Skin Preparation:** Ensure your skin is clean and free of makeup, lotions, or perfumes before treatment. Shave the treatment area prior to use.
- **Eye Protection:** Always wear the provided protective goggles during treatment to shield your eyes from the intense light.
- **Skin Contact Sensor:** The device is equipped with a skin contact sensor to ensure proper application and prevent accidental flashes.
- **Sun Exposure:** Avoid excessive sun exposure and tanning (including sunbeds) for at least one week before and after treatment. Sunbathing can cause skin darkening or sunburn, leading to redness, swelling, or inflammation. Always apply sunscreen to treated areas if unavoidable.
- **Water Contact:** Avoid contact with water for at least 6 hours after hair removal.
- **Post-Treatment Care:** After 6 hours, moisturize the treated area with a hydrating lotion. Avoid using perfumes, face masks, and cosmetics with exfoliating beads on the treated skin.
- **Medical Conditions:** If any discomfort occurs, consult a dermatologist promptly. Consult a doctor if you plan to undergo any plastic surgery procedures in the treated area after hair removal.
- **Contraindications:** The device is not suitable for individuals who are pregnant or breastfeeding.

3. PACKAGE CONTENTS

The INIA Fond IPL Laser Hair Removal Device package includes the following items:

- INIA Fond IPL Hair Removal Device
- Power Adapter
- Protective Goggles
- Razor
- Instruction Manual
- Quick Start Guide
- Warranty Card



Image: The INIA Fond IPL Laser Hair Removal Device, protective goggles, and packaging.

4. DEVICE OVERVIEW

The INIA Fond device features a sleek, ergonomic design for comfortable handling. It includes a flash window, power button, mode button, and indicator lights for intensity levels and treatment modes.

TIKTOK #1 BEST-SELLING

Hair Removal Device



Image: Close-up view of the INIA Fond device, highlighting the control panel with mode and intensity indicators.

5. SKIN TONE & HAIR COLOR SUITABILITY


The effectiveness and safety of IPL treatment depend on your skin tone and hair color. Please refer to the chart below to determine if the INIA Fond device is suitable for you.

INIA

Skin Tone & Body Hair Type Chart

Please review the charts below to ensure your eligibility for safe & effective IPL use.

Suitable Skin Tones

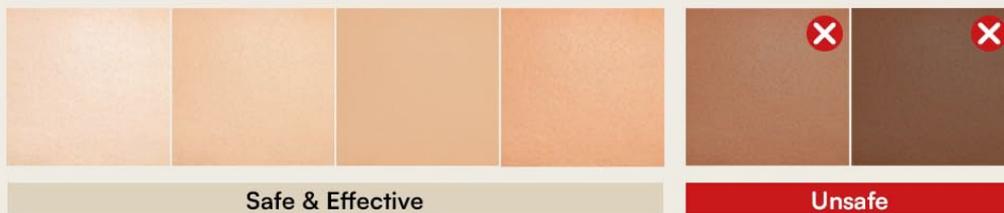

Suitable Body Hair Colors



Image: A chart illustrating suitable skin tones (lighter tones are safe and effective, darker tones require caution or are not applicable) and suitable hair colors (darker hair colors are effective, lighter colors like red, blonde, and gray are ineffective) for IPL treatment.

General Guidelines:

- IPL is most effective on lighter skin tones with darker hair.
- Darker skin tones may absorb too much light energy, leading to skin irritation or burns.
- Very light hair colors (red, gray, white, very light blonde) may not contain enough melanin for the IPL to effectively target the follicle.

6. SETUP

Before your first treatment, follow these steps to prepare the device and your skin:

1. **Prepare and Shave:** Ensure the treatment area is clean and dry. Shave the hair in the area you intend to treat. Do not wax or pluck, as the hair root needs to be present for IPL to be effective.
2. **Wear Goggles:** Put on the provided protective goggles to shield your eyes from the intense light flashes.
3. **Power Connection:** Safely plug the device into a power outlet using the provided adapter. The device is designed for corded use.

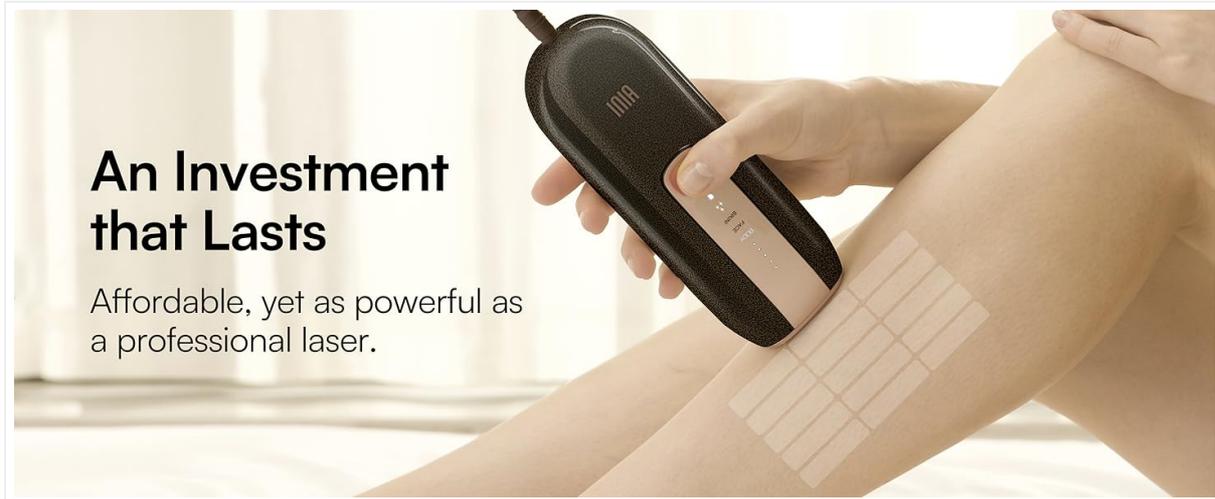


Image: Visual guide for device setup, showing shaving, wearing goggles, plugging in the device, and initiating hair removal.

7. OPERATING INSTRUCTIONS

The INIA Fond device offers customized modes and adjustable intensity levels for various body parts.

1. **Power On:** Long press the power button to turn on the device.
2. **Select Mode:** Briefly press the mode button to cycle through the available treatment modes: **Body**, **Face**, and **Bikini**. Select the appropriate mode for your target area.
3. **Adjust Intensity:** Briefly press the power button to cycle through the 5 adjustable intensity levels. Start with Level 1 for your first use and gradually increase as your skin tolerates it.
4. **Apply to Skin:** Place the flash window firmly against your skin, ensuring full contact at a 90-degree angle. The device will only flash when proper skin contact is detected.
5. **Start Treatment (Manual Mode):** Press the flash button to emit a single flash. Move the device to the next area and repeat.
6. **Start Treatment (Auto-Glide Mode):** For larger areas, hold down the flash button for 2 seconds until the auto mode indicator light illuminates. The device will automatically flash at regular intervals as you glide it across your skin.

3 Custom Modes

Head to Toe Smooth in 12 Minutes

1

BIKINI Mode

Say farewell to stubborn hair in just 2 minutes!



2

BODY Mode

Effortless body hair removal in 8 minutes!



3

FACE Mode

Transform your facial hair with a 2-minute session!



Image: Illustration of the three custom modes (Bikini, Body, Face) and their recommended usage for different body parts.

Your browser does not support the video tag.

Video: An official instructional video demonstrating how to use the INIA Fond device, including setup, mode selection, and application techniques.

8. TREATMENT SCHEDULE

Consistent use is key for optimal results. Follow this recommended treatment schedule:

- **Weeks 1-4:** Three times per week.

- **Weeks 5-8:** Twice per week.
- **Week 9 and after:** Once per fortnight (every two weeks) for maintenance.

An Investment that Lasts

Affordable, yet as powerful
as a professional laser.



Image: Visual comparison of hair growth before and after 4 weeks of treatment, demonstrating visible hair reduction.

9. MAINTENANCE

Proper care ensures the longevity and effectiveness of your device:

- **Cleaning:** After each use, ensure the device is unplugged and cool. Gently wipe the flash window and the device body with a soft, dry cloth. Do not use abrasive cleaners or immerse the device in water.

- **Storage:** Store the device in a cool, dry place, away from direct sunlight and moisture. Keep it in its original packaging or a protective case to prevent damage.

10. TROUBLESHOOTING

If you encounter any issues with your INIA Fond device, please refer to the following general advice:

- **Device Not Turning On:** Ensure the power adapter is securely plugged into both the device and a working power outlet.
- **No Flash:** Check that the flash window is in full contact with your skin. The device has a safety mechanism that prevents flashing without proper skin contact. Ensure the device is not in auto-glide mode if you intend to manually flash.
- **Skin Irritation:** If you experience redness or discomfort, reduce the intensity level for your next treatment. Ensure you are following all pre- and post-treatment care instructions.
- **Lack of Results:** Ensure you are following the recommended treatment schedule and using the appropriate intensity level for your skin and hair type. Consistent use over several weeks is necessary to see significant hair reduction.

For further assistance, please consult the full instruction manual or contact INIA customer support.

11. SPECIFICATIONS

Feature	Detail
Model Number	D-1198
Product Dimensions	1 x 2 x 5 inches
Item Weight	1.96 Pounds
UPC	198715096736
Manufacturer	INIA

12. WARRANTY & SUPPORT

The INIA Fond IPL Laser Hair Removal Device comes with a 1-year warranty. This warranty can be extended to 2 years upon product registration. Free replacements are provided within the warranty period for any manufacturing defects or if the device runs out of flashes.

For warranty claims, technical support, or any product-related inquiries, please contact INIA Official US customer service.