

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

- › [ORDO](#) /
- › [ORDO Lite Sonic Toothbrush User Manual](#)

## ORDO SONICLITE

# ORDO Lite Sonic Toothbrush User Manual

Model: SONICLITE

## INTRODUCTION

---

Thank you for choosing the ORDO Lite Sonic Toothbrush. This manual provides essential information for the safe and effective use of your new electric toothbrush. Please read it thoroughly before first use and retain it for future reference.



Image: The ORDO Lite Sonic Toothbrush in Petal color, shown alongside its retail packaging.

## **SAFETY INFORMATION**

---

- Replace the toothbrush head every 3 months for optimal oral hygiene.
- Keep out of reach of children under 8 years old.
- This appliance contains batteries that are non-replaceable.
- Do not attempt to disassemble or repair the toothbrush.
- Do not use if the product appears damaged.
- Consult a dentist before use if you have severe gum disease or have undergone oral surgery.

## PACKAGE CONTENTS

---

Verify that all items are present in your ORDO Lite Sonic Toothbrush package:

- ORDO Lite Sonic Toothbrush Handle
- 1 x ORDO Sonic Toothbrush Head
- USB-C Charging Cable
- User Manual (this document)

## PRODUCT OVERVIEW AND FEATURES

---

The ORDO Lite Sonic Toothbrush is designed for effective oral hygiene with advanced features.

- **Superior Sonic Technology:** Delivers 35,000 sonic pulses per minute for efficient plaque removal.
- **Long-Lasting Battery:** Lithium-Ion battery provides over 5 weeks of use on a single charge.
- **Dual Brushing Modes:** "Clean" for daily cleaning and "Sensitive" for gentle brushing.
- **Smart Timer:** A 2-minute timer ensures recommended brushing duration.
- **Quad-Pacer:** Guides you to brush all quadrants of your mouth evenly.
- **Waterproof Design:** Suitable for use in wet environments.
- **USB-C Charging:** Convenient and universal charging.
- **Anti-Roll Design:** Prevents the toothbrush from rolling on flat surfaces.



Image: The ORDO Lite Sonic Toothbrush highlighting its advanced sonic technology, two brushing modes, 5+ weeks battery life, and USB-C charging.

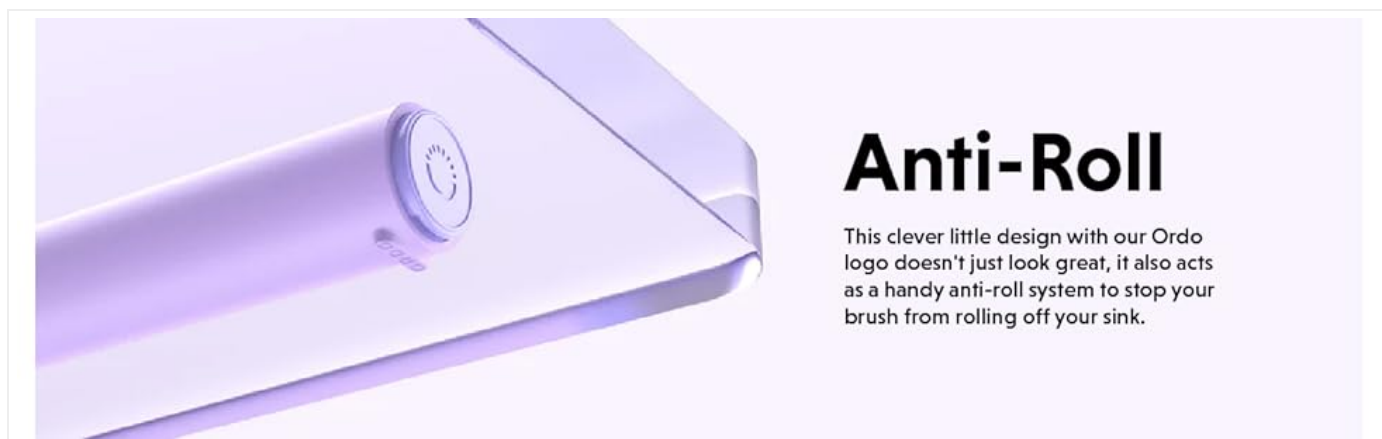


Image: A close-up view of the ORDO Lite Sonic Toothbrush handle, illustrating its anti-roll design feature.

## SETUP

---

### 1. Charging the Toothbrush

Before first use, fully charge your ORDO Lite Sonic Toothbrush. A full charge typically takes several hours.

1. Open the protective cover at the bottom of the toothbrush handle to reveal the USB-C port.
2. Insert the provided USB-C charging cable into the port.
3. Connect the other end of the USB-C cable to a compatible USB power adapter (not included) or a USB port.
4. The LED indicator lights will illuminate to show charging status. When fully charged, the lights will indicate a full battery.
5. Once charging is complete, remove the cable and securely close the protective cover to maintain water resistance.



Image: A close-up of the ORDO Lite Sonic Toothbrush's USB-C charging port with the protective cover open and a USB-C cable inserted.

# Effortless Charging

The Sonic Lite toothbrush can be charged at home or on the go with our USB-C to USB-A cable.



Image: An illustration demonstrating the effortless USB-C charging process for the ORDO Lite Sonic Toothbrush.

## 2. Attaching the Brush Head

Align the brush head with the metal shaft on the toothbrush handle and push it down firmly until it clicks into place.

## OPERATING INSTRUCTIONS

---

### 1. Brushing Your Teeth

1. Wet the bristles and apply a small amount of toothpaste.
2. Place the brush head against your teeth at a slight angle towards the gum line.
3. Press the power button to turn on the toothbrush. It will start in the default "Clean" mode.
4. Gently move the brush head slowly across your teeth, allowing the sonic vibrations to do the work. Do not scrub.
5. The built-in 2-minute timer will automatically turn off the toothbrush after the recommended brushing time. The quad-pacer will pause briefly every 30 seconds to signal you to move to the next quadrant of your mouth.
6. After brushing, rinse your mouth and the brush head thoroughly.

### 2. Changing Brushing Modes

The ORDO Lite Sonic Toothbrush features two brushing modes:

- **Clean Mode:** For a strong, daily clean.
- **Sensitive Mode:** For a gentle approach, suitable for sensitive teeth or new users of electric toothbrushes.


To switch between modes, press the mode button (located below the power button) while the toothbrush is operating. The indicator light next to the selected mode will illuminate.



## 2 brushing modes

clean, sensitive

Image: A detailed view of the ORDO Lite Sonic Toothbrush's control panel, showing the power button and indicators for the "Clean" and "Sensitive" brushing modes.

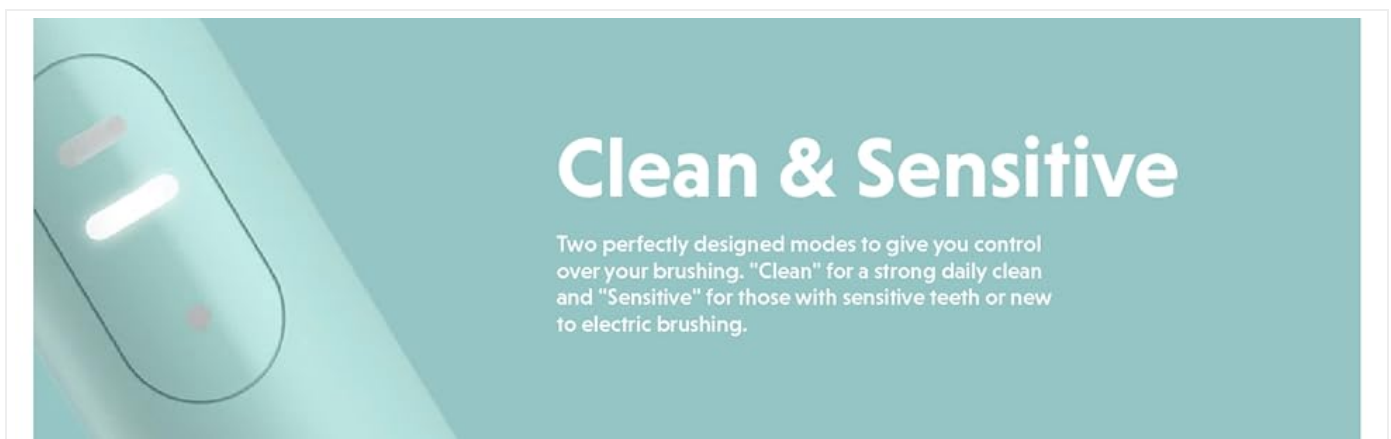

## Clean & Sensitive

Two perfectly designed modes to give you control over your brushing. "Clean" for a strong daily clean and "Sensitive" for those with sensitive teeth or new to electric brushing.

Image: An illustration explaining the two brushing modes, "Clean" for daily use and "Sensitive" for gentle cleaning.

# MAINTENANCE

## 1. Cleaning the Toothbrush

- After each use, rinse the brush head and bristles thoroughly under running water.
- Wipe the toothbrush handle with a damp cloth. Do not immerse the handle in water for extended periods.
- Ensure no water enters the charging port.

## 2. Replacing the Brush Head

For optimal performance and hygiene, replace the brush head every 3 months, or sooner if the bristles become frayed.

1. Grip the brush head and pull it firmly away from the handle.
2. Align the new brush head with the metal shaft and push it down until it clicks into place.

# TROUBLESHOOTING

Problem	Possible Cause	Solution
Toothbrush does not turn on.	Battery is depleted.	Charge the toothbrush fully using the provided USB-C cable.
Weak vibrations or reduced performance.	Low battery or worn brush head.	Recharge the toothbrush. Replace the brush head if it's old or frayed.
Charging indicator lights not working.	Cable or power source issue.	Ensure the USB-C cable is securely connected. Try a different USB power adapter or port.
Brush head does not stay attached.	Incorrect attachment or damaged brush head/handle.	Ensure the brush head is pushed down firmly until it clicks. If the issue persists, replace the brush head. If the handle is damaged, contact support.

# SPECIFICATIONS

Feature	Detail
Model Number	SONICLITE
Brand	ORDO
Sonic Pulses	35,000 per minute
Brushing Modes	Clean, Sensitive
Battery Type	Lithium Ion
Battery Life	Over 5 weeks (based on 2 minutes, twice daily)
Charging	USB-C
Timer	2-minute smart timer with 30-second quad-pacer



Feature	Detail
Waterproof Rating	Waterproof design
Age Range	Adult
Item Firmness	Soft
Package Dimensions	9.61 x 3.66 x 2.05 inches
Item Weight	3.42 ounces

## WARRANTY AND SUPPORT

ORDO products are manufactured to high-quality standards. For information regarding warranty coverage, please refer to the warranty card included with your purchase or visit the official ORDO website.



If you encounter any issues or have questions about your ORDO Lite Sonic Toothbrush, please contact ORDO customer support:

- **Online Support:** Visit the [ORDO Store on Amazon](#) for FAQs and contact options.
- **Manufacturer:** Ordo



© 2025 ORDO. All rights reserved.

## Related Documents - SONICLITE

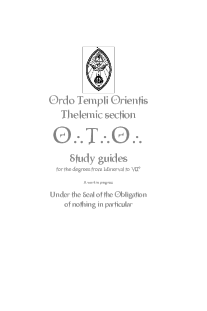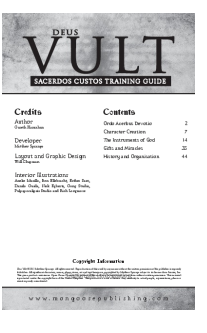	<p><a href="#">Ordo Templi Orientis Thelemic Section Study Guides</a></p> <p>A comprehensive guide to the study of the Ordo Templi Orientis (O.T.O.) Thelemic section, covering degrees from Minerval to VII. This document outlines recommended readings, memorization, and practices for initiates.</p>
	<p><a href="#">Deus Vult: Sacerdos Custos Training Guide</a></p> <p>A comprehensive training guide for Sacerdos Custos characters within the secretive Ordo Acerbus Devotio in the Deus Vult role-playing game. Explore the order's history, structure, character creation, combat, skills, and lore.</p>

