

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [Sejoy](#) /
- › [Sejoy FC5091 Cordless Water Flosser User Manual](#)

Sejoy FC5091

Sejoy FC5091 Cordless Water Flosser User Manual

Model: FC5091 | Brand: Sejoy

INTRODUCTION

Thank you for choosing the Sejoy FC5091 Cordless Water Flosser. This portable oral irrigator is designed to provide effective plaque removal and improve gum health through powerful water pulsation. Please read this manual thoroughly before use to ensure proper operation and maintenance.

IMPORTANT SAFETY INFORMATION

To ensure safe operation, please observe the following precautions:

- Avoid ingesting water used with the flosser.
- Exercise caution when using the device for the first time, especially if you have sensitive gums.
- Do not submerge the entire device in water for extended periods.
- Ensure the charging port is dry and covered before exposing the device to water.
- Keep out of reach of children.


PRODUCT OVERVIEW

The Sejoy FC5091 Cordless Water Flosser features a compact design with advanced functionalities for comprehensive oral care.

Key Features:

- **Multiple Nozzles:** Includes 6 nozzles (3 standard, 1 periodontal, 1 orthodontic, 1 tongue cleaner) for various oral care needs.
- **360° Rotating Nozzle:** Allows for easy access to all areas of the mouth.
- **4 Cleaning Modes:** Soft (25-65psi), Normal (70-110psi), Pulse (60-110psi), and DIY (20-110psi).
- **300ML Detachable Water Tank:** Large capacity for uninterrupted flossing and easy cleaning.
- **IPX7 Waterproof:** Safe for use in the shower and bathroom.
- **USB Rechargeable Battery:** 2000mAh battery provides up to 30 days of continuous use on a 4-hour charge.

6 NOZZLES with Different Function

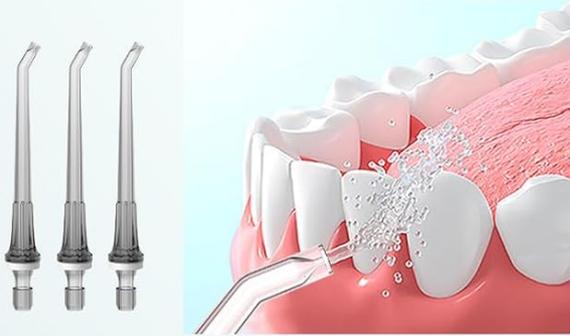


3×Standard Nozzles

Daily Cleaning




1×Periodontal Tip

Gum Cleaning

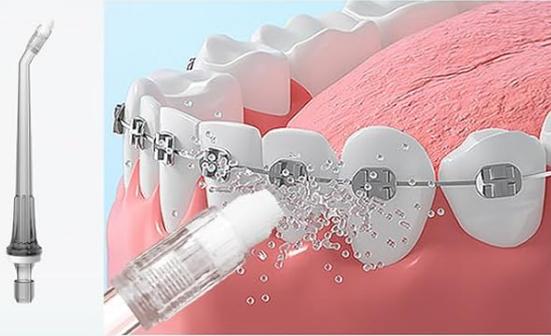

1×Orthodontic Tip

Remove food debris in the brace

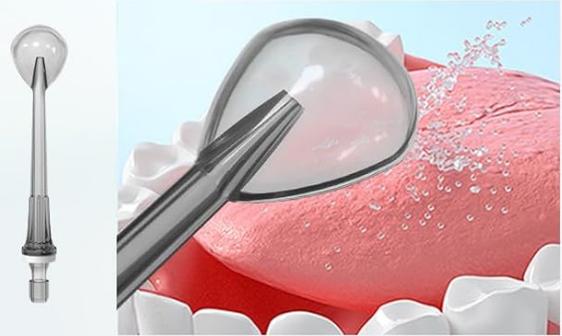

1×Tongue Cleaner

Remove bacteria build up for fresher breath

Image: The Sejoy FC5091 Cordless Water Flosser, showcasing its sleek design.

IPX7 waterproof



Safe use in the bathroom
anytime!

NOTE: Make sure the silicone plug is in place

Image: Illustration of the six included nozzles: 3 standard, 1 periodontal tip, 1 orthodontic tip, and 1 tongue cleaner, each designed for specific oral care needs.

How to adjust the Mode

- 1 Short press , the indicator lights up, enter Normal mode. Normal
- 2 Short press , the indicator lights up, enter Soft mode. Soft
- 3 Short press , the indicator lights up, enter DIY mode. (How to adjust the DIY intensity, long press 3s , the DIY indicator is flashing, the gear begins to change from weak to strong, gradually to the appropriate intensity, short press , you can fix the DIY intensity.) DIY
- 4 Press again and the three indicators will flash automatically in turn to enter Pulse mode

After selecting the appropriate mode, press and hold the  2S to start using it



Image: Diagram illustrating the components of the water flosser, including the nozzle, uninstal button, switch button, mode selection button, indicator lights for modes (High-level, Low-level, DIY), and water tank. Dimensions are also shown.

SETUP GUIDE

1. Initial Charging

Before first use, fully charge the device. The Sejoy FC5091 is USB rechargeable. Connect the provided USB cable to the charging port on the device and a suitable USB power source. A full charge takes approximately 4 hours and provides up to 30 days of continuous use.

4 Professional Modes



25-65psi  Soft mode Suitable for first time users or sensitive teeth	70-110psi  Normal mode Deep cleaning of teeth and debris in between teeth
20-110psi  DIY mode Water pressure can be freely adjusted according to your dental needs.	60-110psi  Pulse mode Water pressure is cycled from low to high to clean and massage the teeth gums and tongue

Image: The water flosser showing its 2000mAh powerful battery and 300ml large water tank, with an illustration of the USB charging port.

2. Attaching the Nozzle

Select the desired nozzle (standard, periodontal, orthodontic, or tongue cleaner). Insert the nozzle into the hole on the top of the controller until you hear a "click" sound, indicating successful installation. Rotate the nozzle to adjust the direction of water outflow as needed.



Image: Step-by-step visual guide showing how to add water, insert the nozzle tip, and press the mode switching button.

3. Filling the Water Tank

The 300ML water tank is detachable for easy filling and cleaning. Twist the bottom part of the device in the direction of the arrow to detach the tank. Fill the tank with lukewarm water. Reattach the tank by aligning it and twisting until secure. Ensure the silicone plug for the charging port is in place to maintain IPX7 waterproofing.

2000mAh

Powerful Battery

Fully charged in
4 hours

Up to **30 days**
of use according to 2
minutes a day



300ml

Large Water Tank

Large diameter, easy to fill

Translucent water tank with ml and oz scales, water level is clearly visible.



Image: The water flosser submerged in water, highlighting its IPX7 waterproof rating, with a close-up of the charging port cover.

OPERATING INSTRUCTIONS

1. Selecting a Mode

Press the mode button to cycle through the available cleaning modes. The indicator lights will show the selected mode:

- **Normal Mode (70-110psi):** For effective daily cleaning.
- **Soft Mode (25-65psi):** Recommended for first-time users or those with sensitive teeth/gums.
- **DIY Mode (20-110psi):** Allows for custom adjustment of water pressure.
- **Pulse Mode (60-110psi):** For gum massage and improved circulation.

How to use DIY mode?



- 1 Long press the power button
- 2 Short press the mode button to switch to DIY mode
- 3 Long press mode button 3 seconds after letting go, at this time the DIY indicator light flashes, the water pressure will be from low to high to low cycle, when the water pressure to your satisfaction gear, short press mode button. DIY mode setup is complete!

Unlock nozzle button

Power button

Mode button



Image: Visual representation of the four cleaning modes: Soft, Normal, DIY, and Pulse, with their respective pressure ranges and recommended uses.

How to remove the tank



- 1 As shown in the picture, twist in the direction of the arrow



- 2 Pull down firmly

Image: Detailed instructions on how to adjust the cleaning modes using the mode selection button, including how to set the DIY intensity.

2. Using DIY Mode

To adjust the intensity in DIY mode (20-110psi):

1. Long press the power button to turn on the device.
2. Short press the mode button to switch to DIY mode.
3. Long press the mode button for 3 seconds. The DIY indicator light will flash, and the water pressure will cycle from low to high.
4. Short press the mode button again when the desired pressure is reached to fix the intensity.



Image: Visual guide explaining the steps to activate and adjust the water pressure in DIY mode, showing the power button, mode button, and unlock nozzle button.

3. Flossing Procedure

After selecting the mode, press and hold the power button for 2 seconds to start the device. Place the nozzle tip into your mouth, directing the water jet at a 90-degree angle to your teeth. Keep your lips slightly closed to prevent splashing, but allow water to flow into the sink. Move the nozzle along your gumline and between teeth, pausing briefly between each tooth. Clean all areas of the mouth, including the back teeth and difficult-to-reach areas. The device has a 120-second auto-timer for optimal cleaning duration.



Image: A sequence of images demonstrating the flossing procedure: adding water, inserting the tip, selecting a mode, positioning the nozzle in the mouth, and turning on the power.

MAINTENANCE

1. Cleaning the Water Tank and Nozzle

After each use, empty any remaining water from the tank. Detach the water tank by twisting it off. Rinse the tank thoroughly with clean water. For the nozzle, press the uninstall button and remove it. Rinse the nozzle under running water. It is recommended to replace nozzles every 3 months for hygiene.

How to replace the replacement head



Insert the nozzle into the hole on the top of the controller. You should hear a "click" sound on successful installation. Rotate the nozzle to adjust the direction of water outflow.



If you need to replace the nozzle, press and hold the uninstall button, then remove the nozzle

Image: Two images demonstrating how to remove the water tank: twisting it in the direction of the arrow and then pulling it down firmly.

How to use

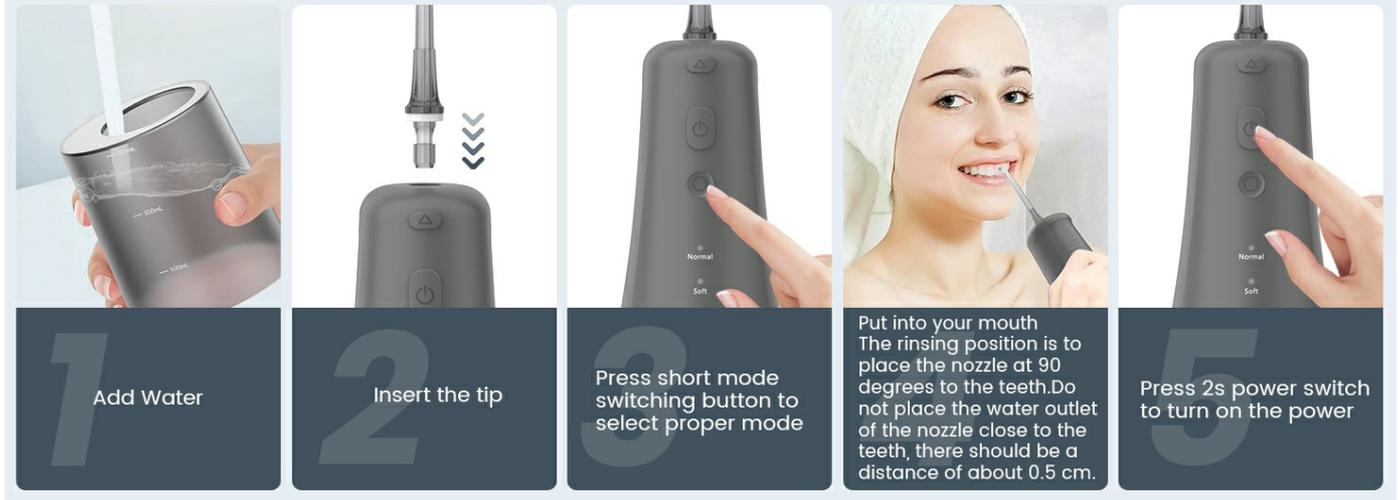


Image: Two images showing how to insert a new nozzle into the controller until it clicks, and how to remove a nozzle by pressing the uninstall button.

2. Device Storage

Store the water flosser in a cool, dry place away from direct sunlight when not in use. Ensure the water tank is empty and clean before storage.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery or not charged.	Charge the device fully (approx. 4 hours).
Weak water pressure.	Low battery, insufficient water in tank, or clogged nozzle.	Recharge the device, refill water tank, or clean/replace nozzle. Ensure DIY mode is set to desired pressure.
Water leakage.	Water tank not properly attached or silicone plug for charging port is open.	Ensure water tank is securely twisted into place. Close the silicone plug firmly.
Nozzle does not stay in place.	Nozzle not fully inserted.	Push the nozzle firmly until you hear a "click".

SPECIFICATIONS

- **Model Number:** FC5091
- **Brand:** Sejoy
- **Water Tank Capacity:** 300ML
- **Pressure Range:** 20-110 PSI (varies by mode)
- **Battery:** 2000mAh Lithium Battery
- **Charging Time:** Approximately 4 hours
- **Usage Time:** Up to 30 days (based on 2 minutes/day)
- **Waterproof Rating:** IPX7
- **Package Dimensions:** 9.53 x 3.43 x 3.03 inches
- **Item Weight:** 14.11 ounces

WARRANTY AND SUPPORT

For warranty information or technical support, please refer to the product packaging or contact Sejoy customer service directly. Keep your purchase receipt as proof of purchase.

Legal Disclaimer: Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition. This product is intended for oral hygiene purposes only.

