

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [ADIUPUL](#) /

› [ADIUPUL Back Massager Belt K212 User Manual](#)

ADIUPUL K212

ADIUPUL Back Massager Belt K212 User Manual

Model: K212

INTRODUCTION

Thank you for choosing the ADIUPUL Back Massager Belt. This device combines red light therapy, soothing heat, and vibration massage to provide comfort and relief. Please read this manual thoroughly before use to ensure safe and effective operation.

Important Safety Instructions

- Do not use if pregnant, have a pacemaker, or any serious medical conditions without consulting a doctor.
- Avoid direct contact with water. The device is not water-resistant.
- Do not use on broken or irritated skin.
- Keep out of reach of children.
- The device features an automatic shut-off after 20 minutes for safety.
- Only use the provided charging cable.

Package Contents

- ADIUPUL Back Massager Belt
- Magnetic Controller (5000mAh)
- Extension Belt
- Charging Cable
- User Manual

Product Overview

The ADIUPUL Back Massager Belt is designed for targeted relief and relaxation. It integrates three core functions:

- **Red Light Therapy:** Features 110 LED lights (660nm visible red light and 850nm invisible infrared light) for deep tissue penetration and surface skin benefits.
- **Soothing Heat:** Offers 5 adjustable heating levels (38°C/100.4°F to 48°C/118.4°F) for comfortable warmth.
- **Vibration Massage:** Equipped with 4 powerful vibrating motors and 3 intensity levels to stimulate muscles.



Image: The ADIUPUL Back Massager Belt showing the magnetic controller and key features like red light therapy, heat, vibration, and auto shut-off.



5 Heating Levels
110pcs Red Light Therapy

38°C/100.4°F 40°C/104°F

42°C/107.6°F 45°C/113°F 48°C/118.4°F

Image: A woman wearing the massager belt, highlighting the 5 heating levels and 110 red light therapy LEDs with temperature ranges.

4 Powerful Vibrating Motors **3 Vibration Modes & Intensity Levels**



Image: A man relaxing on a couch with the massager belt, illustrating the 4 powerful vibrating motors and 3 vibration modes/intensity levels.



5000mAh Magnetic Controller



5 Strong
Magnetic Buttons



Touch
Screen Control

Image: A close-up of the magnetic controller, showing its 5000mAh capacity, 5 strong magnetic buttons, and touch screen control with battery and temperature display.

Setup Instructions

1. **Charge the Controller:** Before first use, fully charge the magnetic controller using the provided charging cable. The controller detaches from the belt for convenient charging.
2. **Attach the Controller:** Align the magnetic controller with the designated magnetic points on the belt. It will securely snap into place. The device will not activate until the controller is properly attached.
3. **Wear the Belt:** Wrap the massager belt around your waist, ensuring the red light and vibration panels are positioned over the desired area of your back. Secure the belt using the hook-and-loop fasteners.
4. **Use Extension Belt (if needed):** If the standard belt size (25-48 inches) is too small, attach the extension belt to achieve a comfortable fit.

Video: Demonstration of attaching the magnetic controller and wearing the cordless back massager belt.

The Perfect Gift for Loved Ones

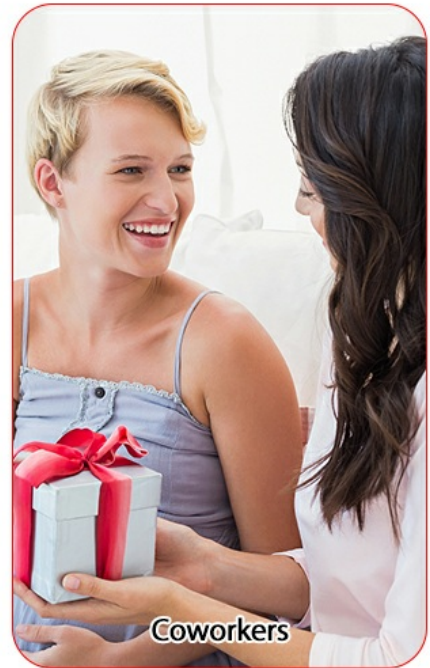


Image: Visual guide on how to attach the extension belt for a larger fit.

Operating Instructions

1. **Power On:** Press and hold the power button on the magnetic controller to turn on the device. The display will light up.
2. **Adjust Heat:** Press the 'Heat' button (represented by wavy lines) to cycle through 5 heating levels. The current temperature will be displayed.
3. **Adjust Vibration:** Press the 'Vibration' button (represented by zigzag lines) to select from 3 vibration modes and intensity levels.
4. **Combine Functions:** You can use the heat and vibration functions independently or simultaneously.
5. **Automatic Shut-off:** For safety, the device will automatically turn off after 20 minutes of continuous use. You can restart it if further use is desired.
6. **Power Off:** Press and hold the power button to turn off the device.

Video: Demonstration of operating the magnetic controller to adjust heat and vibration settings.

Care and Maintenance

- **Cleaning:** Disconnect the magnetic controller before cleaning. Wipe the belt with a soft, damp cloth. Do not submerge the belt or controller in water.

- **Storage:** Store the device in a cool, dry place away from direct sunlight when not in use.
- **Battery Care:** To prolong battery life, avoid fully discharging the controller frequently. Charge it regularly.

Troubleshooting

Problem	Possible Cause	Solution
Device does not turn on.	Battery is low or dead. Controller not properly attached.	Charge the magnetic controller. Ensure the controller is securely attached to the magnetic points on the belt.
Heat or vibration not working.	Function not activated. Device automatically shut off.	Press the respective heat or vibration button on the controller. Restart the device if 20 minutes have passed.
Belt feels too tight/loose.	Incorrect sizing.	Adjust the hook-and-loop fasteners. Use the included extension belt if a larger size is needed.

Specifications

- **Model Number:** K212
- **Power Source:** Battery Powered (5000mAh)
- **Heating Levels:** 5 (38°C/100.4°F to 48°C/118.4°F)
- **Vibration Modes:** 3
- **Vibration Intensities:** 3
- **LED Lights:** 110 Pcs (660nm Red Light, 850nm Invisible Infrared Light)
- **Automatic Timer:** 20 minutes
- **Product Dimensions:** 13.7 x 9 x 2.5 inches
- **Product Weight:** 2.16 Pounds
- **Belt Fit:** 25-48 inches (with extension belt for larger sizes)
- **Water Resistance:** Not Water Resistant

Warranty and Customer Support

For warranty information or customer support, please refer to the contact details provided with your purchase or visit the official ADIUPUL website. Keep your purchase receipt for warranty claims.