

Teglu FHC-2222

Teglu FHC-2222 Paraffin Bath Instruction Manual

Model: FHC-2222

INTRODUCTION

This manual provides essential instructions for the safe and effective use of your Teglu FHC-2222 Paraffin Bath. Please read this manual thoroughly before operating the device and retain it for future reference. The Teglu FHC-2222 is designed for paraffin wax treatments for hands, feet, and elbows, offering a convenient and efficient way to enjoy the benefits of warm paraffin therapy.

IMPORTANT SAFETY INSTRUCTIONS

To reduce the risk of burns, electrocution, fire, or injury, always follow basic safety precautions when using electrical appliances. Read all instructions before use.

- Do not reach for an appliance that has fallen into water. Unplug it immediately.
- Do not use while bathing or in a shower.
- Do not place or store appliance where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or other liquid.
- Always unplug this appliance from the electrical outlet immediately after using.
- Unplug this appliance before cleaning.
- Close supervision is necessary when this appliance is used by, on, or near children or invalids.
- Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair.
- Keep the cord away from heated surfaces.
- Never use while sleeping or drowsy.

- Never drop or insert any object into any opening.
- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not use on sensitive skin, open wounds, or areas with skin conditions. Consult a physician if you have any concerns.
- Always test the wax temperature on a small area of your skin before full immersion to ensure it is comfortable and safe.
- Ensure the unit is placed on a stable, heat-resistant surface.
- Do not overfill the paraffin bath.

PACKAGE CONTENTS

Verify that all items are present in the package:

- 1 x Teglu FHC-2222 Paraffin Bath (4L capacity)
- 3 x 450 g Paraffin Wax (Lavender scent)
- 100 x Plastic Gloves Liner for Hands
- 100 x Plastic Gloves Liner for Feet
- 2 x Thermal Gloves for Hands
- 1 x Silicone Scraper/Brush
- 1 x Thermal Insulation Grid


PRODUCT FEATURES

- **Large Capacity:** 4-liter capacity with a large inner tank (32 x 17 cm) suitable for hands, elbows, and feet (up to UK shoe size 12 / EU 46.5).
- **Fast Melting:** 200W high-performance heating with an aluminum liner melts 3 x 450g paraffin wax in approximately 45 minutes.
- **Digital Control:** LCD display for temperature and time, with touchable buttons. Adjustable temperature range from 30-70°C (86-158°F).
- **Automatic Lid Opening:** Intelligent system automatically opens the lid and initiates cooling based on melting progress.
- **Constant Temperature & Overheating Protection:** Maintains a constant temperature and stops heating if the temperature exceeds 70°C (158°F) for safety.

SETUP INSTRUCTIONS

1. **Placement:** Place the paraffin bath on a flat, stable, and heat-resistant surface. Ensure it is away from water sources and out of reach of children.
2. **Add Wax:** Open the lid and place the paraffin wax blocks into the inner pot. Ensure the wax is evenly distributed.

Easy to Use

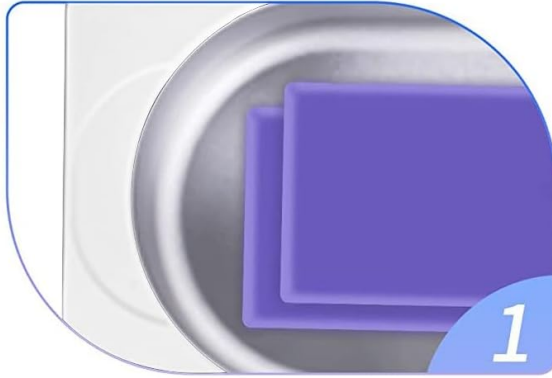

1. Put the paraffin wax into the pot and start heating.



2. Place the heating insulating grid on the pot after the wax cools to proper temperature.



3. Slowly dip your hand or feet into the melted wax, keep in the wax for 20 seconds, then keep them out till the wax is dried. Repeat it 3 to 5 times until it forms a thick layer.



4. Put on the gloves and remove the wax after 10-15 minutes.

Image: Step 1 - Place paraffin wax into the pot.

- Power On & Melt:** Plug the unit into a suitable electrical outlet. The display will light up. Set the desired melting temperature (e.g., 70°C / 158°F) using the touch controls. The 200W heating element and aluminum liner will melt 3 x 450g of wax in approximately 45 minutes.

200W Fast & Gentle Heat

The use of ALUMINUM LINER enables the melting time of paraffin to be reduced to 30 mins(2 pack of paraffin)-45mins(3 pack of paraffin), which is 4 TIMES FASTER than traditional plastic liner (2-3 hours)



Image: The 200W heating system melts paraffin wax quickly.


4. **Automatic Lid Opening:** The intelligent system will automatically open the lid once the wax has melted and begins to cool to the appropriate treatment temperature.

Intelligent Automatic Lid Opening

The system will automatically open the lid and start cooling according to the melting progress, so that the paraffin wax can quickly reach the appropriate temperature.



Wax melts

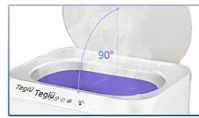

Automatic opening




Cool to fine temperature

Image: The automatic lid opening feature.

5. **Insert Insulation Grid:** Once the wax is fully melted and has cooled to a safe temperature (typically around 50-55°C / 122-131°F), carefully place the thermal insulation grid into the pot. This grid prevents direct contact with the heating element and provides a stable surface for immersion.

Easy to Use

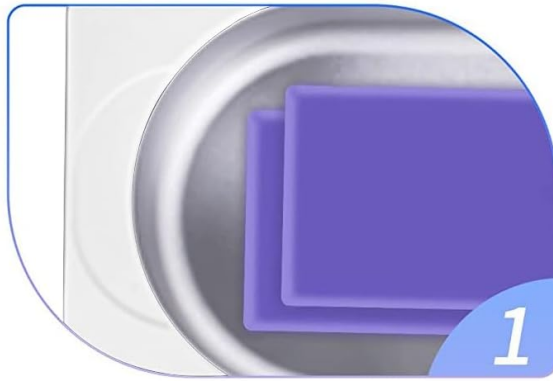

1. Put the paraffin wax into the pot and start heating.



2. Place the heating insulating grid on the pot after the wax cools to proper temperature.



3. Slowly dip your hand or feet into the melted wax, keep in the wax for 20 seconds, then keep them out till the wax is dried. Repeat it 3 to 5 times until it forms a thick layer.



4. Put on the gloves and remove the wax after 10-15 minutes.

Image: Step 2 - Place the heating insulating grid on the pot.

OPERATING INSTRUCTIONS (PARAFFIN TREATMENT)

1. **Prepare Skin:** Ensure your hands, feet, or elbows are clean and dry. Remove any jewelry.
2. **Test Temperature:** Before immersing, test a small amount of wax on your wrist or inner arm to ensure the temperature is comfortable and not too hot. The recommended treatment temperature is typically between 50-55°C (122-131°F).
3. **Immerse:** Slowly dip your hand, foot, or elbow into the melted wax. Keep it submerged for approximately 20 seconds, then slowly remove it. Allow the wax to dry slightly on your skin. Repeat this process 3 to 5 times until a thick layer of wax forms.


4L Large Capacity & Large Opening Design

The large liner size (32x17cm/12.6x6.7inch) is suitable for your hands, elbows and feet, and it can fit feet up to 12(UK)/46.5(EU).



Image: The 4L capacity accommodates feet up to EU size 46.5.

Easy to Use

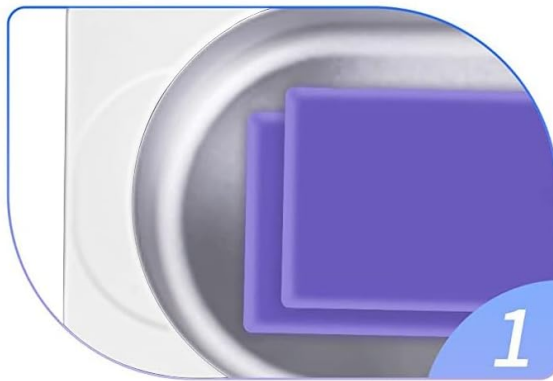

1. Put the paraffin wax into the pot and start heating.



2. Place the heating insulating grid on the pot after the wax cools to proper temperature.



3. Slowly dip your hand or feet into the melted wax, keep in the wax for 20 seconds, then keep them out till the wax is dried. Repeat it 3 to 5 times until it forms a thick layer.



4. Put on the gloves and remove the wax after 10-15 minutes.

Image: Step 3 - Slowly dip your hand or foot into the melted wax.

4. **Cover:** After applying the desired layers of wax, immediately cover the treated area with a plastic glove liner (provided). For enhanced heat retention, place a thermal glove over the plastic liner.
5. **Relax:** Allow the wax to remain on your skin for 10-15 minutes. This allows the heat to penetrate and moisturize the skin.
6. **Remove Wax:** Gently peel off the hardened wax, starting from the wrist or ankle. The wax should come off easily, taking dead skin cells with it. Dispose of used wax; do not reuse it for hygienic reasons.

Easy to Use



1. Put the paraffin wax into the pot and start heating.

2. Place the heating insulating grid on the pot after the wax cools to proper temperature.



3. Slowly dip your hand or feet into the melted wax, keep in the wax for 20 seconds, then keep them out till the wax is dried. Repeat it 3 to 5 times until it forms a thick layer.



4. Put on the gloves and remove the wax after 10-15 minutes.

Image: Step 4 - Put on gloves and remove the wax after 10-15 minutes.

7. **Moisturize:** After removing the wax, apply a moisturizer to your skin to lock in the hydration.

MAINTENANCE

• Cleaning:

- Always unplug the unit and allow it to cool completely before cleaning.
- To clean the interior, allow the remaining wax to solidify. Once solid, gently push on one end of the wax block to lift it out. Use the provided silicone scraper if needed.
- Wipe the interior and exterior surfaces with a soft, damp cloth. Do not use abrasive cleaners or immerse the unit in water.
- Ensure the unit is completely dry before storing or next use.

• **Wax Replacement:** It is recommended to replace the paraffin wax after approximately 20-30 uses, or if it appears dirty or contaminated.

• **Storage:** Store the paraffin bath in a cool, dry place, away from direct sunlight and moisture. Keep the lid closed to prevent dust accumulation.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Unit does not turn on.	Not plugged in; power outlet issue; damaged cord/plug.	Ensure the unit is securely plugged into a working outlet. Check the power cord for damage.
Wax not melting or melting slowly.	Incorrect temperature setting; too much wax; unit malfunction.	Increase the temperature setting. Ensure you are using the recommended amount of wax. If issues persist, contact support.
Wax is too hot for treatment.	Temperature set too high; insufficient cooling time.	Lower the temperature setting. Allow more time for the wax to cool to the recommended treatment temperature (50-55°C). Always test before immersion.
Lid not opening automatically.	Sensor obstruction; malfunction.	Check for any obstructions around the lid mechanism. If the issue persists, contact support.

SPECIFICATIONS

Feature	Detail
Model Number	FHC-2222
Brand	Teglu
Capacity	4 Liters
Inner Tank Dimensions	32 x 17 cm (12.6 x 6.7 inches)
Power	200 W
Temperature Range	30-70°C (86-158°F)
Melting Time (3x450g wax)	Approx. 45 minutes
Overheating Protection	Automatic shut-off at 70°C (158°F)
Country of Origin	China
Shipping Dimensions	41 x 29.8 x 17.7 cm
Shipping Weight	3.96 kg

WARRANTY AND SUPPORT

For warranty information or technical support, please refer to the documentation included with your purchase or contact Teglu customer service. Keep your purchase receipt as proof of purchase. For further assistance, please visit the official Teglu website or contact their customer support channels.