

## AMIRO SD-1186

# AMIRO Opal IPL Hair Removal Device SD-1186 User Manual

Model: SD-1186

## 1. INTRODUCTION

The AMIRO Opal IPL Hair Removal Device (Model SD-1186) utilizes Intense Pulsed Light (IPL) technology for effective and convenient at-home hair removal. This device is designed to provide a comfortable experience with its integrated ice cooling mode and offers multiple energy levels and operating modes for various body areas.



Figure 1: AMIRO Opal IPL Hair Removal Device

## 2. SAFETY INFORMATION

- Always read the user manual thoroughly before first use.
- Wear the provided protective goggles during operation to shield your eyes from the pulsed light.
- Ensure your skin tone and hair color are suitable for IPL treatment. Refer to the skin tone and hair color chart in the manual.
- Do NOT use the device on eyebrows or tattoos.
- Perform a patch test on a small skin area 24 hours before full treatment to check for adverse reactions.
- If any irritation or discomfort occurs, discontinue use and consult a medical professional.
- The device is designed for use on naturally black or brown hair. It may not be as effective on white, gray, blonde, or red hair.

### 3. PACKAGE CONTENTS

The AMIRO Opal IPL Hair Removal Device package includes:

- AMIRO Opal IPL Hair Removal Device (Model SD-1186)
- Power Adapter (60W)
- Protective Goggles
- Shaver
- User Manual

### 4. SETUP

1. Unpack the device and all accessories from the packaging.
2. Inspect the device for any damage. Do not use if damaged.
3. Connect the power adapter to the device's power input port.
4. Plug the power adapter into a suitable electrical outlet.

The device is now ready for operation.

### 5. OPERATING INSTRUCTIONS

#### 5.1 Power On/Off

Press and hold the power button (usually located on the device) for a few seconds to turn the device on or off. The device will emit a sound and the display will light up.

#### 5.2 Mode Selection (HR/RA/SC)

The device features 3-in-1 functions: HR (Hair Removal), RA (Skin Rejuvenation), and SC (Acne Clearance). Use the mode selection button to cycle through these options. The display will indicate the currently selected mode.

#### 5.3 Energy Levels

The device offers 8 energy levels. Start with level 1 and gradually increase the intensity as your skin tolerates it. Higher levels generally yield better results. Adjust the level using the designated buttons on the device.

## 8 ENERGY LEVELS 13J



Figure 2: 8 Energy Levels and Treatment Progression

### 5.4 Auto/Manual Modes

The IPL hair remover supports two operating modes:

- **Manual Mode (Single Flash):** Ideal for smaller, sensitive, or curved areas such as bikini lines, upper lip, chin, face, and armpits. Press the flash button for each pulse.
- **Auto Mode (Continuous Flash):** Suitable for larger areas like the back, chest, stomach, arms, and legs. Hold the flash button to activate continuous flashing as you glide the device over the skin.

## 2MODES



Manual Mode



Auto Mode

Figure 3: Manual and Auto Modes for Targeted Treatment

### 5.5 Ice Cooling Mode

Activate the ice cooling mode to cool the skin to approximately 50°F (10°C) during treatment. This feature helps reduce pain and discomfort, providing a more comfortable experience, especially at higher energy levels.

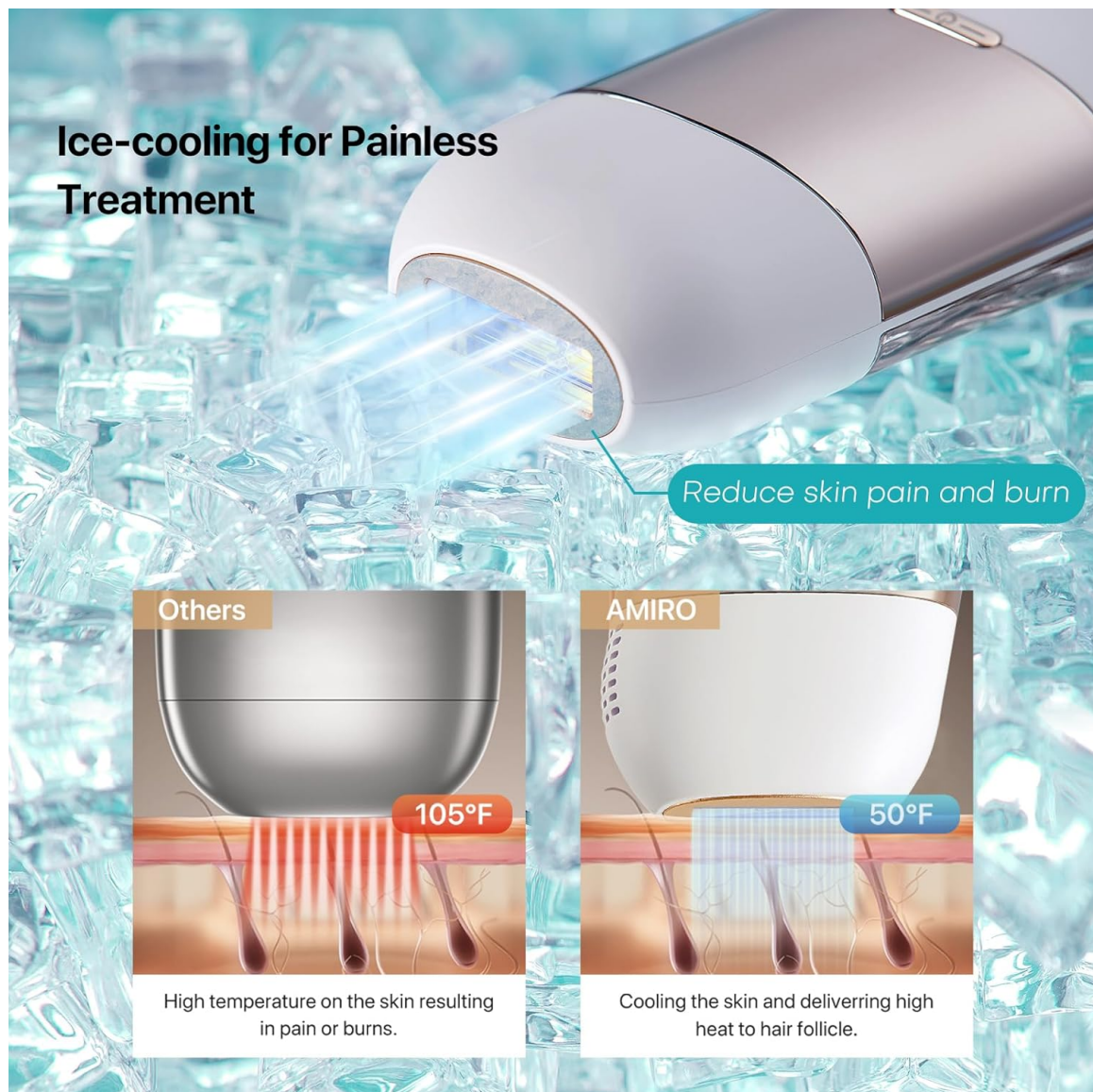


Figure 4: Ice Cooling for Painless Hair Removal

## 5.6 Skin Tone and Hair Color Compatibility

The effectiveness and safety of IPL treatment depend on your skin tone and hair color. Refer to the charts below to determine if the device is suitable for you.





Figure 5: Skin Tone and Hair Color Compatibility Chart

## 5.7 General Usage Steps

1. Shave the treatment area thoroughly before each session.
2. Cleanse the skin and ensure it is dry and free of lotions or oils.
3. Wear protective goggles.
4. Select the desired mode (HR/RA/SC) and energy level.
5. Place the device's flash window firmly against the skin.
6. Press the flash button (Manual Mode) or glide the device (Auto Mode) to deliver pulses.
7. Move the device to the next area, ensuring slight overlap to cover all desired regions.

Video 1: Demonstration of Ice Cool Painless Hair Removal and Auto/Manual Modes.

Video 2: Quick overview of the hair removal device in use.

Video 3: Detailed demonstration of the AMIRO Opal Hair Removal Device features and usage.

## 6. TREATMENT SCHEDULE

For optimal results, follow a consistent treatment schedule:

- **Weeks 1-4:** Use 2 times a week.
- **Weeks 5-8:** Use 1 time a week.
- **After 8 Weeks:** Use 1-2 times each month for maintenance.

Visible hair reduction typically begins after 3-4 uses, with up to 94% hair reduction after 7-8 uses.

## 7. MAINTENANCE

- Clean the flash window after each use with a soft, dry cloth to ensure optimal performance and hygiene.
- Store the device in a cool, dry place away from direct sunlight.
- Do not immerse the device in water or any other liquids.

## 8. TROUBLESHOOTING

- **Device not turning on:** Ensure the power adapter is securely connected to both the device and the power outlet.
- **No flash:** Make sure the flash window is fully in contact with your skin. The device has a skin contact sensor that prevents accidental flashes.
- **Discomfort/Pain:** Activate the ice cooling mode. Reduce the energy level if discomfort persists.
- **No visible results:** Ensure consistent use according to the treatment schedule. IPL results vary by individual and may take several weeks to become apparent.

## 9. SPECIFICATIONS

Feature	Detail
Model Number	SD-1186
Product Dimensions	6.89 x 2.56 x 1.57 inches
Item Weight	10.09 ounces
Energy Levels	8 Levels (up to 13J)
Functions	HR (Hair Removal), RA (Skin Rejuvenation), SC (Acne Clearance)
Cooling Technology	Ice Cooling Mode (cools skin to 50°F/10°C)



Modes	Auto/Manual
Manufacturer	Shenzhen Beauty Every Moment intelligent electric Co.,Ltd.

## 10. WARRANTY AND SUPPORT

AMIRO provides a one-year service coverage for this device. If you have any questions or require support, please contact customer service through the retailer where the product was purchased.