

## Operan OR989

# Operan Electric Toothbrush OR989 User Manual

Model: OR989

## INTRODUCTION

---

Thank you for choosing the Operan Electric Toothbrush OR989. This manual provides essential information for the safe and effective use of your new electric toothbrush. Please read it thoroughly before first use and keep it for future reference.

## WHAT'S IN THE BOX

---

Your Operan Electric Toothbrush package includes the following items:

- 1 x Operan Sonic Electric Toothbrush Handle
- 8 x Dupont Brush Heads (replace every 2-3 months)
- 1 x Travel Case
- 1 x Original Type-C Charging Cable (Wall adapter not included)
- 1 x Instruction Manual

# PROFESSIONAL level care

Recommendation by Dental  
practitioner



*Dr. Aaron*



Image: The package contents of the Operan Electric Toothbrush, showing the handle, eight brush heads, a travel case, and a charging cable.

## PRODUCT FEATURES

- **5 Brushing Modes & 3 Intensity Levels:** Customize your brushing experience with White, Clean, Sensitive, Polish, and Massage modes. Each mode offers three intensity levels (light, medium, heavy) for personalized oral care.
- **Smart Features:** A built-in 2-minute smart timer ensures optimal brushing duration, while a 30-second interval reminder prompts you to move to a different quadrant of your mouth, promoting thorough cleaning habits.
- **40,000 VPM Sonic Motor:** The high-vibration sonic motor delivers up to 40,000 vibrations per minute, effectively removing up to 99.99% of plaque and stains along the gumline and between teeth.

- **100-Day Battery Life & IPX7 Waterproof:** A full 6-hour charge provides approximately 100 days of usage (based on two uses per day). The toothbrush is IPX7 waterproof, allowing safe use in the bath or shower. Ensure the toothbrush is dry before charging.
- **W-Shaped Polished Bristles:** Designed to remove plaque effectively without damaging gums, reaching every corner for a deeper clean.

# Smart Toothbrush with Advanced Technology



Image: An exploded view diagram of the toothbrush highlighting its advanced technology and features.

# W-SHAPED Polished Bristles

Remove plaque without damaging your gums



Image: Detailed view of the W-shaped polished bristles designed for effective and gentle cleaning.

## SETUP

### 1. Initial Charging

Before first use, fully charge your Operan Electric Toothbrush. A full charge takes approximately 6 hours.

1. Connect the Type-C charging cable to the charging port at the bottom of the toothbrush handle.
2. Plug the USB end of the cable into a compatible wall adapter (not included) or a USB port.
3. The charging indicator light will illuminate to show charging status. It will turn solid when fully charged.



Image: The toothbrush handle connected to its Type-C charging cable, highlighting its long battery life.

## 2. Attaching the Brush Head

Align the brush head with the metal shaft on the toothbrush handle and push it down firmly until it clicks into place.

## OPERATING INSTRUCTIONS

### 1. Power On/Off

Press the power button once to turn the toothbrush on. Press and hold the power button for 2 seconds to turn it off.

## 2. Changing Modes and Intensity

After powering on, press the mode button repeatedly to cycle through the 5 brushing modes: White, Clean, Sensitive, Polish, and Massage. The indicator light next to the selected mode will illuminate.

To adjust intensity, press the intensity button (usually located below the power button) to select between 3 levels (low, medium, high). The corresponding indicator lights will show the selected intensity.

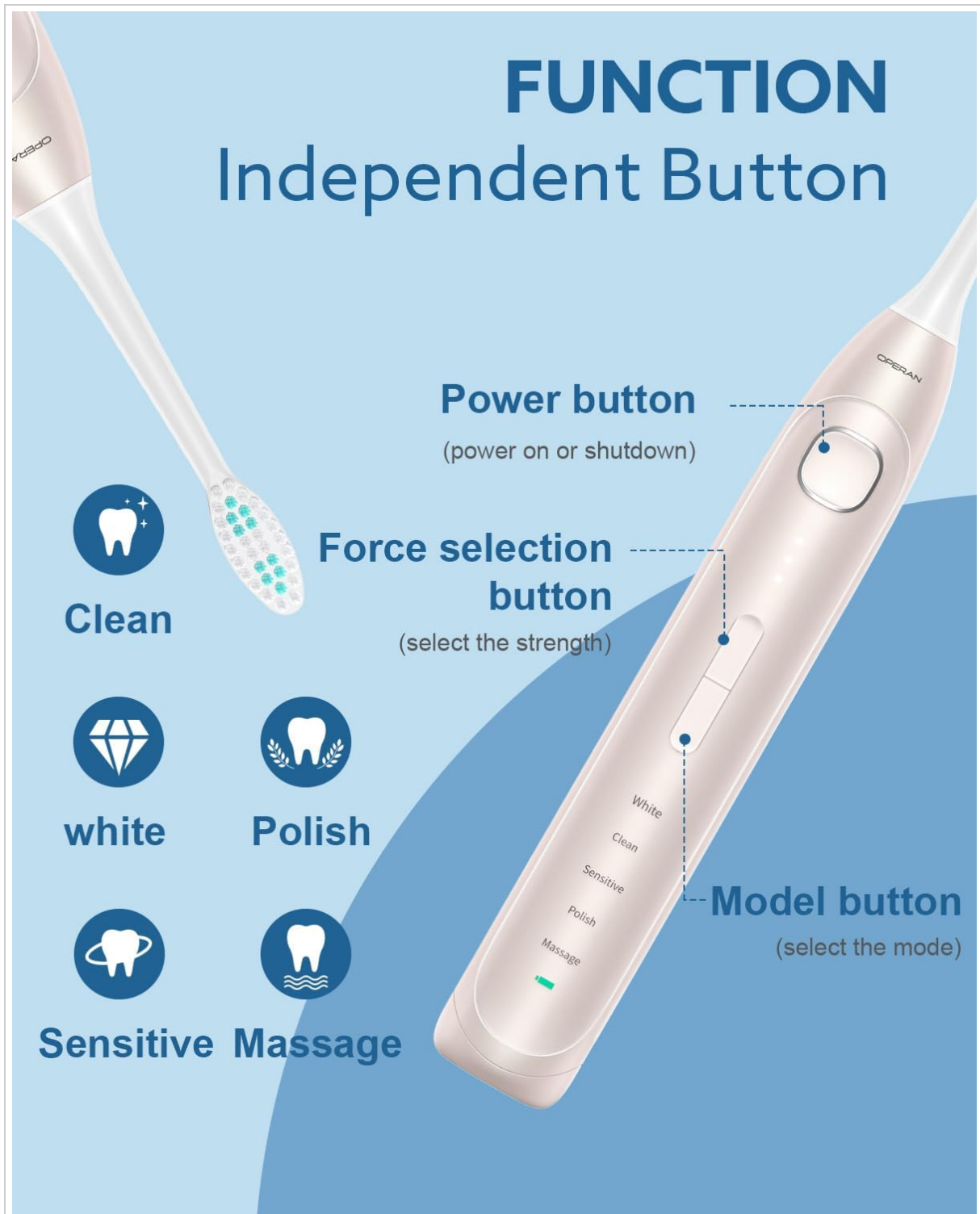


Image: The toothbrush handle with clearly labeled buttons for power, intensity, and mode selection.

## 3. Brushing Technique

Apply a small amount of toothpaste to the brush head. Place the brush head against your teeth at a slight angle towards the gum line. Gently move the brush across your teeth, allowing the sonic vibrations to do the work. The smart timer will guide you through a 2-minute brushing cycle with 30-second intervals for each quadrant of your mouth.

## MAINTENANCE

### 1. Cleaning the Toothbrush

- After each use, rinse the brush head thoroughly under running water.
- Remove the brush head and rinse the metal shaft area.
- Wipe the handle with a damp cloth and dry it. Do not immerse the handle in water for prolonged periods, especially the charging port area.

### 2. Replacing Brush Heads

It is recommended to replace the brush head every 2-3 months, or sooner if the bristles become frayed. To replace, simply pull the old brush head off and push a new one onto the shaft until it clicks.

# 8X

Eight toothbrush heads  
Available for up to **2.5+** years



## High Quality Travel Box

Your best travel companion



Image: Eight replacement brush heads and the included travel case.

## TROUBLESHOOTING

---

- **Toothbrush not turning on:** Ensure the toothbrush is fully charged.
- **Weak vibrations:** The battery might be low. Recharge the toothbrush.
- **Brush head not vibrating:** Ensure the brush head is securely attached to the handle.
- **Charging issues:** Check if the charging cable is properly connected and the wall adapter is functioning. Ensure the charging port is dry.

## SPECIFICATIONS

---

<b>Model Number</b>	OR989
<b>Power Source</b>	Battery Powered (1 Lithium Ion battery included)
<b>Vibrations Per Minute (VPM)</b>	Up to 40,000
<b>Waterproof Rating</b>	IPX7
<b>Battery Life</b>	Approx. 100 days (2 uses/day)
<b>Charging Time</b>	Approx. 6 hours for full charge
<b>Item Firmness Description</b>	Medium

## SAFETY INFORMATION

---

- Avoid prolonged water exposure.
- Always dry the toothbrush before charging.
- Keep out of reach of small children.
- Do not attempt to disassemble or repair the device yourself.

## WARRANTY AND SUPPORT

---

Operan provides a standard manufacturer's warranty for this product. For detailed warranty information, product support, or to purchase replacement parts, please visit the official Operan website or contact customer service through your retailer.