

ANLAN 1

ANLAN Ultrasonic Skin Scrubber User Manual

Brand: ANLAN | Model: 1

1. INTRODUCTION

Thank you for choosing the ANLAN Ultrasonic Skin Scrubber. This device is designed to provide professional-grade facial cleansing, blackhead extraction, and anti-aging benefits through advanced ultrasonic vibrations, LED light therapy, and EMS technology. Please read this manual carefully before use to ensure proper operation and maintenance.

2. SAFETY INFORMATION

- Do not use the device on damaged, inflamed, or broken skin.
- Keep out of reach of children.
- Ensure the device is fully charged before first use.
- Do not immerse the charging base in water. The device itself is IPX7 waterproof, but the charging base is not.
- Consult a dermatologist if you have sensitive skin conditions or concerns.
- Discontinue use if irritation or discomfort occurs.

3. PRODUCT OVERVIEW

3.1 Components

- ANLAN Ultrasonic Skin Scrubber Device
- Magnetic Charging Base
- USB Charging Cable
- User Manual

3.2 Key Features

- **Deep Cleaning:** Ultrasonic vibrations (up to 24,000/s) effectively remove blackheads, dirt, and exfoliate the skin.
- **LED Therapy:** Red light stimulates collagen, reduces pores and wrinkles; Blue light has anti-inflammatory effects, helps with acne.
- **EMS Anti-Aging:** Activates facial muscles, firms skin, and reduces signs of aging.
- **Enhanced Nutrient Absorption:** Negative ions facilitate deeper penetration of skincare products.
- **IPX7 Waterproof:** Safe for use in wet environments.
- **Magnetic Charging Base:** Convenient and space-saving charging.



Figure 1: ANLAN Ultrasonic Skin Scrubber and its accessories.

4. SETUP

4.1 Charging the Device

1. Connect the USB charging cable to the magnetic charging base.
2. Plug the other end of the USB cable into a suitable power adapter (not included) or a USB port.
3. Place the ANLAN Ultrasonic Skin Scrubber onto the magnetic charging base. Ensure the device is properly aligned for charging.
4. The indicator light on the device will show charging status. Charge fully before first use (approximately 2-3 hours).

LA BASE DE CARGA MAGNÉTICA AHORRA ESPACIO



Figure 2: The magnetic charging base provides convenient and space-saving power.

5. OPERATING INSTRUCTIONS

5.1 General Usage Steps

1. **Preparation:** Cleanse your face thoroughly. For best results, use a warm towel or facial steamer for 5 minutes to open pores.
2. **Apply Product:** Apply a facial cleanser, toner, serum, or cream to your face, ensuring the skin remains moist throughout the treatment.
3. **Power On:** Press and hold the power button to turn on the device.
4. **Select Mode:** Press the power button repeatedly to cycle through the different modes.
5. **Treatment:** Gently glide the spatula over your skin in an upward and outward motion, following the contours of your face. Do not press too hard.
6. **Post-Treatment:** After use, rinse your face with lukewarm water and gently pat dry. Clean the device as per maintenance instructions.



Figure 3: Illustrated steps for using the ANLAN Ultrasonic Skin Scrubber.

5.2 Modes of Operation

The ANLAN Ultrasonic Skin Scrubber features multiple modes to address various skin concerns:

Mode	Function	Recommended Use
Ultrasonic Cleaning (Peeling)	High-frequency vibrations atomize water to deeply cleanse pores, remove blackheads, dead skin, and impurities.	Use on wet skin with a cleanser. Glide the spatula at a 45-degree angle.
Ion+ (Ionic Drainage)	Utilizes positive ions to draw out stubborn dirt, oil, and makeup residue from pores.	After cleansing, use with toner or cleansing water.
Ion- (Nutrient Introduction)	Utilizes negative ions to help skincare products (serums, essences, moisturizers) penetrate deeper into the skin for better absorption.	Apply serum or cream, then glide the spatula to massage.
EMS Lifting (Microcurrent)	Microcurrent stimulation activates facial muscles, promotes collagen production, and helps to firm and lift the skin, reducing fine lines and wrinkles.	Use with a conductive gel or moisturizing cream. Glide upwards from chin to forehead.

4 FUNCTIONAL ULTRASONIC

Restore Radiant Skin

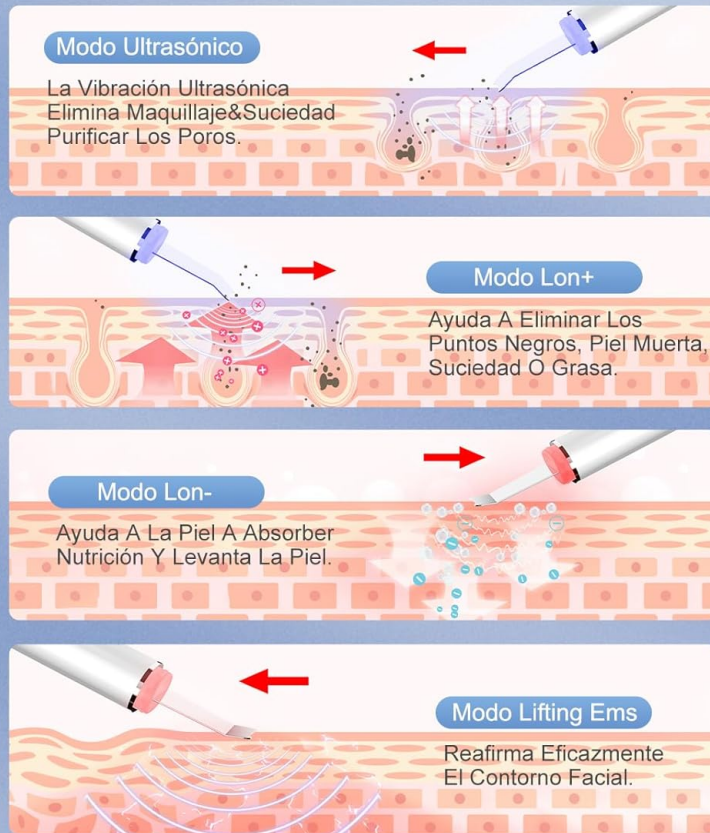


Figure 4: Overview of the four functional modes for radiant skin.



Figure 5: The ultrasonic process for deep cleaning and blackhead removal.

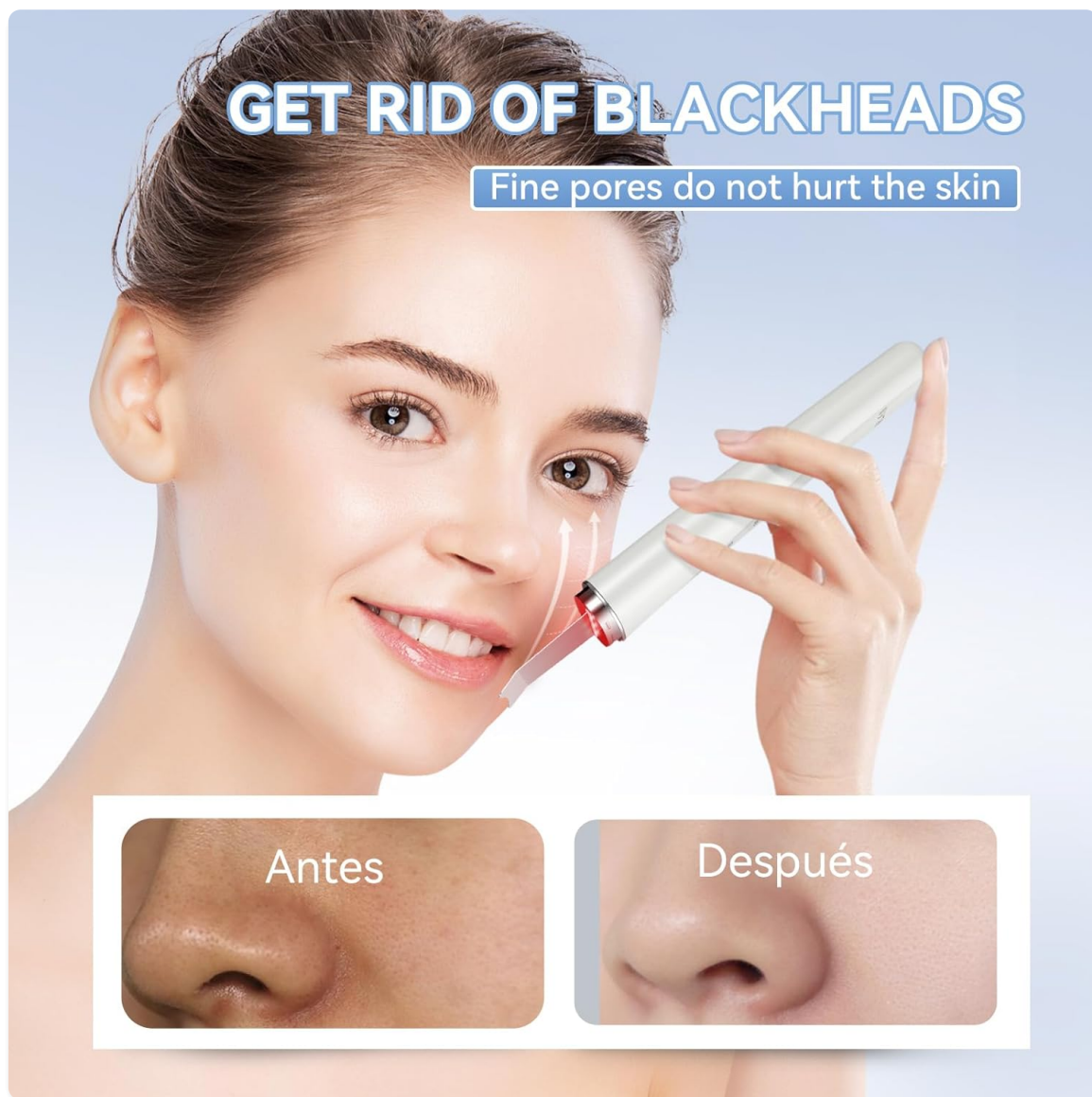


Figure 6: Visible improvement in blackheads after using the skin scrubber.



Figure 7: The device acts as a dermatology assistant, showing improvements in various skin conditions.

Massage nettoyant pour le visage



Trouble et dommageable pour les pores de la peau



Figure 8: The device provides an effective facial cleansing massage.

6. MAINTENANCE AND CARE

- **Cleaning the Spatula:** After each use, wipe the stainless steel spatula with a soft cloth or cotton pad moistened with water or a mild disinfectant.
- **Cleaning the Device Body:** Wipe the main body of the device with a dry or slightly damp cloth. Do not use abrasive cleaners or solvents.
- **Waterproof Feature:** The device is IPX7 waterproof, meaning it can withstand immersion in water up to 1 meter for 30 minutes. This allows for easy rinsing under running water.
- **Storage:** Store the device in a cool, dry place, away from direct sunlight and extreme temperatures.
- **Charging Base:** Keep the charging base dry. Do not rinse or immerse it in water.



Figure 9: The IPX7 waterproof design allows for easy cleaning and use without fear of water damage.



Figure 10: The device's IPX7 rating ensures it can be used and cleaned safely in the bathroom.

7. TROUBLESHOOTING

If you encounter any issues with your ANLAN Ultrasonic Skin Scrubber, please refer to the table below for common problems and solutions:

Problem	Possible Cause	Solution
Device does not turn on.	Low battery or not charged.	Charge the device fully using the provided charging base and cable.
No ultrasonic vibration.	Incorrect mode selected or insufficient contact with skin.	Ensure 'Ultrasonic Cleaning' mode is active. Make sure the spatula is in firm contact with moist skin.
Weak EMS sensation.	Skin is not moist enough or incorrect mode.	Ensure skin is adequately moist with a conductive gel or cream. Select 'EMS Lifting' mode.
Device stops working during use.	Battery is low.	Recharge the device.
Charging indicator light not working.	Improper connection or faulty cable/base.	Ensure the device is correctly placed on the magnetic base and the cable is securely connected. Try a different USB port or power adapter.

If the problem persists after trying these solutions, please contact ANLAN customer support for further assistance.

8. SPECIFICATIONS

Feature	Specification
Model Number	1
Brand	ANLAN
Material	Plastic
Color	White
Waterproof Rating	IPX7
Ultrasonic Frequency	Up to 24,000 Hz
Power Source	Rechargeable Battery
Charging Method	Magnetic Charging Base

Feature	Specification
Dimensions (Package)	19 x 12.9 x 3.2 cm
Weight	250 grams
Country of Origin	China




9. WARRANTY AND SUPPORT



ANLAN products are manufactured to the highest quality standards. For information regarding warranty coverage, returns, or technical support, please refer to the warranty card included with your product or visit the official ANLAN website. You may also contact ANLAN customer service directly for assistance with any product-related inquiries.

We are committed to providing a satisfactory solution for any confusion you may have during use.

© 2024 ANLAN. All rights reserved.

Related Documents - 1

 <p>ANLAN Ultrasonic Skin Scrubber User Manual</p> <p><small>If this product is not a medical device for diagnosis and treatment, but an auxiliary product for general household use. Thank you for purchasing ANLAN Ultrasonic Skin Scrubber. Before using this product, please read this instruction carefully. In addition, this user manual also served as a warranty card, so please keep it in a safe place for future use.</small></p>	<p>ANLAN Ultrasonic Skin Scrubber User Manual</p> <p>User manual for the ANLAN Ultrasonic Skin Scrubber (Model: JD-CP002), detailing its features, safety precautions, usage instructions, maintenance, and FAQs.</p>
 <p>ANLAN NECK BEAUTY DEVICE USER MANUAL</p>	<p>ANLAN Neck Beauty Device User Manual</p> <p>User manual for the ANLAN Neck Beauty Device, detailing its features, functions, usage instructions, and warranty information. The device utilizes positive and negative ions, vibration, warm massage, and red, blue, and green light therapy for skin cleansing, EMS care, and moisturizing.</p>
 <p>ANLAN NECK BEAUTY DEVICE USER MANUAL</p>	<p>ANLAN Neck Beauty Device User Manual: Features, Functions, and How to Use</p> <p>Comprehensive user manual for the ANLAN Neck Beauty Device (Model ES-081). Learn about its features, functions (EMS, warm, red/blue/green light), safety precautions, and how to use it for skin rejuvenation and care.</p>

 <p>ANLAN Body Beauty Device Model OFY-2102</p>	<p>ANLAN Body Beauty Device User Manual - Operation, Safety, and Specifications</p> <p>Comprehensive user manual for the ANLAN Body Beauty Device (Model OFY-2102), covering safety precautions, product details, usage instructions, troubleshooting, and warranty information.</p>
 <p>ANLAN Multifunctional Skincare Device Model 01-ADRY14-10A</p>	<p>ANLAN Multifunctional Skincare Device User Manual</p> <p>User manual for the ANLAN Multifunctional Skincare Device (Model 01-ADRY14-10A), providing comprehensive instructions on usage, features, safety, and maintenance for effective home skincare.</p>