

Comfytemp B0CLV5S4ZK

Comfytemp Wireless TENS Unit Machine User Manual

Model: B0CLV5S4ZK | Brand: Comfytemp

1. INTRODUCTION

This user manual provides comprehensive instructions for the safe and effective operation of your Comfytemp Wireless TENS Unit Machine. Please read this manual thoroughly before using the device and retain it for future reference.

The Comfytemp Wireless TENS Unit is designed for pain relief and muscle stimulation, offering a portable and wire-free solution for various body areas. It features both app control and on-device button control for ease of use.

2. IMPORTANT SAFETY INFORMATION

WARNING: Consult your physician before using this device, especially if you have pre-existing medical conditions, are pregnant, have a pacemaker, or any other implanted medical device. Do not use on open wounds, irritated skin, or near the heart or head.

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

- Do not use while driving or operating machinery.
- Keep out of reach of children.
- Do not immerse the device in water.
- Use only the electrode pads provided or recommended by Comfytemp.
- Discontinue use if you experience any discomfort or adverse reactions.

3. PACKAGE CONTENTS

Your Comfytemp Wireless TENS Unit package should include the following items:

- 2 x TENS Unit Hosts
- 2 x Butterfly Shape Electrode Pads (1-piece, 5.2 x 2.4 inches)
- 2 x Line Shape Electrode Pads (2-piece, 3.5 x 2.8 inches)
- 1 x Dual-head Charging Cable (USB-C)
- 1 x Storage Bag

- 1 x User Manual (this document)



Figure 3.1: Included components in the Comfytemp Wireless TENS Unit package.

4. SETUP AND FIRST USE

4.1 Charging the Device

Before first use, fully charge both TENS unit hosts. Use the provided dual-head charging cable.

1. Connect the dual-head charging cable to a USB power adapter (not included) or a computer USB port.
2. Attach the magnetic ends of the cable to the charging ports on the TENS unit hosts.
3. The indicator light on the unit will show charging status (e.g., red for charging, blue/green for fully charged).

Stable Connection

Strong magnetic buckle, not easy to fall off



Figure 4.1: Charging the Comfytemp Wireless TENS Units.

4.2 Attaching Electrode Pads

Ensure your skin is clean and dry before applying the electrode pads for optimal adhesion and conductivity.

1. Carefully peel the electrode pad from its protective liner.
2. Attach the TENS unit host to the designated magnetic connector on the electrode pad.
3. Apply the electrode pad with the attached TENS unit to the desired area of your body. Ensure full contact with the skin.

Wires Gone, Comfort On

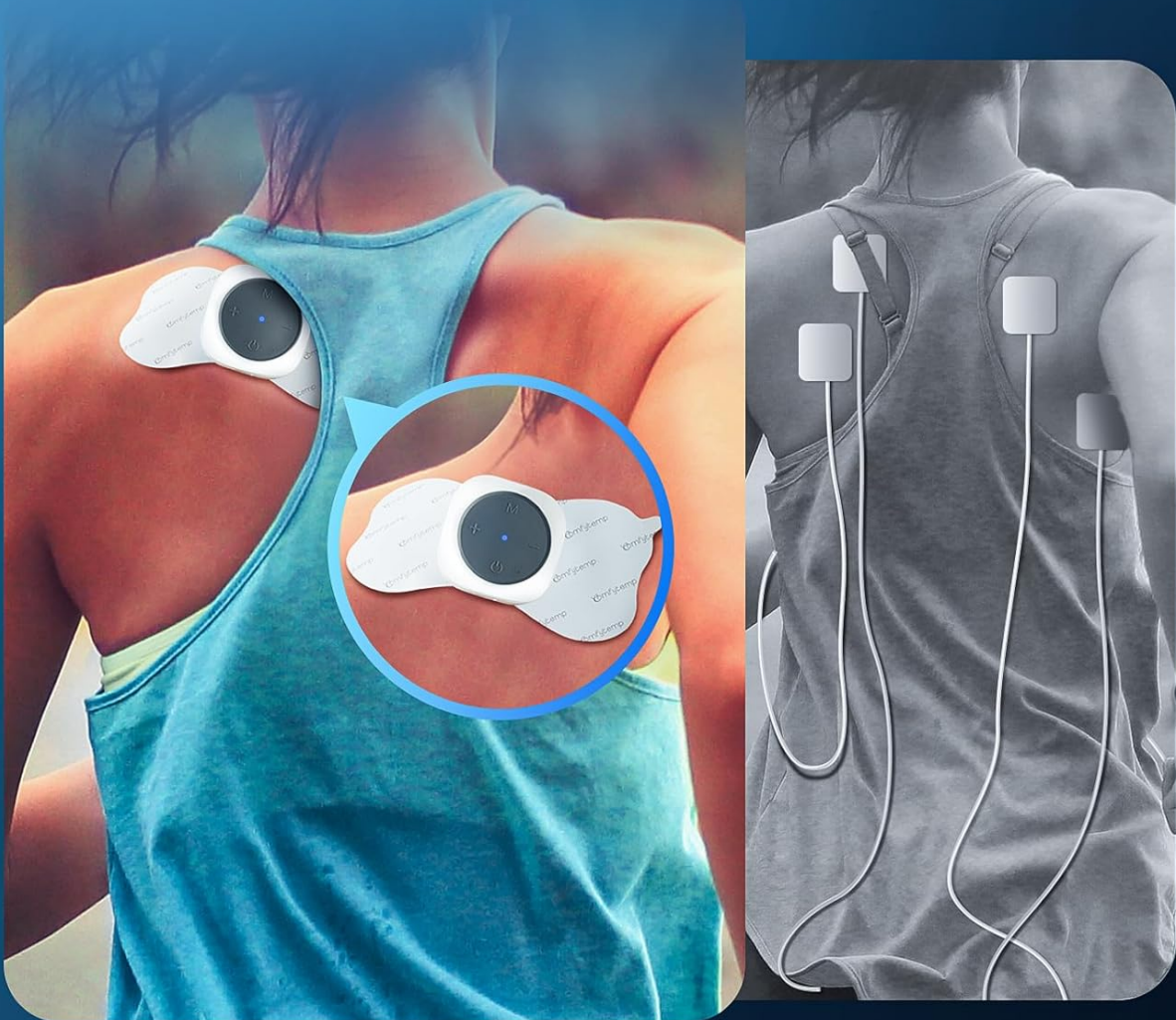


Figure 4.2: Steps for applying the TENS unit and electrode pads.

5. OPERATING INSTRUCTIONS

5.1 Using with the Comfytemp App (Recommended)

The Comfytemp app provides full control over 30 modes and 20 intensity levels, offering a personalized therapy experience.

1. **Download the App:** Scan the QR code in the physical manual or search for "Comfytemp" in your mobile device's app store (no registration required).
2. **Connect the Device:** Ensure Bluetooth (BLE) is enabled on your phone. Short press the power button on the TENS unit host to turn it on. The LED will flash blue, indicating it's ready to connect. Open the app and follow the on-screen prompts to pair.
3. **Select Mode and Intensity:** Once connected, use the app interface to select from 30 available modes and adjust the intensity level (0-20).
4. **Start/Stop:** Initiate or stop the therapy session directly from the app.

APP Control

Pain Relief at Your Fingertips



30 Modes



20 Intensities



One-click Stop



**Please Download APP.
No include remote!**

Figure 5.1: App control interface for the Comfytemp Wireless TENS Unit.



Figure 5.2: Step-by-step guide for app usage.

Your browser does not support the video tag.

Video: Demonstrates the use of the Wireless TENS Unit with both App Control and Button Control. This video is provided by the seller and shows the product in action.

Your browser does not support the video tag.

Video: Provides instructions on how to correctly use the TENS unit. This video is provided by the seller and offers practical guidance.

Your browser does not support the video tag.

Video: Answers frequently asked questions about the Wireless TENS Unit. This video is provided by the seller and addresses common user queries.

5.2 Using without the App (Button Control)

For quick and direct control, the TENS unit hosts also feature on-device buttons with 5 modes and 20 intensity levels.

1. **Turn On:** Short press the power button () on the unit.
2. **Select Mode:** Short press the 'M' button to cycle through the 5 available modes.
3. **Adjust Intensity:** Use the '+' and '-' buttons to increase or decrease the intensity level (0-20).
4. **Turn Off:** Long press the power button () to turn off the unit.



Figure 5.3: Button controls on the TENS unit host.



Figure 5.4: Step-by-step guide for button usage.

6. MAINTENANCE AND CARE

6.1 Electrode Pad Care

Proper care of electrode pads extends their lifespan and ensures optimal performance.

- After each use, place the pads back on their protective liners.
- If the electrode pads lose adhesion, gently rinse the adhesive surface with a small amount of water. Scrape off any surface residue.
- Allow the pads to air-dry completely (do not wipe). Adhesion should renew, enabling reuse for 60+ cycles.
- Store pads in a cool, dry place, away from direct sunlight.

6.2 Device Cleaning

Clean the TENS unit hosts with a soft, damp cloth. Do not use abrasive cleaners or immerse the units in water.

7. TROUBLESHOOTING

Problem	Possible Cause	Solution
No sensation or weak sensation	Unit not turned on. Intensity too low. Pads not making full contact with skin. Pads have lost adhesion or are worn out. Battery low.	Press power button to turn on. Increase intensity using '+' button or app. Ensure pads are firmly pressed onto clean, dry skin. Clean pads or replace with new ones. Recharge the unit.
App cannot connect to device	Bluetooth is off on phone. Device is off or not in pairing mode. Interference.	Enable Bluetooth on your phone. Short press power button on unit to turn on and enter pairing mode (flashing blue light). Move closer to the device; restart app and device.
Pads losing adhesion quickly	Skin not clean/dry. Excessive sweat/oil. Pads are old/worn.	Clean and dry skin thoroughly before application. Clean pads as per Section 6.1. Replace pads if cleaning does not restore adhesion.

8. SPECIFICATIONS




- **Model:** B0CLV5S4ZK
- **Modes:** Up to 30 (via App), 5 (via Button Control)
- **Intensity Levels:** 20
- **Battery:** Lithium Ion (included)
- **Charging:** USB-C Dual-head Cable
- **Dimensions (Package):** 6.77 x 4.45 x 2.32 inches
- **Weight (Package):** 8.15 ounces
- **Manufacturer:** Shenzhen Yicaiwangluo Technology Co.,Ltd
- **Country of Origin:** China

9. WARRANTY AND SUPPORT

Comfytemp products are designed for quality and reliability. For warranty information, technical support, or any inquiries regarding your TENS unit, please contact Comfytemp customer service directly. Refer to the product packaging or the official Comfytemp website for contact details.

Please retain your proof of purchase for warranty claims.

Related Documents - B0CLV5S4ZK

	Comfytemp K6145 TENS Unit User Manual and Specifications Comprehensive user manual for the Comfytemp K6145 TENS unit, covering device description, usage instructions, technical specifications, safety precautions, and troubleshooting for effective pain relief.
	Comfytemp K9224 Massage Heating Pad User Manual and Safety Guide Comprehensive user manual for the Comfytemp K9224 Massage Heating Pad, detailing product specifications, operating instructions, cleaning, troubleshooting, safety warnings, power bank usage, and warranty information.
	Comfytemp K4005 Red Light Pad: User Manual & Therapy Guide Explore the Comfytemp K4005 Red Light Pad for neck and shoulder pain relief and muscle recovery. This user manual details its red and near-infrared light therapy functions, specifications, usage, and safety precautions.
	Comfytemp K4044 XXL Flexible Red Light Pad User Manual and Guide User manual for the Comfytemp K4044 XXL Flexible Red Light Pad, detailing its functions, specifications, usage instructions, maintenance, troubleshooting, and safety precautions for red light therapy.
	Comfytemp K9274 Heating Belt User Manual Comprehensive user manual for the Comfytemp K9274 Heating Belt, detailing its features, operation, safety instructions, troubleshooting, and specifications. Learn how to use the far-infrared heat and support functions for temporary pain relief.



[Comfytemp K4009 Heating Pad User Manual](#)

User guide for the Comfytemp K4009 heating pad, detailing setup, app control, specifications, troubleshooting, and safety precautions for optimal use.