

[Manuals.plus](#) /

› [Redsack](#) /

› Redsack Electric Portable Bidet Instruction Manual

Redsack PB-E3

Redsack Electric Portable Bidet Instruction Manual

Model: PB-E3

PRODUCT OVERVIEW

The Redsack Electric Portable Bidet (Model PB-E3) is a compact and rechargeable personal cleaning device designed for travel and on-the-go hygiene. It features a high-pressure sprayer with adjustable water flow, an integrated water tank, and a foldable design for easy portability. This manual provides essential information for the safe and effective use, maintenance, and troubleshooting of your device.



Figure 1: RedSack Electric Portable Bidet, showing the device, its protective carrying case, and the USB-C charging cable.

SETUP AND CHARGING

Before first use, ensure the bidet is fully charged and properly filled with water.

Charging the Device

1. Locate the USB-C charging port at the base of the bidet. It is protected by a waterproof cover.

2. Open the waterproof cover.
3. Connect the provided USB-C cable to the bidet and the other end to a compatible USB power adapter (not included) or a computer USB port.
4. The charging indicator light will illuminate. A full charge provides approximately 200 uses.
5. Once charging is complete, disconnect the cable and securely close the waterproof cover to maintain water resistance.

Comes with a storage box

Portable travel, Protection product

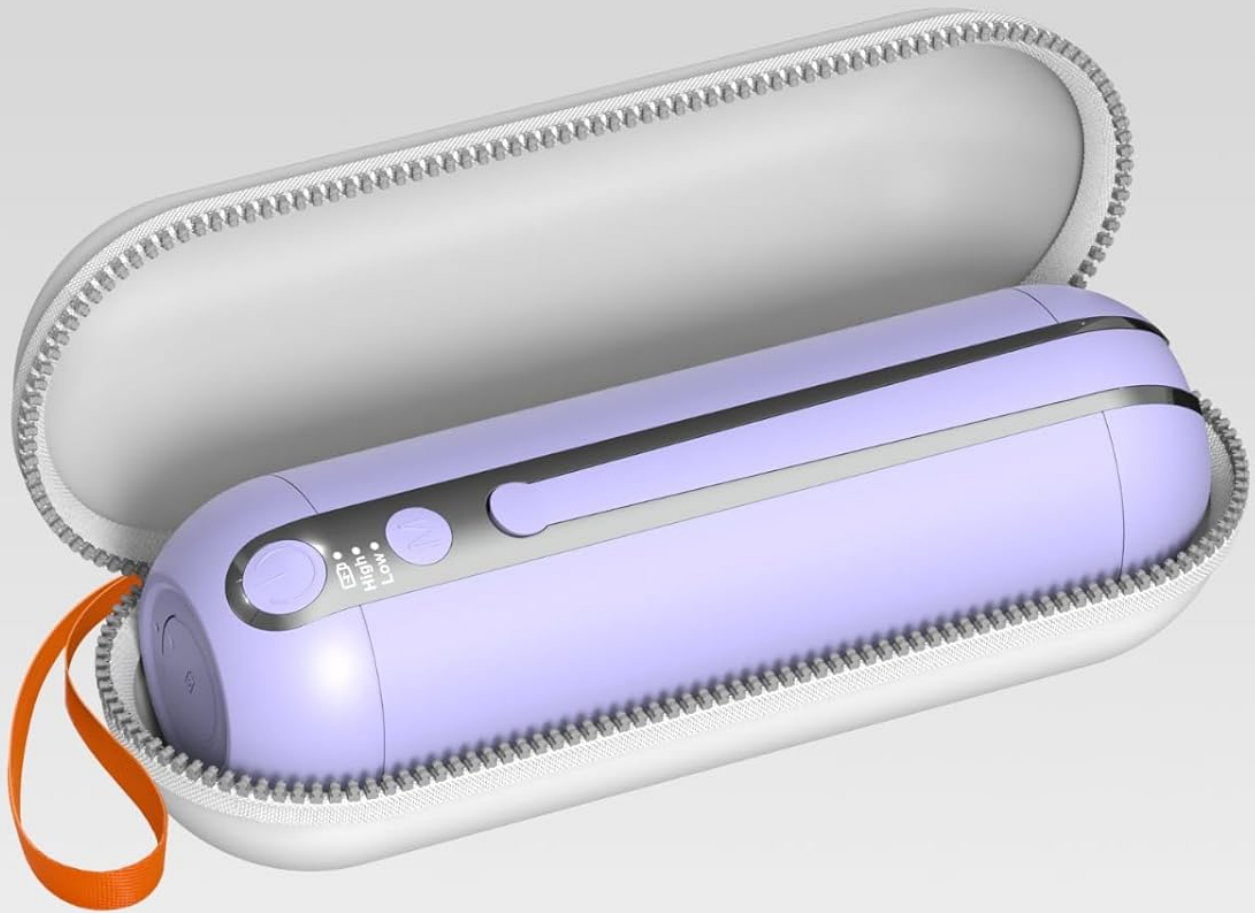


Figure 2: USB-C charging port and battery indicators on the bidet.

Filling the Water Tank

1. Ensure the bidet is in its folded, compact state.
2. Locate the flip-top lid on the top of the bidet.
3. Open the lid and fill the integrated 200ml (6.8oz) water tank with clean water. Warm water (35-40 degrees Celsius) can be used for increased comfort.
4. Close the lid securely to prevent leaks.

Precision cleaning



Hold the threaded part



Open to 180 degrees



Select gear



Avoid misuse



Figure 3: Demonstrating the quick water filling process for the bidet.

OPERATING INSTRUCTIONS

Follow these steps for effective and safe operation of your portable bidet.

1. **Prepare for Use:** Hold the bidet by the anti-slip threaded part at the end of the handle.
2. **Activate Power:** Fold open the spray bar to an angle between 150 and 180 degrees. The device is designed with a safety feature that prevents accidental activation; power will only connect when the spray bar is opened to the correct angle.
3. **Select Water Pressure:** Use the control buttons to select your desired water pressure:
 - **High:** For powerful cleaning.
 - **Low:** For a gentler rinse or massage-like experience.
4. **Initiate Spray:** Press the power button to start the water spray. The front-facing design of the controls allows for clear operation and accurate rinsing.
5. **Cease Spray:** Press the power button again to stop the water flow.
6. **Fold and Store:** After use, fold the spray bar back into its compact position. This action automatically disconnects the power, preventing accidental sprays.

Quick water filling

Integrated flip-top water tank, no need to rotate



Figure 4: Operational steps for precision cleaning, including holding, opening, and pressure selection.

Avoid accidental activation

No power



Trigger angle not reached

No power



Power connected

Start using

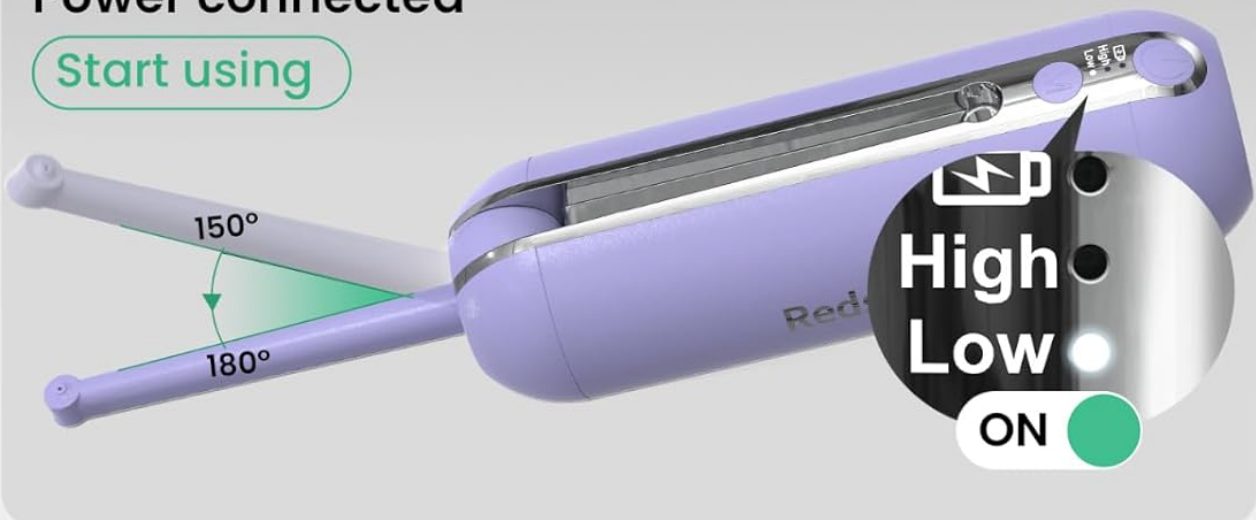


Figure 5: Understanding the folding switch design for power activation and accidental spray prevention.

MAINTENANCE

Proper maintenance ensures the longevity and hygiene of your Redsack Portable Bidet.

Cleaning

- After each use, empty any remaining water from the tank.

- Wipe down the exterior of the bidet with a clean, damp cloth.
- The bidet has an IPX6 waterproof certification, allowing for easy rinsing under running water. Ensure the USB-C port cover is securely closed before rinsing.
- Allow the bidet to air dry completely in a shaded area before storing.
- Do not use abrasive cleaners or harsh chemicals, as these may damage the device.

Storage

- Store the bidet in its provided protective storage case when not in use. This prevents external contamination and protects the device from bumps and damage.
- Keep the bidet in a cool, dry place, away from direct sunlight and extreme temperatures.

Rechargeable Portable Bidet

A full charge can spray \approx **200**
bottles of **200ML/6.8OZ** water



Figure 6: The bidet stored in its protective case, highlighting its portability and protection.

TROUBLESHOOTING

If you encounter issues with your Redsack Portable Bidet, refer to the following common problems and solutions.

Problem	Possible Cause	Solution
---------	----------------	----------

Problem	Possible Cause	Solution
Device does not spray water.	<ul style="list-style-type: none"> • Low battery. • Spray bar not fully opened. • Water tank is empty. 	<ul style="list-style-type: none"> • Charge the bidet fully. • Ensure the spray bar is opened to 150-180 degrees until power activates. • Refill the water tank.
Water pressure is weak.	<ul style="list-style-type: none"> • Low battery. • Set to 'Low' pressure mode. 	<ul style="list-style-type: none"> • Charge the bidet. • Switch to 'High' pressure mode.
Accidental spray when handling.	<ul style="list-style-type: none"> • Spray bar not fully closed. 	<ul style="list-style-type: none"> • Ensure the spray bar is fully folded into its compact position to disengage power.

SPECIFICATIONS

Feature	Detail
Model Name	portable bidet
Brand	Redsack
Color	Purple
Material	Plastic
Product Dimensions	2.6 x 2.6 x 8.3 inches
Item Weight	11.64 ounces
Item Model Number	PB-E3
Battery	1 Lithium Ion battery (included), 2000mAh
Water Tank Capacity	200ml (6.8oz)
Water Pressure Levels	High, Low
Waterproof Rating	IPX6

WARRANTY AND SUPPORT

Specific warranty details for the Redsack Electric Portable Bidet (Model PB-E3) are not provided within this manual. For warranty information, technical support, or customer service inquiries, please refer to the product packaging or contact Redsack customer support directly through their official website or the retailer where the product was purchased.

Legal Disclaimer: Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.



