

[Manuals.plus](#) /

> [Bitvae](#) /

> Bitvae R1 Rotating Electric Toothbrush Instruction Manual

Bitvae R1

Bitvae R1 Rotating Electric Toothbrush Instruction Manual

Model: R1
Brand: Bitvae

INTRODUCTION

This manual provides detailed instructions for the safe and effective use of your Bitvae R1 Rotating Electric Toothbrush. Please read thoroughly before first use to ensure optimal performance and longevity of your device.

Key Features

- **Superior Clean at a Dentist Level:** Equipped with a powerful motor, its rotating movements effectively clean each tooth and reach areas manual brushing might miss.
- **8 Replacement Heads for 2 Years of Use:** Includes eight round-end bristle heads, providing a long-lasting supply for consistent oral hygiene.
- **Gentle yet Effective Cleaning with 5 Modes:** Offers 'Clean', 'Sensitive', 'White', 'Gum Care', and 'Tongue' modes to cater to various oral care needs.
- **Single Charge Lasts Up to 60 Days:** Features a 3-hour fast-charge battery, providing up to 60 days of use, ideal for daily routines and travel.
- **2-Minute Smart Timer & IPX7 Waterproof:** A built-in timer pulses every 30 seconds to guide brushing across all mouth quadrants, and its waterproof design allows for safe use in the bathroom.

PACKAGE CONTENTS

Please verify that all items are included in your package:

- 1 Bitvae R1 Toothbrush Handle
- 8 Toothbrush Heads
- 1 USB Cable (Adapter not included—please charge the toothbrush using a 5V-1A adapter)
- 1 Instruction Manual

SETUP

Charging the Toothbrush

Before first use, fully charge the toothbrush.

1. Connect the USB cable to the charging port at the bottom of the toothbrush handle.
2. Plug the USB cable into a 5V-1A adapter (not included).
3. The charging indicator light will illuminate. A full charge takes approximately 3 hours and provides up to 60 days of use.



Image: The Bitvae R1 electric toothbrush handle with text indicating '3 Hours charge = 60 Days uses', highlighting its long battery life.

Attaching a Brush Head

Align the brush head with the metal shaft on the toothbrush handle and push down until it clicks into place.

8 Toothbrush Heads

*provide 2 years of use



Image: A hand is shown attaching a blue-bristled brush head to the top of the Bitvae R1 electric toothbrush handle, demonstrating the attachment process.

OPERATING INSTRUCTIONS

Power On/Off and Mode Selection

1. Press the power button to turn the toothbrush on.
2. Press the mode button (below the power button) repeatedly to cycle through the 5 cleaning modes: Clean, Sensitive, White, Gum Care, Tongue.
3. The indicator light next to the selected mode will illuminate.
4. To turn off the toothbrush, press and hold the power button for 3 seconds, or cycle through all modes until it turns off.

Customize Your Clean

with 5 modes



Clean



Sensitive



White



Tongue



Gum Care

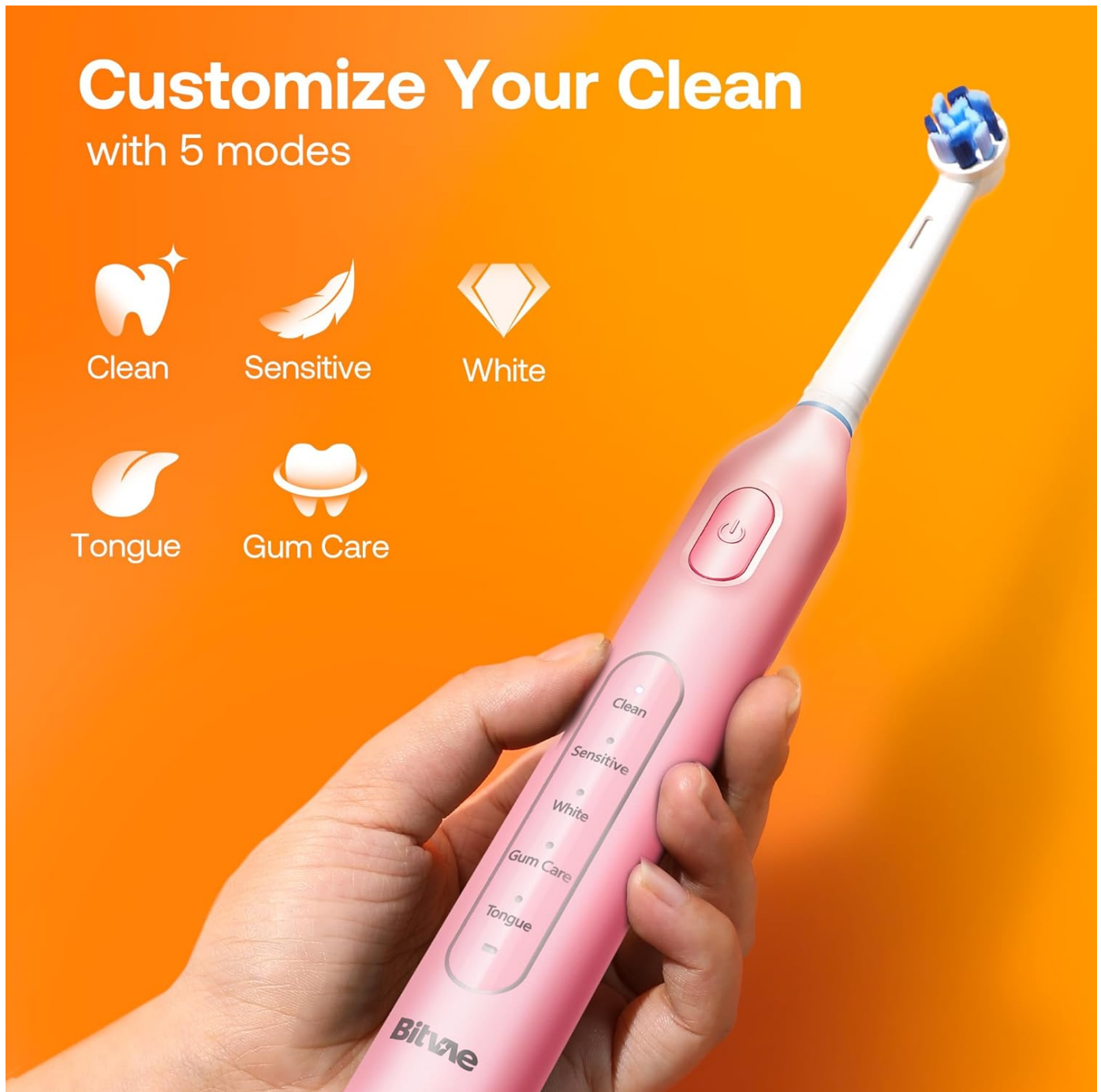


Image: A pink Bitvae R1 electric toothbrush handle is displayed, highlighting the five cleaning modes: Clean, Sensitive, White, Gum Care, and Tongue, each with an indicator light.

Brushing Technique

Apply a small amount of toothpaste to the brush head.

Place the brush head at a 45-degree angle to your teeth, gently moving it across each tooth. Do not scrub aggressively.

Smart Timer

The built-in 2-minute smart timer ensures you brush for the dentist-recommended duration. The toothbrush will pause briefly every 30 seconds to signal you to move to a different quadrant of your mouth.

2-Minute Smart Timer



Image: A visual representation of a mouth divided into four quadrants, with a timer indicating '2 MINS' total brushing time and '30S change brushing area', guiding users on proper brushing duration per section.

MAINTENANCE

Cleaning the Toothbrush

After each use, rinse the brush head thoroughly under running water.

Wipe the handle with a damp cloth. Do not immerse the handle in water for extended periods, despite its IPX7 rating.



Long Battery Life

3
HOURS
charge

=

60
DAYS
uses

Clean
Sensitive
White
Gum Care
Tongue

Bitvae

Image: A pink Bitvae R1 electric toothbrush is shown under running water, demonstrating its IPX7 waterproof capability for easy cleaning.

Replacing Brush Heads

Dentists recommend replacing brush heads every 3 months for optimal cleaning performance.

To replace, gently pull the old brush head off the shaft and snap a new one into place.

Provide an Ultimate Clean



Image: A collection of eight individually packaged replacement brush heads with blue and white bristles, designed for the Bitvae R1 electric toothbrush, emphasizing the 2 years of use provided.

TROUBLESHOOTING

Toothbrush Not Turning On

- Ensure the toothbrush is fully charged. Connect it to a 5V-1A adapter and charge for 3 hours.
- Check if the power button is being pressed firmly.

Reduced Performance

- Replace the brush head if it has been used for more than 3 months or if bristles appear worn.
- Ensure the toothbrush is fully charged.

SPECIFICATIONS

Feature	Detail
Model	R1
Power Source	Battery Powered (1 Lithium Ion battery)
Battery Life	Up to 60 days on a 3-hour charge
Modes	5 (Clean, Sensitive, White, Gum Care, Tongue)
Waterproof Rating	IPX7
Timer	2-minute smart timer with 30-second quadrant reminders
Brush Head Type	Rotating, Round
Bristle Material	Nylon
Handle Material	Plastic
Dimensions	8.89 x 1.06 x 1.1 inches
Weight	9.14 ounces

WARRANTY AND SUPPORT

The Bitvae R1 Electric Toothbrush comes with a 2-year warranty.

For customer support or warranty claims, please contact Bitvae directly. Refer to the contact information provided in your original product packaging or visit the official Bitvae website.

PRODUCT VIDEOS

Bitvae R1 Electric Toothbrush Overview

Your browser does not support the video tag.

Video: A seller-created video providing an overview of the Bitvae R1 Electric Toothbrush, showcasing its features and benefits for adults and kids.

Bitvae R1 Battery Life and Memory Function

Your browser does not support the video tag.

Video: A seller-created video discussing the battery life and memory function of the Bitvae R1 Electric Toothbrush, addressing common user questions about its performance.