

Osaki OS-Pro 4D DuoMax

Osaki OS-Pro 4D DuoMax SL-Track Massage Chair User Manual

INTRODUCTION

The Osaki OS-Pro 4D DuoMax is an advanced 4D massage chair engineered to deliver a comprehensive relaxation experience. It integrates intelligent sensors to scan various body metrics, customizing the massage to individual preferences and targeting specific areas for optimal therapeutic benefit.

SAFETY INFORMATION

Always consult the user manual for specific safety instructions and warnings provided by the manufacturer. Ensure proper body positioning during use. Avoid using the massager if you have certain medical conditions or implants; consult a physician if unsure. Adhere to all electrical safety guidelines for powered massagers, including proper grounding and avoiding use near water. Do not attempt to repair or modify the chair yourself.

SETUP

Unboxing and Placement

Carefully remove all packaging materials. Ensure the chair is placed on a level, stable surface. Due to its weight (354 Pounds), it is recommended to have assistance during placement. Allow adequate space around the chair for recline functions and air circulation. The chair dimensions are 68 inches (length) x 31.5 inches (width) x 50 inches (height).



This image illustrates the physical dimensions of the Osaki OS-Pro 4D DuoMax massage chair. The chair measures 68 inches in length, 31.5 inches in width, and 50 inches in height, providing a clear understanding of its footprint.

Power Connection

Connect the power cord securely to the chair's power input and then to a grounded electrical outlet. Ensure the power source matches the chair's voltage requirements. The chair operates via a corded electric power source.

OPERATING INSTRUCTIONS

Initial Use and Settings

Before first use, ensure the chair is properly assembled and connected to power. The chair features an intelligent health detection system that scans various metrics to tailor the massage. Utilize the 4D massage mechanism and the 8-roller track massage system for a comprehensive experience.



A detailed view of the massage chair's digital control panel, which features a touchscreen interface. The screen shows various massage settings, including shoulder, foot, 4D intensity, speed, and air pressure adjustments for different body areas. Power and pause buttons are also visible.

Key Features and Modes

- **Zero Gravity:** Reclines the chair to a position that elevates your legs above your heart, reducing pressure on the spine.
- **Dual Roller 4D Mechanism:** Provides advanced massage movements that extend in and out, up and down, and side to side, mimicking human hands.
- **Intelligent HealthPro:** Utilizes sensors to detect body metrics and customize massage programs.
- **8 Roller-Track Massage:** Offers extensive coverage along the SL-Track for a full-body massage.
- **2 Heating Cores:** Provides soothing heat therapy to enhance muscle relaxation.
- **Piano Key Massage:** A unique massage technique integrated into the chair's functionality.

Using the Control Panel

The chair is operated via its integrated control panel. Select desired settings based on target body areas and preferred massage intensity. Utilize the various techniques available, including the 4D massage mechanism and roller-track massage. Ensure adequate space around the massager during operation, especially when reclining.



This image displays the Osaki OS-Pro 4D DuoMax SL-Track Massage Chair, highlighting its sleek black design and ergonomic structure. The chair features a modern aesthetic with integrated armrests and a comfortable seating area.

MAINTENANCE

Cleaning

To clean the faux leather material, use a soft, damp cloth with a mild, non-abrasive cleaner. Avoid harsh chemicals or solvents that could damage the finish. Wipe dry immediately. Regularly dust the chair to prevent buildup.

Care and Storage

Keep the chair away from direct sunlight and extreme temperatures to prevent material degradation. When not in use for extended periods, it is recommended to cover the chair to protect it from dust. Ensure the power cord is unplugged when cleaning or moving the chair.

TROUBLESHOOTING

Common Issues

- **Chair not powering on:** Check if the power cord is securely plugged into both the chair and the electrical outlet. Verify the outlet is functional.
- **Massage rollers not moving:** Ensure the chair is not in a paused state. Check the control panel for any error messages.
- **Unusual noises:** If you hear abnormal sounds, immediately stop using the chair and contact customer support.

For issues not listed here, or if troubleshooting steps do not resolve the problem, please contact Osaki customer support.



SPECIFICATIONS





| Feature | Detail |
|--------------------|-----------------------|
| Product Dimensions | 68 x 31.5 x 50 inches |
| Item Model Number | OS-Pro 4D DuoMax |
| Manufacturer | Osaki |
| Material | Faux Leather |
| Item Weight | 354 Pounds |
| Power Source | Corded Electric |
| Use For | Back |

WARRANTY AND SUPPORT


For detailed warranty information, please refer to the documentation included with your purchase or visit the official Osaki website. For technical support, service, or parts inquiries, please contact Osaki customer service directly. Contact information can typically be found on the manufacturer's website or in your product registration materials.

Related Documents - OS-Pro 4D DuoMax

| | |
|---|---|
|  | <p>Osaki OS-Pro 4D DuoMax Massage Chair: Operation Instructions and User Manual</p> <p>This document provides comprehensive operation instructions and a user manual for the Osaki OS-Pro 4D DuoMax Massage Chair. It details safety precautions, installation steps, various massage functions and settings (including 4D intensity, air pressure, voice control, and health detection), maintenance, troubleshooting, and product specifications. Learn how to set up and operate your Osaki massage chair for optimal relaxation and therapeutic benefits.</p> |
|  | <p>Osaki OS-Pro 4D DuoMax Massage Chair Operation Instructions</p> <p>Comprehensive operation manual for the Osaki OS-Pro 4D DuoMax Massage Chair, covering safety precautions, installation, usage methods, maintenance, troubleshooting, and product specifications.</p> |

| | |
|---|---|
|  | <p>Osaki OS-Pro 4D DuoMax Massage Chair Operation Instructions</p> <p>Comprehensive operation instructions and safety guide for the Osaki OS-Pro 4D DuoMax massage chair, covering installation, usage, maintenance, troubleshooting, and product specifications.</p> |
|  | <p>Osaki OS-4D Pro Ekon+ Deluxe Multi-function Massage Chair User Manual</p> <p>User manual for the Osaki OS-4D Pro Ekon+ Deluxe Multi-function Massage Chair. Covers installation, operation, features, safety precautions, maintenance, and troubleshooting.</p> |
|  | <p>Osaki OS-PRO MAESTRO Massage Chair Instruction Manual</p> <p>Comprehensive instruction manual for the Osaki OS-PRO MAESTRO massage chair, covering safety guidelines, product parts, before-use procedures, controller instructions, operation guide, cleaning and maintenance, troubleshooting, and specifications.</p> |
|  | <p>Osaki OS-Honor Massage Chair User Manual and Safety Instructions</p> <p>Comprehensive guide to the Osaki OS-Honor massage chair, covering important safety instructions, features, operation, and specifications. Learn how to use your massage chair safely and effectively.</p> |


Documents - Osaki – OS-Pro 4D DuoMax

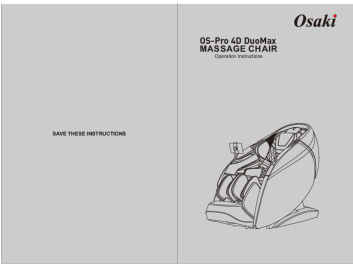

[\[pdf\]](#) User Manual Instructions Troubleshooting Guide Specifications Guide

1 Administrator User Guide Osaki OS Pro 4D DuoMax Massage Chair Bedplanet Manual cdn shopify s files 1835 3687 v 1692641793 |||


OS-Pro 4D DuoMax MASSAGE CHAIR Operation Instructions SAVE THESE INSTRUCTIONS Thank you for purchasing this product. Please read this manual instruction carefully in order to know how to use this product correctly before operation. Please pay more attention to some important safety warning and kee... lang:en score:50 filesize: 2.67 M page_count: 30 document date: 2023-08-21

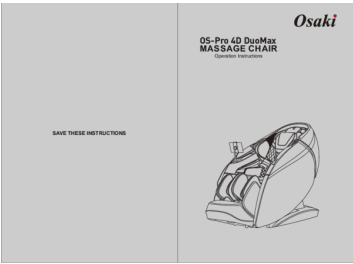

[Osaki OS-Pro 4D DuoMax Massage Chair: Operation Instructions and User Manual](#)

This document provides comprehensive operation instructions and a user manual for the Osaki OS-Pro 4D DuoMax Massage Chair. It details safety precautions, installation steps, various massage functions and settings (including 4D intensity, air pressure, voice control, and health detection), maintenance, troubleshooting, and product specifications. Learn how to set up and operate your Osaki massage chair for optimal relaxation and therapeutic benefits.


lang:en score:47 filesize: 2.63 M page_count: 30 document date: 2024-07-25

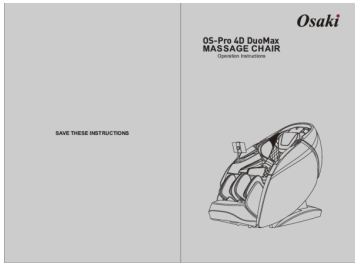

[Osaki OS-Pro 4D DuoMax Massage Chair Operation Instructions](#)

Comprehensive operation manual for the Osaki OS-Pro 4D DuoMax Massage Chair, covering safety precautions, installation, usage methods, maintenance, troubleshooting, and product specifications.

lang:en score:39 filesize: 2.73 M page_count: 30 document date: 2023-05-25

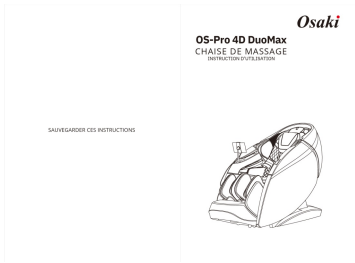


[\[pdf\] User Manual Instructions](#)

Osaki OS Pro 4D DuoMax Manual french DAGHY 4hSwl BAELpfuGUN8 Alia Khan DownloadOS Canadacdn shopify s files 1 0624 9634 2154 v 1718834001 |||

OS-Pro 4D DuoMax CHAISE DE MASSAGE INSTRUCTION D UTILISATION SAUVEGARDER CES INSTRUCTIONS Merci d avoir achet ce produit. Veuillez lire attentivement ce manuel afin de savoir comment utiliser correctement ce produit avant de l utiliser. Veuillez prter plus d attention certains avertissements de s...


lang:fr score:36 filesize: 18.28 M page_count: 30 document date: 2024-06-19

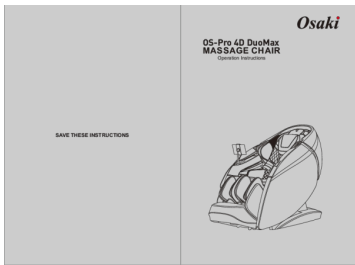

[Osaki OS-Pro 4D DuoMax Massage Chair Operation Instructions](#)

Comprehensive operation instructions and safety guide for the Osaki OS-Pro 4D DuoMax massage chair, covering installation, usage, maintenance, troubleshooting, and product specifications.

lang:en score:29 filesize: 2.63 M page_count: 30 document date: 2024-07-25

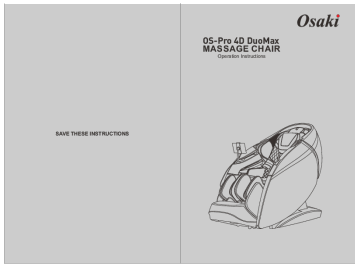

[\[pdf\] User Manual Instructions Troubleshooting Guide Specifications](#)

1 Administrator Osaki DuoMax Manual OS Pro 4D Massage Chair — Prime Chairs cdn shopify s files 2444 6637 v 1693256549 |||

OS-Pro 4D DuoMax MASSAGE CHAIR Operation Instructions SAVE THESE INSTRUCTIONS Thank you for purchasing this product. Please read this manual instruction carefully in order to know how to use this product correctly before operation. Please pay more attention to some important safety warning and kee...

lang:en score:24 filesize: 2.58 M page_count: 30 document date: 2023-08-28