



[Manuals.plus](#) /

> [comrelax](#) /

> Comrelax Shiatsu Back Massager with Heat (Model HY-6122) User Manual

comrelax HY-6122

Comrelax Shiatsu Back Massager with Heat

Model: HY-6122 User Manual

IMPORTANT SAFETY INFORMATION

WARNING: Always consult a healthcare professional before using this massager if you have any pre-existing medical conditions or concerns. Do not use the massager while driving or operating heavy machinery. This product is not intended for children.

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

PRODUCT OVERVIEW

The Comrelax Shiatsu Back Massager with Heat (Model HY-6122) is designed to provide a relaxing massage experience. It features deep kneading shiatsu nodes, an optional heating function, and vibration massage for the seat area. This versatile massager can be used on various chairs at home, in the office, or during travel.



Figure 1: Comrelax Shiatsu Back Massager with Heat, including power adapter and car adapter.

Key Features:

- **Deep Kneading Shiatsu Massage:** Four rotating nodes move up and down, simulating a professional massage therapist's hands. Features 3 adjustable intensity levels.
- **Soothing Heat Function:** Optional heat delivered through the Shiatsu nodes to help relieve muscle tension.
- **Vibration Seat Massage:** Three adjustable intensity levels for the hips and thighs. *Note: The seat area provides vibration only, no heat.*
- **Versatile Design:** Securely fits most chairs including office chairs, recliners, and sofas.
- **Smart Timer:** Automatic shut-off after 5, 10, or 15 minutes to prevent overheating.
- **Comfortable Materials:** Constructed with luxurious PU leather and soft mesh.



Figure 2: Illustration of 3D Rolling & Shiatsu massage action on the back.



Back Heating Nodes

Note: The Heat Is On The Roller Nodes, It's Not A Heating Pad.



Figure 3: Visual representation of the back heating nodes providing warmth.

PACKAGE CONTENTS

- Comrelax Shiatsu Back Massager with Heat (Model HY-6122)
- AC Power Adapter
- Car Adapter
- User Manual

SETUP

1. **Choose a Suitable Chair:** Place the massager on a sturdy chair such as an office chair, recliner, sofa, or dining chair. Ensure the chair provides adequate back support.
2. **Secure the Massager:** Use the integrated strap system to securely fasten the massager to the chair. This prevents slipping during use.

3. Connect Power:

- For home or office use, connect the AC power adapter to the massager's power input and then plug it into a standard electrical outlet.
- For car use, connect the car adapter to the massager's power input and then plug it into your vehicle's 12V power outlet.

4. **Position Yourself:** Sit comfortably on the massager, ensuring your back is aligned with the massage nodes. The massager is designed to comfortably accommodate users between 5'1" and 6'0" (155-183 cm).



Figure 4: The massager is designed to fit most common chair types.

Fit Most Chairs



Wooden Chair



Lazy Couch



Office Chair



Folding Recliner



Rocking Chair



Home Sofa

Figure 5: The massager can be used in a car with the included car adapter.


OPERATING INSTRUCTIONS

The massager is controlled by a wired remote. Familiarize yourself with the buttons and their functions before use.

Notes Before Purchase!



① Massage strength varies for everyone. If it feels too tough, throw a small blanket on the massager—it'll feel softer and cozier!

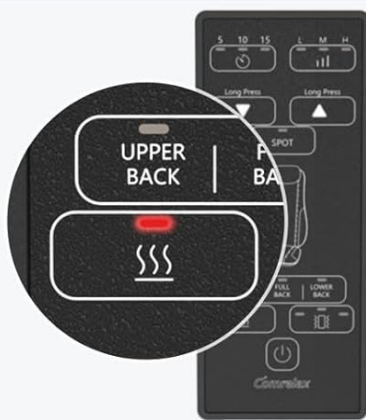

② The heating of the massager is  on the massage balls.

Click the heating button and the heating effect will be quite obvious after waiting for few minutes.

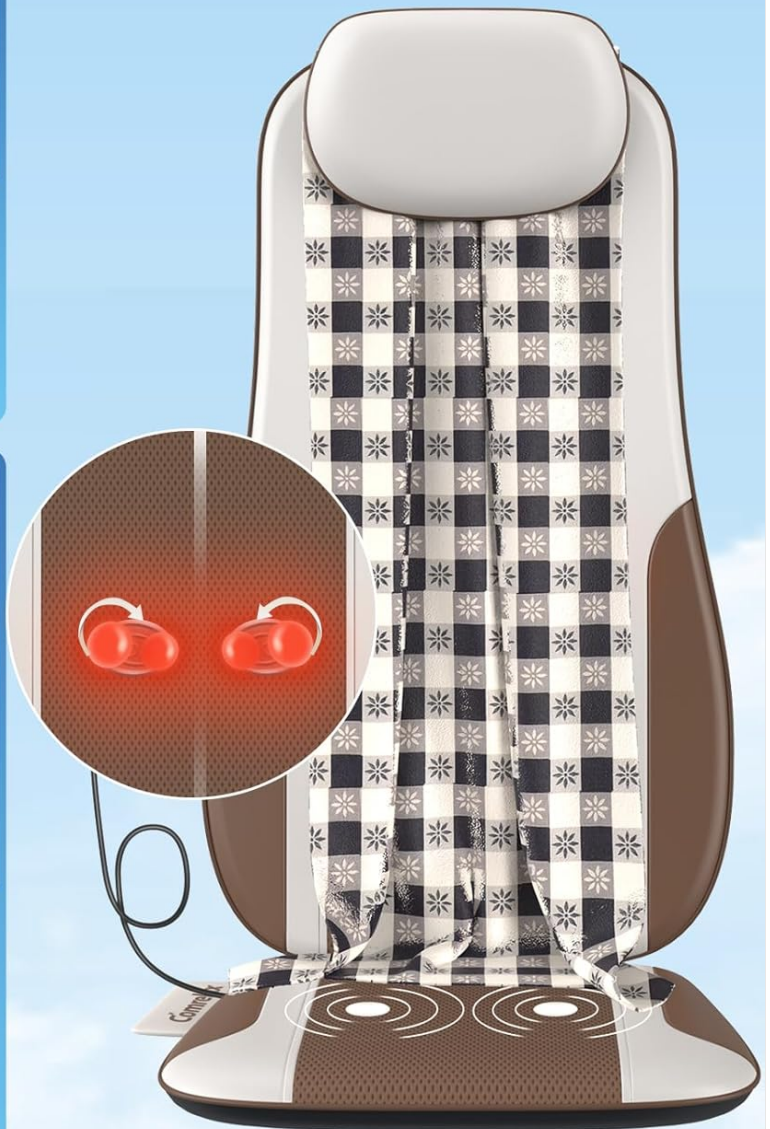


Figure 6: Comrelax Wired Controller Functions.

Controller Functions:

- **Power Button:** Press to turn the massager ON or OFF.
- **Time Control:** Press to set the massage timer to 5, 10, or 15 minutes. The massager will automatically shut off after the selected time.
- **Intensity Button:** Press to adjust the Shiatsu massage intensity across 3 levels (Low, Medium, High).
- **Shiatsu Area Button:** Select the desired massage area:
 - **FULL BACK:** Activates Shiatsu massage across the entire back.
 - **UPPER BACK:** Focuses Shiatsu massage on the upper back region.
 - **LOWER BACK:** Focuses Shiatsu massage on the lower back region.
- **Spot Button:** Press to pause the Shiatsu function at a specific spot. Use the up/down arrow buttons to adjust the node position.
- **Heat Button:** Press to turn the heating function ON or OFF. When heat is active, the lights on the roller nodes will

turn red.

- **Vibration Button:** Press to adjust the vibration intensity for the seat cushion across 3 levels, or to turn it ON/OFF.



Figure 7: The massager targets various areas including shoulders, waist, back, and thighs.

Comrelax Controller

Time Control:

Press to set the Massage Timer at 5, 10, 15 minutes.

Spot Button:

Press to keep shiatsu function in one area. Press arrow buttons to adjust the shiatsu nodes position up or down.

Heat Button:

Press this button to turn on/off the heating function, when the heating is on, the lights on the roller are red.



Intensity Button:

Press this button to adjust the shiatsu intensity

Shiatsu Area Button:

Select Full, Upper or Lower back to designate the area covered by the shiatsu massage.

Vibration Button:

Press this button to adjust the vibration modes, or turn on/off it

Figure 8: Tips for use, including using a blanket for softer massage and noting heat is on massage balls.

MAINTENANCE AND CARE

- **Cleaning:** Disconnect the massager from the power source before cleaning. Wipe the surface with a soft, damp cloth. Do not use abrasive cleaners or immerse the unit in water.
- **Storage:** When not in use, store the massager in a cool, dry place, away from direct sunlight and moisture. Avoid placing heavy objects on the massager.
- **Avoid Sharp Objects:** Keep sharp objects away from the massager to prevent damage to the fabric and internal components.

TROUBLESHOOTING

| Problem | Possible Cause | Solution |
|---------|----------------|----------|
|---------|----------------|----------|

| Problem | Possible Cause | Solution |
|----------------------------------|--|--|
| Massager does not turn on. | Not properly plugged in; power outlet issue; power button not pressed. | Ensure power adapter is securely connected to both the massager and a working electrical outlet. Press the power button firmly. |
| Heat function is not noticeable. | Heat function not activated; insufficient warm-up time. | Press the heat button on the controller. Allow a few minutes for the heating nodes to warm up. Note that heat is on the massage balls, not a separate heating pad. |
| Massage feels too intense. | Intensity setting is too high; direct contact with nodes. | Adjust the intensity to a lower setting using the controller. Place a thin blanket or towel between your back and the massager for a softer experience. |
| Vibration is not working. | Vibration function not activated or set to lowest intensity. | Press the vibration button on the controller to cycle through intensity levels or turn it on. |

SPECIFICATIONS

| | |
|-----------------------------|--|
| Brand | comrelax |
| Model Number | HY-6122 |
| Product Dimensions | 10 x 6 x 5 inches |
| Item Weight | 1.92 ounces |
| Material | Leather, Polyester |
| Power Source | Corded Electric |
| Use For | Back, Hip, Thigh |
| Manufacturer | Xiamen High top Electronic Technology Co.,Ltd. |
| Date First Available | November 15, 2023 |



Figure 9: Product dimensions and material details.

WARRANTY AND CUSTOMER SUPPORT

For information regarding product warranty, returns, or technical support, please refer to the purchase documentation or contact Comrelax customer service directly. Contact details can typically be found on the product packaging or the official Comrelax website.

Protection plans are available for purchase separately. Please refer to your retailer for details on available protection plans.

