

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [WELLFINE](#) /

› [WELLFINE Portable Cordless Water Dental Flosser F5030B User Manual](#)

WELLFINE F5030B

WELLFINE Portable Cordless Water Dental Flosser User Manual

MODEL: F5030B

For optimal oral hygiene and product longevity

Introduction

This manual provides instructions for the proper use and maintenance of your WELLFINE Portable Cordless Water Dental Flosser, Model F5030B. Please read this manual thoroughly before first use to ensure safe and effective operation.

Safety Information

Important Safety Precautions:

- If you have braces, implants, or other dental work, consult your dentist to ensure safe use.
- Do not immerse the main unit in water beyond the IPX7 waterproof rating.
- Keep out of reach of children.
- Use only with water or dentist-recommended solutions.

Package Contents

Verify that all items are present in the package:

- WELLFINE Portable Cordless Water Dental Flosser (Main Unit)
- 5 Specialized Nozzles (including 2 Standard Jet Tips, 1 Orthodontic Tip, 1 Curved Nozzle, 1 Tongue Scraper)
- USB Charging Cable
- User Manual



Image: Main unit, packaging, and included nozzles.

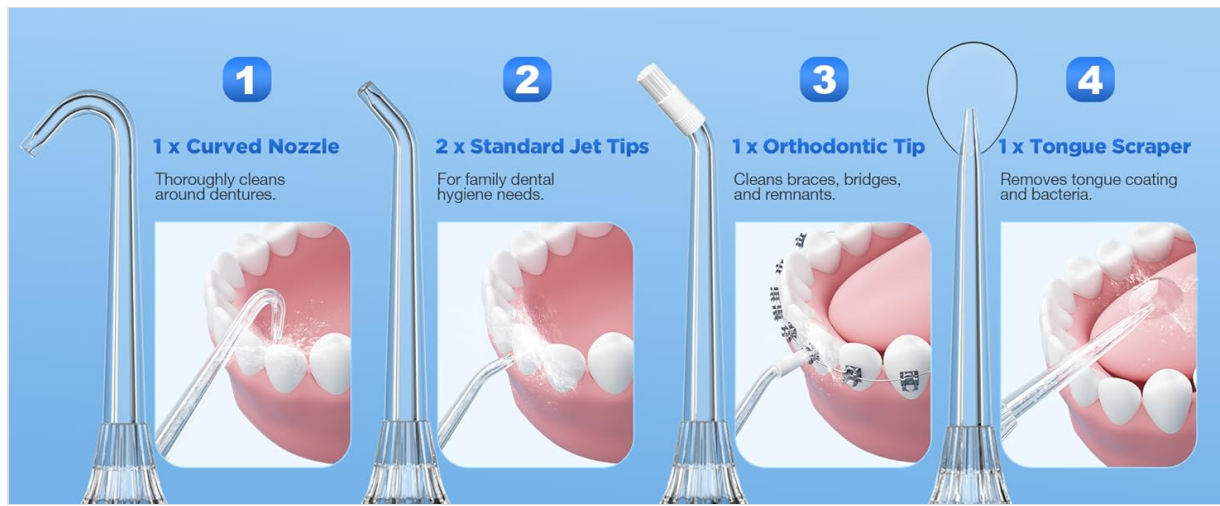


Image: Detailed view of the five specialized nozzles.

Product Features

- **Deep Cleaning Technology:** The flosser delivers 1,400-1,800 high-frequency pulses per minute with adjustable water pressure from 30-120 PSI. This action effectively removes food debris and plaque from between teeth and along the gumline, and gently massages gums.

EFFECTIVELY



REMOVES UP PLAQUE

Wellfine



0.3mm+0.3mm
Water column



Up to 1800
Pulsations per minute



Up to 120 psi
Strong water pressure



Up to 11
Cleaning modes

OTHERS

0.66mm
Water column

Up to 1500
Pulsations per minute

Up to 110 psi
Strong water pressure

Up to 3
Cleaning modes



Image: Illustration of deep cleaning effectiveness and specifications.

- **4 Customizable Modes:** Includes Normal, Soft, Pulse, and a Custom mode with 8 adjustable pressure settings (totaling 11 pressure settings). A memory function recalls the last used setting.

3+8 CLEANING MODES



SUITABLE FOR ALL AGES



NORMAL

Powerful
Clean Teeth

SOFT

For beginners
& sensitive teeth

PULSE

For Massage
& Protect Gums

CUSTOM

8 Water Pressures
adjustable by you



SOLVE MULTIPLE ORAL PROBLEMS

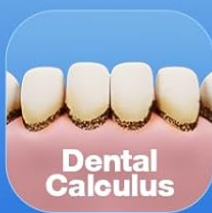


Image: Control panel displaying the four operating modes.

- **Long Battery Life & Portability:** A single charge provides approximately 30 days of use. The compact and lightweight design facilitates travel.

1 CHARGE USE 30 DAYS



TRIPS, TRAVEL, OFFICE, CAMPING

It is the Best Companion for Business Trips, Travel, Office and Camping.



Image: Flosser highlighting its long battery life and portability for travel.

- **Ergonomic Design & 300ML Water Tank:** Features an ergonomic grip for comfortable handling. The 300ML water tank allows for a complete cleaning session without frequent refills. The tank is detachable for easy filling and cleaning.

300ML LARGER



REMOVABLE WATER TANK

It's super clean!

✓ OURS

✗ OTHERS
Small water inlet
harbors dirt



More Larger
Double capacity
One water fill is enough

Easy
to Fill



Image: Illustration of the 300ML removable water tank and its features.

- **IPX7 Waterproof Rating:** The device is rated IPX7 waterproof, allowing for safe use in the shower.

IPX7 WATERPROOF



USE IN THE SHOWER

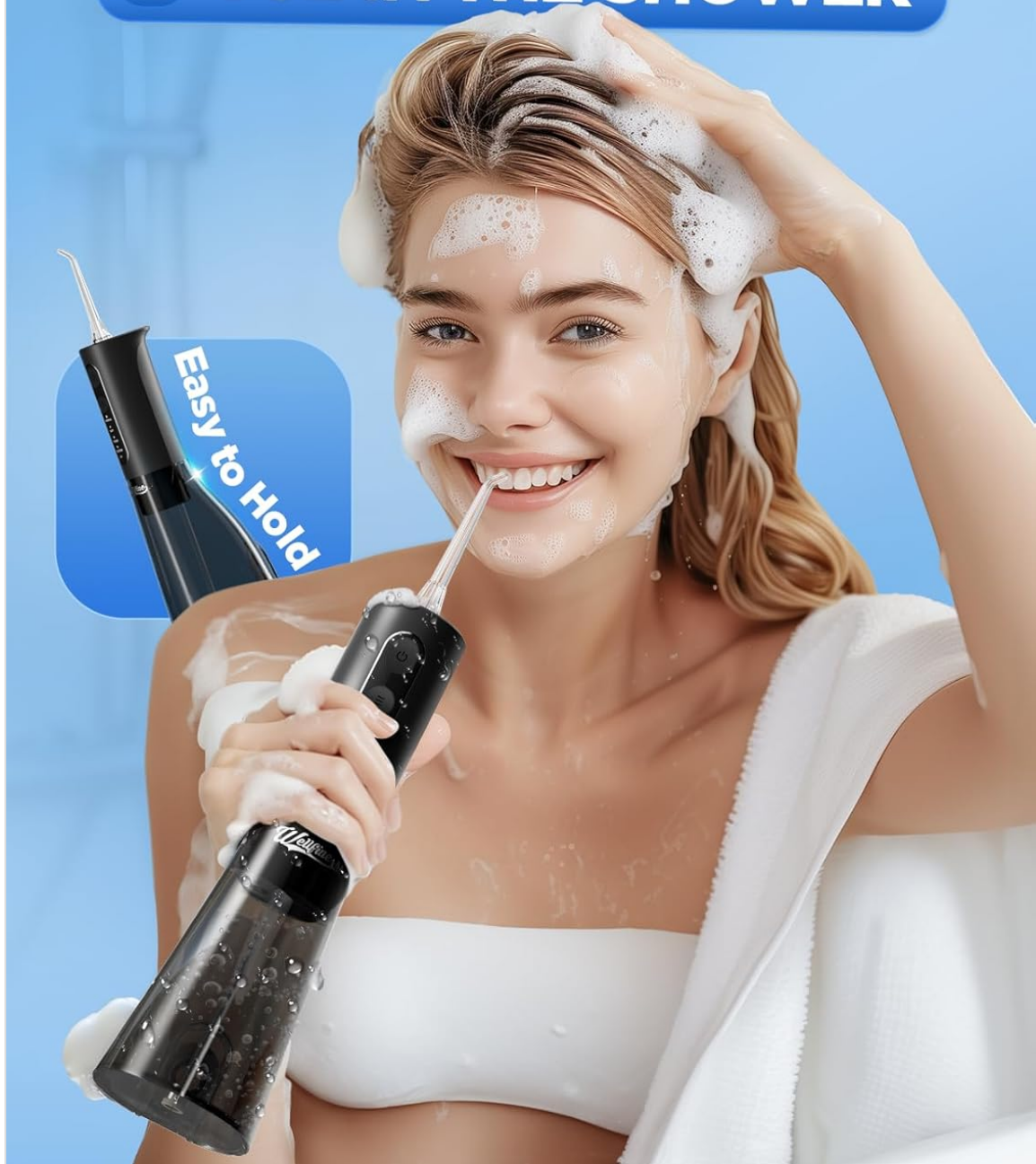


Image: Product shown in a shower environment, emphasizing IPX7 waterproofing.

Setup

1. **Initial Charge:** Fully charge the device before using it for the first time. Connect the USB charging cable to the flosser and a compatible USB power adapter (not included). The indicator lights will show charging status.
2. **Attach Nozzle:** Select the desired nozzle and insert it into the hole on the top of the flosser until it clicks into place.
3. **Fill Water Tank:** Open the water tank cover and fill the 300ML tank with warm water. Close the cover securely.

Operating Instructions

Follow these steps for effective use:

1. **Add Water:** Ensure the 300ML water tank is filled with warm water.
2. **Insert Nozzle:** Attach the desired nozzle firmly to the flosser.
3. **Select Mode:** Press the 'MODE' button to cycle through Normal, Soft, Pulse, and Custom modes. The corresponding indicator light will illuminate.
4. **Position Nozzle:** Lean over a sink. Place the nozzle tip into your mouth, directing it at a 90-degree angle to your gumline.
5. **Start Flossing:** Press the power button to turn on the device. Move the nozzle along your gum line and between your teeth, pausing briefly between each tooth.
6. **Turn Off:** After use, press the power button again to turn off the device.



Image: Visual guide for operating the water flosser.

Maintenance

Regular cleaning ensures optimal performance and hygiene:

- **After Each Use:** Empty any remaining water from the tank. Run the flosser for a few seconds to clear water from the internal components.
- **Clean Water Tank:** The detachable water tank can be removed and rinsed thoroughly with warm water. For deeper cleaning, use a mild soap and water solution, then rinse completely.
- **Clean Nozzles:** Remove the nozzle and rinse it under running water. Replace nozzles every 3-6 months, or sooner if they show signs of wear.
- **Clean Main Unit:** Wipe the exterior of the main unit with a soft, damp cloth. Do not use abrasive cleaners.

Specifications

Brand	WELLFINE
Model Number	F5030B
Power Source	Battery Powered

Water Tank Capacity	300ML
Water Pressure Range	30-120 PSI
Pulse Frequency	1,400-1,800 pulses/minute
Waterproof Rating	IPX7
Operating Modes	Normal, Soft, Pulse, Custom (with 8 pressure settings)
Package Dimensions	9.02 x 3.82 x 3.03 inches
Item Weight	13.44 ounces