

[manuals.plus](#) /

› [LUSTRE](#) /

› [LUSTRE ClearSkin Renew Pro Facewear LED Light Therapy Mask User Manual](#)

LUSTRE PR4001US

LUSTRE ClearSkin Renew Pro Facewear LED Light Therapy Mask User Manual

Model: PR4001US

1. INTRODUCTION

The LUSTRE ClearSkin Renew Pro Facewear LED Light Therapy Mask is an advanced skincare device designed to improve skin health using specific light wavelengths. This wireless, app-controlled mask integrates Blue (415nm), Red (630nm), and Near-Infrared (830nm) light therapy to target various skin concerns, including acne, fine lines, wrinkles, and overall skin tone and texture. This manual provides essential information for the safe and effective use of your device.



Image: LUSTRE ClearSkin Renew Pro Facewear LED Light Therapy Mask and its accompanying smartphone application.

LUSTRE
RENEW
LED light therapy face mask



Image: Two users demonstrating the LUSTRE ClearSkin Renew Pro Facewear, emphasizing key features.

Tested. Proven. Trusted.



FDA
cleared




Dermatologist
tested



Clinically proven
wavelengths



Award
winning

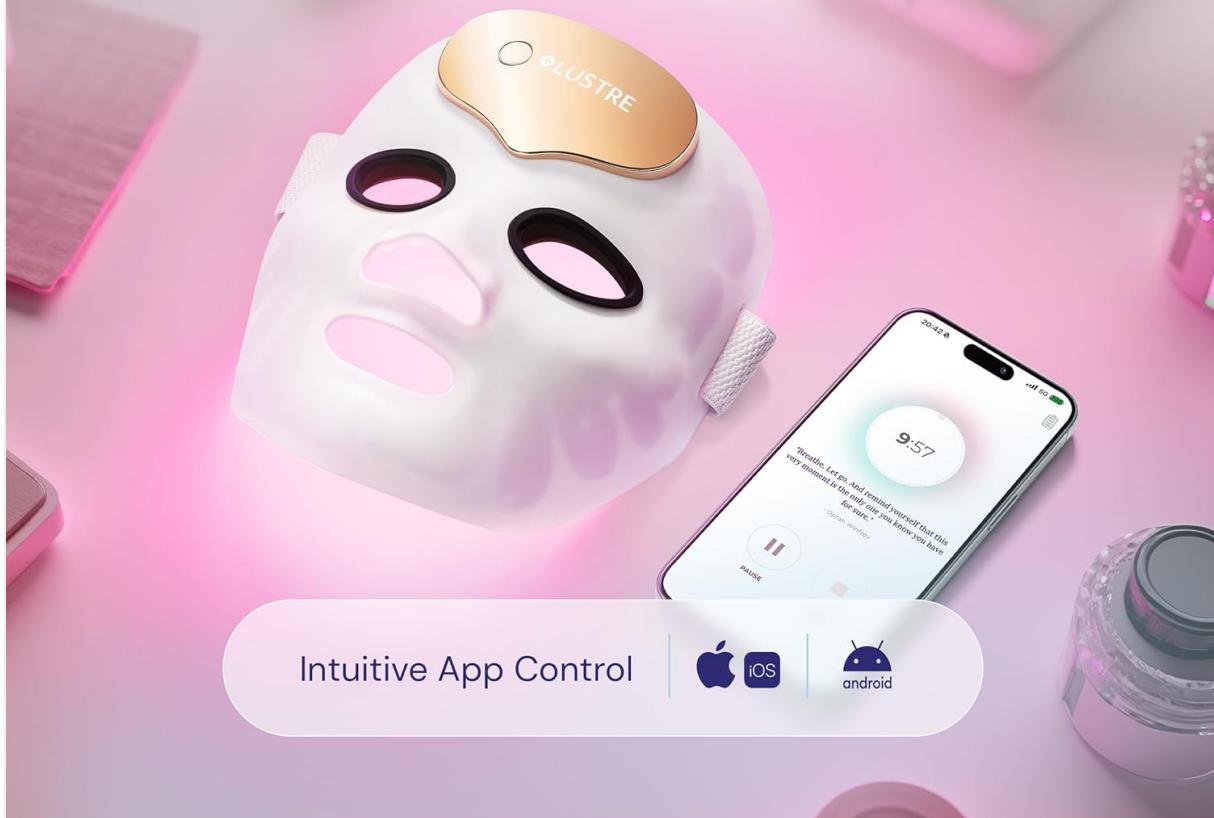

Intuitive App Control



Image: The LUSTRE ClearSkin Renew Pro Facewear highlighting its tested, proven, and trusted attributes.

Modern skincare

No chemicals, no needles, no irritation



Acne &
breakouts



Large
pores



Smile
lines



Fine lines
& wrinkles



Skin
elasticity



Uneven
tone

Endorsed by dermatologists in the USA and UK

Image: User wearing the mask, illustrating the range of skin benefits provided by the device.



Image: Product packaging showcasing the scientific backing and key benefits of the device.

2. IMPORTANT SAFETY INFORMATION

Please read all safety instructions carefully before using the LUSTRE ClearSkin Renew Pro Facewear. Failure to follow these instructions may result in injury or damage to the device.

- Do not look directly into the light emitted by the device.
- Do not place the device near water or on wet surfaces.
- Perform a light sensitivity test before the first use of the device.

- This device is designed for use by adults 21 years or older.
- This product shall not be used by adolescents under 21 years of age.
- Consult the full instruction manual for contraindications. Do not use the product if any contraindications apply to you.
- This device is suitable for Fitzpatrick Skin Types I-IV only.
- Refrain from using this device if you have severe acne.
- Distinguish your skin type or severity of acne according to the information provided in the full instruction manual before use.



Image: Product packaging displaying important safety warnings and usage guidelines.

3. WHAT'S IN THE Box

Verify that all components are present in the packaging:

- LUSTRE ClearSkin Renew Pro Facewear Light Therapy Mask
- Adjustable Velcro Head Strap
- Built-in Eye Protection (Removable)
- USB-C Charging Cable
- Luxury Pouch
- Instructions for Use Manual

What's in the Box?



Light Therapy Mask



Adjustable Velcro Head Strap



Built-in Eye Protection (Removable)



USB-C Cable



Luxury Pouch




Instructions for Use Manual

LUSTRE®ClearSkin

RENEW Pro Facewear Mask



5 060282 690394

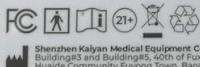

Shenzhen Kalyan Medical Equipment Co., Ltd
Building#3 and Building#5, 40th of Huxin Street,
Hutong Community Fuxing Subdistrict, Bao'an District,
Shenzhen, Guangdong 518103, China
Alba Innovation Centre, Alba Campus, Livingston
Scotland, UK, EH54 7GA. Infous@lustreskin.com
www.lustreclearskin.us

Quantity: 1 pcs
Input: 5V=1A
The Wavelength: 630nm, 830nm, 415nm
Battery: 3.7V, 1500 mAh, 5.55Wh



Risk group 1
CAUTION: Do not stare at the light source.
DO NOT STARE AT THE LIGHT SOURCE
EYES CAN BE DAMAGED BY THIS DEVICE



2024-06-15 900041
Product Number: PRA001US
(0)5060282690394
(1)20240615
(0)P210008KV003
(2)PRA00100004

Designed in the EU
Assembled in China
900012.2

Image: Diagram illustrating the components included in the LUSTRE ClearSkin Renew Pro Facewear package.



Image: Front packaging of the LUSTRE ClearSkin Renew Pro Facewear.

4. SETUP GUIDE

4.1 Initial Charging

Before first use, fully charge the device. Connect the provided USB-C cable to the charging port on the mask and to a suitable power adapter (not included). A full charge typically takes 3-5 hours. Ensure the battery is fully charged for optimal performance.

4.2 App Installation and Connection

Download the LUSTRE ClearSkin App from your device's app store (available for iOS and Android). Follow the in-app instructions to connect your mask via Bluetooth. The app allows for personalized session control and tracking.

4.3 Adjusting the Mask

Place the mask on your face and adjust the Velcro head strap for a secure and comfortable fit. The built-in eye protection should align with your eyes. Ensure the mask sits flush against your skin for effective light delivery.

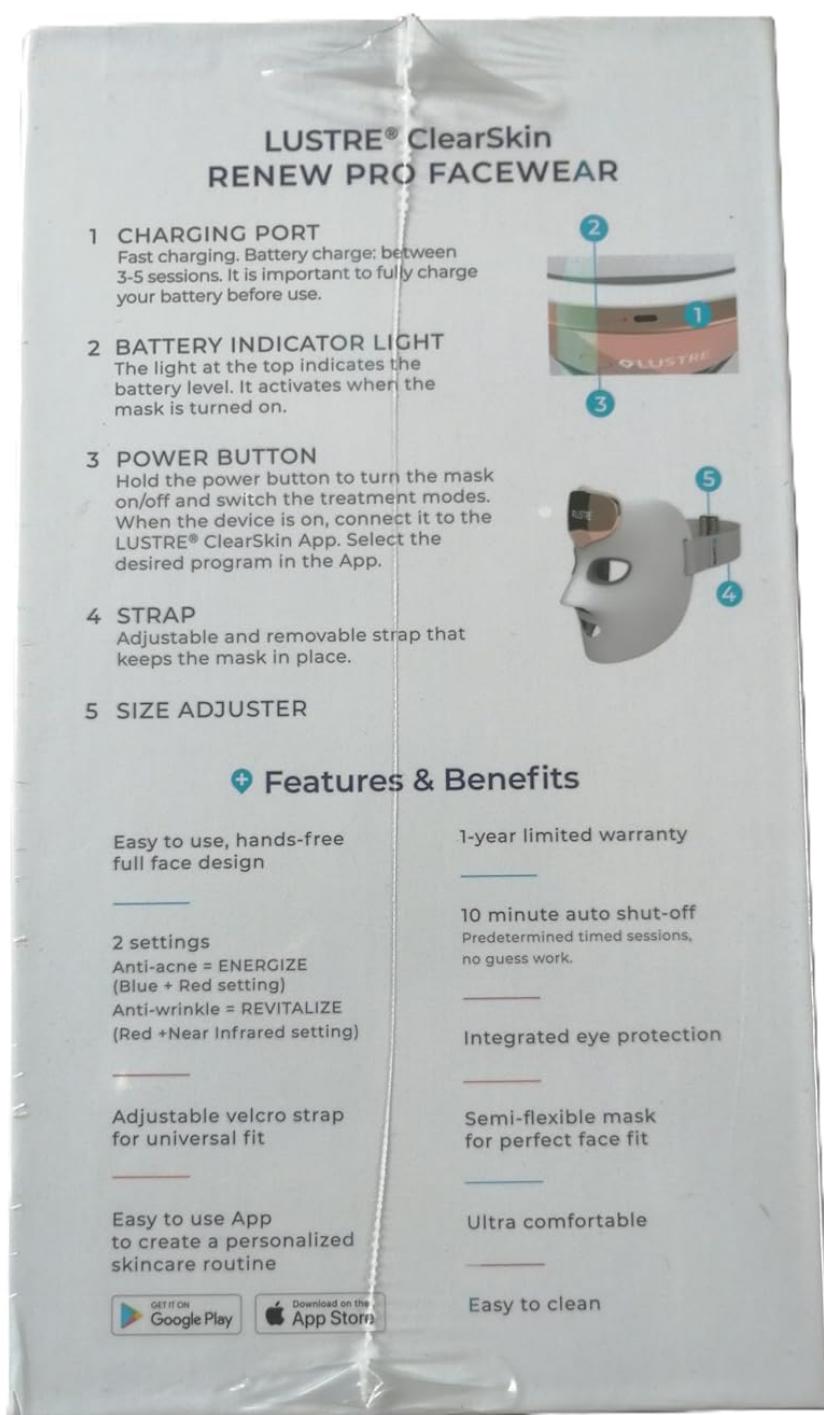


Image: Product packaging highlighting key features and components for setup and use.

5. OPERATING INSTRUCTIONS

5.1 Powering On and Mode Selection

Press the power button on the mask to turn it on. You can cycle through the pre-set treatment modes by

pressing the power button again, or connect to the LUSTRE ClearSkin App for advanced control.

5.2 Treatment Modes

- **ENERGIZE Mode (Anti-Acne):** Combines Blue (415nm) and Red (630nm) light. Blue light targets acne-causing bacteria and prevents breakouts. Red light helps reduce inflammation and redness, promoting skin renewal.
- **REVITALIZE Mode (Anti-Wrinkle):** Combines Red (630nm) and Near-Infrared (830nm) light. Red light boosts collagen production and reduces fine lines. Near-infrared light supports cellular recovery and skin repair.
- **Personalized Mode (via App):** Use the LUSTRE ClearSkin App to customize light therapy for specific facial areas (forehead, eyes, nose, cheeks, chin) based on your skin's needs.

5.3 Session Duration

Each treatment session is 10 minutes. The device features an automatic shut-off timer. Consistent daily use is recommended for optimal results.



Image: Product packaging detailing usage steps and the scientific benefits of different light wavelengths.

The science of beauty

Blue light

Helps reduce acne, prevent breakouts, and soothe irritation

Red light

Boosts collagen for anti-aging and anti-wrinkle effect

Near infrared light

Supports skin renewal, calms inflammation, and accelerates skin repair

CUSTOM
MODE

Personalize your facial

Save time by receiving both anti-aging and anti-acne treatment during one session: pick personalized treatment for specific areas of your face

Image: Scientific diagram explaining the penetration and effects of different light wavelengths on skin layers.

Swipe right for better skin



Customize light therapy modes and zones



Treat acne and wrinkles at the same time



Track progress and sessions effortlessly

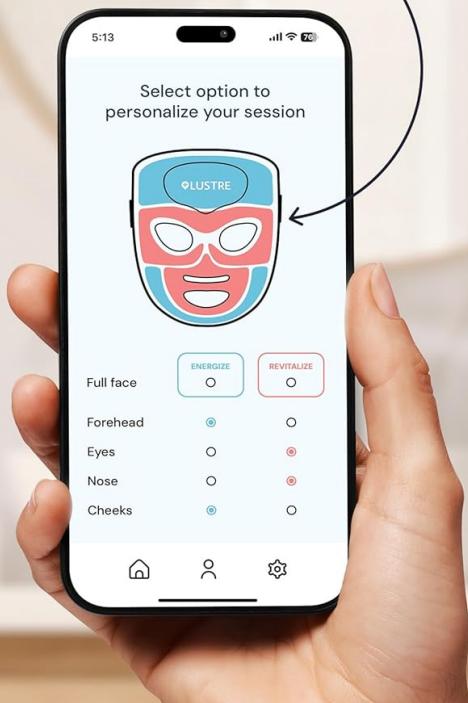


Image: The LUSTRE ClearSkin app interface demonstrating personalized facial treatment zone selection.

Pro tech. Home comfort.

Salon-grade skin transformation
available 24/7 in your personal space

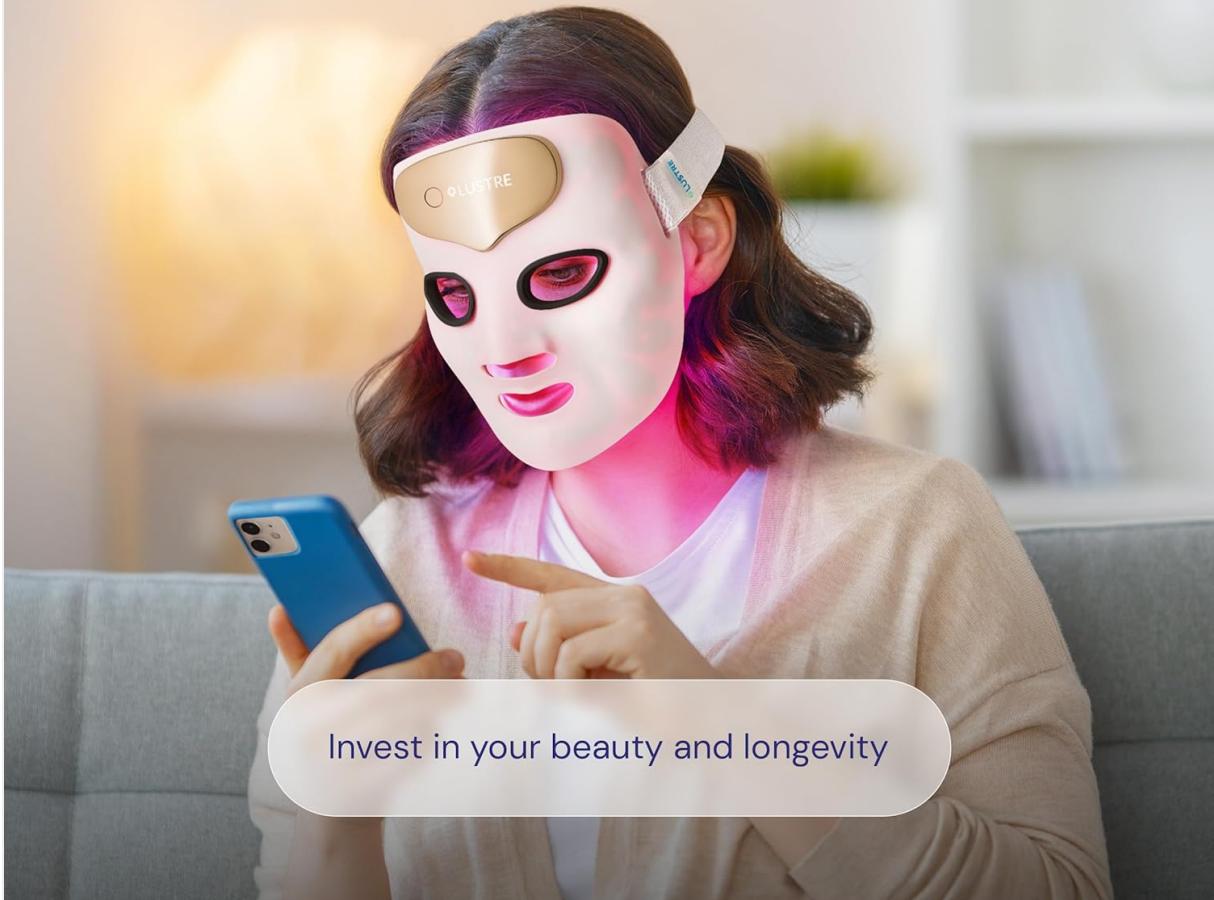


Image: User enjoying a light therapy session at home, highlighting the device's comfort and app control.

6. CARE AND MAINTENANCE

To ensure the longevity and hygiene of your LUSTRE ClearSkin Renew Pro Facewear:

- Clean the mask regularly with a soft, damp cloth. Do not use abrasive cleaners or immerse the device in water.
- Store the mask in its luxury pouch in a cool, dry place away from direct sunlight and extreme temperatures.
- Ensure the charging port is free from dust and debris.

7. PRODUCT SPECIFICATIONS

| Feature | Detail |
|-------------------|--|
| Model Number | PR4001US |
| Light Wavelengths | Blue (415nm), Red (630nm), Near-Infrared (830nm) |

| | |
|------------------|--------------------------------------|
| Battery Type | 1 Lithium Ion battery (included) |
| Battery Capacity | 3.7V, 1500 mAh, 5.55Wh |
| Charging Time | 3-5 hours |
| Session Duration | 10 minutes (automatic shut-off) |
| Connectivity | Wireless, App-Controlled (Bluetooth) |
| Dimensions | 8.82 x 7.95 x 5.04 inches |
| Weight | 2.2 Pounds |
| Manufacturer | Lustre |
| ASIN | B0CPMFKWF2 |

8. WARRANTY AND SUPPORT

8.1 Limited Warranty

The LUSTRE ClearSkin Renew Pro Facewear comes with a 1-year limited warranty from the date of purchase. This warranty covers manufacturing defects under normal use. Please retain your proof of purchase for warranty claims.

8.2 Customer Support

For further assistance, technical support, or warranty inquiries, please refer to the contact information provided in the full instruction manual or visit the official LUSTRE ClearSkin website. You can also find additional resources and FAQs within the LUSTRE ClearSkin App.