

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [ARTOLF](#) /

› [ARTOLF IPL Laser Hair Removal Device User Manual](#)

ARTOLF FZ-200

ARTOLF IPL Laser Hair Removal Device User Manual

Model: FZ-200

1. INTRODUCTION TO YOUR DEVICE

The ARTOLF IPL Laser Hair Removal Device (Model FZ-200) is an advanced at-home hair removal solution utilizing 3.0 IPL technology. It is designed to effectively and thoroughly remove unwanted hair by targeting hair follicles and disrupting the hair growth cycle. This device is suitable for use on various body areas including the face, legs, and arms.

What's Included:

- ARTOLF IPL Laser Hair Removal Device
- Power Adapter
- Protective Eyewear
- Razor
- Cleaning Cloth
- User Manual (this document)



Figure 1: The ARTOLF IPL Laser Hair Removal Device with included protective eyewear and razor.

For detailed instructions and safety guidelines, please refer to the complete user manual provided with your device. A digital version of the user manual can be found [here](#).

2. IMPORTANT SAFETY INFORMATION

Before using the ARTOLF IPL Laser Hair Removal Device, it is crucial to read and understand all safety instructions. Failure to do so may result in injury or device malfunction.

- **Eye Protection:** Always wear the provided protective eyewear during operation to shield your eyes from the intense light pulses.
- **Skin and Hair Type Compatibility:** This device is most effective on certain skin tones and hair colors. Refer to the compatibility chart below to determine if this device is suitable for you. Using the device on incompatible skin or hair types may lead to adverse reactions or ineffective results.
- **Skin Sensitivity:** For first-time users, it is recommended to start with the lowest energy level (Level 1) and gradually increase it based on your skin's adaptability.
- **Temperature:** The device may become warm during use. This is normal.

- **Ingredients:** The device uses IPL (Intense Pulsed Light) technology for hair removal. No topical ingredients are involved in the process.

Check if IPL is right for you

SKIN TONE / HAIR COLOUR		SKIN TONE						
								
White Grey		⊗	⊗	⊗	⊗	⊗	⊗	⊗
Light Blond		⊗	⊗	⊗	⊗	⊗	⊗	⊗
Red		⊗	⊗	⊗	⊗	⊗	⊗	⊗
Dark Blond		✔	✔	✔	✔	⊗	⊗	⊗
Light Brown		✔	✔	✔	✔	⊗	⊗	⊗
Dark		✔	✔	✔	✔	⊗	⊗	⊗



It's match!



Try a different hair removal method

Figure 2: Skin Tone and Hair Color Compatibility Chart. Check if IPL is right for you.

Note: The device is not suitable for white/grey, light blond, or red hair, nor for very dark skin tones. Consult the chart for detailed compatibility.

3. SETUP AND PREPARATION

1. **Cleanse Skin:** Ensure the treatment area is clean, dry, and free of any lotions, gels, or deodorants.
2. **Shave Hair:** Shave the hair in the treatment area. Do not wax, pluck, or use depilatory creams, as the hair follicle needs to be present for IPL to be effective.
3. **Connect Power:** Plug the power adapter into the device and then into a power outlet.
4. **Wear Protective Eyewear:** Put on the provided protective glasses before turning on the device.

4. OPERATING INSTRUCTIONS

The ARTOLF IPL device offers multiple energy levels and operating modes for a customized hair removal experience.

Power On/Off:

Press and hold the power button to turn the device on or off. The display will show the current energy level and remaining flashes.

Energy Levels:

The device features 9 adjustable energy levels (1-9), with 18J energy distributed more evenly. Higher energy levels provide more effective hair removal. For your first treatment, begin with Level 1 and gradually increase as your skin tolerates it. Use the appropriate buttons on the device to adjust the energy level.



Figure 3: The device offers 9 energy levels for varied intensity.

Operating Modes:

The device supports two operating modes:

- **Manual Mode:** Ideal for small, precise areas such as the bikini line, armpits, and upper lip. Press the flash button for each pulse.
- **Auto Mode:** Suitable for larger areas like arms, legs, and back. Hold the flash button down for continuous pulses.



Figure 4: Manual mode for small areas, Auto mode for large areas.

Application:

Place the device's light window firmly against your skin in the desired treatment area. Ensure full contact for optimal results. Press the flash button (Manual Mode) or hold it down (Auto Mode) to emit light pulses. Move the device slowly across the skin, ensuring each area is treated once per session.

600-1200nm Red Color Wave

Directly hit the hair follicles for precise hair removal

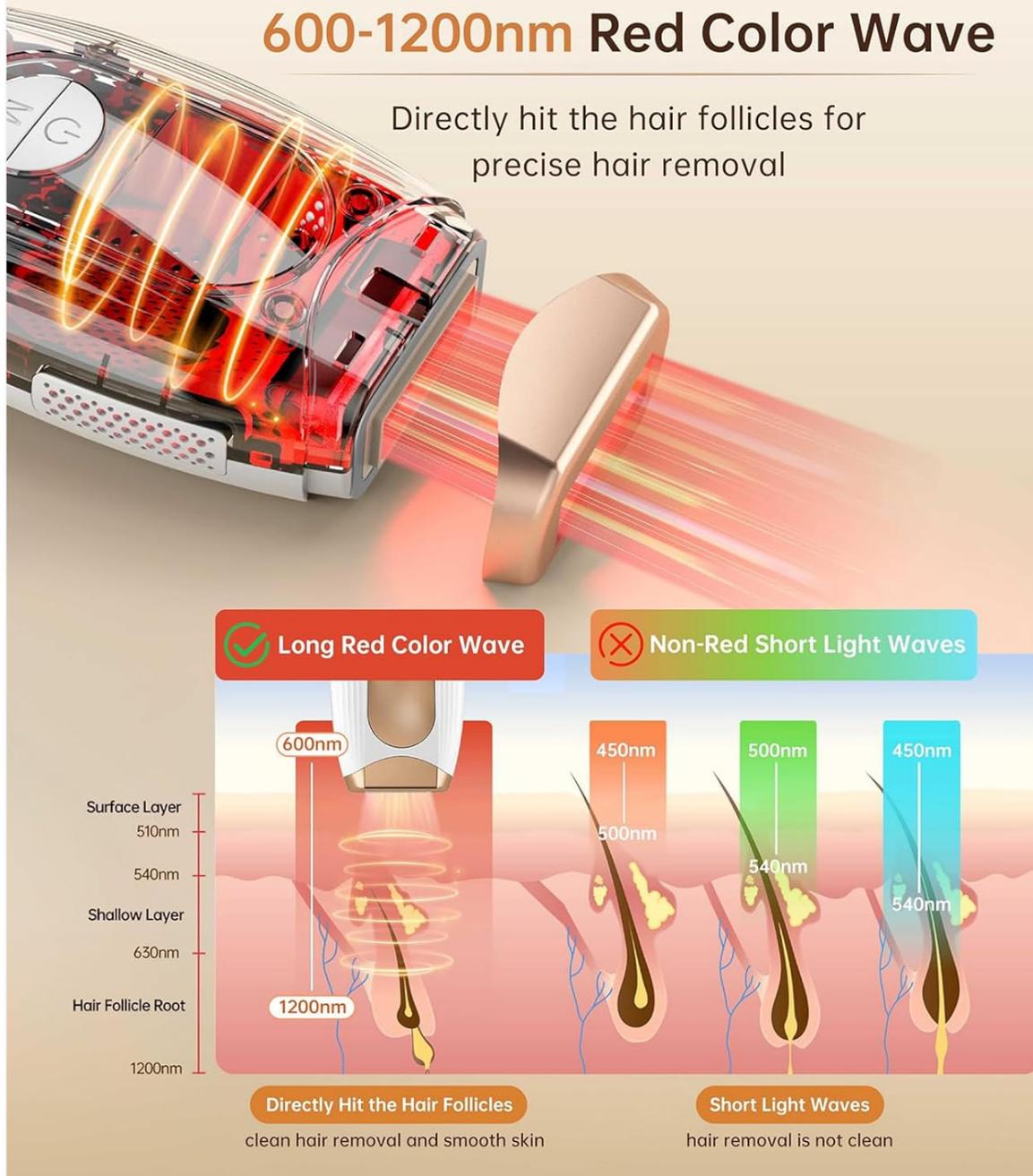


Figure 5: The 600-1200nm red color wave directly targets hair follicles for precise hair removal.


Hair Removal Cycle:

Consistent use is key for effective hair reduction. The typical hair removal cycle is as follows:

- **First Month:** Use 2-3 times per week.
- **Second Month:** Use 1-2 times per week.
- **Third Month Onwards:** Use 2 times per month for maintenance.

An efficient hair removal rate of up to 94% can be achieved after 7-8 weeks of consistent treatments.

Hair Removal Cycle

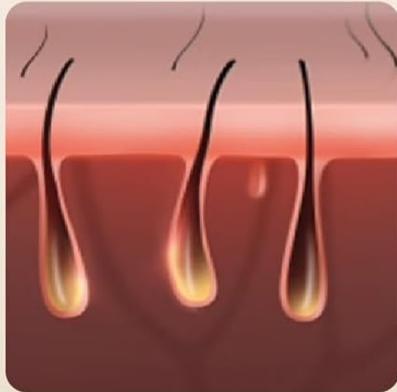

First month
(use 2-3 times per week)



Second month
(use 1-2 times per week)



Third month
(use 2 times per month)



Before After 8 weeks (95% hair reduction)

Figure 6: Visual representation of hair reduction over time with consistent use.

Additional Functions:

The upgraded device includes 3 functions: HR (Hair Removal), SC (Skin Care), and RA (Rejuvenation). These functions allow for skin care benefits alongside hair removal.

Safe, Fast and Effective IPL Hair Removal



Figure 7: Overview of device features including the 3-in-1 functions.

5. MAINTENANCE AND CARE

- **Cleaning:** After each use, ensure the device is powered off and unplugged. Use the provided cleaning cloth to gently wipe the light window and the device surface. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the device in a cool, dry place, away from direct sunlight and moisture. Keep it out of reach of children.

6. TROUBLESHOOTING

If you encounter any issues with your ARTOLF IPL device, please refer to the following general advice. For more specific problems, consult the full user manual.

- **Device Not Turning On:** Ensure the power adapter is securely connected to both the device and a working power outlet.
- **No Light Pulse:** Check if the light window is in full contact with your skin. The device has a safety mechanism that prevents flashing if contact is not complete. Ensure the device is not in standby mode.
- **Skin Discomfort:** If you experience discomfort, reduce the energy level. Ensure your skin is clean and dry before treatment.

- **Ineffective Results:** Ensure you are following the recommended treatment schedule and using an appropriate energy level for your skin and hair type. Consistency is key.

7. PRODUCT SPECIFICATIONS

Brand	ARTOLF
Model Number	FZ-200
Product Dimensions	10 x 8.11 x 3.23 inches
Item Weight	1.72 Pounds
Flashes	999,900+
Energy Levels	9
Wavelength Range	600 nm to 1200 nm

8. WARRANTY AND SUPPORT

The ARTOLF IPL Laser Hair Removal Device is designed for long-term use, with over 999,900 flashes, ensuring it is sufficient for lifetime use. This provides a cost-effective alternative to salon treatments.

For any product inquiries, technical support, or warranty claims, please refer to the contact information provided in your product packaging or visit the official ARTOLF website.

