

MEGELIN 05080508

MEGELIN 7 Colors LED Light Therapy Facial Mask

Model: 05080508 - User Manual

INTRODUCTION

This user manual provides comprehensive instructions for the safe and effective operation, maintenance, and troubleshooting of your MEGELIN 7 Colors LED Light Therapy Facial Mask. Please read this manual thoroughly before first use and retain it for future reference.

The MEGELIN LED Light Therapy Facial Mask is designed to provide various skin benefits through the application of specific light wavelengths. Its flexible silicone design and multiple light settings offer a customizable skincare experience.

SAFETY INFORMATION

WARNING: Read all safety warnings and instructions before using this product. Failure to follow the warnings and instructions may result in electric shock, fire, and/or serious injury.

- Do not use the device if you are pregnant, breastfeeding, or have a history of seizures triggered by light.
- Consult a healthcare professional before use if you have any pre-existing skin conditions, light sensitivity, or are taking medications that cause photosensitivity.
- Always use the provided eye protection during treatment to shield your eyes from direct light exposure.
- Do not use the device while it is charging.
- Keep the device away from water and do not immerse it in any liquid.
- Do not use the device if it is damaged or appears to be malfunctioning.
- Keep out of reach of children.
- This device is for cosmetic use only and is not intended to diagnose, treat, cure, or prevent any medical condition.

PACKAGE CONTENTS

Verify that all items listed below are included in your package:

- 1 x LED Light Face Mask

- 1 x Instruction Manual
- 1 x Charging Cable (Type-C)
- 1 x Eye Protection (Eye Pads)
- 1 x Fixing Strap

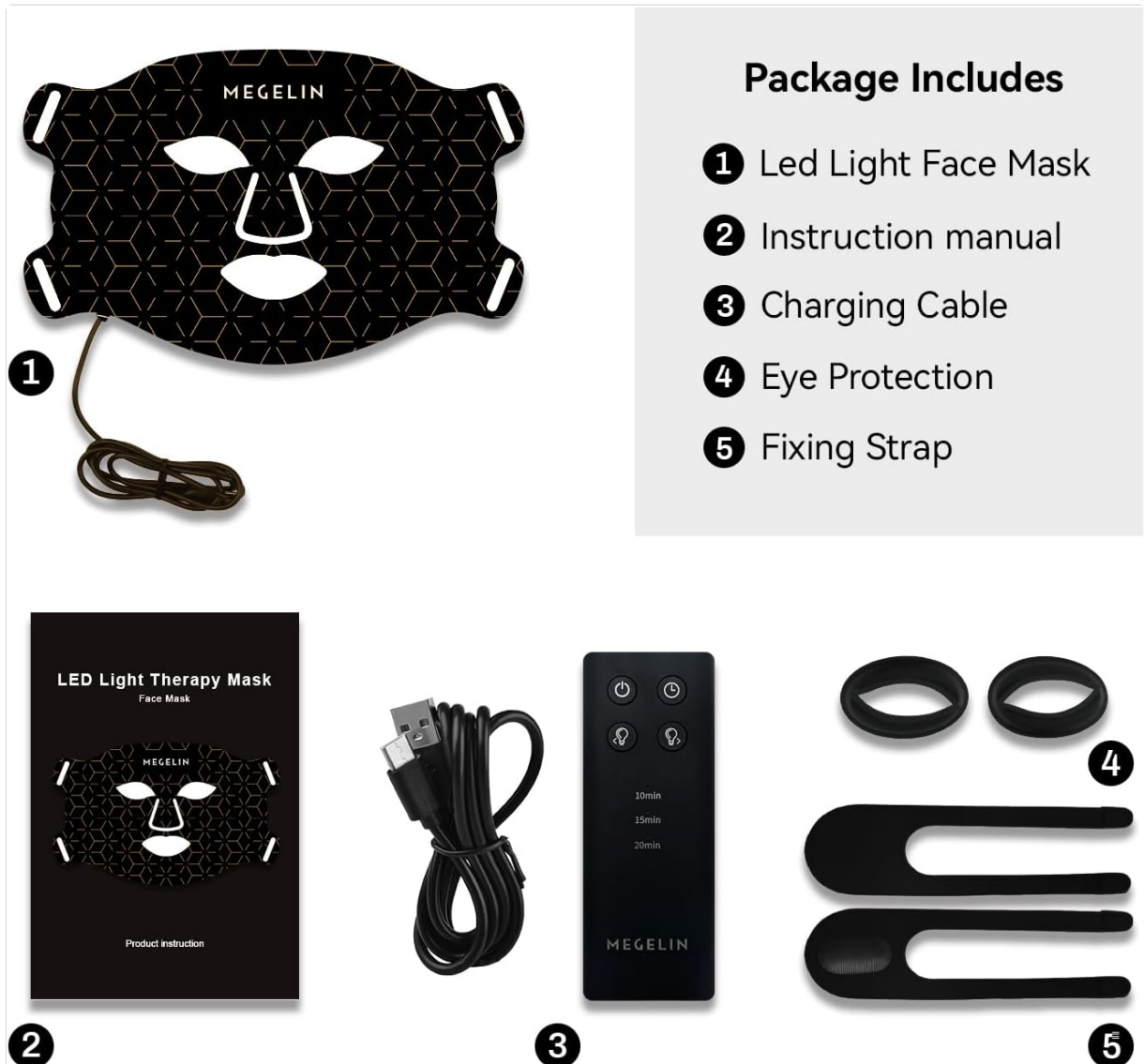


Image: The package contents laid out, including the LED mask, instruction manual, charging cable, eye protection, and fixing strap.

SETUP

Charging the Device

Before first use, fully charge the device. A full charge typically takes approximately [Insert Charging Time if available, otherwise omit].

1. Connect the Type-C charging cable to the charging port on the mask's controller.
2. Connect the other end of the cable to a suitable USB power adapter (not included) and plug it into a power outlet.
3. The battery indicator light on the controller will flash during charging and remain solid when fully charged.

How to recharge the battery?



1

Battery and Type c



2

Connect the type c to plug



3

Choose one of the type c sockets



4

The battery will be flashing

Image: Visual guide showing the Type-C cable, controller, and power adapter connection for charging the mask.

Preparing for Use

1. Ensure your face is clean and dry before applying the mask. Remove all makeup and skincare products.
2. Install the provided eye pads into the designated areas on the mask to protect your eyes from the LED light.
3. Adjust the fixing strap to ensure a comfortable and secure fit on your face. The mask should sit snugly without causing discomfort.

Easy-To-Use



Image: A sequence of images demonstrating how to clean the face, install eye pads, and adjust the mask for wear.

OPERATING INSTRUCTIONS

Using the Mask

1. After preparing your face and fitting the mask, press and hold the power button on the controller for 2 seconds to turn the mask on.
2. Press the power button briefly to cycle through the 7 available light colors. Select the desired light mode for your treatment.
3. Use the timer buttons (10min, 15min, 20min) on the controller to set your desired treatment duration.
4. Relax and enjoy the light therapy session for the selected duration. The mask will automatically turn off when the time is complete.
5. After the session, remove the mask and proceed with your regular skincare routine.

Easy-To-Use





1 Keep face clean



2 Install the eye pads and adjust straps to wear



3 Press the button 2 seconds
To turn on and off,
Press once to switch to a color



4 Lie down or sit back
And enjoy the skin
For 10-20 minutes



5 Follow your normal
Skin care routine

Image: A sequence of images showing the user turning on the mask, selecting settings, and relaxing during treatment.

Light Modes and Their Benefits

The MEGELIN LED Light Therapy Facial Mask features 7 distinct light wavelengths, each targeting specific skin concerns:

Light Color	Wavelength	Primary Benefit
Red Light	630nm	Anti-aging, increases collagen production.
Yellow Light	590nm	Fades spots, reduces skin sensitivity.
Green Light	550nm	Balances pigmentation, fades scars.
Cyan Light	490nm	Boosts cell metabolism.
Purple Light	420nm	Improves dull skin, tightens skin.
Blue Light	390nm	Calms sensitive skin, reduces creases.
Infrared Light	850nm	Penetrates superficial and deep tissues for enhanced rejuvenation.

7 Colors LED Light

Meet Different Skin Needs



Infrared Light

850nm

Penetrates superficial,
Deep tissues



Red Light

630nm

Anti-aging,
Increase collagen



Yellow Light

590nm

Fade spots,
Reduce skin sensitivity



Green Light

550nm

Balance pigmentation,
Fades scars



Cyan Light

490nm

Boost cell,
Metabolism



Purple Light

420nm

Improve dull skin,
Tightens skin



Blue Light

390nm

Clams sensitive skin,
Crease reduction

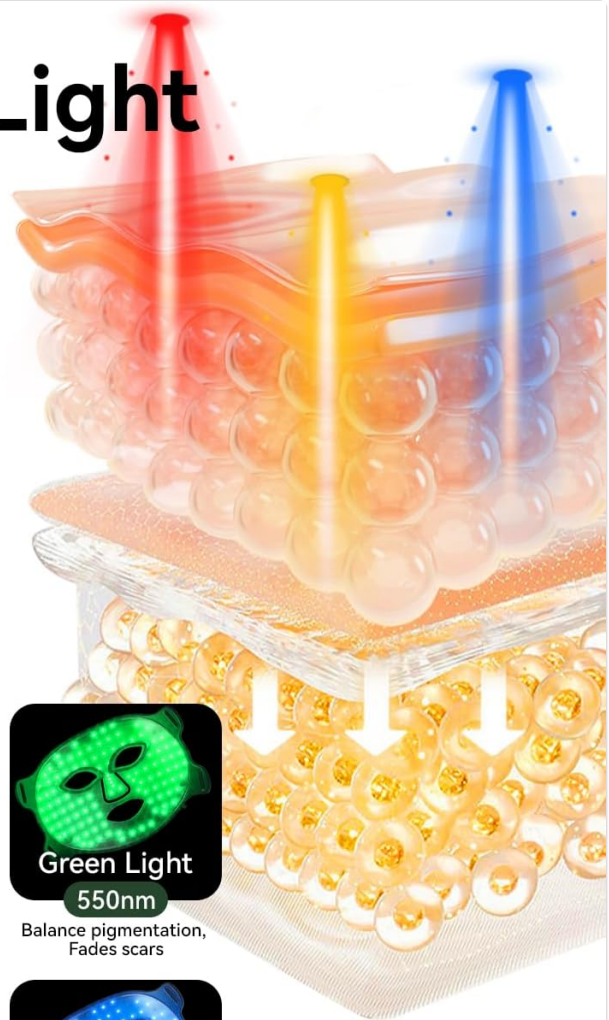


Image: A visual representation of the different LED light colors, their wavelengths, and how they penetrate the skin to provide various benefits.

Recommended Usage Frequency: It is recommended to use the Facial Light Therapy 3-4 times a week, with each session lasting 10-20 minutes. Consistent use for 4-6 weeks may help improve skin condition. Individual results may vary.

MAINTENANCE AND CARE

- **Cleaning:** After each use, gently wipe the mask with a soft, damp cloth. Do not use abrasive cleaners, alcohol, or harsh chemicals. Ensure the mask is completely dry before storing.
- **Storage:** Store the mask in a cool, dry place, away from direct sunlight and extreme temperatures. Keep it in its original packaging or a protective case to prevent damage.
- **Battery Care:** To prolong battery life, avoid fully discharging the battery frequently. Charge the device regularly, even if not in frequent use.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery; improper power button press.	Charge the device fully. Ensure you press and hold the power button for 2 seconds.
Lights are dim or flickering.	Low battery.	Recharge the device.
Cannot change light color.	Incorrect button press.	Ensure you are pressing the power button briefly to cycle through colors, not holding it.
Mask feels uncomfortable.	Strap too tight; improper positioning.	Adjust the fixing strap for a looser fit. Reposition the mask on your face.
Device stops working after a short period.	Battery depleted; timer expired.	Recharge the device. Check the timer setting; the device automatically turns off when the set time expires.

If the problem persists after attempting these solutions, please contact customer support for further assistance.

SPECIFICATIONS

- **Model Number:** 05080508
- **Dimensions:** 14.96 x 10.62 x 2.75 inches
- **Weight:** 3.15 Pounds
- **LED Count:** 288 LEDs (3-in-1 Chip)
- **Light Wavelengths:** 7 Colors (Red 630nm, Yellow 590nm, Green 550nm, Cyan 490nm, Purple 420nm, Blue 390nm, Infrared 850nm)
- **Power Density:** 35mW/cm²
- **Battery:** 1 Lithium Ion battery (included)
- **Charging Interface:** Type-C
- **Manufacturer:** MEGELIN

WARRANTY AND SUPPORT

MEGELIN offers a 365-day warranty for this product. If you encounter any issues or have questions regarding your MEGELIN 7 Colors LED Light Therapy Facial Mask, please do not hesitate to contact our customer support team.

For support, please visit the official MEGELIN store on Amazon or refer to the contact information provided with your purchase documentation.

Online Resources:

- Official MEGELIN Store: [Visit the MEGELIN Store on Amazon](#)
- Product Page: [MEGELIN 7 Colors LED Light Therapy Facial Mask](#)

Related Documents

	<p>MEGELIN LED Light Therapy Machine Spray Mist Type User Manual</p> <p>Comprehensive user manual for the MEGELIN LED Light Therapy Machine Spray Mist Type, covering safety precautions, product composition, installation, usage methods, light modes, specifications, maintenance, and warranty information.</p>
	<p>Megelin LED Light Therapy Machine User Manual</p> <p>User manual for the Megelin LED Light Therapy Machine, covering safety precautions, product composition, installation, usage methods, product specifications, maintenance, and warranty information.</p>