

## Braun PL5262

# Braun Silk-expert Pro 5 IPL Hair Removal System User Manual (Model PL5262)

[Product Overview](#)   [Safety Warnings](#)   [Setup](#)   [Operating](#)  
[Instructions](#)   [Maintenance](#)   [Specifications](#)   [Warranty & Support](#)

## 1. PRODUCT OVERVIEW

The Braun Silk-expert Pro 5 is an advanced Intense Pulsed Light (IPL) hair removal system designed for effective, semi-permanent hair reduction. It is not a laser hair removal device but utilizes IPL technology to deliver smooth skin with visible results in as little as 3 weeks, when following the recommended treatment plan. Individual results may vary.

Key features include:

- **SkinPro 2.0 (SensoAdapt™) Technology:** Automatically and continuously adapts light intensity to your skin tone for optimal balance between fast results and skin comfort.
- **Fast and Efficient:** Delivers up to 125 light flashes per minute, ensuring fewer missed spots and thorough coverage for effortless application.
- **Gentle Modes:** Features 2 sensitive modes to reduce light intensity for a softer application, ideal for sensitive areas.
- **Versatile Attachments:** Includes standard, precision, and wide heads for effective treatment across various body parts, from head to toe.
- **Complete Kit:** Comes with a storage pouch and a Venus razor for pre-treatment preparation.

# BRAUN

Silk-expert Pro 5

## 1 Jahr seidig glatte Haut\*

Auch für den Intimbereich geeignet

\*Unter Einhaltung des Behandlungsplans.  
Individuelle Ergebnisse können variieren



Image 1.1: The Braun Silk-expert Pro 5 in use, demonstrating smooth skin results. This device helps achieve smooth skin for up to one year when used according to the treatment plan.



### SkinPro 2.0 (SensoAdapt™)

Passt die Leistung jedes Lichtimpulses automatisch an den Hautton an.



### Weniger ausgelassene Bereiche

Lichtimpulse alle 0,5 Sekunden für ein müheloses Gleiten und eine gründliche Abdeckung.



Image 1.2: Visible results can be achieved in as little as three weeks with consistent use of the Braun Silk-expert Pro 5.




## 2 Sensitiv-Modi

Für eine angenehme Behandlung, auch in empfindlichen Bereichen.

### Ist IPL das Richtige für Sie?

Geeignete Hauttöne

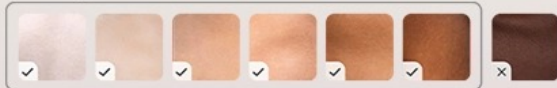

Geeignete Körperhaartöne



Image 1.3: The device automatically adapts its power output to match your specific skin tone, ensuring both safety and effectiveness.

## 2. SAFETY WARNINGS

Please read all safety instructions carefully before using the Braun Silk-expert Pro 5 IPL device. Failure to follow these warnings may result in injury or damage to the device.

- **Water Exposure:** Never use this device in the shower or bath. Do not allow the device to come into contact with water. It is designed for use on dry skin only.
- **Charging:** Never use the device while it is charging.
- **Skin Conditions:** Do not use the device on sensitive or irritated skin.
- **Application Areas:**
  - Never use the device on facial hair.
  - Never use the device on genital areas.
- **Post-Use Care:** After each use, clean the device thoroughly and store it in a dry place.

## 3. SETUP

Before your first treatment, ensure the device is properly set up and your skin is prepared.

1. **Skin Preparation:** Cleanse your skin thoroughly to remove any oils, lotions, or creams. Shave the area to be treated to ensure the light energy reaches the hair follicle effectively.
2. **Attachment Selection:** Select the appropriate attachment head for the body area you intend to treat.
  - *Standard Head:* For larger body areas like legs and arms.
  - *Precision Head:* For smaller, more sensitive areas like the bikini line and underarms.
  - *Wide Head:* For efficient treatment of large areas.
3. **Power Connection:** Connect the device to a power outlet using the provided power cord. The device is

designed for corded use only.



Image 3.1: The Braun Silk-expert Pro 5 with its included accessories, including different treatment heads and a Venus razor.

## 4. OPERATING INSTRUCTIONS

---

Follow these steps for effective and safe use of your IPL device.

1. **Power On:** Plug in the device. The indicator lights will illuminate.
2. **Skin Tone Adaptation:** Place the treatment window firmly against your skin. The SkinPro 2.0 (SensoAdapt™) sensor will automatically detect your skin tone and adjust the light intensity. The indicator lights on the side of the device will show the detected intensity level.
3. **Select Intensity (Optional):** If desired, you can select one of the 2 sensitive modes for a gentler treatment, especially in sensitive areas. Refer to the user manual for specific button functions.
4. **Apply Treatment:**
  - For larger areas, use the *gliding mode*: Hold down the activation button and slowly glide the device over your skin. The device will flash continuously.
  - For smaller or curved areas, use the *stamp mode*: Place the device on the skin, press and release the activation button, then lift and move to the next spot.

Ensure the treatment window is in full contact with the skin to activate the flashes. Keep your skin firm and guide the device in the direction of hair growth.

5. **Post-Treatment:** After completing the treatment, you may apply a soothing lotion to the treated skin if desired.
6. **Treatment Schedule:** For best results, follow the recommended treatment schedule provided in the full user manual. Typically, initial treatments are performed weekly for 4-12 weeks, followed by maintenance treatments as needed.



Image 4.1: The Braun Silk-expert Pro 5 is suitable for various body areas, including underarms, legs, and bikini line. Specific attachments are provided for different areas.

## 5. MAINTENANCE

Proper maintenance ensures the longevity and optimal performance of your device.

- **Cleaning:** After each use, ensure the device is unplugged and cool. Wipe the treatment window and the outer casing with a dry, soft cloth. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the device and its attachments in the provided storage pouch in a dry, cool place, away from direct sunlight and moisture.
- **Inspection:** Periodically inspect the device and power cord for any signs of damage. If any damage is found, discontinue use and contact customer support.

## 6. TROUBLESHOOTING

If you encounter any issues with your Braun Silk-expert Pro 5, refer to the following common solutions:

Problem	Possible Cause	Solution
Device does not turn on.	Not properly plugged in; power outlet issue.	Ensure the power cord is securely plugged into both the device and a working electrical outlet.
Device does not flash.	Treatment window not in full contact with skin; skin tone too dark; device malfunction.	Ensure the treatment window is fully pressed against the skin. Check if your skin tone is suitable for IPL treatment (refer to the skin tone chart in the full manual). If the issue persists, contact customer support.
Skin feels warm or slightly irritated after use.	Normal reaction; intensity setting too high for sensitive skin.	This is a normal sensation. Apply a soothing lotion. For future treatments, consider using a sensitive mode or reducing the intensity if available.

For more detailed troubleshooting, please consult the comprehensive user manual or contact Braun customer support.

## 7. SPECIFICATIONS

Feature	Detail
Model Name	Silk-expert Pro 5
Model Number	PL5262
Brand	Braun
Power Source	Corded Electric
Voltage	230 Volts
Item Weight	1218 Grams
Water Resistance Level	Not Water Resistant
Recommended Uses	Body Hair Removal
Manufacturer	Procter & Gamble
Country of Origin	United Kingdom

## 8. WARRANTY AND SUPPORT

Braun products are manufactured to meet high quality standards. In the unlikely event that your product requires service, please refer to the warranty information provided in the original packaging.

- **Warranty Period:** Please check your product's specific warranty card or the official Braun website for details on the warranty period applicable to your region.
- **Customer Support:** For technical assistance, troubleshooting beyond this manual, or warranty claims, please contact Braun customer service. Contact details can typically be found on the Braun official website or in your product's documentation.
- **Online Resources:** Visit the official Braun website for additional product information, FAQs, and support videos.