

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

> [LUBEX](#) /

> [LUBEX IPL Hair Removal Device with Sapphire Ice-Cooling \(Model D-1196C\) User Manual](#)

LUBEX D-1196C

LUBEX IPL Hair Removal Device with Sapphire Ice-Cooling

MODEL D-1196C USER MANUAL

1. Introduction

Thank you for choosing the LUBEX IPL Hair Removal Device. This device utilizes advanced Intense Pulsed Light (IPL) technology combined with Sapphire Ice-Cooling to provide a comfortable and effective hair removal experience at home. Please read this manual thoroughly before first use to ensure safe and optimal operation.



Image: The LUBEX IPL Hair Removal Device, including the main unit, safety goggles, two razors, and product packaging.


2. Safety Information

Your safety is paramount. Adhere to all warnings and precautions to prevent injury or damage to the device.

- **Eye Protection:** Always wear the provided safety goggles during operation to protect your eyes from the intense light pulses.
- **Skin Preparation:** Ensure skin is clean, dry, and free of lotions, gels, or deodorants before treatment.
- **Shaving:** Shave the treatment area prior to using the device to ensure the light energy targets the hair follicle effectively.
- **Skin Contact:** The device must be in full contact with the skin for the light pulse to activate.
- **Avoid Overlapping:** Do not treat the same area multiple times in one session to prevent skin irritation.
- **Not for All Skin/Hair Types:** IPL technology is most effective on darker hair colors and lighter skin tones. It is not suitable for very dark skin or very light (blonde, red, white, grey) hair.

Can I use IPL?

Skin Color

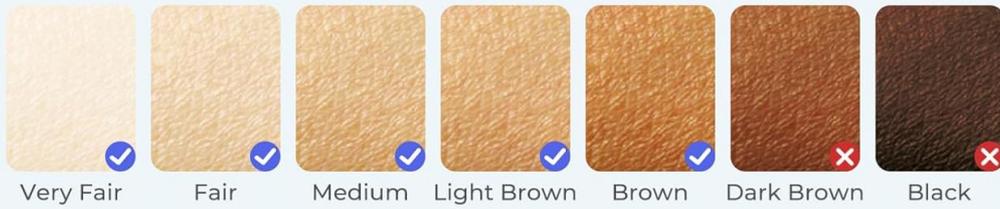

Hair Color



IPL works primarily on melanin, so higher melanin levels in dark skin increase the risk of burns. Light hair color lack enough melanin for IPL to be effective.

Image: A chart illustrating suitable skin tones (Very Fair to Brown) and hair colors (Black, Dark Black, Brown, Light Brown) for IPL treatment, and unsuitable types (Dark Brown, Black skin; Red, Blond, White Grey hair).

3. Package Contents

Verify that all items are present in the package:

- LUBEX IPL Hair Removal Device (Model D-1196C)
- Power Adapter
- Safety Goggles
- Disposable Razors (2 units)
- User Manual



Image: The LUBEX IPL device shown with its power adapter, safety goggles, and disposable razors, indicating the complete package contents.

4. Product Overview and Features

The LUBEX IPL device is designed for efficient and comfortable hair removal with several key features:


- **Sapphire Ice-Cooling Technology:** Provides instant cooling to below 41°F (5°C) upon skin contact, ensuring a painless and comfortable experience during treatment. This cooling effect is maintained for over 30 minutes of continuous use.
- **21J IPL Energy Output:** Delivers powerful light pulses to effectively target hair follicles, leading to significant hair reduction.
- **Multiple Modes:** Features dedicated modes for Face, Bikini, and Body, optimizing treatment for different areas.
- **Adjustable Intensity Levels:** 15 sensitive levels allow for customized treatment based on skin sensitivity and area.
- **Ergonomic Design:** The 4-angle diamond shape and 80-degree even light distribution ensure comfortable handling and thorough coverage without missed spots.

Nearly **Painless** Hair Removal At Home



Image: A close-up of the device's treatment window, highlighting the Sapphire Ice-Cooling technology with a constant temperature of 47°F, ensuring safe and comfortable use.

Sapphire Ice Skin Sensation

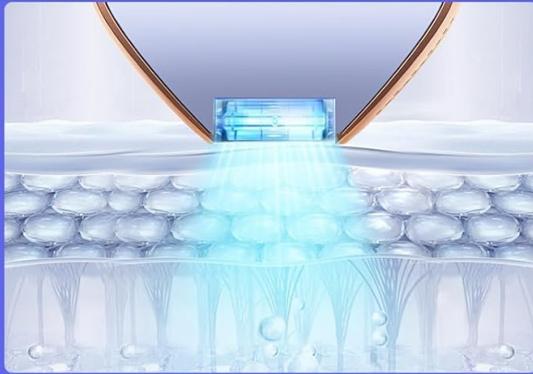

Speedy Cooling:

1ms instant cooling to 41°F



Extreme Cooling:

Lowers to 41°F



Long-lasting Cooling:

Maintains 47°F for 30 minutes of continuous use

Image: An illustration detailing the Sapphire Ice Skin Sensation, showing 1ms instant cooling to 41°F and long-lasting cooling for 30 minutes of continuous use.

Boosted 21J Energy

Efficiency improved to 2X and eliminates stubborn hair comfortably



21J IPL
Energy Output



Auto
Continuous Flash



Image: A man using the LUBEX IPL device on his armpit, illustrating the 21J IPL energy output and Auto Continuous Flash mode for efficient hair removal.

5. Setup

Follow these steps to prepare your device for use:

1. **Shave:** Shave the visible hair from the area you intend to treat.
2. **Connect Power:** Plug the power adapter into the device and then into a power outlet.
3. **Wear Goggles:** Put on the provided safety goggles to protect your eyes.



Image: A visual guide showing four steps: 1. Shave the surface hair, 2. Connect the power plug, 3. Wear safety goggles, 4. Activate the device and place it vertically at 90 degrees against the skin.

6. Operating Instructions

Once set up, operate the device as follows:

1. **Power On:** Press and hold the power button to turn on the device.
2. **Select Mode:** Use the mode button to cycle through the available modes: Face, Bikini, and Body. Each mode is optimized for specific areas.
3. **Adjust Intensity:** Select an appropriate intensity level (1-15) using the controls. Start with a lower level (1-2) for sensitive areas or initial treatments, gradually increasing as your skin tolerates it.
4. **Apply to Skin:** Place the device's treatment window firmly and vertically (90°) against the skin.
5. **Flash:**
 - **Manual Mode (Single Flash):** Press the flash button for a single pulse. Ideal for small, sensitive areas like the upper lip or bikini line.
 - **Auto Mode (Continuous Flash):** Hold the flash button or select auto mode (if available) and glide the device over larger areas like legs or arms. The device will automatically emit pulses when in full contact with the skin.
6. **Move to Next Area:** After each flash or glide, move the device to an adjacent area. Avoid treating the same spot repeatedly.
7. **Power Off:** After your session, press and hold the power button to turn off the device and unplug it.

3 Comfort Mode & Smart Touchscreen

Precision for whole body hair removal



1 Face Mode
Face, Under Lips



2 Bikini Mode
Bikini, Under Armpits



3 Body Mode
Arms, Legs, and Chest

Image: The LUBEX IPL device displaying its smart touchscreen with three comfort modes: Face Mode (for face, upper lips), Bikini Mode (for bikini, underarms), and Body Mode (for arms, legs, chest).

15 Adjustable Sensitive Levels for Painless Treatments

Quick 8 Minutes Treatments
15 Intensity Levels, Less Pain

----- Bikini Mode

Sensitive 1-2: For those with sensitive skin, for sensitive areas like bikini, upper lip, balls.



----- Face Mode

Sensitive 3-4: For those with moderate sensity or the later weeks of treatments on your more sensitive areas.



----- Body Mode

Sensitive 4-5: For those with overall sensitivity, less sensitive areas, like arms, legs, back, and shoulders.



Image: An illustration showing 15 adjustable intensity levels categorized by mode: Bikini Mode (Sensitive 1-2), Face Mode (Sensitive 3-4), and Body Mode (Sensitive 4-5), indicating recommended levels for different body parts and sensitivities.


7. Usage Schedule

For optimal results, follow this recommended treatment schedule:

- **Weeks 1-2:** 3-4 times per week.
- **Weeks 2-3:** 2-3 times per week.
- **After 4 Weeks:** 2 times per month for long-lasting results.

Visible hair reduction can be observed in as little as 3 weeks, with successful hair removal typically achieved within 6-8 weeks of consistent use.

3 Weeks Visible Result

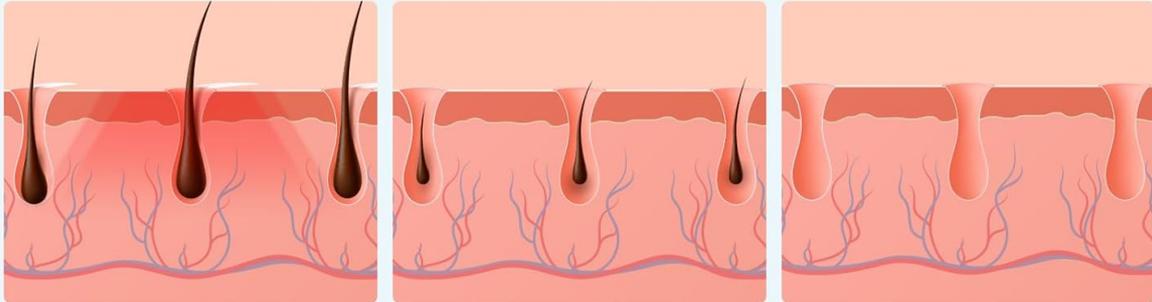

1-2 weeks
(Slow hair growth)

3 weeks
(Visibly reduced hair)

6-8 weeks
(Successful hair removal)



Image: A visual timeline showing hair growth reduction: 1-2 weeks (slow hair growth), 3 weeks (visibly reduced hair), and 6-8 weeks (successful hair removal), with corresponding images of an armpit area.

8. Maintenance

Proper care ensures the longevity and performance of your device:

- **Cleaning:** After each use, ensure the device is powered off and unplugged. Gently wipe the treatment window with a soft, dry cloth. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the device in a cool, dry place, away from direct sunlight and out of reach of children.

9. Troubleshooting

If you encounter any issues, refer to the following common solutions:

- **Device Not Turning On:** Ensure the power adapter is securely connected to both the device and a working power outlet.
- **No Flash:** The device requires full contact with the skin to emit a flash. Ensure the treatment window is pressed firmly against the skin. Also, check if the device is in manual mode and the flash button needs to be pressed.
- **Skin Discomfort:** If you experience discomfort, reduce the intensity level. Ensure your skin is clean and dry before treatment.
- **Poor Results:** Ensure you are following the recommended usage schedule and preparing your skin correctly (shaving before treatment). IPL effectiveness varies by skin and hair type; refer to the suitability chart in Section 2.

If problems persist, please contact customer support.

10. Specifications

Model Number	D-1196C
Product Dimensions	7.4 x 5.48 x 6 inches
Weight	2.25 Pounds
Manufacturer	LUBEX

11. Warranty and Support

The LUBEX IPL Hair Removal Device comes with a **24-month warranty** from the date of purchase. This warranty covers manufacturing defects and malfunctions under normal use.

For warranty claims, technical support, or any questions regarding your device, please contact LUBEX customer service. Refer to the product packaging or the official LUBEX website for current contact information.

LUBEX
Enjoy smooth skin

Image: The LUBEX brand logo, representing the manufacturer of the device.