



Manuals.plus /

› PhoneSoap /

› PhoneSoap UV Toothbrush Sanitizer, Dryer, and Case (Model BGv1G) User Manual

PhoneSoap BGv1G

PhoneSoap UV Toothbrush Sanitizer, Dryer, and Case (Model BGv1G) User Manual

Enhance your dental hygiene routine effortlessly.

INTRODUCTION

The PhoneSoap UV Toothbrush Sanitizer, Dryer, and Case is designed to prioritize your oral health by effectively eliminating bacteria from your toothbrush. This device combines UVC sanitization with a built-in drying fan, ensuring a hygienic brush for every use. It also functions as a convenient storage and travel case.

Your toothbrush encounters various bacteria from your mouth and environment. The PhoneSoap UV Toothbrush Sanitizer's UVC technology is tested and proven to eliminate 99.9%* of E. coli, contributing to a cleaner brushing experience.

*Tested and proven to eliminate 99.9% of E. coli.

SETUP

Charging the Device

Before initial use, ensure your PhoneSoap UV Toothbrush Sanitizer is fully charged. The device features a USB-C charging port. Connect a compatible USB-C cable to the port and a power source. Charging takes approximately 80 minutes for a full charge.



Image: The USB-C charging port and indicator lights on the PhoneSoap UV Toothbrush Sanitizer.

Mounting Options

The PhoneSoap UV Toothbrush Sanitizer includes a magnetic dock for flexible placement. You can attach it to a mirror, wall, or any convenient spot for easy access and storage.



Image: The PhoneSoap UV Toothbrush Sanitizer mounted on a mirror, holding two toothbrushes.

Toothbrush Compatibility

The sanitizing bay is designed to accommodate various toothbrush heads. It fits brush heads up to 12mm in width.

Fits all these toothbrushes and more...



Image: An illustration showing different toothbrush head sizes that fit the sanitizer.

OPERATING INSTRUCTIONS

Sanitization and Drying Cycle

1. Simply place your toothbrush neck inside the rubber holder to secure the brush head within the unit.
2. Ensure your toothbrush bristles are facing towards the LED light for effective sanitization.
3. Close the lid. The device will automatically activate the sanitization and drying process.




Image: A toothbrush being inserted into the open PhoneSoap UV Toothbrush Sanitizer.

Quieter Operation

For quieter operation, you can press and hold the button on the outside of the device. This action will turn off the drying fan while maintaining the sanitization process.

Key Features

The device utilizes UVC light to sanitize and includes a built-in fan for drying.

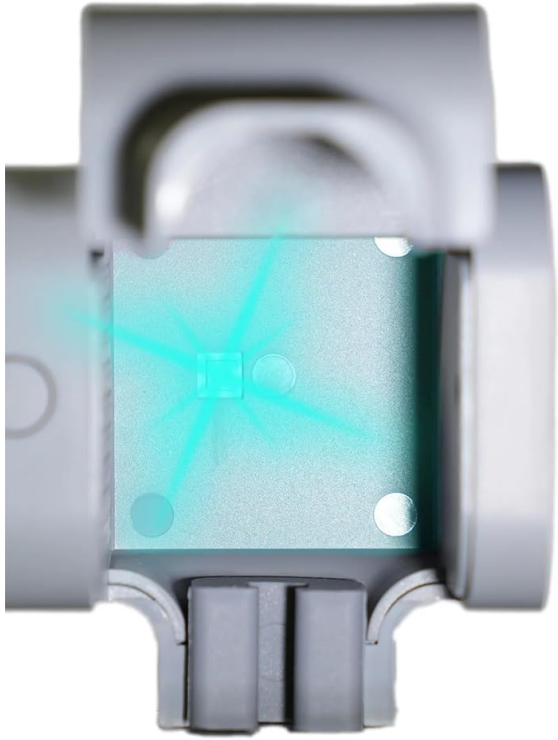

99.9% Sanitizing UV-C

Image: A close-up view of the UV-C sanitization area inside the device.



Built-in Drying Fan

Image: A close-up view of the built-in drying fan within the device.

Battery Life

A full charge provides:

- 28 fan + sanitizing cycles
- 34 fan only cycles
- 40 sanitizing only cycles

Travel Functionality

The device doubles as a travel case, protecting your toothbrush from germs when on the go or stored in a drawer.



Image: The PhoneSoap UV Toothbrush Sanitizer packed in a suitcase, highlighting its travel-ready design.

Product Video

Your browser does not support the video tag.

Video: An official PhoneSoap video demonstrating the features of the Toothbrush UV Sanitizer, including its UVC sanitization and drying capabilities.

MAINTENANCE

To ensure optimal performance and hygiene, regularly clean your PhoneSoap UV Toothbrush Sanitizer. Use a soft, damp cloth to wipe down the exterior. For the interior, ensure the device is powered off and gently clean the sanitizing bay with a dry or slightly damp cloth. Avoid using harsh chemicals or abrasive cleaners, as these can damage the device or its UV-C components.

TROUBLESHOOTING

Device Not Activating

- Ensure the lid is fully closed. The sanitization cycle begins automatically upon lid closure.

- Check the battery level. If the battery is low, charge the device using the provided USB-C port.

Fan Not Operating

- If the fan is not running but the sanitization light is on, you may have activated the quieter operation mode. Press and hold the button on the outside to re-engage the fan.

Toothbrush Not Fitting Securely

- Ensure the toothbrush neck is properly seated within the rubber holder.
- Verify that your toothbrush head does not exceed the maximum compatible width of 12mm.

SPECIFICATIONS

Model Number	BGv1G
UPC	724049646219
Product Dimensions	1.5 x 1.5 x 1.5 inches
Unit Weight	2.4 ounces (7.05 ounces with packaging)
Sanitizing Bay Dimensions	1.5in x 1.5in x 1.5in
Battery Capacity	3.7V 500mAh
Input	DC 5V 1A
Manufacturer	PhoneSoap

WARRANTY AND SUPPORT

For warranty information and customer support, please refer to the official PhoneSoap website or contact their customer service directly. You can also visit the PhoneSoap Store on Amazon for additional resources and contact details.