



Manuals.plus /

› DentJet /

› DentJet FC5240 Cordless Oral Irrigator User Manual

DentJet FC5240

DentJet FC5240 Cordless Oral Irrigator User Manual

Model: FC5240 | Brand: DentJet

INTRODUCTION

Thank you for choosing the DentJet FC5240 Cordless Oral Irrigator. This device is designed to provide effective and convenient oral hygiene, helping to remove plaque and food debris from between teeth and below the gumline. Please read this manual thoroughly before first use to ensure proper operation and maintenance.

IMPORTANT SAFETY INFORMATION

- Do not immerse the charging base in water.
- Keep out of reach of children.
- Do not use if the device or charging cable is damaged.
- Use only with water or dentist-recommended mouthwash.
- Avoid directing the water jet directly into the eyes, nose, or ears.
- If you have sensitive gums or have undergone recent oral surgery, consult your dentist before use.

PACKAGE CONTENTS

Verify that all items are present in your package:

- DentJet FC5240 Oral Irrigator Unit
- Standard Jet Tips (various types for multifunctions)
- USB Charging Cable
- User Manual
- Travel Pouch (if included in your specific kit)



Image: Contents of the DentJet FC5240 package.

PRODUCT OVERVIEW

The DentJet FC5240 features a compact, cordless design with a 300ml water reservoir and multiple pressure modes for personalized cleaning.



Image: Front view of the DentJet FC5240 Oral Irrigator.

Key Components:

- **Nozzle Release Button:** For detaching jet tips.
- **Water Inlet:** For filling the reservoir.
- **Power Button:** To turn the device on/off.
- **Mode Button:** To cycle through pressure modes.
- **Indicator Lights:** Display selected mode and battery status.
- **Water Reservoir:** 300ml capacity.
- **Charging Port:** Located at the back, covered by a waterproof plug.

SETUP

1. **Initial Charge:** Before first use, fully charge the device for approximately 4 hours. Ensure the charging port cover is securely closed after charging.
2. **Attach Nozzle:** Insert the desired jet tip into the opening on the top of the unit until it clicks into place.
3. **Fill Reservoir:** Twist the bottom of the unit counter-clockwise to detach the reservoir. Fill it with lukewarm water or mouthwash. Reattach the reservoir by twisting clockwise until secure.

OPERATING INSTRUCTIONS

1. **Select Mode:** Press the Mode button to cycle through the available pressure settings (Normal, Medium, Soft, Pulse).
2. **Position Nozzle:** Lean over a sink and place the jet tip into your mouth, aiming it at your gumline at a 90-degree angle.
3. **Start Flossing:** Press the Power button to start the water flow. Keep your mouth slightly open to allow water to drain into the sink.
4. **Guide the Jet:** Move the jet tip along your gumline, pausing briefly between teeth. Clean both the inside and outside surfaces of your teeth.
5. **Finish:** Press the Power button again to turn off the device. Empty any remaining water from the reservoir.

MODES AND NOZZLES

Pressure Modes:

- **Normal Mode:** Powerful cleaning for general use.
- **Medium Mode:** Daily cleaning for regular oral care.
- **Soft Mode:** Gentle cleaning for beginners, children, or sensitive teeth.
- **Pulse Mode:** Massages gums and helps improve gum health.



Image: Four distinct pressure modes for varied cleaning needs.

Multifunction Nozzles:

- **Standard Jet Tip:** For general cleaning, removing food and dirt between teeth.
- **Orthodontic Tip:** Designed for cleaning around braces, bridges, and other dental work.
- **Periodontal Tip:** For effectively removing dental plaque and preventing dental calculus.
- **Tongue Cleaner Tip:** For cleaning dental implants, dentures, dental bridges, and crown cleaning.

4 Nozzles for Multifunctions

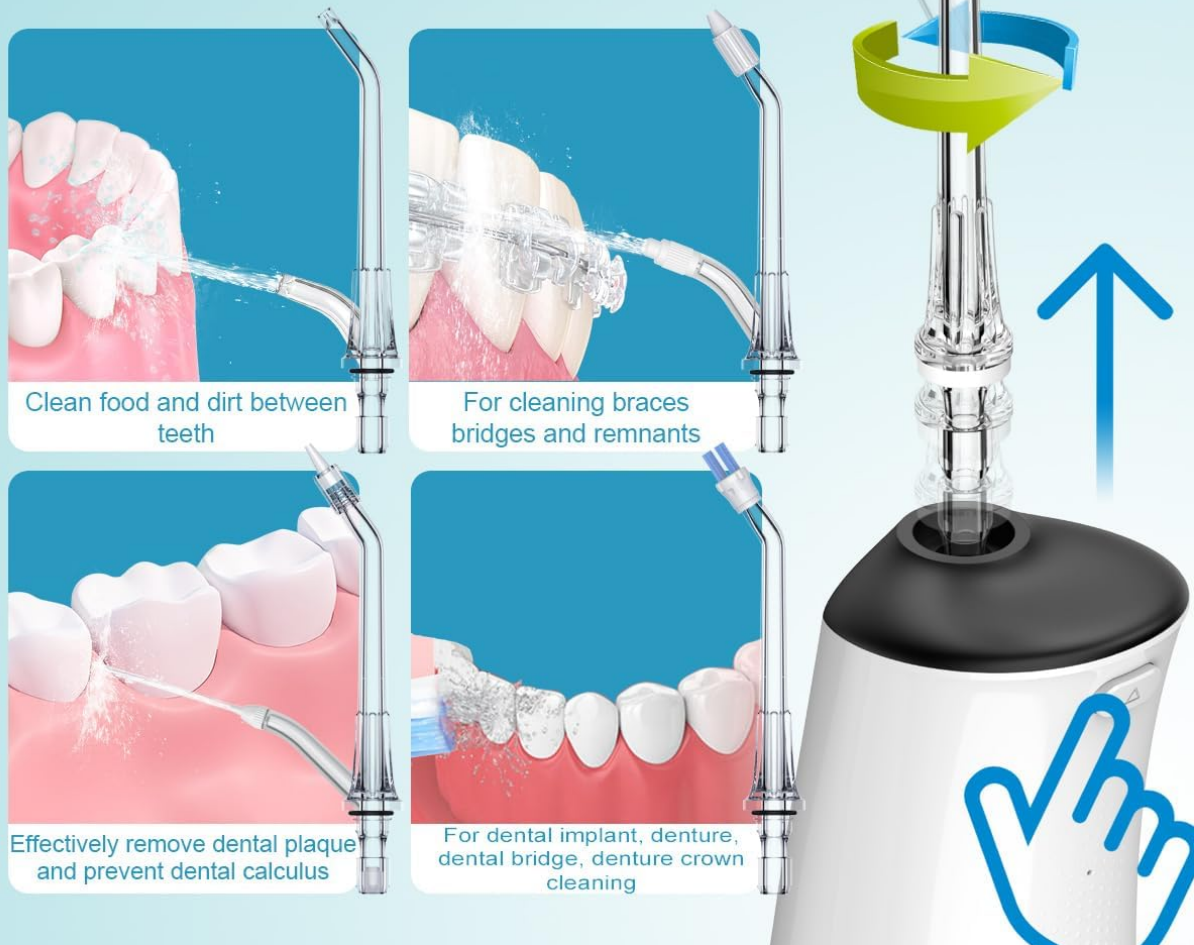


Image: Four specialized nozzles for comprehensive oral care.

CHARGING

The DentJet FC5240 is equipped with a 2000 mAh Li-Battery. A full charge takes approximately 4 hours and provides up to 15 days of use (based on 2 minutes of use, twice a day in Normal Mode).

- Connect the USB charging cable to the charging port on the back of the device.
- Plug the other end into a compatible USB power adapter (not included).
- The indicator light will show charging status.
- Ensure the charging port cover is tightly closed before using the device in water.



Image: Charging time and battery life details.

MAINTENANCE AND CLEANING

Regular cleaning ensures optimal performance and extends the life of your oral irrigator.

- **After Each Use:** Empty any remaining water from the reservoir. Run the device briefly with clean water to flush out the internal components.
- **Cleaning the Reservoir:** The reservoir can be detached and washed with mild soap and water. Rinse thoroughly.
- **Cleaning the Unit:** Wipe the exterior of the unit with a soft, damp cloth. Do not use abrasive cleaners.
- **IPX7 Waterproof:** The device is IPX7 waterproof, meaning it can be safely washed under running water or used in the shower. Always ensure the charging port cover is sealed.
- **Nozzle Cleaning:** Remove the jet tip and rinse it under running water. Replace jet tips every 3-6 months or when they show signs of wear.

IPX7 100% Waterproof

Perfect Work in Bathing or Showering

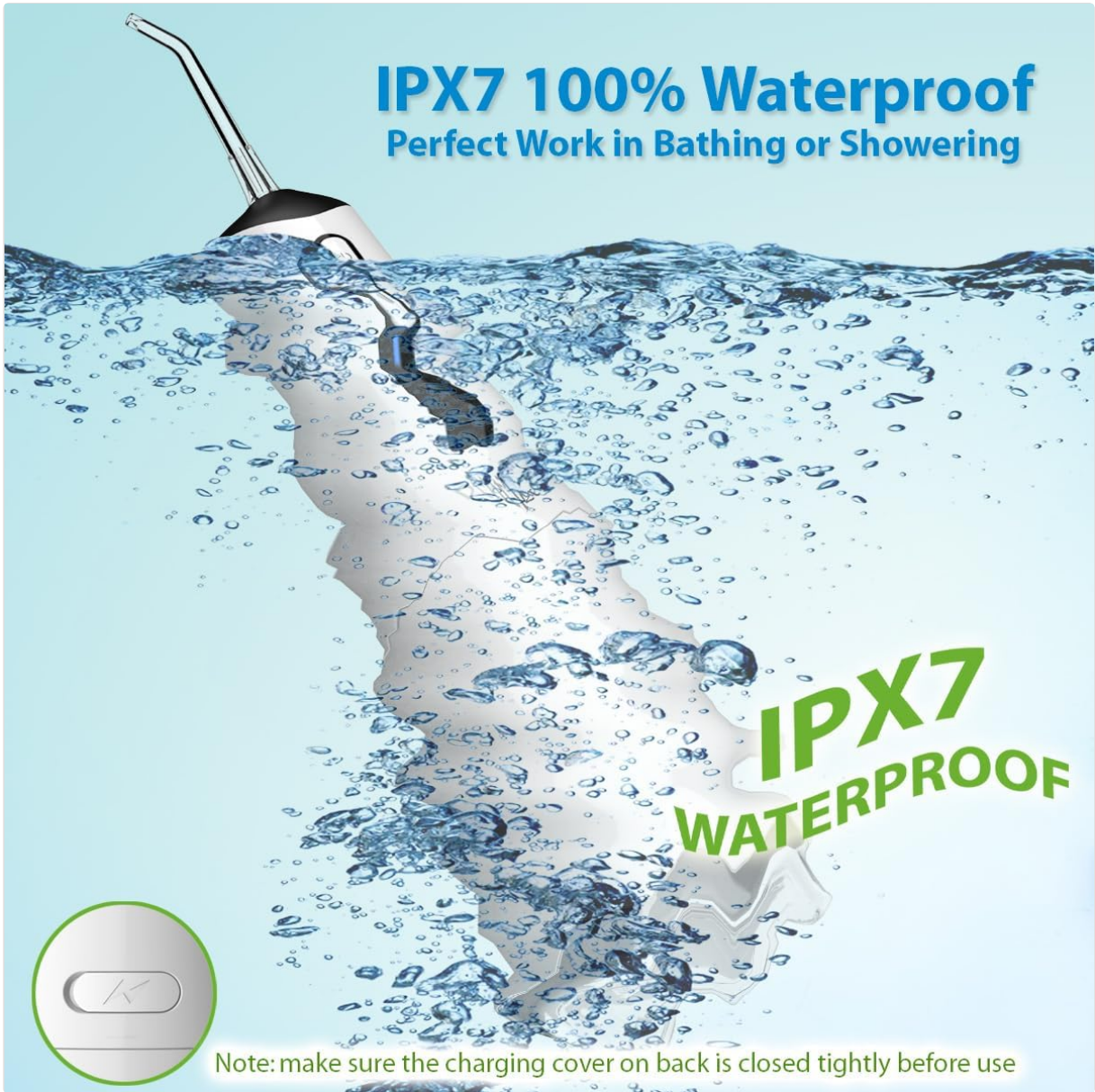


Image: The device is IPX7 waterproof for easy cleaning and safe use in wet environments.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery or no charge.	Charge the device for 4 hours.
Weak water pressure.	Low battery, clogged nozzle, or insufficient water in reservoir.	Recharge the device. Clean or replace the nozzle. Refill the reservoir.
Water leakage.	Reservoir not properly attached or charging port cover not sealed.	Ensure the reservoir is twisted on securely. Check that the charging port cover is tightly closed.
Device stops working unexpectedly.	Battery depleted or internal malfunction.	Recharge the device. If the problem persists, contact customer support.

SPECIFICATIONS

Model Number: FC5240

Brand: DentJet

Water Reservoir Volume: 300 ml

Power Source: Rechargeable Li-Battery (2000 mAh)

Charging Time: Approximately 4 hours

Operating Time: Up to 15 days (based on typical use)

Waterproof Rating: IPX7

Dimensions (L x W x H): 7.9 x 7 x 22.4 cm

Weight: 410 grams

Special Features: Cordless, Rechargeable, Waterproof, Multiple Pressure Modes, Multifunction Nozzles

WARRANTY AND SUPPORT

DentJet products are manufactured to high-quality standards. For warranty information or technical support, please refer to the warranty card included in your package or contact DentJet customer service through the retailer where the product was purchased.

For further assistance, please visit our official website or contact our support team.